

RENPHO

Smart Health Simplified

18888

18888

RENPHO

RENPHO is committed to our environment. To reduce paper usage, please scan the QR Code or enter the following URL to download the full version of the electronic User Manual:



<https://renpho.com/pages/user-manual-body-composition-scale>

01

04

RENPHO

Smart Health Simplified

US Importer: JOICOM CORPORATION

14129 The Merge Street, Building 3 Unit A, Eastvale, CA 92880

Manufacturer: Shenzhen Ruiyi Business Technology Co., Ltd.

Qianhai Complex A201, Qianwan Road 1, Qianhai Shenzhen-Hong Kong Cooperation Zone, Shenzhen, 518000 P.R.China

Made in China

NA2

02

05

IMPORTANT SAFETY INSTRUCTIONS

READ ALL INSTRUCTIONS CAREFULLY BEFORE USE.

1. The scale is for household use only and not for medical or commercial purposes.

2. Children under 8, the elderly, or individuals who are pregnant, have osteoporosis, a pacemaker, artificial joints, portable ECG or other metal implants should use the scale for weighing purposes only and make sure to wear shoes to prevent any electric current when stepping onto the scale.

3. All body composition data provided is for reference and should not be used as a substitute for advice from a licensed healthcare professional.

4. Do not overload the scale. (Maximum Capacity: 400lb / 180kg)

5. Inspect the scale for any damage before each use. Do not use a damaged scale to avoid inaccuracy or injury.

6. In case of any malfunction, check the battery power first. If necessary, replace the batteries.

7. Keep the batteries out of reach of children.

8. To ensure accuracy, use the scale on a hard, flat and dry surface.

9. To avoid interference or danger, do not use the scale in an environment near equipment with strong electromagnetic fields.

10. To prevent injury, do not stand on the edges of the scale, jump on or stomp it.

11. The scale is a delicate and precise measuring device. Handle with care to avoid breakage.

12. The scale is not waterproof. Keep it dry and never submerge it in water. To clean the surface, use a damp cloth or glass cleaner. Do not use soap or chemical cleaning agents.

13. Do not expose the scale to direct sunlight or use it in high temperatures or extremely humid areas. Store the scale in a cool, dry place.

14. Do not attempt to disassemble the scale for any reason, as this may cause damage or void the warranty.

03

06

ABOUT YOUR SCALE

LED Display

Electrodes

Weighing Platform

Foot Pad

Battery Compartment

Icons

CAL

Calibrating

Calibrate the scale before each use.

Err

Overload

Maximum Capacity: 180kg/400lb

Lo

Low Battery

Replace the batteries timely.

0000

Measuring Body Fat

Bluetooth

Bluetooth is disconnected.

Bluetooth

Bluetooth is connected.

BMI

BMI

Body Fat

Body Fat

User Identification

User Identification

Bare Feet

Bare Feet

07

06

Method 1

a. Open the Renpho Health App and step onto the scale's electrodes with bare, dry feet.

b. Tap "Connect" to pair when a window pops up on the app.

08

07

CALIBRATING YOUR SCALE

Calibrate the scale before each use to ensure accuracy.

a. Place the scale on a flat, solid and even surface. Avoid any carpeting or uneven flooring.

b. Step on the scale to activate it and then step off.

c. The scale will show "CAL" and then "0.0" or "0.00" to indicate calibrating is complete.

09

08

USING YOUR SCALE

Before taking your first measurement, it is important to have a user ID set up in the Renpho Health App. This ensures your measurements are accurately recorded and associated with your profile.

a. Ensure the scale is placed on a flat, solid and even surface.

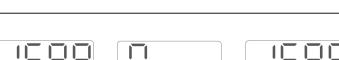
b. After calibrating the scale, step on the scale's electrodes with bare and dry feet and wait for the measurement(s) to be shown on the scale's display.

c. After measuring, the scale will transmit the data to the app via Bluetooth. Open the app to sync and view your detailed measurements, track your progress, and access additional features and insights.

10

09

1. Display Issues

Display	Description/Cause
	Scale has been connected to the app to use, and this measurement was taken with bare feet.
	Weight, BMI, Body Fat
	Scale has been connected to the app to use and this measurement was taken with socks or shoes.
	Weight, BMI
	Scale has not been connected to the app to use, and this measurement was taken with socks or shoes.
	Weight
	Scale has not been connected to the app, and this measurement was taken with bare feet.

11

10

Method 2

a. Tap "+" on the top right corner of the "Device" page to pair manually.

b. Tap "Body Fat Scale" and follow the in-the-app instruction.

c. Tap "+".

12

11

WARRANTY POLICY

Your RENPHO product purchase is covered by a one-year limited manufacturer warranty from the date of delivery.

For warranty terms and conditions, please visit:

<https://renpho.com/pages/warranty-terms-and-conditions>

Note: Product registration is not required for the warranty.

If you choose not to register your product, it will not diminish the product warranty.

CUSTOMER SERVICE

Please feel free to contact us if you have any questions or concerns. RENPHO Customer Service Team guarantees a quick response and hassle-free solutions to any issue you may have within business hours.

TEL: +1(844) 417 0149 (US ONLY)

Monday-Friday 9:00AM-4:30PM

1800 MY RENPHO (1800 69 7367) (AU ONLY)

Monday-Friday 9:00AM-6:00PM (AEST)

Email: support@renpho.com (US&CA)

support-au@renpho.com (AU)

* For defective products or the return of items, please contact us with your order number within the specified warranty period. DO NOT dispose of any product parts as they may be required for inspection/repair.

13

12

FCC Regulatory Compliance Policy

This device complies with part 15 of the FCC Rules. Operation is subject to the following two conditions:
(1) This device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.
Warning: Changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.
Note: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications.
However, there is no guarantee that interference will not occur in a particular installation.
If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:
- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.
RF Exposure Compliance
This equipment complies with FCC radiation exposure limits set forth for an uncontrolled environment. This transmitter must not be co-located or operating in conjunction with any other antenna or transmitter.
ISED Regulatory compliance
This device contains licence-exempt transmitters that comply with Innovation, Science and Economic Development Canada's licence-exempt RSS(s).
Operation is subject to the following two conditions:
(1) This device may not cause interference.
(2) This device must accept any interference, including interference that may cause undesired operation of the device.
Le présent appareil est conforme aux CNR d'Industrie Canada applicables aux appareils radio exempts de licence. L'exploitation est autorisée aux deux conditions suivantes:
(1) L'appareil ne doit pas produire de brouillage, et (2) l'utilisateur de l'appareil doit accepter tout brouillage radioélectrique subi, même si le brouillage est susceptible d'en compromettre le fonctionnement.
This equipment complies with IC RSS-102 radiation exposure limits set forth for an uncontrolled environment.
Cet équipement est conforme aux limites d'exposition aux radiations IC CNR-102 établies pour un environnement non contrôlé.

尺寸: 96x130mm

材质: 108g铜版纸 双面印刷 四色印刷 十字折

虚线边框为裁切框不印刷