



Product Brief Introduction



Thank you for choosing our company's product Acoustic and Dynamic Repair Instrument for Human Body.

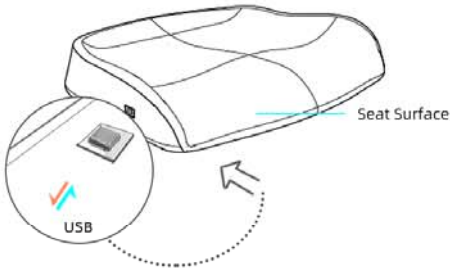

Application Method

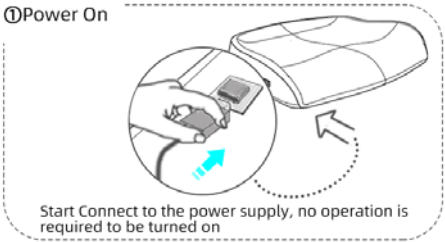


Note: The charging interface should be charged with an adapter charger. The power supply parameters should meet: output voltage 5V, output current 2A.

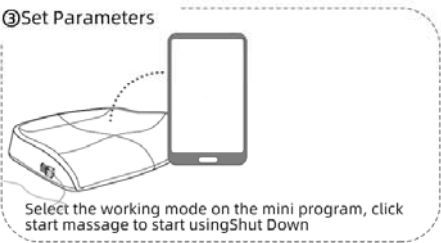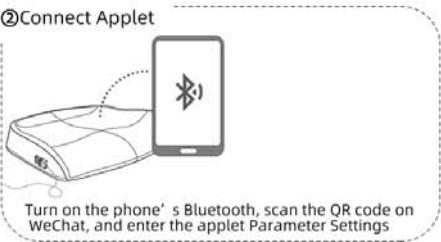

Note: If the applet connection is not successful, please check:

1. Whether the cushion is plugged in;
2. Whether the Bluetooth of the mobile phone is turned on;
3. If your device use android OS, please enable the Location function .

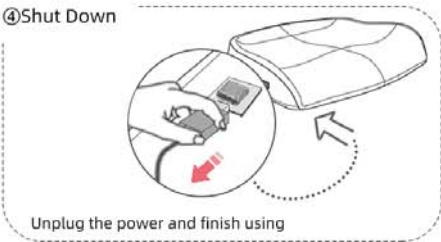

Configuration List


Item	Massage Cushion	Power Cord	User's Manual
Quantity	1	1	1

Precautions

- When not using this product, please be sure to unplug the power plug
- It is recommended to use it for no more than 40 minutes every day. According to your own situation, choose a comfortable massage intensity, and the time can be appropriately extended.
- It is forbidden to approach fire sources, corrosive liquids and various unfavorable environments
- Please keep the ports, plugs and other parts dry, and do not operate the power supply with wet hands.
- Please store in a cool and dry environment, avoid moisture and direct sunlight.

Banned people

- Children, the elderly or the disabled should use this product under supervision and guidance. Ensure that children come into contact with this product under supervision and do not use it as a toy. Do not place the product in a place where children can reach it.
- The surface of the product has heating function, people who are not sensitive to heat should use it with caution.
- Patients with skin sensation impairment or skin abnormalities are not allowed to use.
- If you have hip sensory impairment or hip disease, please consult your doctor before use.
- Prohibited by pregnant women.

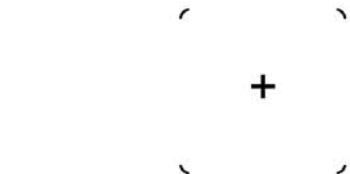

Service Hotline : 400-607-6107
Corporate Name : Taian Kangyu Medical Equipment Co. Ltd
Corporate Add : Kangyu Health Industrial Park, Daiyue District,
Taian, P.R.China

Product Warranty Description

From the date of purchase by the customer, our company is responsible for providing after-sales service. The specific contents are as follows:

- 1.Thank you for your kind purchase of our products.
- 2.We will confirm your Warranty Certificate to provide warranty service according to the warranty terms.
- 3.In order to ensure that you can get technical service support in time, please contact us in time after encountering any malfunctions.
- 4.Our products have lifetime maintenance with this card, and a free warranty for one year from the date of purchase.
5. The vulnerable parts of this instrument are not included in the free warranty, but we will provide customers with long-term preferential prices.

Non-warranty Regulations:

- 1.Unauthorized repairs, misuse, collision, negligence, abuse, fluid intake,man-caused accidents, incorrect use of non-product accessories, or tearing or altering labels or anti-counterfeiting marks;
- 2.It has exceeded the validity period of the three guarantees;
- 3.Damage caused by force majeure;
- 4.This product and its accessories are malfunctioning due to human reasons.

Attention :

In order to facilitate the after-sales service of the product, please keep the inner and outer packaging of the product intact. All costs incurred due to packaging problems shall be borne by the user!

Warranty

Name : _____ Tel : _____

E-mail : _____

Add : _____

Item : _____

Model : _____

Serial No. : _____

Warranty Record

Delivery Date	Cause of Failure	Delivery Acceptance Date	Service Engineer

Certificate Of Quality

Checker:

Inspect Date:

HKANGYU

Manual

Acoustic and Dynamic
Repair Instrument for Human Body

Please read this instruction manual carefully before
using the product and keep it properly!

This device complies with part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) this device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

Any changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

NOTE: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation.

If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

The device has been evaluated to meet general RF exposure requirement. The device can be used in portable exposure condition without restriction

FCC ID: 2A324HFK-4002