




## FY-A300A

### BEAUTY BLACKHEAD INSTRUMENT


Operation Manual

Please read this manual carefully  
Before using and keep this in a safe place for future reference


## Function




**Microcrystalline skin-smooth**  
The probe with microcrystalline particles has a large surface area with high absorption. without hurting the skin, The probe can massage and exfoliate.



**Comedo suction, skin Tightening**  
The large circular probe carries a strong suction. It is best used for oily spots and blackheads. Its unique mechanism uses a vacuum pressure to massage the skin by micro-alift which deeply effects the underlying skin.The device Stimulates micro-circulation operation promotes blood and lymphatic circulation that allows the skin to rejuvenate again.



**Cleansing Pores**  
The small circular probes suction is weaker. Place the probe on the blackhead acne to remove. If the skin starts to hurt or becomes swollen, it is recommended to stop use of the device. If the blackhead is deep in the pores, use the larger size probe The small size probe is more suitable for sensitive skin.



**Smooth Fine Lines**  
The special oval probe is designed to have a micro pulling massage for the canthus, ala nasi, the corner of the mouth and other areas where there are fine lines, areas difficult to care for, or sensitive skin. It stimulates the micro circulation preventing the appearance of wrinkles.

## Tips

❶ The tightening function shapes and pulls.Injury and scarring is possible. It is not recommended to keep the probe in the same place for too long. Generally 3-5 seconds is appropriate.

❷ The probe with microcrystalline particles need to be used according to your skin: once a week for dry skin, twice for oily skin and, no more than 5 minutes every time;

❸ Please don't press the probe firmly on your skin;the machine has strong suction.

❹ After skin care. the skin may feel a little tingly and flush because of the blood circulation. Depending on your skin type, the redness phenomenon will fade away between 0.5-24 hours.

❺ Please pay more attention to moisturizing and sun care after use.

❻ Please do not use scrub cream or exfoliator before or after skin care.

## Connect To App And Usage Method

1. Scan the QR code below to install app (Android Apple general)  
(As shown in Figure 1)



Figure 1

2. The installation is successful  
(As shown in Figure 2)



Figure 2



Figure 3

3. Turn on the key of the blackhead suction instrument, keep it on, and open the app to show that it is not connected.(As shown in Figure 3 and Figure 4)

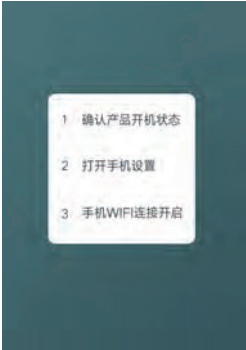


Figure 4

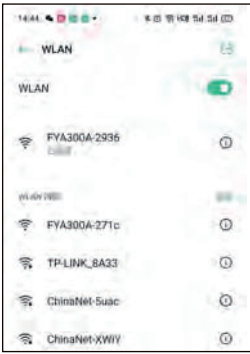


Figure 5

4. Turn on the mobile phone and set up WLAN to search for "fya300a -?".(As shown in Figure 5)

5.Connect, click the start button on app, the camera light is on, normal use.(As shown in Figure 6and Figure 7)

6. Exit app and the camera light goes out.




Figure 6


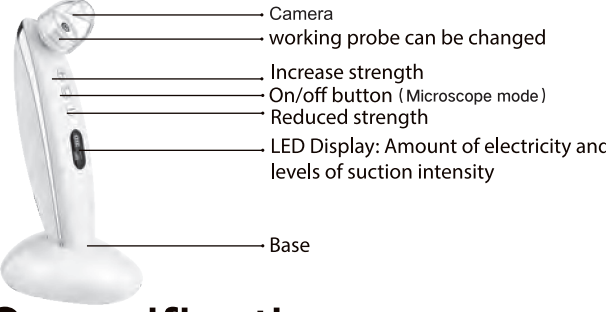


Figure7

## Product Feature



- Camera
- working probe can be changed
- Increase strength
- On/off button (Microscope mode)
- Reduced strength
- LED Display: Amount of electricity and levels of suction intensity
- Base

## Speccification

Name: Comedo Suction Beauty Device  
Model: FY-A300A  
Input voltage: 5V-0.5A  
weight: 410g  
Power: 5W  
Size: 211.5(L)X135.5(W)x52(H)mm

## Using area

**Microcrystalline Probe:**  
Any area on face except around the eyes. It is suitable to use in areas such as elbows, knees and ankles;

**Big size circular from probe:**  
Any area on face except around the eyes.

**Small size circular from probe:**  
Can be used on sensitive areas such as; around the eyes or any other areas of the face.

**Oval probe:**  
Can be used on sensitive areas around the eyes or any other areas of the face.

## Attention

❶Please dry your hands before using the device. The water may get into the device from its air exit hole, which will damage the device;

❷Before using it, please make sure the sponge has been put into the cup. The sponge must be dry, or the liquid will damage the device;

❸ Please don' t wash the device;

❹ This product is not for people who are physically or mentally disabled lack of experience and knowledge(including children) without parental/guardian supervision;

❺ This product is not suitable for body parts that have had cosmetic surgery or rhinoplasty, Pregnant woman and heart problem patients should not use this product.

## Introduction

Microcrystalline Comedo Suction Device is a beauty device used for skin-smoothing acne removing,skin firming, pore cleansing and comedo suction.the microcrystalline probe with 1 million natural mineral particles that work with vacuum suction for removing surface skin blemishes and the dirt from pores,making your skin soft and smooth immediately.

This technology is not invasive and a skin irritant. The probe has different degrees of roughness to control the degree of exfoliation. Each of the four probes have different functions: exfoliation, skin firming, pore cleansing and comedo suction.

## Working Principle

Vacuum Vapour-pressure Type Slimming Face Technology

❶ Massages your skin gently, which can accelerate the circulation of blood promote metabolism, make skin more elastic and make skin come back white and smooth;

❷ Strengthen skin's osmosis by penetrating deeply into the skin, which can shrink pores, moisturize, nourish and firm the skin;

❸ Accelerate refreshing of epidermis cells, strengthen skin's immunity and resistance of ultraviolet radiation, reduces epidermal pigmentation and inhibits melanin production for a flawless, mono-chromatic appearance;

❹ Stimulates the growing of textile fiber metrocyte, produces the massive high quality tension collagenous fibre continuously and helps to tighten the skin.

❺ Improves the micro-circulation of shin, promotes metabolism of melanin,which can control new forming melanin.

## Microcrystalline Probe

The natural mineral particles on the probe can gently exfoliate, deep cleanse skin, and activate the face in an instant. It can eliminate blackheads age spots and suction the dirt from the skin

This allows accelerating circulation of blood to help cells refresh, making skin smooth and rejuvenated.

## Working Principle

**Skin-smooth and tighten function**

Clean your face and dry, use the skin care or steamer to open pores.

❶ Change to the microcrystalline probe or large probe;

❷ Start from your chin and move the device from bottom to top;

❸ Use on the cheek and move it from side to side;

❹ Use on the forehead area start from the middle and work your way to the sides;

❺ Use on T-area and move from bottom to top.

**Clean the pore**

❶ Change to large or small size circular probe(depending on your skin type);



❷ Place the round hole on the blackhead acne to remove, The levels can be adjusted according to your need;

**Method:** When working on a blackhead, please move slowly in the same direction. Proceed to repeat it several times in order to avoid the skin redness marks. Please don't stay too long in one position.

**Smooth fine lines**

Change to the oval probe, use on canthus, alanasi, the corner of the mouth and other areas where it's easy to have fine lines. The levels can be changed according to your need;

**Method:** Suction on the areas of smooth fine line, gently pull away quickly over and over.It stimulates the micro circulation, preventing the appearance of wrinkles.



(Attention: Please don't leave the suction on the same area for long period of time, it will make the skin swell and may leave red marks.If the skin swells, stop of use immediately.)

FCC ID:2A3H4-FY-A300A

This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) this device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

NOTE: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

-- Reorient or relocate the receiving antenna.

-- Increase the separation between the equipment and receiver.

-- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.

-- Consult the dealer or an experienced radio/TV technician for help.

**RF Exposure Information**

The device has been evaluated to meet general RF exposure requirement.

The device can be used in portable exposure condition without restriction.