


5514说明书

日期：2021.01.26

材质：200克双铜纸

装订方式：骑马钉

尺寸：130.18x130.18

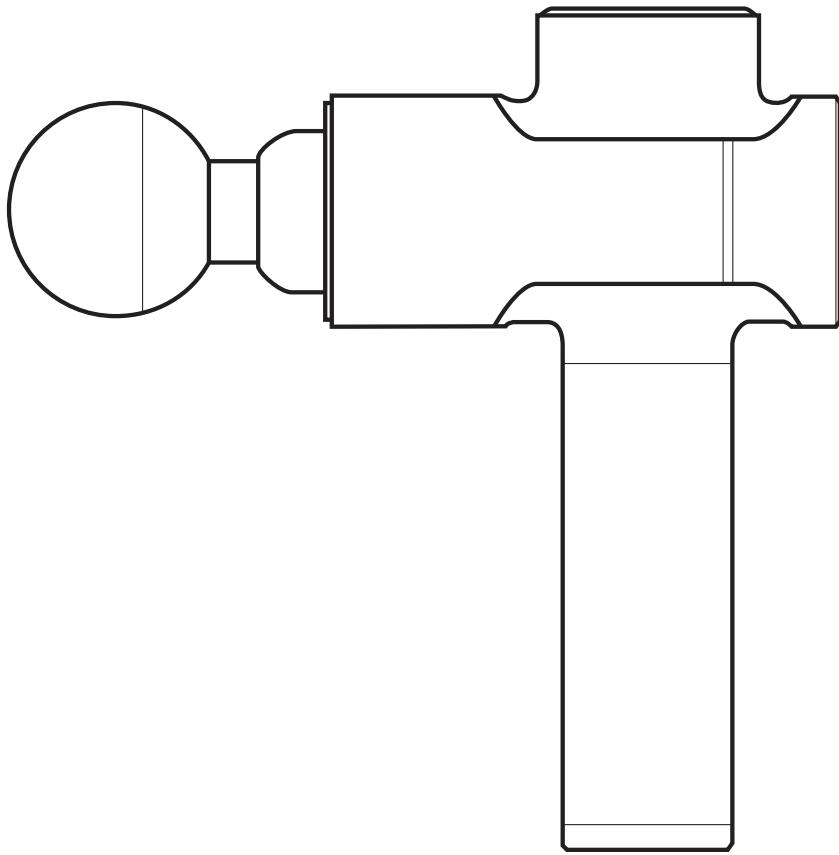

MAX

Operating Instructions | Mode d'emploi | Instrucciones de operación |
Bedienungsanleitung | Instruções de operação | Istruzioni per l'uso

EN

PLEASE READ PRIOR TO USE

The Recovapro Max is a state of the art handheld percussion massage device that helps relieve muscle soreness and stiffness, accelerates warmup and recovery before and after physical activity, optimises mobility and promotes blood circulation.

CHARGING

Before first use, make sure to charge your Recovapro Max for six (6) hours.

To charge, plug the charger to the wall outlet using the correct adapter, then connect the charger to the charging port on the base of the device. The battery indicator will flash to signify the device is charging. When all four battery lights are solid, the battery is 85%-100% charged. Alternatively, plug the charger into the docking station and place the Recovapro Max on the docking station.

OPERATING

Press and hold the Power (R) Button for 2 seconds to Switch your Recovapro Max ON. Once your device is powered ON a White light will illuminate on the first speed setting, plus (up to) four blue lights for battery level on the side of the device. Your Recovapro Max is now ready for use.

Selecting speeds (1-5) - By pressing the Power (R) button, the second speed is selected which is indicated by the 2nd White light above the Power (R) button. Each press of the R button will increase the speed of the device which is indicated by the light above the R button.

To switch off the device, Press and hold the Power (R) button for 2 seconds at any speed level, a full ring Blue Led will appear and after 5 seconds, your device will switch off. The device will stop running after 10 minutes of continuous use. Restart again if needed.

The Recovapro Max is equipped with Intelligent pressure sensor light: The more pressure is applied, the brighter the light becomes. The Recovapro max can deliver up to 60 lbs of pressure before stalling which is more than enough to tackle stubborn areas.

To change head attachments insert the desired head directly into the shaft / opening and then rotate it clockwise (about 90°), until the vertical line on the attachment is aligned with that on the device

CLEANING AND STORAGE

Make sure the device is OFF and the charger is not attached. Use a damp, clean, cloth and gently wipe your Recovapro Max. Store in a clean, cool place when not in use.


SPECIFICATIONS

1. Battery level indicator
2. Interchangeable head attachments (5)
3. Aluminum Shell
4. Handle with precision engineered grip
5. Charging port
6. Speed level indicators (5 LED lights)
7. Power (R) button
8. Intelligent Pressure Sensor Light
9. Vertical line on attachments and device

AC charger: 110-240V 50/60Hz - .5A - 26VDC- 1A

Frequency (Hz): level 1 - 1400/23, level 2 - 1800/30, level 3 - 2200/37, level 4 - 2600/43, level 5 - 3200/53

Weight: 1.2 kg / 2.6 lbs

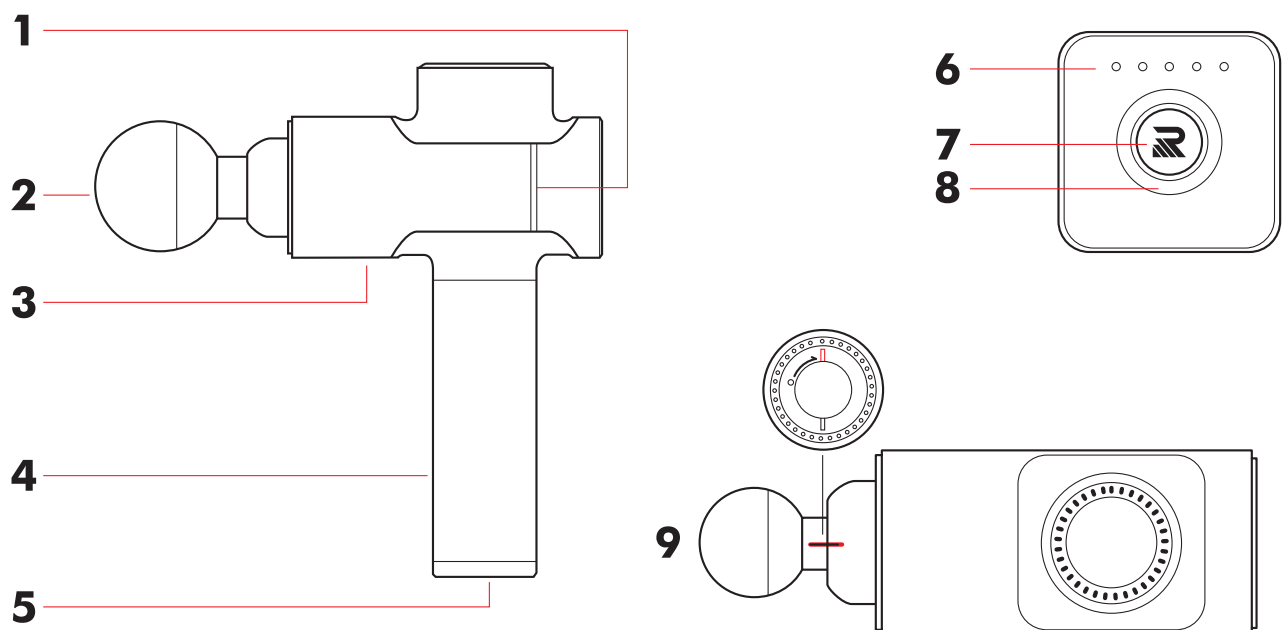



DANGER: TO REDUCE RISKS OF ELECTRIC SHOCK, FIRE, AND PERSONAL INJURY, OR PROPERTY DAMAGE, THIS PRODUCT MUST BE USED IN ACCORDANCE WITH THE FOLLOWING WARNINGS, CAUTIONS AND SAFETY INSTRUCTIONS:

IMPORTANT SAFETY INSTRUCTIONS

The Recovapro Max is for adult use only and should always be used in accordance with the following warnings and instructions.

When using the Recovapro Max use only on dry, clean surfaces of the body by lightly pressing and moving across the body for approximately sixty(60) seconds per region. The Recovapro Max should only be used on soft tissue areas of the body and should not be used on any hard or bony areas of the body, including the head. Please immediately stop use of the Recovapro Max upon experiencing any pain or discomfort, apart from light muscle soreness. Bruising can occur from use of the Recovapro Max, regardless of the pressure setting, and if bruising does occur you should not use the Recovapro Max on the bruised area until the bruise is fully healed. Use of the Recovapro Max should also be avoided on all bruises, contusions, rashes or irritated or injured areas of the skin until they are fully healed. When using the Recovapro Max be sure to keep fingers, toes, hair, and other body parts away from the back of the applicator to avoid pinching or hair entanglement.

The Recovapro Max should be kept dry, should never be immersed in water and should be kept away from liquids and heat sources at all times. The ventilation holes on the Recovapro Max should be kept free of dust, debris and/or any other substance that could penetrate the holes, including but not limited to liquids, dirt, sand or food. You should never operate the Recovapro Max for more than one (1) continuous hour at a time and after each use please allow the Recovapro Max to rest for a minimum of thirty (30) minutes before using it again.

THE RECOVAPRO MAX SHOULD ONLY BE CHARGED WITH THE RECOVAPRO MAX CHARGER THAT WAS SUPPLIED WITH THE UNIT. DO NOT CHARGE THE RECOVAPRO MAX OVERNIGHT OR LEAVE THE RECOVAPRO MAX UNATTENDED WHEN CHARGING OR IN USE.

After charging the Recovapro Max, the charger should be unplugged from the Recovapro Max's charging port prior to use. Prior to disposing of a Recovapro Max unit the battery must be removed. The battery must be complied with safely at an appropriate e-waste disposal or recycling facility.

The Recovapro Max should never be used for any purpose other than its intended use, should not be dropped, and should never be tampered with or altered in any way.

PLEASE DO NOT USE THE RECOVAPRO MAX, OR ANY PERCUSSION OR VIBRATION DEVICE WITHOUT FIRST OBTAINING APPROVAL FROM YOUR PHYSICIAN IF ANY OF THE FOLLOWING APPLY:

Pregnancy, diabetes with complications such as neuropathy or retinal damage, wear of pace-makers, recent surgery or injury, epilepsy or migraines, herniated disks, spondylolisthesis, spondylolysis, or spondylosis, recent joint replacements or IUD's, metal pins or plates or any concerns about your physical health. Frail individuals and children should be accompanied by an adult when using any percussion or vibration device. These contra indications do not mean that you are not able to use a percussion or vibration device but we advise that before doing so you consult your physician. Ongoing research is being done on the effects of percussion massage for specific medical disorders that may result in the shortening of the list of contra indications as shown above.



Caution: The user is cautioned that changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) this device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

NOTE: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation.

If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

FCC Radiation Exposure Statement

This equipment complies with FCC radiation exposure limits set forth for an uncontrolled environment.

This transmitter must not be co-located or operating in conjunction with any other antenna or transmitter.

