Technical overview

In this box you will find:

moonbird

your personal breathing coach

click-on charging cable

connects your moonbird to any USB-A power supply for charging

moonbird protection bag


a safe nest for you moonbird

user manual

instructions on how to use and take care of your moonbird

Your moonbird has the following functions

- ① **thumb sensor**
photoplethysmography (PPG) sensor to measure your heart rate
- ② **charging port**
to charge your moonbird

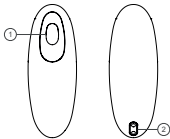Note: For firmware updates go to the Settings section in the app. Click on 'My moonbird' and check for firmware updates.

Charging your moonbird

Make sure your moonbird is fully charged before starting an exercise. When using your moonbird, the app will notify you when it needs charging.

To charge your moonbird:

- ① Plug one end of the micro USB cable in a USB power supply, such as your computer or the charger of your mobile phone.
- ② Use the magnetic end to connect to the charging port of your moonbird.
- ③ It takes approximately 3 hours before the battery is fully charged. Your moonbird is fully charged when the battery icon shown in the home screen in the app is full.

Note: Always use the magnetic USB charging cable provided in the box to charge your moonbird.

Taking care of your moonbird

Cleaning

Clean your moonbird every once in a while. To clean it, use a damp cloth with water, mild soap or alcohol gel (70%). Your moonbird is water resistant but not waterproof: try to avoid it being completely immersed by water.

It is important to completely dry your moonbird before charging it.

Storing

Keep your moonbird dry and safe in its protection bag.

Things to avoid

- Don't sterilize or boil your moonbird
- Don't completely immerse it with water

- Don't use cleaning agents other than soap or alcohol gel
- Don't put it in the microwave
- Don't put it in the dishwasher or washing machine
- Don't use your moonbird in bath or under the shower
- Don't leave it in the sun or in extreme heat for too long
- Don't drop your moonbird
- Don't try to disassemble the device in any way

Health notice

Moonbird is a consumer product that is intended to induce relaxation. Moonbird is not a medical device, nor are its services or products intended as a substitute for a diagnosis or medical treatment by a doctor.

Usually, slow-paced breathing exercises do not cause any health problems. If you have any health issues such as asthma or COPD, or if you don't feel well as a result of the breathing exercises, please ask your doctor for advice.

When you are new to breathing exercises, your body might need time to adjust to this new way of breathing. After a few weeks you should feel more calm and connected with your body.

Certification

The designated product is in conformity with the legislation in the markets in which it is offered:

Regulations	Code
RED	2014/53/EU
RoHS	2011/65/EU
Ecodesign Directive	2009/125/EC
IC	27965-MPROV1
FCC	2A3NN-MPROV1
UK	Radio Regulation 2017

Hereby, Moonbird BV, declares that the radio equipment type MBPROV1.3 is in compliance with Directive 2014/53/EU and UK Radio Regulation 2017. The full text of the EU/UK declaration of conformity is available at the following internet address:

<https://www.legislation.gov.uk/ukSI/2017/1206/regulation/13/made>

FCC and ISSED Declarations

Compliance statement (part 15.19)

This device complies with part 15 of the FCC Rules and to RSS of Industry Canada. Operation is subject to the following two conditions: (1) this device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

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This portable transmitter with its antenna complies with FCC/IC RF exposure limits for general population / uncontrolled exposure.

Warning (part 15.21)

Changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

Information to the User (Part 15.105 (b))

Note: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

This Class B digital apparatus complies with Canadian ICES-003.

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Recycling information

This symbol means the product must not be discarded as household waste, and should be delivered to an appropriate collection facility for recycling. Proper disposal and recycling helps protect natural resources, human health and the environment. For more information on disposal and recycling of this product, contact your local municipality, disposal service, or the shop where you bought this product.

Technical information

To use moonbird, iPhones must have iOS10.0 or newer (recommended devices are iPhone 5 and later). Android phones must have version 8.0.0 Oreo or above.

Feature	Description
Power supply	5VDC 500mA
Battery type	Lithium ion battery
Battery lifetime	6 x 20 min = 2 hours
Battery charging time	3 hours
Battery charging type	Magnetic charging cable to USB A male. Please fit this in a power adaptor with USB A female min 500mA.
Dimensions	122 mm x 54 mm x 44 mm
Weight	186 gr
Operating temperature	5-40°C
Storage temperature and humidity	-5 - 35°C, < 75% RH
Material	Biocompatible, medical grade silicone skin (RoHS and phthalates compliant)
Sensor	Optical reflectance PPG sensor to measure heart rate and derivatives
Bluetooth connectivity	Bluetooth® 5 - Bluetooth® Low Energy (BLE) BLE frequency range: 2402 - 2480MHz BLE max. transmitting output: -11.34 dBm

Warranty

Moonbird is made from biocompatible, medical-grade silicone and is covered by a limited manufacturer's warranty against any defects in materials and/or workmanship for two years from the date of original purchase.

Defects discovered within the warranty period should be reported to Moonbird BV. Moonbird BV reserves the right to exercise its discretion in determining whether or not to replace the defective moonbird free of charge.

Any warranty claims must be supported with reasonable evidence including proof of date of purchase, alongside these warranty conditions. This warranty does not cover cosmetic deterioration or damage caused by general wear and tear, physical or natural destruction, accident, misuse, neglect or other external causes. Any attempt to take moonbird apart will invalidate this warranty.

Disclaimer

Users of moonbird do so at their own risk. Neither moonbird nor its resellers assume any responsibility or liability for any injuries or damages, physical or otherwise, resulting, directly or indirectly, from the use or possession of this product.

Your device is not a toy, keep it out of reach of small children.

Moonbird reserves the right to revise this user guidance and to make changes from time to time to its contents without obligation to notify any person of such revision or changes. The latest user manual can be found at moonbird.life/usermanual.

The product model may be changed for improvements without notice at any time at the discretion of moonbird. For complete information on limited product warranties, disclaimer of warranties, limitation of liability, please see the moonbird website at moonbird.life.

June 2022, version 8.0_English
All rights reserved
Designed and made in Belgium
MODEL: MBPROV1.3

English

For other languages visit moonbird.life/usermanual

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Moonbird LTD