


可视冰热敷黑头仪

产品说明书

产品部件

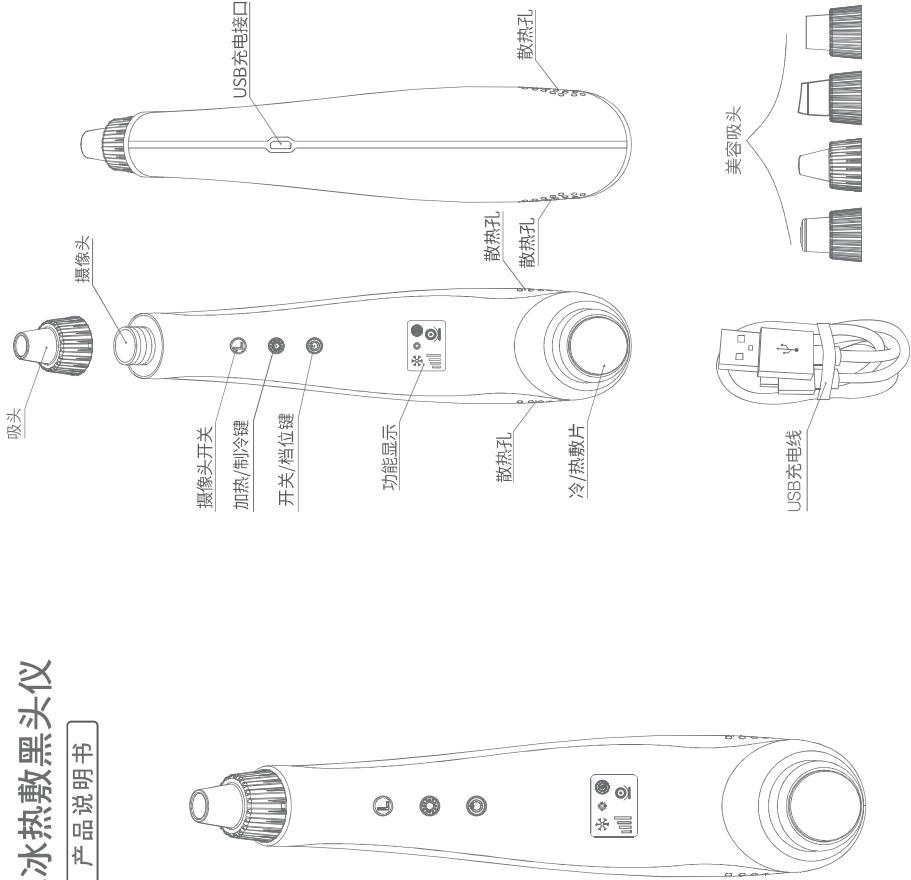

产品简介

本产品是一款集吸黑头与热敷/冷敷加可视化的多功能美容仪，
A:采用真空吸技术对肌肤进行轻柔的拉伸按摩并吸出毛孔中的黑头，污垢，让皮肤更加剔透亮泽。
B:搭配30万像素的高清摄像头加20倍放大镜片能清楚看到毛孔油脂黑头等杂质，边看边吸，清洁毛孔不留盲区。
C:采用热敷/冷敷、令肌肤弹性提升与塑造紧致饱满的肌肤。

产品参数

型号: ZSX-H300-P	产品尺寸: 212.0X44.0MM
充电电压: 5.0V	净重: 155g
额定功率: 6.5W	电池容量: 1000 Mah
续航时间: 2H	像素: 30万高清/放大20倍
热敷温度: 42-45℃	冷敷温度: 10℃以下
产品颜色: 黑白色	吸力: 65 Kpa

操作指示

大圆孔美容吸头

吸力比较强,用于吸黑头,吸黑头时,需要慢慢在皮肤上移动,请勿长时间在一个部位使用,这吸头重点针对黑头,粉刺,油脂,灰尖,化妆品残留等。脆弱、敏感皮肤者及眼部周围慎重使用。



小圆孔美容吸头

吸力小,用于皮肤紧致,敏感肌肤,而不适用太强的吸力,以及吸鼻翼角薄弱部位,就用这吸头。



椭圆孔美容吸头

主要用于祛除细纹,需吸住皮肤后快速拨开反复多次,提拉脸部,可起到促进皮肤微循环,增加皮肤弹性。用于脸颊,眼角,鼻翼等部位。



微晶美容吸头

细微晶可以轻柔去除肌肤老化角质,令肌肤细腻又嫩滑。



使用方法

一：长按开关/档位键3秒，产品进入待机状态，显示屏指示亮

A:开机状态

- 1.按开关/档位键一次，则进入吸黑头模式的1档（低档），显示屏吸力指示亮低档。
- 2.按开关/档位键二次，则进入吸黑头模式的2档（中档），显示屏吸力指示亮中档。
- 3.按开关/档位键三次，则进入吸黑头模式的3档（中高档），显示屏吸力指示亮中高档。
- 4.按开关/档位键四次，则进入吸黑头模式的4档（高档），显示屏吸力指示亮高档。
- 5.按开关/档位键五次，重新回到待机状态。

B:可视功能

- 1.按摄像头开关键一次，摄像头亮起。
- 2.打开手机APP/进入设置栏/打开连接WiFi/找到相应的网络连接。
- 3.退出设置返回到画面，点播按钮进入可视窗口即可看到画面。
- 4.点击“录像”后可将吸黑头仪摄像头传输的图像保存到手机内。
- 5.点击“拍照”后可将吸黑头仪摄像头传输的图像以图片的形式保存到手机内。
- 6.后再按摄像头开关键，关闭可视功能。

C:产品进入开机状态

- 1.按加热/制冷键一次，则进入热敷模式，显示屏太阳指示亮。
- 2.按加热/制冷键二次，则进入冷敷模式，显示屏雪花指示亮，3分钟后自动关机。

二：长按开关/档位键3秒，则产品退出开机状态，即关机

注意事项：

- 1.使用前，请将手上的水或其它液体擦干，以免水或其它液体从按键与气孔进入机体内，从而损坏机体。
- 2.使用前，后，请不要用磨砂膏与其它角质产品。
- 3.请不要用力将吸头长时间按压在皮肤上，会使本机产生过强的吸力，损伤皮肤。
- 4.由于带有微晶头的吸头的磨皮功能，建议干性的皮肤一周使用一次，油性皮肤一周使用二次，每次使用不超过5分钟，具体视皮肤而定。
- 5.请不要在同一位置持续使用，会引起皮肤红紫，出现这种情况，请立即停止使用。
- 6.本产品不能给敏感肌肤、精神紧张与无经验的人使用，特别是不能给儿童使用。
- 7.本产品不适合整容、隆鼻与皮肤美容手术的部位使用。
- 8.孕妇与心脏病患者不能使用本产品。

吸黑头具体操作：

① 热敷面部

使用本机热敷功能敷脸打开毛孔，或用热毛巾热敷至5分钟打开毛孔。

③ 选择档位

调节好适合自己皮肤的档位，可参上页的图表（起润滑作用）。

⑤ 冲凉清除面部

吸出黑头后用清水清洗面部。

② 涂上导出液

把导出液滴到化妆棉上敷3至5分钟，可以利用产品的热敷功能同时使用。

④ 吸附脸部

在脸上有黑头处缓慢移动（每次使用1至2分钟）。

⑥ 护理脸部

涂上收缩水或面膜，打开水雾模式，水雾面部喷雾毛孔或敷补水面膜收缩毛孔。

可视功能具体操作

- 1.手机微信扫描二维码下载APP: (苹果系统和安卓系统)，中英文根据使用设备语言自动变更。
- 2.扫描二维码下载对应的APP:




ISO 二维码 安卓二维码

- 3.按摄像头开关键一次，摄像头亮起。



- 4.APP连接摄像头

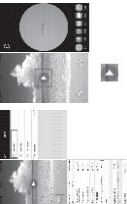

1.打开APP，点击右上角设置
2.用户可根据使用设备语言选择，简体中文或繁体中文或英文。
3.选择语言后，点击右上角Sheng Long，即可开始使用。

备注：因手机系统差异出现的速度也会不同，正常3-5S可以出图。

注意事项

可视吸头取下后用清水冲洗干净，请勿用大量接冲液仪器机身

以下情况，不建议使用本产品，以免发生意外

- 1.正接受医师治疗和身体感觉异常者
- 2.恶性肿瘤者
- 3.孕妇、产妇、儿童
- 4.感染病者
- 5.有皮肤知觉障碍者、或皮肤异常者
- 6.急性病者

保修卡

商品名称	可视黑头仪	型号	ZSX-H300-P
购买方信息	姓名 地址 电话	省 市/县 区	
经销商信息			
故障描述			
维修时请填写详细写保修卡上要求的内容			本产品整机保修1年


合格证

生产日期


质检

VISUAL COLD HOT COMPRESS BLACKHEAD
REMOVAL APPLIANCE

Manual

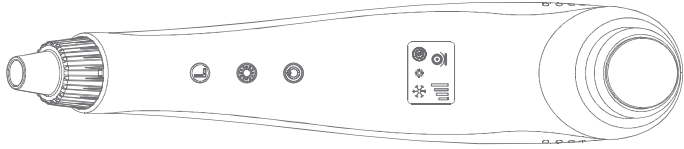

Product Parts

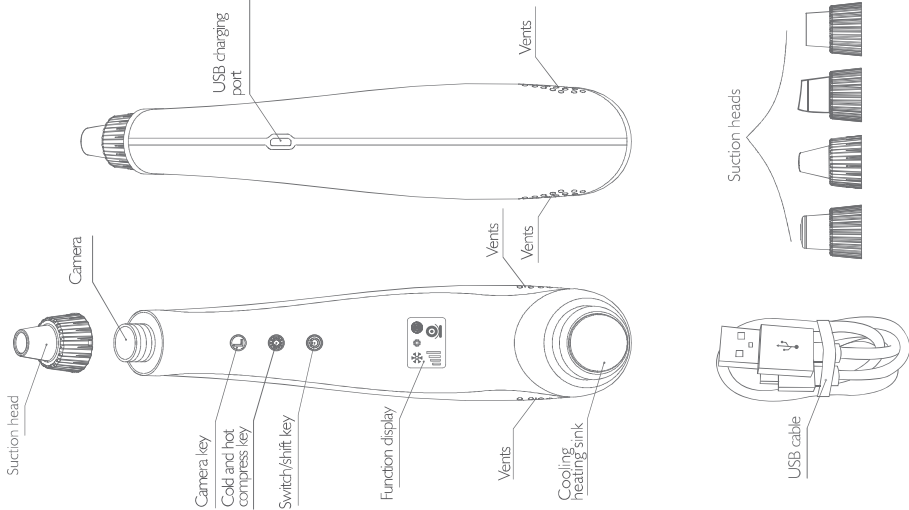

Product Introduction

The product is a multi-functional blackhead removal appliance with cold and hot compress, visualization.

A: Using vacuum suction technology of the product to stretch and massage your skin gently to help suck out the blackheads, dirt deep in the pore, which makes your skin smoother and shinier.

B: The product is equipped with an HD camera with 300,000 Pixel and a 20x magnifier, which can make you see the impurities such as the blackheads, grease in the pores. Watching while you suck out the blackheads, and clean up the pores without leaving a blind spot.

C: Adopting cold and hot compress to improve skin elasticity and shape firm and plump skin.

Product Parameters

Model: ZSX-H300-P
Product size: 212.0X44.0mm
Net weight: 155g
Charging voltage: 5.0V
Battery capacity: 1000 Mah
Rated power: 6.5W
Pixel: 300,000 HD/ 20-times magnifier
Battery life: 2H
Hot compress temperature: 42°C-45°C
Cold compress temperature: <10°C
Color: Black and white
Suction power: <65 Kpa

Operating Instructions

Large round hole suction head

This suction head with a powerful force that is not suitable for fragile, sensitive skin or the skin around your eyes. It is used to remove blackheads, grease, dirt, cosmetic residues on the skin. It needs to be moved slowly on your skin when you operate the product and please do not use it in the same area for a long time.



Small round hole suction head

This suction head with low force, which is suitable for tender, sensitive skin, and the weak parts of the corner of the nose.



Oval suction head

This suction head is suitable for removing blackheads on cheeks, corners of eyes, and wings of the nose. Use the vacuum absorption function of the product to absorb facial skin, lift the skin, and then quickly remove the product. Repeated operations will promote the blood circulation of the skin, increase skin elasticity, and reduce wrinkles and fine lines.



Microcrystalline head

The head can gently exfoliate aging skin, leaving the skin delicate and smooth.



How to use

1: Press the switch/shift key for 3s, the product enters into the standby state, and the display indicator lights up.

A. Power on

1. Press the switch/shift key once to enter the first gear (low gear) of the blackhead suction mode, and the indicator on the display will light up.
2. Press the switch/shift key twice to enter the 2nd gear (medium gear) of the blackhead suction mode, and the indicator on the display will light up.
3. Press the switch/shift key three times to enter the 3rd gear (high and medium gear) of the blackhead suction mode, and the indicator on the display will light up.
4. Press the switch/shift key four times to enter the 4th gear (high gear) of the blackhead suction mode, and the indicator on the display will light up.
5. Press the switch/shift key five times to return to the standby state.

B. Visualization

1. Press the camera key for once, and the indicator will light up.
2. Open the App on your mobile phone/enter into the setting bar//turn on WiFi /Use the corresponding network connection.
3. Exit the Settings to return to the desktop, click the Play button to enter the Visual window to see the picture.
4. Click "record" to save blackhead removal Video transmitted by the product camera in the phone.
5. After clicking "photograph", the images transmitted by the product camera can be saved in the mobile phone in the format of a photo.
6. Then press the camera switch button to turn off the visual function.

C. Power-on State

1. Press the cold and hot compress key for once to enter into the hot compress mode, and the sun indicator will light up.
2. Press the cold and hot compress key twice to enter into the cold compress mode, and the snow indicator will light up.

1: Press the switch/shift key and hold it for 3s to power off.

Precautions :

1. Before using, please dry up your hands to prevent water or other liquids from entering the products from the vents or buttons, thereby damaging the product.
2. Please do not use exfoliating creams and other exfoliating products at the same time.
3. Do not press hard on the skin with the suction head for a long time, which will create a strong pulling force, thus damage your skin.
4. Because of the exfoliating function of the microcrystalline heads, it is recommended to use it once a week for dry skin and twice a week for oily skin, and should not exceed 5 minutes for each use. It depends on the skin.
5. Please do not use the product continuously in the same position, or it will make the skin red and injured. If this happens, please stop using it immediately.
6. This product is not recommended to be used by people who lack five senses, people who have a mental disease, and inexperienced people. Especially not for children.
7. The product is not suitable for people who get cosmetic surgery, such as eyelids surgery/ face lift, nose reshaping, or other cosmetic treatments.
8. The product is not suitable for pregnant women and patients with heart disease.

Operation of blackhead removal



① Face hot compress

Using the hot compress function of the product, or put a warm towel to your face for 3 to 5 minutes to help open the pores of the skin.



③ Suction level choice

Choose a proper suction level for your skin, can apply a little lotion before using (lubrication function)



④ Adsorption

Put the product on the part where there are blackheads, and then slowly to suck them out (1-2 minutes each time)



⑤ Clean up your face

Wash up your face with water after sucking out the blackheads.



⑥ Face care

Put on some toner or a mask, or use the cold compress of the product to help close up the pores.



② Apply the blackhead derived liquid

Drip the liquid onto the makeup cotton and then put it on where the blackheads are for 3-5 minutes. You can use the hot compress function of the product at the same time

Visual function operations

1. Scan the QR code on WeChat to download the corresponding APP




2. Scan the QR code to download the corresponding APP



QR code for ISO QR code for Android system

3. Press the camera key for once, and the indicator will light up.

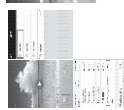

4. Camera connection

Remarks:
The speed of drawing pictures will be different according to the different mobile phone systems, and pictures can be generated within 3-5s normally.

Precautions

These suction heads can be removed and rinsed with clean water, but do not make the body of the product directly will water.

To avoid any accidents this product is not recommended to be used by the following people:

1. Those who are being treated by a doctor and have abnormalities in their body.
2. Patients with malignant tumors.
3. Pregnant women, puerpera, children.
4. Infectious patients.
5. People with skin sensation impairment or skin abnormalities.
6. Acutely ill patients.

Warranty Card

Product name	Visual cold hot compress blackhead removal appliance	Model	ZSX-H300-P
Buyer Information	Name Address Telephone	Province City/ County District	
Dealer information			
Fault description			
Please fill in the content required on the warranty card in detail when repairs are required.			This product is guaranteed for 1 year

Certification

Certification _____

Quality inspection _____

This device complies with part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) This device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

Any Changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

Note: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.