




MASSAGE CHAIR

USER MANUAL

MODEL: MF-600

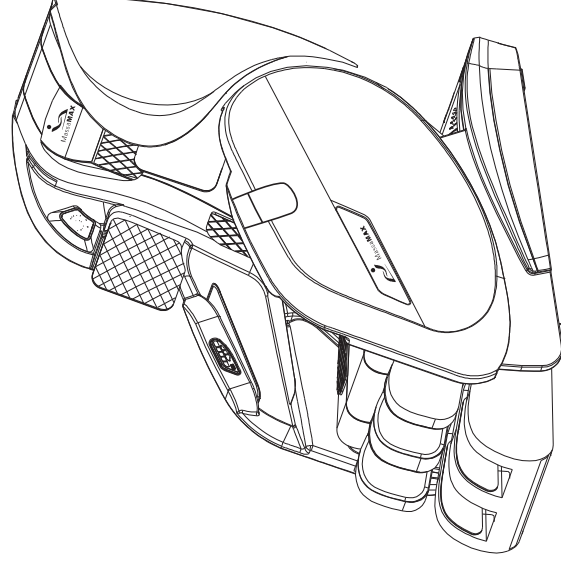

TROUBLESHOOTING

PROBLEM	SOLUTION
The rollers and the backrest pad make frictional sound. Motor sound is clearly hearable during operation. Loud sound is given when tapping and kneading massage is performed	This is the sound produced by mechanical structure instead of failure of internal components. This may occur during normal operation.
Message is suddenly stopped and the massage rollers return to their original position.	*The preset auto-timing has been reached. *In order to avoid the massage armchair from causing accident due to overhear, this chair is equipped with overhear reboot protection equipment. If the chair stops and fails to restart as a result of overhear, turn off the main power switch for at least 30 minutes in order to prevent it from suddenly starting after components are cooled down. *The power button was pressed accidentally.
The remote controller does not work.	*Check the power plug and ensure that it is indeed plugged into the socket. *Check whether the lines on the back of the massage chair are properly connected. *Confirm the main power switch is turned on.
The power cord and plug are excessively hot.	Stop using immediately and contact the service centre.
The massage stops in the process.	Turn off the main switch on the back of the chair and turn on 5 seconds later. If the same problem is occurred again, please turn off the main switch, unplug the power cord and contact the customer service.

CONTENTS

Safety Instructions	02-06
Operating precautions	07-08
Installation Guide	09-13
Names and Functions of Parts	14-15
Before starting massage	15
How to use the chair	16-17
Maintenance and Storage	18
Specifications	19
Troubleshootings	20
Thanks for choosing our massage chair. This operating manual contains important information about how to use the chair safely and correctly. Be sure that you have read all safety instructions thoroughly before using your chair. Keep this manual nearby so that you can refer to it any time.	

IF YOU HAVE ANY PROBLEM PLEASE FEEL FREE TO CONTACT US

SAFETY INSTRUCTIONS

CAUTION

TO AVOID THE RISK OF INJURY

- 1 Persons with any of the conditions listed below are not suitable to use the chair :
 - Osteoporosis patient.
 - People with heart disease including embedded pacemaker or other medical electronic device inside the body.
 - For pregnant women or people with health issues, please consult doctor before use.
 - Patient with a high fever forbid to use this machine.
 - Perceptual disabilities patients forbid to use this machine.
 - People over 80 years old forbid to use this machine.
 - Kids and the drunk people are not allowed to use this machine.
 - Please don't sleep in the chair.
 - Don't use the chair when your body is wet
- 2 This unit should not be used by children. Turn off the power when not in use and store out of reach of children.
- 3 Do not insert your hands, fingers, feet or head into the following spaces:
 - Between the backrest and seat or armrest
 - Between the leg rest and seat or armrest
 - Between the seat and armrest
 - Between the back cover and leg cover
 - Back of the leg rest
- 4 Do not use the unit against bare skin. While thin clothing may increase effectiveness, exposing the skin directly to the massager may irritate the skin.
Do not use the unit while wearing anything hard on your head like a hair accessory, etc. Do not use the unit for massaging the head, abdomen, front neck, elbows or knees.

SPECIFICATIONS

Product Name	Massage Chair
Model No.	MF-600
Rated Voltage	AC120V
Rated Frequency	60Hz
Rated Power Consumption	100W
Cover Cloth	PU synthetic leather
Timer	15/25/30 minutes
Weight Capacity	6.25 FT/242.5 LBS
Dimension	Main Unit:43.30×23.22×18.50" Backrest :38.97×17.51×24.21" Armrest :43.11×13.77×21.85"
Gross Weight	Main carton :92.59LBS Backrest carton:55.66LBS Armrest carton :44.53LBS
Net Weight	Main carton :87.41LBS Backrest carton :49.16LBS Armrest carton :37.91LBS
Length of the wire	Controller wire: 70.8" Power supply wire: 70.8"
Usage condition	Environment temperature: 10°C-40°C Contrasting humidity: 30-85RH
Safety feature	Equipped with overheated and power surge safe protection

※ ※ It is not a medical equipment and not for medical purpose. If you feel uncomfortable during use, please stop using it immediately.

MAINTENANCE AND STORAGE

CLEANING

- ✗ Prior to the maintenance work, make sure to disconnect the power cord. Be careful that hands must not be wet to prevent electric shock or injury.
- ✗ Do not use benzene, thinner or alcohol to clean the chair, nor spray any liquid containing any chemical agents to prevent the risk of electric shock or catching fire.
- ✗ Do not iron to prevent damage to the fabric.

For main body, areas covered with synthetic leather, back pad, pillow and seat cushion

Use a dry cloth moistened with diluted neutral detergent and thoroughly wrung out to clean the chair, and then wipe with a dry cloth to completely remove the residual detergent. Thereafter, allow the parts to dry naturally.

For Remote controller

Use a dry cloth to clean the remote controller.
Never use a wet cloth for cleaning to avoid malfunction.

STORAGE

- ✗ Do not place the unit in areas exposed to direct sunlight, or near heaters, hot air blowers or any other heating devices to avoid discoloration or deterioration on the surface of the holster.

1. Wipe off dust and dirt on the chair each time after use.
2. Keep the chair clean and away from places of high humidity.
3. If the chair will be laid aside for a long period of time, cover it with a dust cover and unplug its power cord.

- 5 Do not use excessively strong massage action on the back of the neck.
Do not stretch the back muscles for more than approximately 16 minutes at a time.
To avoid overstretching the muscles which may result in discomfort, do not use the unit for more than 5 minutes when first using it. You can gradually increase the duration and strength of the massage as you become used to it.
- 6 Stop using the unit immediately should you experience any discomfort and consult the physician.
- 7 Do not go to sleep while using the unit. Do not use the unit after drinking alcohol.
- 8 When the seat back is in a reclined position, do not sit on it or sit with your legs resting on the headrest.
Do not sit on the leg rest as the unit may fall over.
- 9 Do not allow children or pets to play on or around the unit, especially during operation.
- 10 Do not unplug the unit or turn it off during operation. Unplug the unit immediately if there is a power outage.
- 11 Be sure to start on the gentle massage program if using the massager for the first time.
- 12 Do not allow children to play on this mechanized furniture or operate the mechanism. Leg rest folds down on closing so that a child could possibly be injured.
Always leave in an upright and closed position. Keep hands and feet clear of mechanism. Only the occupant should operate it.
- 13 Remove hard objects, such as false fingernails, wristwatches or rings, when massaging hands or arms. Empty pockets of hard objects when massaging the area of the seat or thighs.
- 14 People with any of the conditions listed below should not use the Hip stretch:
 - People suffering from acute lumbago, spinal disc herniation, or misalignment of the vertebrae
 - People with hip joint problems such as hip osteoarthritis
 - People with osteoporosis
 - People who felt numbness or listlessness in the legs It may worsen the symptoms.
- 15 Always return the leg rest to its original position before moving the unit.
Failure to observe this precaution may result in injury due to the leg rest returning to its original position when the unit is being moved.

- 16 Store the attachment screws out of reach of children. Accidental ingestion may occur.
- 17 Do not stand, jump on, or place objects on the seat until the air is completely extracted for the air massage.
- 18 Do not hold the sole massage section when moving the unit.
Doing so may result in the section sliding as you move the unit, causing injuries.
- 19 Do not drag or push the unit in an installed state.
Protect your floor (place a mat or similar) when moving the unit on its castors to avoid damage to your floor.
- 20 If you have not used the unit for a while, carefully read this operating instructions again and check that the unit is operating normally before use.
Failure to do so may result in accident or injury.
- 21 Do not use at the same time as another medical device. Doing so may cause an accident or deteriorated health.
Do not allow children to use the controller holder fixing screw.
- 22 Failure to do so may result in an accident or the child accidentally swallowing the key or screw.

DANGER

TO REDUCE THE RISK OF ELECTRIC SHOCK:

- Always unplug this unit from the electrical outlet immediately after using.
- Always unplug this unit from the electrical outlet immediately before cleaning.

WARNING

TO REDUCE THE RISK OF BURNS, FIRE, ELECTRIC SHOCK, OR INJURY TO PERSONS:

- 1 Fully insert the power plug to reduce the risk of short-circuit and fire.
- 2 The unit should never be left unattended when plugged in. Always unplug the unit when not in use.
- 3 The unit is not to be used by people who cannot move or communicate on their own.
- 4 Only use the unit for the purpose described in these instructions.
- 5 Do not stand on the unit.



Vibrating

This is button for vibrating ON/OFF , press to start vibrating massage, press again to stop



Footroller

This is button for foot roller ON/OFF , press to start foot roller massage, press again to stop



Backrest up and footrest down

Adjust the angle of backrest when press the button for long.

Backrest down and footrest up

Adjust the angle of backrest when press the button for long.



- Power button, the power is on when press the button the first time
- when massage function is going on ,short press will pause all massage function, the power will be off when press again for 3 seconds.



Zero gravity

One key to reach zero gravity position



Heating

This button is for heating on/off

HOW TO USE THE CHAIR

BUTTON AND FUNCTIONS

Music



Bule Tooth ON/OFF Button,press once to turn on,press again to close .

Back kneading



This is button for back kneading on/off

Back extending



This is button for back extending on/off

Airbag intensity




This button is for adjusting the strength of airbags intensity

Message strength



This button is for adjusting the strength of message

Auto mode

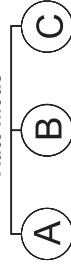

Ease mode Relax mode Sleep mode



Press this button will add Air intensity /Time/Message strength




This button is for adjusting the length of the massage session,the default working time is 15 minutes ,5-30 minutes adjustable ,5 minutes as a adjustable unit.



Press this button will reduce Air intensity /Time/Message strength



Airbags on-off button, press to start airbags massage, press again to stop

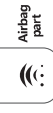

The button allows you turn on/off the air pressure function at certain areas.
Press it to select air massage position in the following order:
Legs → Arms → seats → shoulders → OFF

- 6 Do not use any accessories other than those recommended by the manufacturer.
- 7 Always return the seat to the upright position and make sure the leg rest is completely retracted. Be careful that nothing is obstructing the leg rest while retracting. Keep children away from the leg rest.
- 8 Never operate the unit if it has a damaged power cord or plug. Return the unit to the nearest authorized service center if it is not working properly, if it has been dropped or damaged, or if it has been immersed in water.
- 9 Keep power cords away from heated surfaces. Do not carry the unit by the power cord.
- 10 Never block the air openings during operation and make sure that all air openings are kept free of lint and hair,etc. Do not drop or insert anything into any of the openings of the unit.
- 11 Do not use the unit on top of heating appliances, such as electric carpets, etc. Do not operate the unit under blanket or pillow. Excessive heating can occur and cause fire, electric shock, or injury to persons.
- 12 Do not use the unit in bathrooms or other damp or humid places as this may cause electric shock or cause the unit to malfunction. Do not spill water, etc. onto the controller.
- 13 Do not use the unit outdoors.
- 14 Do not use or operate where aerosol spray products are being used or where oxygen is being administered.
- 15 When unplugging the unit, turn off all controls before removing the plug from power outlet. Do not use the unit with a transformer as this may lead to malfunction or electric shock.
- 16 Connect this unit to a properly grounded outlet only. See Grounding instructions.
- 17 Use heated surfaces carefully. May cause serious burns. Do not use over sensitive skin areas or in the presence of poor circulation. The unattended use of heat by children or incapacitated persons may be dangerous.
- 18 Always raise the comfort cushion to check that the fabric covering the movement area of the massage heads in the backrest has not been ripped before using the unit. Please also check other areas to ensure that the fabric has not been ripped. No matter how small the rip, if you find one then immediately stop using the unit,

unplug the power cord and have the unit repaired at an authorized service center.

- 19 Do not push hands, elbows, arms, feet or legs into the fabric covering the movement area of the massage heads in the backrest while using the unit.
- 20 Do not use the unit in close proximity to loose clothing or jewelry.
- 21 Keep long hair away from the unit while in use.
- 22 Always use sockets and wiring devices at the correct rated value. Use at the specified power supply 120 V, AC.
- 23 Failure to do so may result in electric shock or fire due to overheating.
- 24 When Heat massage heads are used, do not let the massage heads touch the same place for a long time. Do not use the Foot Heat massage for a long time. Low-temperature burn may occur.
- 25 Do not damage the Heat (heater) section. Do not stick pins or needles into the unit. Do not damage with a sharp object. Doing so may damage the internal heater, causing fire or electric shock.
- 26 Sit on the unit after confirming there is no foreign matter between the parts of the unit. (Confirm that there are no foreign objects stuck in the comfort cushion, leg rest or hand & arm massage section.) Sit down and ensure that the leg rest is properly locked in place. Failure to observe this precaution may result in accident or injury.
- 27 Periodically remove dust, etc. from the power plug. Failure to do so may result in insulation failure due to humidity, etc., which may cause fire. (Remove dust using a dry cloth.
- 28 Do not insert or unplug the power plug with wet hands. Doing so may cause electric shock.
- 29 Do not massage any one point for more than 5 minutes at a time. Excessive massaging can overstimulate the muscles and nerves and result in an adverse effect.
- 30 After each massage, slide the power switch, which is located on the right side of the chair, to the off position and turn its lock switch to the lock position and remove the key and power plug.
- 31 Do not attempt to open or disassemble any part of the unit. Only use an authorized service to fix your unit.

REMOTE CONTROLLER



1. TFT Display
2. Tooth Button ON/OFF
3. Air intensity/Time/Massage strength Increase
4. Back kneading ON/OFF Button
5. Air intensity Adjust Button
6. Waist kneading ON/OFF Button
7. Air intensity/Time/Massage strength Decrease
8. Auto Mode A/B/C
9. Vibrating ON/OFF Button
10. Airbag auto ON/OFF Button
11. Airbag part ON/OFF Button
12. Foot roller ON/OFF Button
13. Pause/Power ON/OFF Button
14. Zero-G ON/OFF Button
15. Heating ON/OFF Button
16. Backrest down And Footrest up
17. Backrest up And Footrest down
18. Waist extending ON/OFF Button
19. Massage strength Adjust Button
20. Time Adjust Button
21. Back extending ON/OFF Button

*When the chair is turned off, the LCD display of remote controller shows "OF".

BEFORE STARTING USING THE CHAIR

1.CONNECT THE POWER CORD PLUG WITH THE RECEPTACLE.

Check if the power cord, remote controller cord or power plug is damaged and check whether there is any hair or dust on the plug.

4.MAKE SURE TO CHECK THE SURROUNDING CONDITIONS.


Make sure that there is no person, pet or object behind the chair, in front of and below the footrest.

3.CHECK THE CONDITION OF THE MESSAGE CHAIR.


- Hold up the back pad and check for any tear in the backrest cloth as well as in any other parts.
- Should there be any tear, no matter how small it is, stop using the unit immediately, unplug the power cord and have the unit repaired.
- Check whether the power cord, remote controller cord or any other object is pinched in the main body.

NAMES AND FUNCTIONS OF PARTS

MAIN UNIT

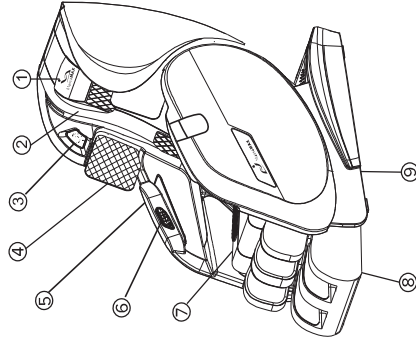

1. Cushion for Head
2. Backpad
3. Speaker
4. Shoulder Airbag
5. Armrest
6. Keyboard shortcut start with USB charging port
7. Seat Cushion

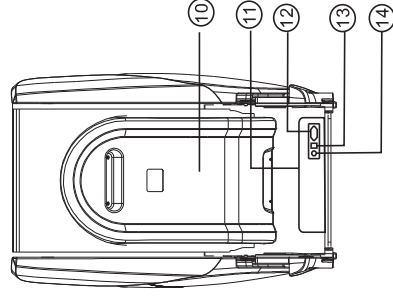

8. Footrest
9. Pedestal
10. Back cover
11. Electric Control Cabinet
12. Power Supply
13. Power Switch
14. Fuse

- 32 If the massager functions abnormally, immediately turn off the power and have the massager checked by an authorized service center.

OPERATING PRECAUTIONS

CAUTION

TO AVOID DAMAGING THE CHAIR, FOLLOW THESE PRECAUTIONS

- 1 Use the chair only when it is placed on an even, flat surface. Failure to observe this may result in upset.
- 2 Do not sit on or drop the controller.
- 3 Do not pull the controller cord or the power cord. Be careful not to catch the cord of the controller between the armrest and the chair back.
- 4 Do not treat the chair roughly, such as turning the chair on its side, turning the chair over or standing on it.
- 5 Do not treat the covering roughly.
- 6 Keep sharp or pointed objects away from the covering of the chair.
- 7 Be careful not to drop lit ashes, lit cigarettes or matches on the chair.
- 8 Exposure to direct sunlight can cause fading or a change in color of the covering.
- 9 Do not repeatedly operate any switch at short intervals. Such action may cause the switch to malfunction.
- 10 Some noise may be heard coming from the massager while in use. This is due to the structure of the massager and is normal.
- 11 Make sure there are no obstacles behind the chair before reclining it. If the back hits a wall or pillar, the chair may malfunction. Be sure to have ample space behind the chair. Recline the chair slowly in order to avoid contact with an obstacle.
- 12 When excessive pressure is applied to the massage heads, they may stop moving for safety.
- 13 When moving the heads up from the lowered position, ones body weight may stop the heads' movement. In such a case, slightly raise your body and allow the massage heads to move up.
- 14 Do not sit on the chair with wet body or hair.
- 15 The unit must not be used "wet rooms" (sauna, swimming pool) or outdoors.

- 16 Always hold the power plug, not the power cord, when unplugging it from a socket. Failure to observe this precaution may result in electric shock or short-circuit.

GROUNDING INSTRUCTIONS

This product must be grounded. If it should malfunction or break down, grounding provides a path of least resistance for electric current to reduce the risk of electric shock. This product is equipped with a cord having an equipment-grounding conductor and a grounding plug. The plug must be plugged into an appropriate outlet that is properly installed and grounded in accordance with all local codes and ordinances.


DANGER

Improper grounding may lead to electric shock. Check with an authorized service center if you are in doubt as to whether the unit is properly grounded. In case the supplied plug does not match the outlet, do not attempt to modify the plug. Please contact a qualified electrician to have a proper outlet installed.


This unit has a grounding plug. Make sure that the unit is connected to a power to a power outlet with the same configuration as the plug. Do not use an adaptor with this unit.

MOVING THE CHAIR


As the following illustration, push the back rest downwards and tilt chair on to its wheels at 45° angle. Slowly push and maneuver the massage chair to desired location. Lower massage chair slowly and smoothly until it is fully flat on the ground.

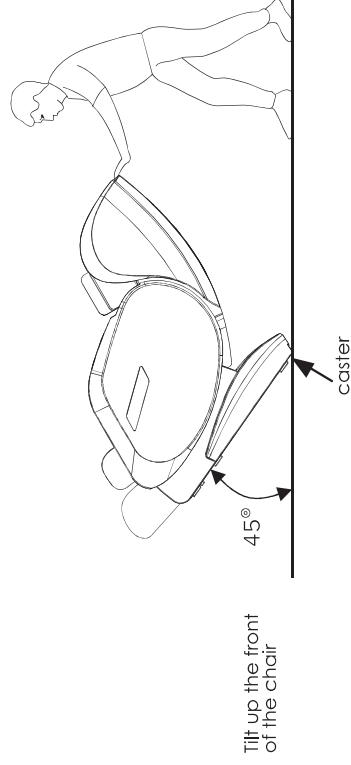

12. Attach the chair pad to the backrest with the zipper and Velcro tapes.

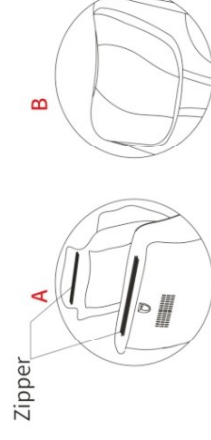


13. Insert the controller into the seat. (Pls notice the gap position);
• Tighten the screw thread.

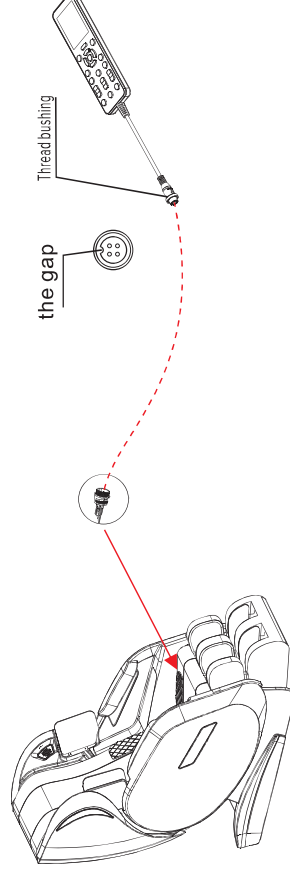

14. Insert the plug to the power socket, the machine on power.


- The machine is power on when the button turn " " and power off when the button turn "O"

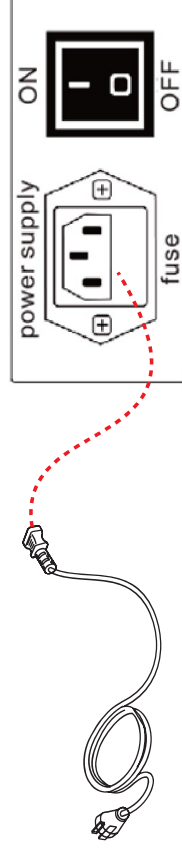

Caution: Please connect the controller first; then turn on the power, in order to avoid the failure and damage of the controller.


9.Connection of left Armrest

- Connect the jack cable from armrest with the female cable from backrest,
- Firmly insert the air plug of main unit into the air pipe of armrest.
- Insert the armrest along the holder O the shaft P; Put Q up the S.
- Use the Allen Key to fasten the screw.

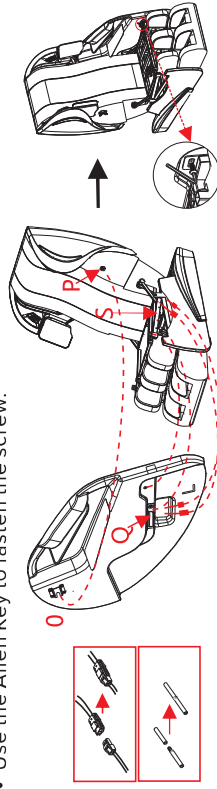


10.Connection of right Armrest

- Connect the jack cable from armrest with the female cable from backrest,
- Firmly insert the air plug of main unit into the air pipe of armrest.
- Insert the armrest along the holder T the shaft U; Put V up the W.
- Use the Allen Key to fasten the screw.

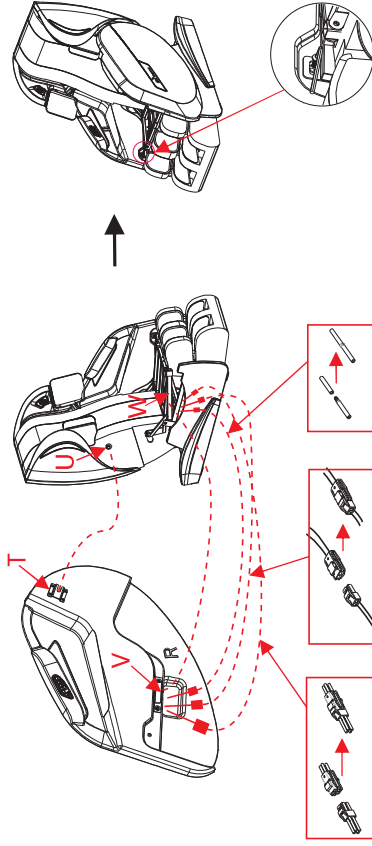

11.Connection of seat cushion


- Take out the seat cushion and place on the seat frame,
- Connect the air plug of the air cushion, and put the cushion on the seat.

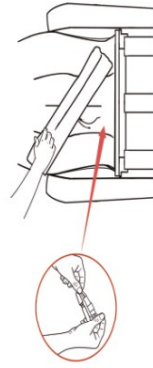

INSTALLATION GUIDE


SETTING UP THE MESSAGE CHAIR

- Ensure adequate space around the unit and place it on a flat surface. Avoid placing the unit exposed to direct sunlight or hear any heating device or in other locations with a high temperature to avoid discoloration or deterioration on the surface of the holster or back pad.
- It is recommended to place a carpet or soft mat underneath the unit to avoid damaging the floor.
- The size of such carpet or mat should be large enough to cover the surface where the unit and the footrest may com in contact with the floor.
- With the backrest reclined backward, the footrest goes up. To avoid collision with any surrounding objects, a clearance of at least 19.5" at the front and at least 19.5" at the back should be reserved around the unit.

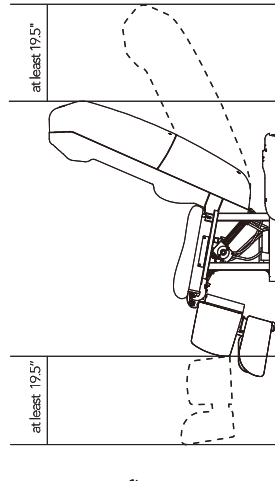

ASSEMBLY

1.Take out all the following parts from the 3 boxes for installation

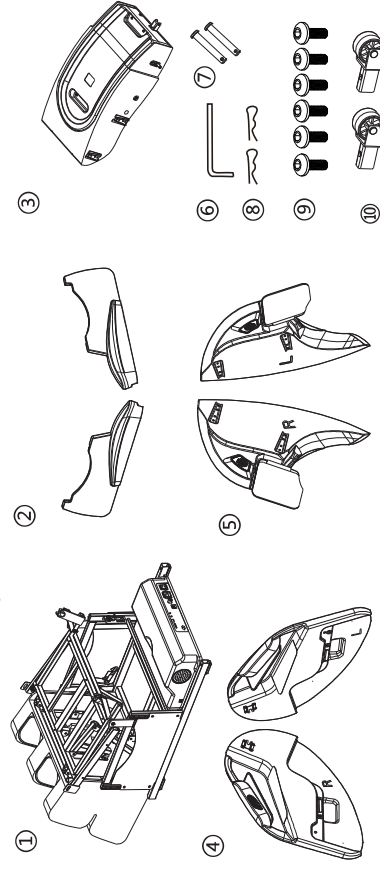


- ⑪ Main unit ⑫ Pedestal(L/R) ⑬ Backrest ⑭ Armrest(L/R)
 ⑮ Shoulder(L/R) ⑯ 5mm-Allen driver ⑰ pins×2 ⑱ R type pin
 ⑲ Screw×4 ⑳ Whells×2 ㉑ Chair pad ㉒ Seat cushion
 ㉓ Power line ㉔ controller

2. Install the wheel;

- Fix the wheel into the square pipe ;
- Use the Allen Key to fasten the screw


3. Fix the pedestal;

- Insert A to B to the correct place ;
- Use the Allen Key to fasten the screw.

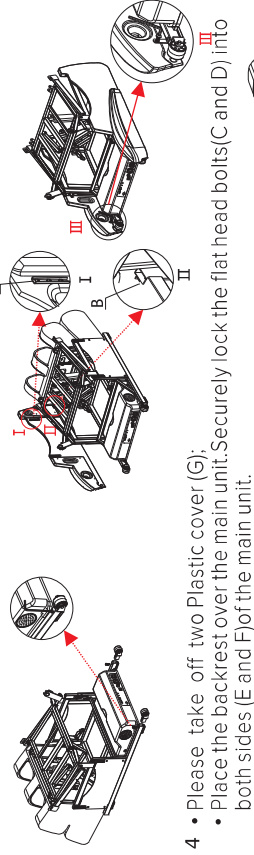


4. Please take off two Plastic cover (G);


- Place the backrest over the main unit. Securely lock the flat head bolts(C and D) into both sides (E and F) of the main unit.

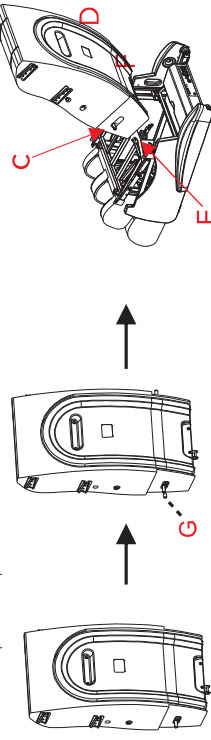

5. Connect the rod flat head rivet(H) with the U stay fork(I), and secure them with a slit pin (J)

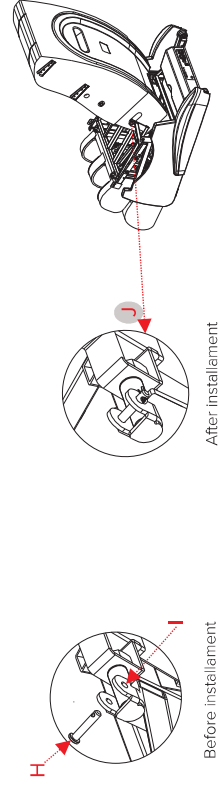


6. Press down the backrest of the main unit, K and L will be locked when a reminding a sound "Ka" is heard.

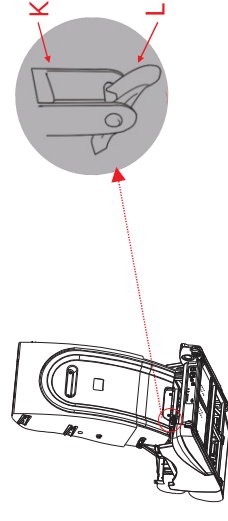


7. Connect the wires on bottom of chair:

- Connect the jack cable from main unit with the female cable from backrest,
- Firmly insert the air plug of main unit into the air pipe of armrest.

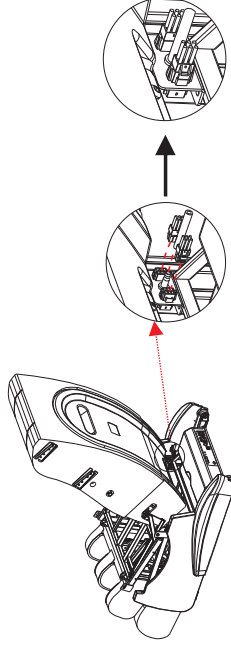

8. Connection of Shoulder:

- Connect the jack cable from Shoulder with the female cable from backrest,
- Firmly insert the air plug of Shoulder into the air pipe of Backrest.
- Insert the shoulder along the hook M the holder N
- Repeat the step to install other sider shoulder unit

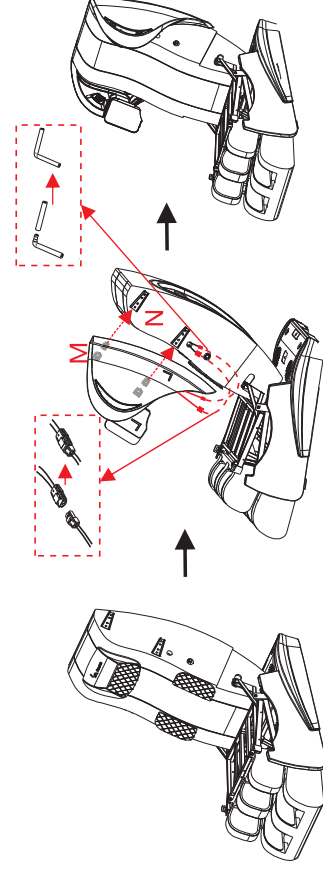

FCC Caution: Any changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate this equipment.

This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions:

- (1) This device may not cause harmful interference, and
- (2) this device must accept any interference received, including interference that may cause undesired operation.

This device and its antenna(s) must not be co-located or operating in conjunction with any other antenna or transmitter.

NOTE:

This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules.

These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications.

However, there is no guarantee that interference will not occur in a particular installation.

If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.