


Valid warranty 1. The warranty card must be engraved or stamped by our sales staff to be valid 2. The user must complete the customer information in order to provide you with effective after-sales service. 3. It is necessary to carry out the warranty card, shopping invoice and product to maintenance. Please check the contents of warranty card when purchasing, whether it is true and the product sequence is complete. If you have any questions, please call our after-sales service hotline.

Thank you for choosing our products Please read this manual carefully before use. and keep it for future reference.

2-2 外面

Product Features After sales and service Storage When it is stored for more than 3 months, it should be unpacked and inspected before recharging the machine Auto mode: Airpressure + Hot pack + Vibration + Music before use. It should not be close to the ground, wall and roof. Ventilation in the air should be good to avoid strong sunlight and corrosive gas erosion. One-click control mode cycle: Transport The product is equipped with simple anti-vibration Air pressure + Hot pack + Vibration + Music facilities, suitable for aviation, railway, road and ship transportation. Rain, snow splash, upside-down and collision should be avoided. Air pressure + Vibration + Music Check the power supply, if its voltage exceeds the Air pressure + Hot pack + Music specified range of the instrument, do not use it. Air pressure + Music Maintenance Take anti-moisture, anti-rat, anti-fouling measures. When cleaning, the device should be turned off. After use, the instrument should be placed in a safe place. Do not place heavy objects on it. When there is a fault, please call for professional personnel for maintenance, non-professionals, do not disassemble the instrument by yourself.

3-3 外面

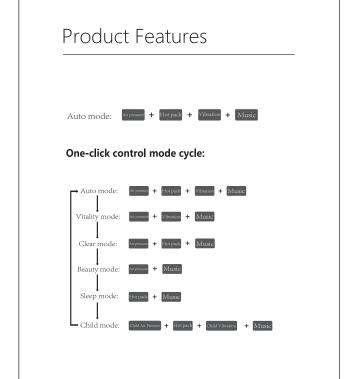

Warning

Before using the product, please be sure to read the operating instructions carefully (User Manual). The company is not responsible for the consequences of incorrect use or failure to follow the corresponding instructions of this product.

- 1. Do not use this massager while driving.
- 2. Do not use it in excessively high temperature and excessively wet places.
- 3. For those who have symptoms of illness (functioning body parts), please use it with caution.
- 4. Do not use the massager under (operating body parts)
- or having a history of surgery.
- 5. If you are sick, please consult a professional physician and use it with your loved ones.
- 6. Do not use it for the person who are with stroke,
- heart disease, hypertension or hypoglycemia. 7. Avoid using after drinking.
- 8. When using this product, do not wear other medical care, beauty, health care items and accessories at the body parts to avoid conflicts resulting in bodily harm
- 9. Please use a power adapter that conforms to national standards, and pay attention to electricity safety.

Product Function instruction 6 7 8 1. Switch / mode (long press for two seconds to turn on / off, short press for mode switching) 2. Volume control key (default medium volume at startup) 3. Air pressure mode 4. Hot pack mode 5. Vibration mode 6. Bluetooth indicator 7. Child air pressure mode 8. Charging indicator

1-1 外面

Warranty

Warranty Conditions

- 1 Within 1 year from the date of purchase, free parts replacement and repair services can be in accordance with a valid warranty card or invoice. (The replacement of product shell leather sheath and Leather belts are not covered by the warranty. The company promises free repairs including maintenance cost and material costs, excluding transportation costs, and the company does not provide on-site repair services). The company is willing to be responsible for recycling, processing, and reuse of the product after it is discarded.
- 2 Any damage or malfunction caused by improper use of the customer, such as catching the wrong voltage, causing an accident, or repairing the product without the maintenance personnel authorized by the company, is not covered by
- 3 The company will charge a corresponding fee for any additional services or transportation services beyond the standard maintenance scope.
- 4 Damage caused by improper use, storage and maintenance of customers; damage caused by disassembled action by maintenance personnel not authorized by the company: no product warranty card or invoice; damage caused by force majeure, in view of the above situations, our company can provide with repair service, but the corresponding accessories and transportation costs will be charged.

Note: Due to the continuous upgrading of services, the company reserves the right to modify the rules without notice

Table of Contents

Product Features

Product configuration

Product function

Use operation

Bluetooth connection and music control

Warning

Fail-safe

After sales and warranty

warranty card

Troubleshooting

If a fault occurs during operation, please perform the following checks. If the fault cannot be eliminated or the product still fails to operate normally, please disconnect the power supply, and contact the company's after-sales center directly to apply for repair service.

Do not disassemble or attempt repairs by yourself

Basic failure:		Basic troubleshooting method:
unable to boot	Is the battery low?	If yes, please charge the machine

Technical Parameters

Product Name:	Eye Massager	
Item No.:	EM-1901	
Power:	5W	
Power Adapter:	DC5V === 1A	
Size:	about 220mmX80mm(L*H)	
N.W.:	about 300g	

Product configuration: Your massager includes: host, user manual, certificate Certificate USB data line User manual

10. This product is a home health product, and is not suitable for diagnostic and therapeutic purposes. 11. This product has a heating function, please pay attention when using it insensitive to heat.

- 12. If you feel unwell during use, please stop using it immediately.
- 13. Do not use it in a humid environment such as a bathroom. 14. Children, pregnant women, menstrual women are
- prohibited from using it. 15. If the massage site has injuries, redness, inflammation,
- or physical abnormality, it is forbidden to use it.
- 16. For cleaning, please wipe it with a soft cleaning cloth. Do not allow the machine to get wet.

Use operation

Preparation before use

Use the power adapter configured by our company to fully charge the product

Operation method

Keep pressing your fingers on " (10) " key for 2 seconds to turn on the product. After power on, the indicator light is always on and the host starts to work. Working time is 10 minutes, playing music.

Switch mode and mode order

Press " (response than one second, being switch modes. Press keys and the modes will cycle in sequence: Auto mode, vitality mode, clear mode, beautiful pupil mode, sleep mode, and child mode.

When the product indicates that the battery is low, please use the charger supplied with the product to charge it.

Note: It is not recommended to use the product when charging.

Bluetooth connection and music control

1. The default is to turn on the music and wait for the Bluetooth connection at startup. Open the Bluetooth to search the name

EM-1901 to indicate that the link is successful.

2. Shortly press the music control button to adjust the volume, or you can adjust the music volume through your phone after

connecting to Bluetooth

This device complies with part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) This device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

Any Changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

Note: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- -Reorient or relocate the receiving antenna.
- -Increase the separation between the equipment and receiver.
- -Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- -Consult the dealer or an experienced radio/TV technician for help.