Fairywill

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V2

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2022/07/14

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
Bitvae

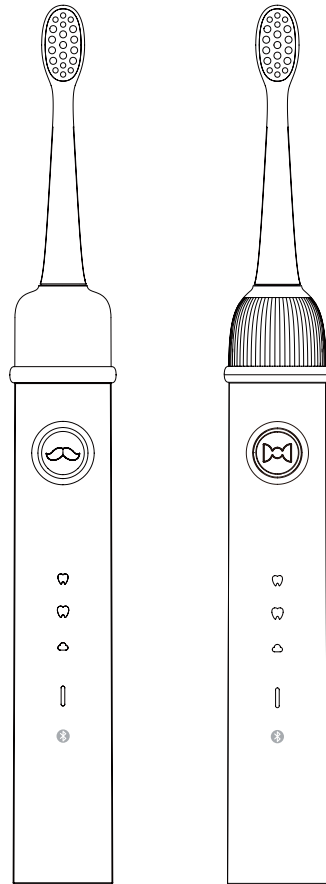Kiddo K7S

Sonic Electric Toothbrush
Instruction Manual

Thank you for choosing the K7S sonic toothbrush. Please read the instructions carefully before using the product, and keep the user manual in a safe place for future reference.

WARNING:

Any Changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

This device complies with part 15 of the FCC Rules. Operation is subject to the following two conditions:

(1) This device may not cause harmful interference.

(2) this device must accept any interference received, including interference that may cause undesired operation.

Note: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

-Reorient or relocate the receiving antenna.

-Increase the separation between the equipment and receiver.

-Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.

-Consult the dealer or an experienced radio/TV technician for help. The device has been evaluated to meet general RF exposure requirement. The device can be used in portable exposure condition without restriction.

This device contains licence-exempt transmitter(s)/receiver(s) that comply with Innovation, Science and Economic Development Canada's licence-exempt RSS(s). Operation is subject to the following two conditions:

(1) This device may not cause interference.

(2) This device must accept any interference, including interference that may cause undesired operation of the device.

This transmitter must not be co-located or operating in conjunction with any other antenna or transmitter.

The device has been evaluated to meet general RF exposure requirement. The device can be used in portable exposure condition without restriction.

1.Features:

- A Waterproof and non-slip design.
- Replaceable brush heads with unique crescent shape bristles design to help keep your teeth whiter and healthier.
- Audio indicator encourages kids to complete their brushing task every time.
- With the pressure indicator, users are able to know when the pressure is too hard to their teeth.
- Over-voltage protection, when the input voltage is over 6.1V, the toothbrush will be auto-off to protect itself and users.
- The Bluetooth technology is built into the toothbrush for application connection, users can learn more about their teeth and gum conditions through the toothbrush connecting to the application.
- The pushed-weekly and monthly generated brushing reports will help parents to learn their children's brushing situation.
- USB charging, travel friendly and compatible with global voltage.

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
☎ 972-679-9068

MADE IN CHINA

- A low-power reminding function, when the voltage of the battery is low on power, the charging indicator light will keep flashing quickly in red or the toothbrush is unable to turn on during the operating process, as a reminder to charge the toothbrush.
- A 2 minute smart timer with interval pausing, that gives a reminder to change the brushing location every 30 seconds.
- There are 3 modes: Mild, Soft and Gum settings built into the toothbrush. There is also a mode memory function built into the brush, that means when you turn on the toothbrush, the default starting mode stays the same as the mode that was set at the time the toothbrush was last shutoff.

2.Designs: (Ribbon for Girls, Mustache for Boys)

- ① Brush Head
- ② Pressure Indicator
- ③ Power/Mode Button
- ④ Mode Indicator (Mild, Soft, Gum)
- ⑤ Charging Indicator / Bluetooth Indicator
- ⑥ Type-C Charging Port

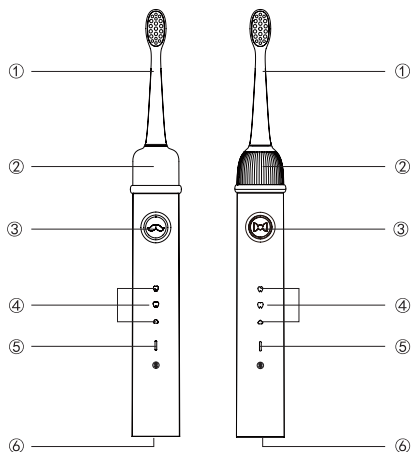3.Operations:

☆ Using the toothbrush

- Slide the toothbrush head down to the handle and make sure it is locked in the correct position.
- This product can be used by connecting the application. Scan the QR code below to download and install the application. And add the devices to the phone according to the hints.

- Turn on the toothbrush by pressing the "Power" button.
- The Bluetooth indicator will light up in White if the device and toothbrush is in connection, please press the "Power" button for 3 seconds to connect to the application that download in phone.
- Shift between the 3 modes ("Mild, Soft and Gum") by continuously pressing the "Power" button.

- Wet the bristles and apply an appropriate amount of toothpaste onto the brush.
- Place the bristles with toothpaste against the teeth at a slight angle towards the gum line.
- Gently move the brush head across the teeth on the front, back, and all chewing surfaces. The brush should be allowed to glide across the teeth with minimal pressure. Do not apply too much pressure to the teeth/gums, or apply the bristles to one area of the mouth for an extended period of time.
- If the pressure is too hard against users' teeth, the pressure indicator light will turn on in Red to remind the users not to force the toothbrush against the teeth too hard. Otherwise, their teeth and gums will be hurt.
- The brush has a 2 minute smart timer as this is recommended by dentists as the optimal time for brushing the teeth and there is interval pausing after every 30 seconds as a reminder to move to another quadrant of the mouth.
- To avoid splashing toothpaste, only turn on the toothbrush after you have placed it against your teeth and only remove it from the mouth after it is turned off.
- The APP guides children to brush their teeth through the audio and real-time cartoon characters, to stimulate children's interest and awareness of brushing to protect gum health.

☆ Changing the toothbrush head

- Install: Align the bottom of the brush head with the driveshaft, push it in firmly until the handle is locked in the position.
- Detach: Firmly hold the handle then, pull the toothbrush head out away from the handle.

☆ Charging the toothbrush

- If the toothbrush is in low battery condition when turning it on for the first time, please ensure it is fully charged before usage.
- If the toothbrush is not in use for more than 3 months, please fully charge it according to the charging instructions before using it.
- If the charging indicator light flashes in red quickly during the operation process, it indicates the toothbrush is in a low battery condition, please charge the toothbrush as soon as possible.
- When you are going to charge your toothbrush, please insert the USB cable into the charging port in the bottom of the toothbrush and ensure the brush is connected to the electric power safely.
- The charging indicator stays on in red during the charging process. When the battery is fully charged, the charging indicator stays on in white for 2 minutes, after that, it will turn off automatically. After the toothbrush is fully charged, it is possible to use it in the "Mild" mode for two minutes per session, twice a day for 60 days.
- There is a waterproof function built into the toothbrush, so it can be used in the bathroom safely.

☆ Modes Introduction

• The first mode: Mild

In this mode, the white indicator light turns on in the "Mild" mode position, the brush swing is mild and the vibration intensity is the softest of the three modes, which is suitable for kids over 6 years old with sensitive teeth or use electric toothbrush for the first time.

• The second mode: Soft

In this mode, the white indicator light turns on in the "Soft" mode position, the swing of the brush is larger and the vibration intensity is the strongest. It is suitable for kids who are 9 to 12 years old and used to electric toothbrushes, to gain the strongest cleaning effect and maintain healthy oral hygiene.

• The third mode: Gum

In this mode, the white indicator light turns on in the "Gum" mode position, the bristle swing and the vibration intensity changes alternatively between the "Mild" and the "Soft" modes. This is for gentle stimulation, designed to improve blood circulation in gums and keep gums firm and healthy.

4.Specifications:

Working temperature: 0-50°C
Working voltage: DC3.7V
Bristles material: Nylon
Input voltage: DC5.0V
Rated power: 2W
Waterproof rating: IPX7
Type of Charging: Type-C Direct Current

5.Precautions:

- For the first usage, it is recommended to start the toothbrush at the "Mild" mode position. If users experience gum discomfort or gum bleeding for more than 2 weeks, stop using the product at once and consult your dentist for help.
- Do not put the toothbrush into a fire or charge it in a high-temperature place.
- Do not place the product in a humid environment or other liquids for a long time to avoid functional defects of the product.
- Do not disassemble, repair, or alter the product or apply other accessories or power supplies without authorization.
- Please avoid the product from falling or strong impacting.
- This product is only for oral cleaning, do not use it for other purposes.
- This product is only designed for personal oral care, not suitable for multi-users.
- Do not use this product while it is being charged for safety purposes.
- Please ensure the charging port or the power adapter are dry enough before charging the toothbrush or connecting it to the power supply.
- Check whether the power cord is damaged regularly. If the cord is damaged or the toothbrush is failed to operate, please stop using it.
- In case of battery leakage, please stop using the product immediately to avoid accidental contact with the skin.
- Do not place the product in the mouth or brush your teeth in a supine position. Otherwise, it may easily lead to accidental swallowing.
- This product should not be played with or disassembled by children.
- This product is not suitable for children under 3 years old. Children between 3 and 8 years old need parental or adult supervision and guidance when brushing their teeth.
- Kids who have severe periodontal disease or oral surgery should consult a doctor before use.
- Users who have problems with pacemakers or other implanted devices should consult a medical professional before using this product.

6.Environmental Protection:

There is a rechargeable lithium battery built in the product. When the toothbrush is being scrapped, please do not mix it with general household garbage. It should be sent to the battery recycling center for disposal in accordance with the local disposal management regulations.

7.After-sales service:

- A 1-year warranty starts on the day of purchase, please contact us if it is necessary.
- Within the warranty period, any defect arising from materials or a product malfunctioning, we will replace the parts or the whole product depending on the situation and its condition.
- When you apply for a warranty based repair, you must submit the toothbrush and proof of the date of purchase (purchasing invoice) to us.
- The following conditions are not included in the warranty:
 - Damages caused by improper use, for example, improper use of the power supply, plugged into improper socket, artificial damage, etc.
 - Ageing and wearing-out caused by normal use.
 - Damages due to product being used for other purposes, or use of components not supplied by us.

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