

The Overview

1. Precautions for use
2. How to use device
3. How to use app
4. FAQ
5. Specification

1. Precautions


Precautions for use

- ▶ Close the lid of device firmly and keep it.
(Preventing dead battery)
- ▶ Do not disassemble device by force.
- ▶ Do not shake the device when you take pills.
- ▶ If the pills are not counted, don't force yourself to put the pills in the bottle.
- ▶ Don't put the physical impact on device.
(The cause of the device failure)
- ▶ Do not drop or throw the device.
(Risk of damage when falling)
- ▶ Don't use next to flammable substances. .
(Combustible substances such as gasoline and thinner, candlelight, fireplace)
- ▶ Do not put it in water or wash it with water.
- ▶ Use a wet tissue or dry cloth for cleaning.
- ▶ Check the remaining battery level and contact the customer service center.

2. How to use device




In place with a lot of Bluetooth devices, there can be errors in connecting to Bluetooth. You need to separate it from other wireless devices.

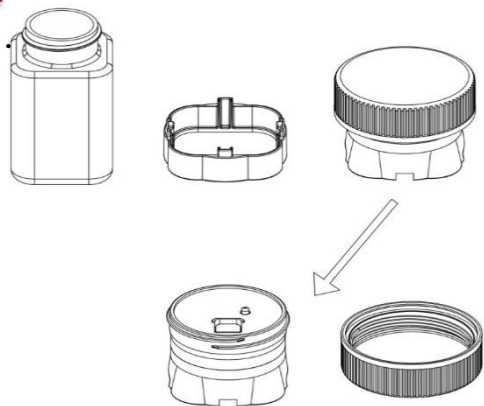


1

Prepare the device in a clean place and open the lid.

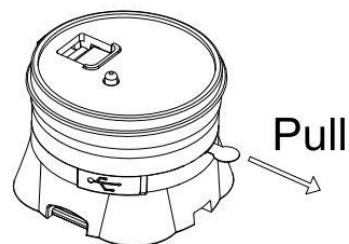

2

Remove the battery contact prevention film from the device.

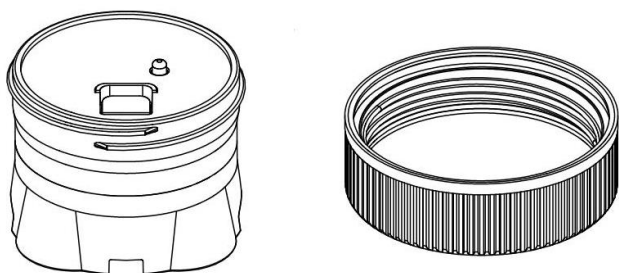


3

Open the lid of device and get ready to connect with the app.

Connect the device with app.


(If Bluetooth number is not recognized, shake the device once and try again.)

4

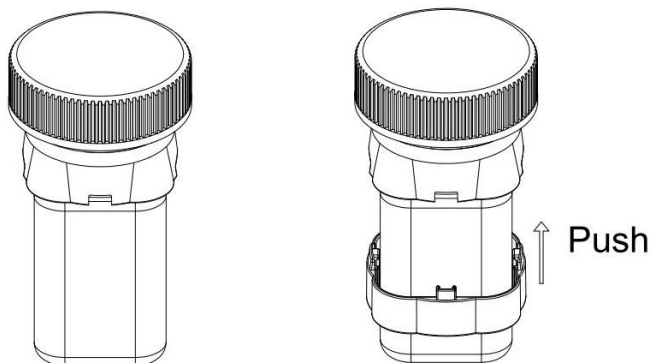

5

1. Combine device and medicine bottle.
2. Fit the safety lock to the protruding shape and pull up.

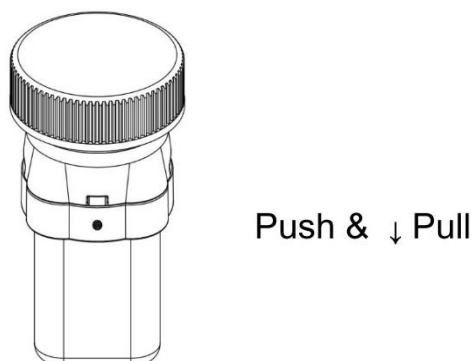


6

Press marked both sides of the safety lock and put it down..

- USB Port is not used for test and does not support the charging mode of the AC Adapter, Only for administrator.

2. How to use device

► How to take out pills

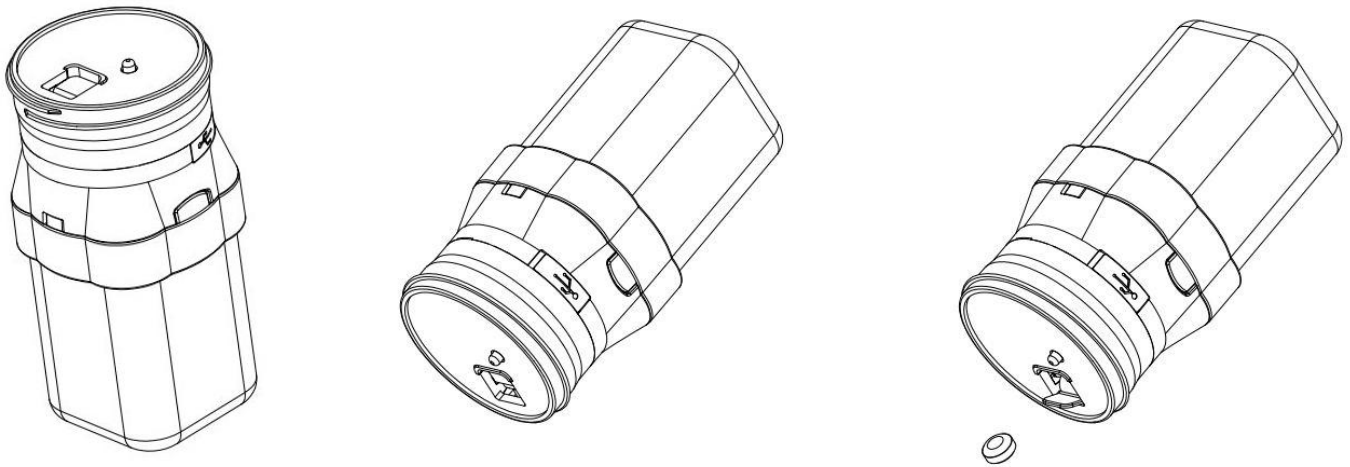

7

Open the lid, turn on the app, set the direction of the bottle to the bottom and tilt it to dispense the pills. If the pills are not dispensed, try again.

► You can check if the device is working :

- 1) If you open the lid, the green light will flash once.
- 2) The Bluetooth number will be searched if you scan the device on the app.

3. How to use App

1. Press registration button and enter personal information including email and password.
2. Log in with the account you signed up.
3. Allow app to access your location.
4. Exclude battery optimization to your app, battery, optimization could delay your medication reminders.
5. Press the setting button to scan device to connect.
6. Click the scanned device, confirm the connectivity after paring the device and app.

4. FAQ

Q1. The app doesn't detect the device operation.

- ▶ If the green LED light does not flash when the lid is opened, it means battery went out. Contact the customer service center.
- ▶ If the green LED light flashes, check the connection status of app and device on the app.

Q2. How do I install the app on my smartphone?

- ▶ Search for Coledy app from the app store, download and install it.

Q3. Is there a problem if the device and smartphone are far away?

- ▶ The device is a Bluetooth wireless device, use it in the same space next to your smartphone.

Q4. Is it possible to exchange battery?

- ▶ Battery replacement is not possible and it can only be used until the battery runs out.

Q5. What if I didn't take the medicine after taking it out?

- ▶ Keep the pills you took out in a sanitary bag and take them.

5. Specification

Smartphone Compatibility

- ▶ Compatible from Android 6.0 version
- ▶ Compatible from iOS 13 version
- ▶ Compatible from Bluetooth 4.2(BLE)

This device complies with part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) This device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

Any changes or modifications (including the antennas) to this device that are not expressly approved by the manufacturer may void the user's authority to operate the equipment.

Note : This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

Modifications not expressly approved by the manufacturer could void your authority to operate the equipment under FCC rules.

RF Exposure Statement

FCC RF Radiation Exposure Statement: This equipment complies with FCC RF Radiation exposure limits set forth for an uncontrolled environment. This device and its antenna must not be co-located or operating in conjunction with any other antenna or transmitter.