



Early-Alert Bed Alarm System
Chair Pad (Optional)
Call Button (Optional)
INSTRUCTION MANUAL

IMPORTANT: Read this manual before use

If you have any issue or question, contact us directly at questions@lundberg.com

Important + customer service

This manual explains how to use the Early-Alert Bed Alarm System and optional Lunderg accessories such as the chair pad and call button. These accessories can be bought separately, as an add-on to your order or as a set with your order on Lunderg.com or at Amazon.com/Lunderg

It is very important that you read all the following pages.

Before you unlock the peace of mind that the Early-Alert Bed Alarm System provides, please make sure all of the parts are included in the box.

If anything is missing, please contact us at questions@lunderg.com and we will make things right. We pride ourselves on our 24 hour customer service and work hard to provide all our customers with the VIP treatment they deserve.

CAUTION: Technology alone cannot prevent falls or injury. In addition to professional caregivers this device is meant to supplement a facility or home patient management system.

If you have any questions at all please contact us at questions@lunderg.com or request a call back at (800) 652-4905.

**REGISTER YOUR WARRANTY WITHIN 30 DAYS OF
RECEIVING THE PRODUCT AT: www.Lunderg.com/warranty**

**REGISTRA TU GARANTÍA NO MAS TARDE DE 30 DIAS
DE RECIBIR EL PRODUCTO en www.Lunderg.com/garantia**

Index

What's in the box	4
Video Guides	5
Identifying the parts of the set	6
Optional Lunderg accessories	6
The Early-Alert bed sensor pad (The Bed Pad)	7
The Early-Alert System	8
Selecting the Pad Settings	9
Where to place the bed pad & how to secure the pad	10-11
The Monitor (The Pager)	12-13
Battery Cover & Batteries	14
How to pair the bed pad with the monitor	15
How to pair the optional chair pad with the monitor	16
How to pair the optional call button with the monitor	17
How to clear the memory	18
Adjusting the Volume	19
LED Light status	19
Tips from the Caregiver	20-21
Care & Maintenance instructions	22
Storage	22
Troubleshooting	23-24
Help + Customer Service	25

What's in the box

- Early-Alert bed sensor pad.
- Monitor with unscrewed battery cover (comes in 2 parts).
- Mini screwdriver (to screw the monitor's battery cover in place).
- 4 x rubber bands and 4 x clips (to secure the bed pad).
- 2 x AA batteries.
- Wall holder for monitor unit.
- 4 x Non-Slip Stickers.
- User manual (this manual you are reading now)

IMPORTANT: If you purchased the Early-Alert bed alarm system together with the chair pad and/or call button as a set or as an add-on to your order, the box will also include a **Chair sensor pad and/or Call button**

Video Guides

For a more dynamic introduction to your products,
watch these video guides

Early-Alert Bed Alarm System

www.Lunderg.com/video guide




Chair Alarm System

www.Lunderg.com/chairvideo




Identifying the parts of the set

The 2 main components of the Early-Alert Bed Alarm System are the Early-Alert Bed Sensor Pad (the bed pad for short) and the Monitor (some folks call it the pager).

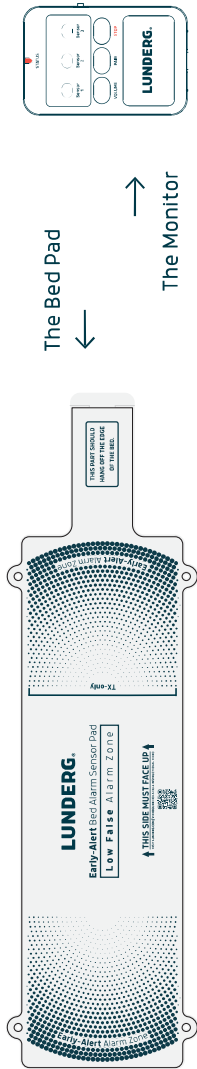


Optional Lunderg Accessories

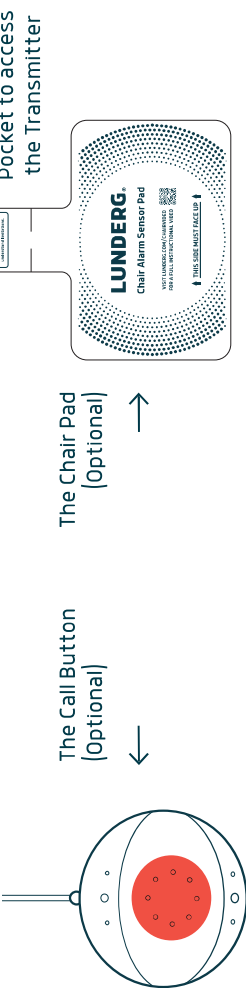

These optional add-ons can be purchased separately at [Lunderg.com](https://lunderg.com) or [Amazon.com](https://amazon.com/lunderg) or [Lunderg.com](https://lunderg.com).

The Early-Alert bed sensor pad (The Bed Pad)

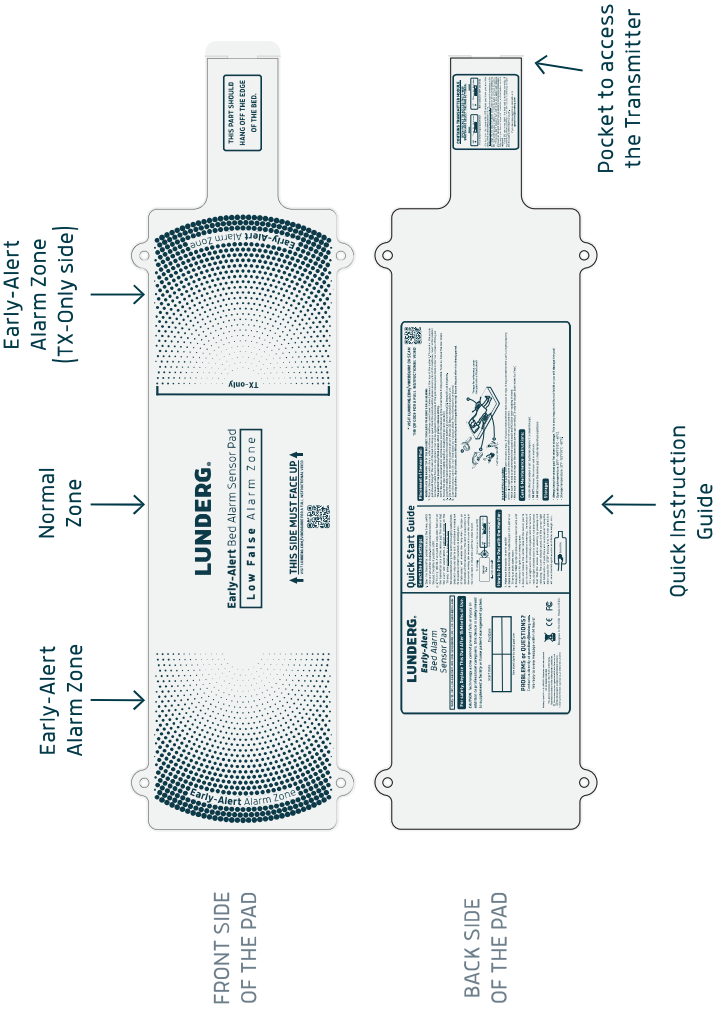



The Early-Alert system


Unlike all the other bed sensor pads in the market, we have developed an innovative Early-Alert system which reduces the false negatives and is able to detect when a patient is about to leave the bed. Traditional systems only detect when there is no pressure directly applied to the bed sensor pad, so it is possible for a patient to move towards the end of the bed pad and even put their feet on the floor, but as they still have their butts on the pad the alarm won't sound until they have totally stood up. That sounds very dangerous, doesn't it?



Good thing that was the past! Our bed sensor pad has two different sensor zones: The Low False Alarm Zone, located at the center of the pad, and the Early-Alert zones, located at the sides of the pad. For those using the pad in a queen sized, king sized or beds with only one side to exit, or even for persons who are restless sleepers that move a lot, it is possible to program it so only the exit side (where the transmitter is) will trigger the alarm using the Early-Alert function, while the other side, remains in a normal / traditional setting.

Selecting the Pad Settings

1. Flip the bed pad to the instructions side and open the transmitter pocket to access the 3-way switch on the transmitter module. It is placed at the neck/end of the pad. The default position is OFF. We do that to maximize the battery life of the transmitter!
2. a) "Tx-Only" will activate the Early-Alert function only on the transmitter side of the bed sensor pad. This means that the alarm will go off when a patient is sitting on the "TX" side edge of the sensor pad and about to leave their bed. The other side remains in the normal setting (which means the alarm will sound when the patient leaves their bed OR there is no longer pressure on the bed pad). We recommend using this setting in large (double) beds, beds with only one side to exit, when there are many low false alarms, or the patient is moving a lot.
b) "Both-sides" will enable the Early-Alert function on both sides of the bed sensor pad. We recommend using this setting in twin-size beds or small beds with two sides to exit.
c) "OFF" is how we ship it so we don't shorten the battery life as its batteries can't be replaced. The only time you should put it in this setting is when you are not going to use the bed pad for long periods of time and want to save on battery life!

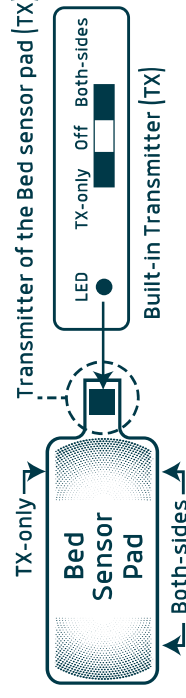

Where to place the bed pad & how to secure the pad

IMPORTANT! NEVER PLACE THE PAD ON TOP OF THE BED SHEET. THIS LEADS TO SLIDING & FALSE ALARMS .




We have included four non-slip stickers that can be applied to the backside of the pad. Please place them in the four corners of the pad.

To secure the pad:

1. Position the sensor pad on top of a flat mattress or mattress cover (never below!), beneath the hips of the patient. The neck of the sensor pad should hang off the side of the bed. In a king or queen bed, the pad should be placed on one side of the mattress, never in the middle!
2. Secure the clips to the rubber bands as shown in the picture to the right (2)
3. Loop the rubber bands with the clips through the eyelets as shown to the right (3). Repeat for all 4 eyelets.
4. Clip to the mattress or mattress cover to secure the pad (4). Bands should be pulled tight!
5. Cover the sensor pad with a fitted sheet so it stays in place. This will reduce false alarms!

Non-slip stickers, clips, bands and a fitted sheet will prevent the pad from moving!

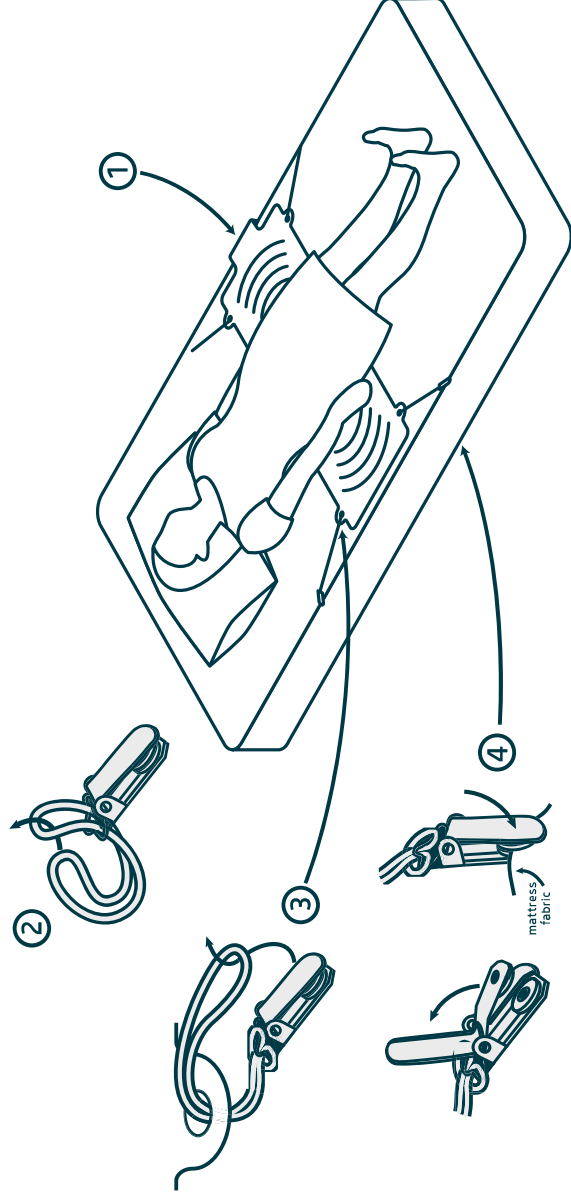


In case of false alarms:

Make sure the pad is beneath the user's hips. If it's placed beneath the lower back or legs, it may not detect the user's weight properly and could trigger an alarm when the patient shifts her weight in the bed.


Make sure the pad is secured and in the correct position EACH DAY when making the bed.

Make sure the pad settings on the transmitter are set up correctly (TX-only for double, Both-sides for Twin).

The Monitor (The Pager)

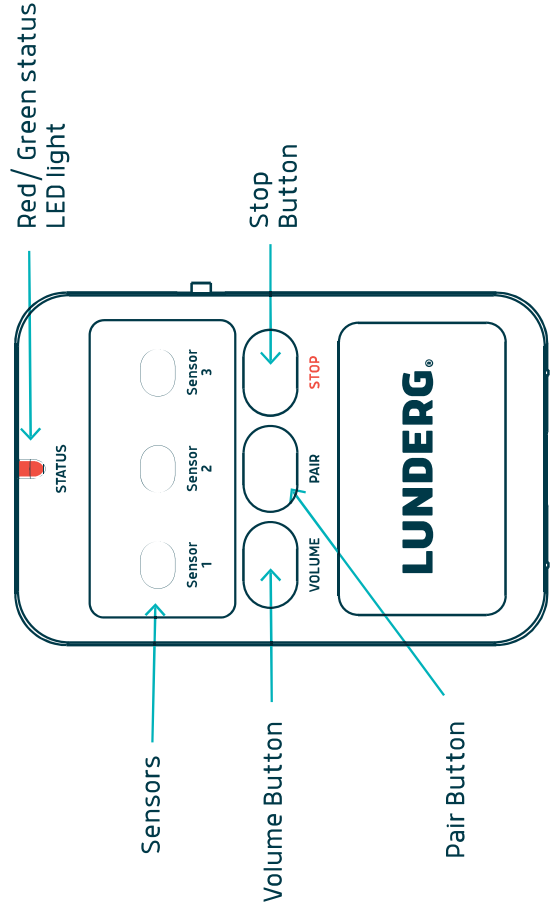

FRONT

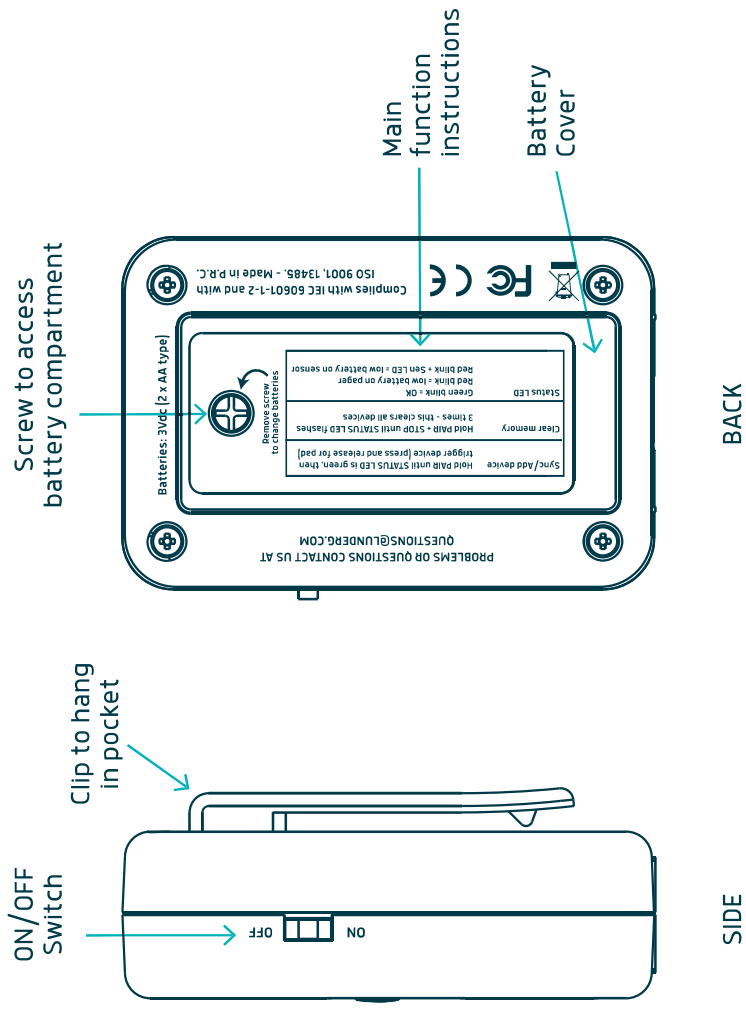

Battery Cover & Batteries

To comply with shipping and warehouse regulations we have to ship the monitor (pager) without the batteries inside of the pager. We have included a small screwdriver and shipped it with the battery cover unscrewed for your convenience!

Simply insert the 2 x AA batteries (the most common ones) in the battery compartment following the + and - symbols. Cover the battery compartment with the battery cover and screw it with the screwdriver provided. The life of the battery will vary based on usage. You are ready to start using your pager!

NOTE: unlike the batteries in the built-in transmitter module of the bed sensor pad and chair sensor pad, which can't be replaced, the 2 x AA batteries of the monitor CAN be replaced.