



LUNDERG™

Early-Alert Bed Alarm System
Chair Pad (Optional)
Call Button (Optional)
INSTRUCTION MANUAL

IMPORTANT: Read this manual before use

If you have any issue or question, contact us directly at questions@lunderg.com

Important + customer service

This manual explains how to use the Early-Alert Bed Alarm System and optional Lunderg accessories such as the chair pad and call button. These accessories can be bought separately, as an add-on to your order or as a set with your order on Lunderg.com or at Amazon.com/Lunderg

It is very important that you read all the following pages.

Before you unlock the peace of mind that the Early-Alert Bed Alarm System provides, please make sure all of the parts are included in the box.

If anything is missing, please contact us at **questions@lunderg.com** and we will make things right. We pride ourselves on our 24 hour customer service and work hard to provide all our customers with the VIP treatment they deserve.

CAUTION: Technology alone cannot prevent falls or injury. In addition to professional caregivers this device is meant to supplement a facility or home patient management system.

If you have any questions at all please contact us at **questions@lunderg.com** or request a call back at (800) 652-4905.

**REGISTER YOUR WARRANTY WITHIN 30 DAYS OF
RECEIVING THE PRODUCT AT: www.Lunderg.com/warranty**

**REGISTRA TU GARANTÍA NO MAS TARDE DE 30 DIAS
DE RECIBIR EL PRODUCTO en www.Lunderg.com/garantia**

Index

What's in the box	4
Video Guides	5
Identifying the parts of the set	6
Optional Lunderg accessories	6
The Early-Alert bed sensor pad (The Bed Pad)	7
The Early-Alert System	8
Selecting the Pad Settings	9
Where to place the bed pad & how to secure the pad	10-11
The Monitor (The Pager)	12-13
Battery Cover & Batteries	14
How to pair the bed pad with the monitor	15
How to pair the optional chair pad with the monitor	16
How to pair the optional call button with the monitor	17
How to clear the memory	18
Adjusting the Volume	19
LED Light status	19
Tips from the Caregiver	20-21
Care & Maintenance instructions	22
Storage	22
Troubleshooting	23-24
Help + Customer Service	25

What's in the box

- Early-Alert bed sensor pad.
- Monitor with unscrewed battery cover (comes in 2 parts).
- Mini screwdriver (to screw the monitor's battery cover in place).
- 4 x rubber bands and 4 x clips (to secure the bed pad).
- 2 x AA batteries.
- Wall holder for monitor unit.
- 4 x Non-Slip Stickers.
- User manual (this manual you are reading now)

IMPORTANT: If you purchased the Early-Alert bed alarm system together with the chair pad and/or call button as a set or as an add-on to your order, the box will also include a **Chair sensor pad and/or Call button**

Video Guides

For a more dynamic introduction to your products,
watch these video guides

Early-Alert Bed Alarm System

www.Lunderg.com/videoguide




Chair Alarm System


www.Lunderg.com/chairvideo

Identifying the parts of the set

The 2 main components of the Early-Alert Bed Alarm System are the Early-Alert Bed Sensor Pad (the bed pad for short) and the Monitor (some folks call it the pager).

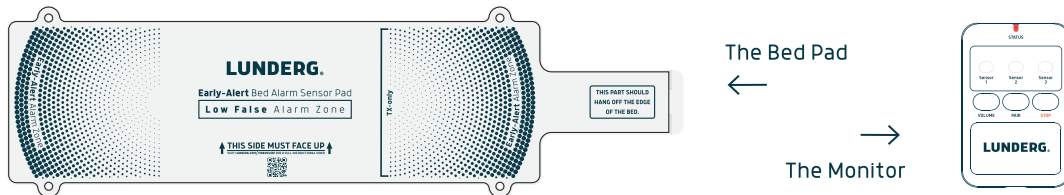


Optional Lunderg Accessories

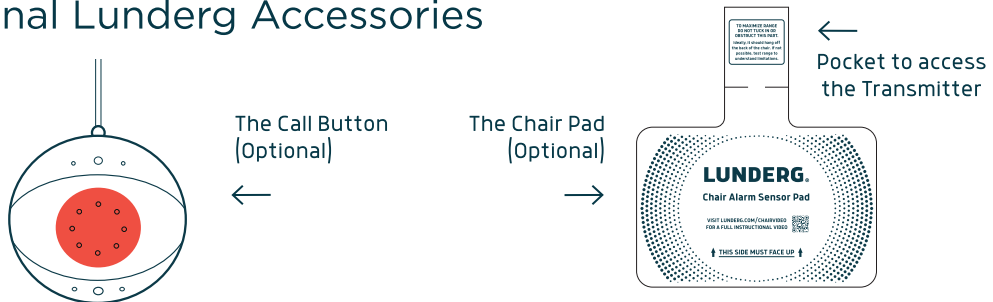

These optional add-ons can be purchased separately at Lunderg.com or Amazon.com/Lunderg.

The Early-Alert bed sensor pad (The Bed Pad)

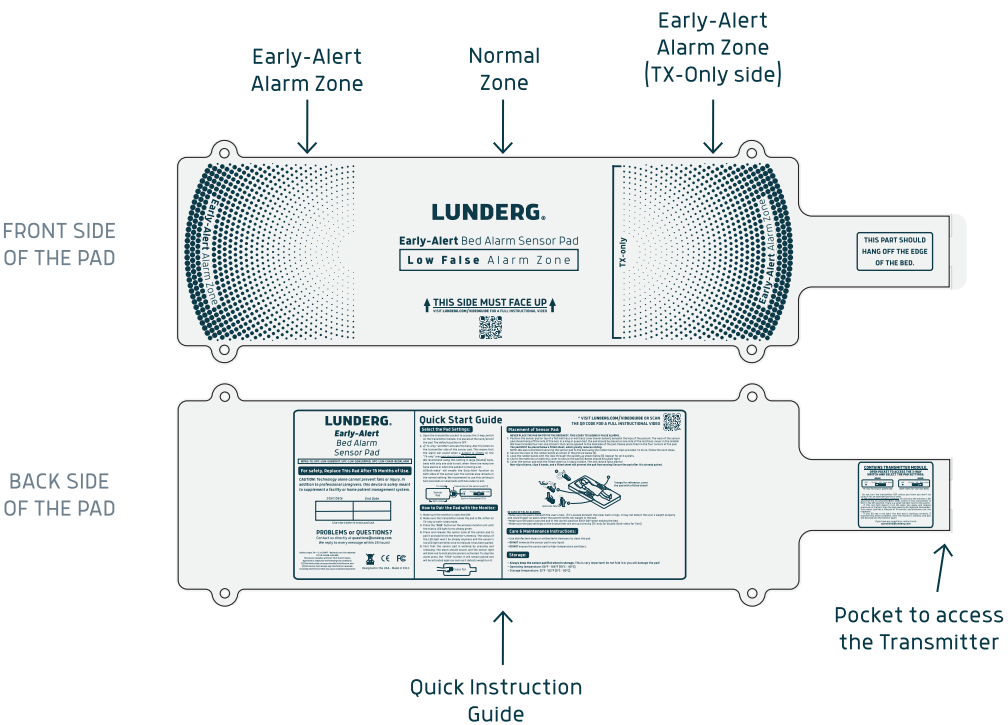

The Early-Alert system

Unlike all the other bed sensor pads in the market, we have developed an innovative Early-Alert system which reduces the false negatives and is able to detect when a patient is about to leave the bed. Traditional systems only detect when there is no pressure directly applied to the bed sensor pad, so it is possible for a patient to move towards the end of the bed pad and even put their feet on the floor, but as they still have their butts on the pad the alarm won't sound until they have totally stood up. That sounds very dangerous, doesn't it?

Good thing that was the past! Our bed sensor pad has two different sensor zones: The Low False Alarm Zone, located at the center of the pad, and the Early-Alert zones, located at the sides of the pad. For those using the pad in a queen sized, king sized or beds with only one side to exit, or even for persons who are restless sleepers that move a lot, it is possible to program it so only the exit side (where the transmitter is) will trigger the alarm using the Early-Alert function, while the other side, remains in a normal / traditional setting.