



Smart eye Massager

User Manual

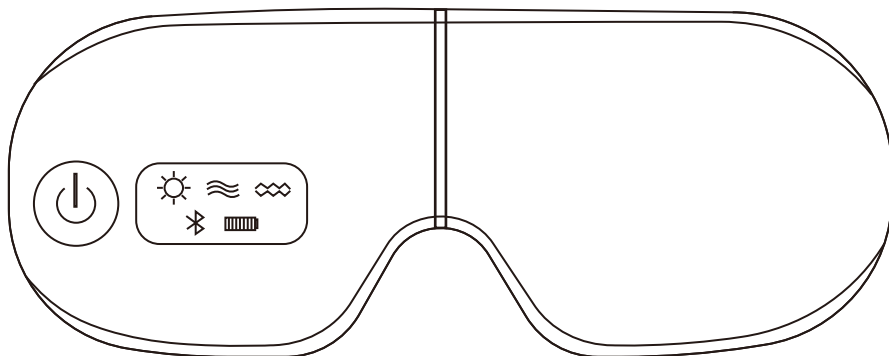


Product figure

Host

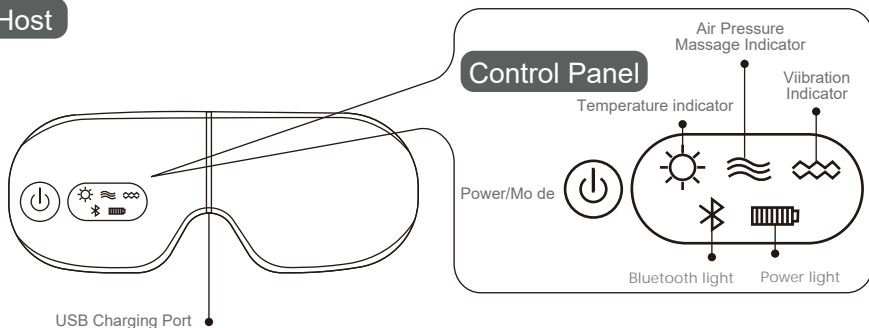

Accessories

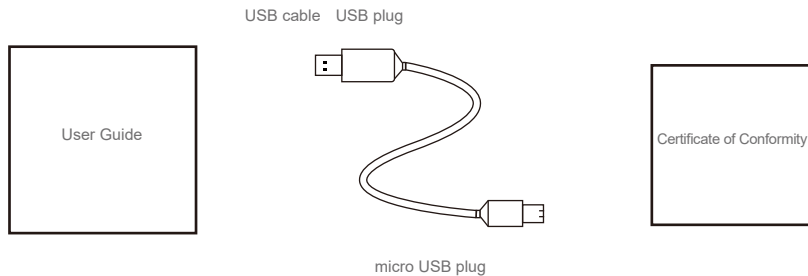

Charging method

- Fold the product first, and use the usb cable to chage
- Please use 5V/1A-2A regular adapter to charge the product.
DO NOT USE ANY ADAPTER WITH OVER 5V
- Please check the indicator light. Please do not use the product when you are using the product. When the product make "Du Du" sound. please charge the machine.
Long time use product with low battery may cause charging problem

Operation

Long term press 3-4 seconds, Turn on/off.

The product has 4 different mode

- When you turn on the machine, it is in Mode 1 (default when turn on the machine)
Mode 1: Low-temperature warm compress + Low-pressure massage+ Low-frequency vibration. One time click the button base on the mode 1, it will turn to mode 1 plus
Mode 1 Plus: High-temperature warm compress + High-pressure massage+ High-frequency vibration
- Double click power button to turn to mode 2
Mode 2: Low-temperature warm compress + Low-pressure massage
One time click the button base on the mode 2, it will turn to mode 2 plus
Mode 2 Plus: High-temperature warm compress + High-pressure massage
- Double click the power button to turn to mode 3
Mode 3: Low-temperature warm compress + Low-frequency vibration
One time click the button base on the mode 3, it will turn to mode 3 plus
Mode 3 Plus: High-temperature warm compress + High-frequency vibration
- Double click the power button to turn to mode 4
Mode 4: Low-temperature warm compress
One time click the button base on the mode 4, it will turn to mode 4 plus
Mode 4 Plus: High-temperature warm compress
- Double Click the power button to turn back to mode 1
- Bluetooth is available for this product. You can link your Bluetooth to listen to music.

Used for

Student, Myopia, Long sighted, Astigmatism, Eye fatigue

People who have fine lines around the eye, dark circles, stay by all night

People who have eyeache, headache, Psychasthenia, Insomnia

Precautions

Using time 15 once is recommended. Overtime use may cause uncomfortable.

Please do not use the product immediately after use eye drops

Please do not use with facial mask together

If you feel your eye over hot, please turn off the product

Please do not use the product when charging

Care Instruction

Do not directly use water to wash it. DO NOT USE alcohol, Sodium Phenyl, or other chemical materials.

USE dry and soft cloth or paper to clean the product

Do not clean it when it charging

Failure solution

Common Problem	Reason	Solution
No charging or all the indicator light, not shiny when charging	The product with low battery due to did not use for a long time	Full charge the product
All the indicator light shiny or There is "Da Da" sound appears.	The product with low battery due to did not use for a long time	Full charge the product
There is a Siren when charging	The product with low battery due to did not use for a long time	Full charge the product
When you feel insufficient temperature when under warm mode	Low battery	Full charge the product
When air pressure mode and warm mode unable to start	Caused by the status of power supply	Full charge the product
Barometric intensity too strong or weak	Adjust the length of the headband	Adjust the headband to the length you need
Murmur sound appear when under air pressure mode	When the airbag start, there would be murmur sound appear	No need to take any action
How to judge charging situation	Indicator light	The indicator light keeps shiny when charging. When the product full charged, the indicator turns off

Specification

Product Name	Eye Massager
Mode	8S
Power Supply	Charger (Built-in lithium battery 3.7V 850mAh)
Charging time	Around 3 hours
Continuous use time	Around 60 min
Rated adaptation time	Around 15min
Temperature in Warm Mode	38℃-40℃
Working Temperature	0℃-40℃
Size	320x80x30mm
Weight	285±5g

Federal Communications Commission (FCC) Statement. This device complies with part 15 of the FCC Rules. Operation is subject to the following two conditions:

- (1) This device may not cause harmful interference, and
- (2) this device must accept any interference received,

including interference that may cause undesired operation. Note: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide Reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

Warning: Changes or modifications made to this device not expressly approved by **Dongguan Shusen Industrial Co., Ltd** may void the FCC authorization to operate this device. Note: The manufacturer is not responsible for any radio or TV interference caused by unauthorized modifications to this equipment. Such modifications could void the user's authority to operate the equipment.

RF exposure statement:

The device complies with RF exposure requirements and can be installed and used without restriction.