


PESSIDO

TCB1 Electric Toothbrush

User Manual



Parameter

Model	TCB1 Electric Toothbrush	Input	5V ---
Gear	High/Medium/Low	Unit weight	93g
Waterproof	IPX7 (Toothbrush Handle)	Waterproof	IPX4 (Charger)
Unit size	236*29.5*29.5mm (including brush)		
Battery	Built-in rechargeable battery (lithium battery)		
Modes	5 brushing modes with 3 intensity settings		
Includes	1handle, 2brush head, 2hygienic cap, 1charger, 1travel case		

※Warning


- Do not share brush heads with others.
- Persons (including young children) who do not have full civil capacity or lack the relevant experience and knowledge must use, clean and maintain the product under the supervision and guidance of a guardian.
- Users of orthodontics, dental restorations, dental implants, etc. should use them under the supervision of a physician.
- This product complies with safety standards for electromagnetic devices and any questions regarding pacemakers or other implantable devices should be directed to a physician or the manufacturer of the implantable device prior to use.
- Use under a doctor's supervision for the following conditions: cardiac patients, seriously ill patients, acutely ill patients, infected patients, hemophiliacs, malignant tumor patients, allergy patients, etc.
- Keep away from sources of fire and heat.
- Do not alter or repair.
- Do not drop the toothbrush hard enough to strongly impact it.
- Unplug the power supply if not used for a long time.

※Kind reminder

Gums may bleed at first use; if bleeding still occurs after 1-2 weeks of use, stop immediately and consult your doctor.

- 1 -

Accessories

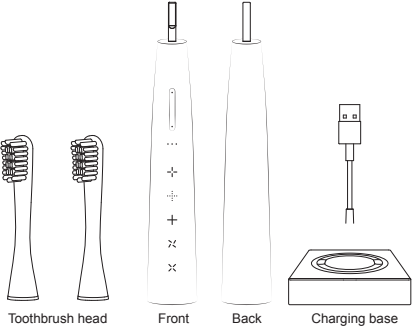

Charging

1. Plug in


Display during charging - gear light display: the 3 gear lights flash in sequence while charging.

- When the 1st gear light flashes, indicating that the charge is not 1/3 full.
- When the 1st gear light is always on and the 2nd gear light is flashing, indicating that the charge is not 2/3 full.
- When the 1st and 2nd gear lights are always on and the 3rd gear light is flashing, indicating that the charge is less than 3/3 full.
- When the 3-stop light is always on, it means the charge is full.

2. Place the charging base on a level table.

3. Toothbrush inserted upright into charging base (as shown).

4. When charging is complete, unplug the power cord.

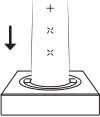

- 2 -

Usage and mode

※ There will be slight discomfort at first use, but you will get used to it after many uses.

1. First attach the brush head to the toothbrush body.

\* Be sure to insert the toothbrush body vertically, do not twist the brush head.

2. Apply toothpaste to the brush head.

3. Press the button "on/off" to start the toothbrush.

4. Press the button "on/off" to switch the mode you want (5 modes).

5. Press the button to switch the gear you want.


\* 3 levels, 1 level is recommended for first use, gradually upgrade the level after getting used to it.

6. Shut down:

① Within 5 seconds of the selected mode, press and hold "Power On" to turn off the toothbrush;

② If the selected mode is more than 5 seconds, short press "Power On" to turn off the toothbrush;

③ The toothbrush will automatically shut down after using for more than 2 minutes.

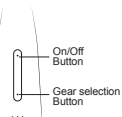

Usage mode:

+

Clean mode

Deep clean, solve the daily residual dirt.

+

Sensitive mode

Gentle touch clean, no burden on sensitive teeth.

+

De-plaque mode

High frequency and thorough cleaning, effectively remove stubborn tartar.

×


Massage mode

Relieve gum discomfort and maintain health.

×

Polishing mode


Polished teeth, radiant.




- 3 -

Recommended Pasteur method of tooth brushing


A smart shift pause reminder is provided every 30s during use.

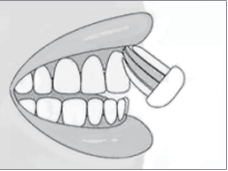


Place the bristles at a 45° angle along the gum line.

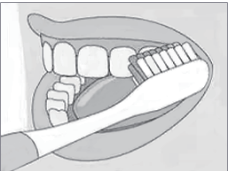

The bristles should touch the surface of the teeth and gums, in groups of 2~3 teeth, vibrate horizontally.

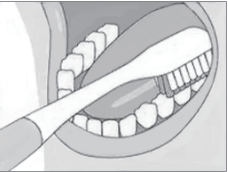

Brush the surface of the teeth, and then move the toothbrush to the next group of 2 to 3 teeth and repeat.

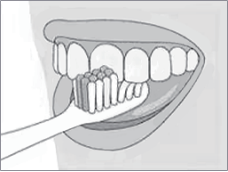


Rest the toothbrush on the occlusal surface of the posterior wisdom teeth and perform gentle back and forth strokes.

- 4 -

Cleaning method

1. Remove the brush head(as shown in the illustration at right) Remove the brush head vertically from the mounting head.

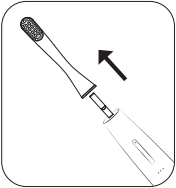


2. Wash the brush head.

Rinse the inside and bottom of your toothbrush well (as shown).

※ Do not pull hard on the bristles.


※ Do not use hot water above 80°C.

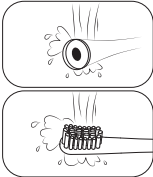

3. Wipe the toothbrush and water from the body and store them away.

※ Do not place the body in water.

※ Do not use hot water above 40°C.



- 5 -

Q&A

Q: The vibration of the toothbrush seems to be getting weaker?

A: The brush head may be connected to the host too tightly.

Please remove the brush head and install it.

Q: The electric toothbrush can't start?

A: Please consider whether the battery is insufficient; if the problem is still not resolved, it may be that the battery has reached the end of its life. Please contact the PESSIDO after-sales for replacement.

Q: I have charged, but the battery life is short?

A: Please make sure whether the toothbrush is placed vertically and steadily during charging; if the problem is still not solved, it may be that the battery has reached the end of its life.

Please contact the PESSIDO after-sales for replacement.

- 6 -

After-sales service

2 Year Warranty

1. The warranty period is two years from the date of purchase. Any quality problems caused by product manufacturing or material defects during normal use belong to the scope of the warranty. he scope of warranty covers the quality problems of parts that affect the function of the product, and does not include the damage caused by the appearance of wear due to long-term use or due to incorrect operation, negligence, or accident.

2. Any damage caused by trying to open or disassemble the product(including its accessories) is not covered by the warranty.

3. Quality assurance claims need to provide reasonable evidence that the product is still within the warranty period. In order to make your product warranty valid, please keep the original proof of purchase properly during the warranty period.

- 7 -

NOTE: This device complies with part 18 of the FCC Rules.

Information to the user.

1. The device has potential interference, but the interference is very small and meets the requirements of the FCC rule

2. Equipment system maintenance is simple, please refer to the manual

3. The user can take simple measures to correct the interference.

Such as staying away from interference sources, turning off interference sources, etc.

\*RF warning for Mobile device:

This equipment complies with FCC radiation exposure limits set forth for an uncontrolled environment. This equipment should be installed and operated with minimum distance 20cm between the radiator & your body.

- 8 -

Warranty Card

Item :	Model :	
Purchasing date :	Purchase Channel :	
Purchaser Information	Name :	Purchase Account :
	Telephone :	
	Address :	
Fault description		

- 9 -