

User Manual

Electronic Moxibustion



This manual is applicable to electronic moxibustion (model / specification: CAJ-9), this is a home electronic heating tools, and it can not be used alone;should be used with other appropriate apparatus (materials) to keep health.

I. Product functional components




On/off key + indicator



Type-C charging port



Heat dissipation hole

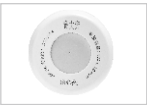

Heat generating piece



Type-C charging cable

II . The Operation Method

1.About Charging

Connect the charging head to the charging cable, then connect charging port of the moxibustion head to start charging. The yellow light turned on which means the moxibustion head is charging. When the light turn to yellow ,it means finished charging.



The charging yellow light goes on



Full yellow lights go out

Note:The power adapters on the market that meet the following requirements can charge the moxibustion head, input: AC100-240V50/60Hz; output: DC5V1A

2. Switching on and off

2.1. Power on

In the off state, press and hold the on/off button for 2 seconds to turn on the power

2.2. Power off

Press and hold the on/off button for 2 seconds to turn off the power, and the working light will be off.

3.Three levels of heating

Press "button" to switch the temperature, there are low, medium and high three temperature levels can be regulated. When the indicator light turn to green and start working, at this time in the low temperature state: press the button again indicator light turn to blue, the product is in the middle temperature state; press again, light turns to yellow, the product is in the high temperature state: and so on cycle. (Note that the indicator color varies slightly depending on the model, the specific object shall prevail)



Indicator light "green"Indicates the current "low temperature"



Indicator light "blue" current "medium temperature"



Indicator "yellow" indicates current "high temperature"

Attention:

a, temperature control shall be based on personal acceptance and comfort, not the blind pursuit of high temperature;
b, if you need accurate temperature regulation and control time, you can connect the product through the APP Bluetooth precision control. (APP only support Chinese)

Warranty Regulations for Follow the Heart Moxibustion Products

Dear Customer:

Thank you for purchasing this product! In order to maintain your legitimate rights and interests, please understand the following in detail:

I. Warranty policy

1) non-human factors lead to the failure of the product, with a valid warranty card can enjoy the following services:
1 from the date of sale of the product (subject to the actual date of receipt) within 7 days, enjoy the return of the same type of product or free repair services;

2) If there is a quality problem, enjoy the replacement of the same type of product or free repair service within 15 days; 3) Within one year from the date of purchase, enjoy the free repair service of Zhongshan Jiaxun Medical Equipment Co.Ltd.

II. The following cases will not be returned or exchanged

1) Damage to the tamper-proof label;
2) the outer packaging of goods, accessories, gifts, instructions are incomplete, missing or altered sales vouchers;
3) the warranty card contains the body number and will be returned to the machine body number does not match;

III. the following cases are not covered by the warranty (can enjoy the manufacturer's paid repair)

1) tamper-proof label is damaged;
2) more than the free warranty period (one year);
3) warranty card contained in the body number and the machine to be warranty body number does not match;
4) damage caused by poor storage, such as falling, folding, into the liquid;
5) damage caused by force majeure (such as natural disasters, earthquakes, lightning strikes);
6) the use of the process, the machine is infected with viruses caused by damage;
7) human factors caused by the failure.

IV. Warranty coverage

Does not include the body shell.

Second, the warranty way

1. warranty, please present a valid warranty slip;
2. products need to be sent to the designated repair point or customer service center.

Note: The final interpretation of the above warranty regulations belong to Zhongshan Jiaxun Medical Equipment Co.

III. Normal working conditions

1. Ambient temperature 5℃~40℃;
2. Relative humidity ≤80%;
3. Atmospheric pressure 86kPa~106kPa;
- 4.Charging input DC5V, 1A;
- 5.Working voltage DC37V±0.5V.

IV.the main performance indicators

Internal power supply:37±0.5V
Battery capacity: 500±50mAh/
heating piece maximum load: 60W
Moxibustion head external size: 22.5mmX45mm
moxibustion head weight: 30g/pc (±05g)
Adjustable time range: within 0 to 70 minutes continuously adjustable

V. Transportation and Storage

- 1.The product should not be heavily pressed during transportation, and it is prohibited to be exposed to sun and rain.
- 2.Pay attention to dustproof and moisture-proof, the product should be kept for a long time without use, and stored in the ambient temperature -5℃ ~ 55℃ relative humidity is not more than 93% no corrosive gases, well-ventilated indoor.
- 3.The product should be stored in a place out of the reach of children, the mentally ill and the consciousness-impaired.
- 4.Special note: The product must be charged after use and then stored, and must be charged once a month to maintain the normal life of the battery when it is rarely used or not used; otherwise, it will lead to battery power depletion, reduce the battery life or even damage such as failure to turn on normally.

VI. Environmental protection

When this product and its accessories exceed the use period or other circumstances need to be scrapped, please do local environmental protection related regulations for disposal.

VII. Special instructions for safe use

- 1.This product is an electronic product, only provide heating use, can not replace drugs and medical devices, prohibited to use alone!
- 2.The body is equipped with pacemakers, stents and metal in the body to use this product with caution.
- 3.The first use of this product, the temperature setting should be adjusted from low to high, and according to different people, season, environmental temperature, etc. timely adjustment, and close observation of the use of the condition adjustment, otherwise the temperature is high prone to burns and other accidents.
- 4.Product work period, prohibit the body, limbs, etc. pressed on the product, the product can not be covered with heavy objects, to prevent contact too tight, resulting in burns.
- 5.After the work of the product, the product should be stored away from flammable and explosive substances to avoid accidents.
- 6.The product can not be charged immediately after work, need to wait for the product to cool down before charging.
7. Any time the product is strictly prohibited to immerse in water, so as not to short-circuit the accident.

VIII.Eight, Emergency handling methods

1. When charging abnormalities occur,please stop powering;
2. Forbidden to disassemble the moxibustion head privately. In the process of using any dysfunction, strictly prohibited Disassemble and repair by yourself to avoid danger. Should promptly contact with our after-sales service department, contact number: 400-966-8187.

Warranty

Product Name	electronic moxibustion		Model / specification		CAJ-9		
Body code			Production date		See body/package		
User Name			User Phone				
User's address	Province	city	district (county)		Zip code		
Purchase date	year	month	day		sales unit name and seal		
Sales address	Province	city (county)					
Sales phone							
Fill in the matter	Please fill in your information correctly and completely. Warranty service will not be implemented if the product model, body number, user name, telephone number, address, date of purchase, sales address or no seal of the seller is not filled in or altered on the warranty form.						

Note:Other contents of product information are detailed in the manual.

Warranty service records					
Delivery date of repair	fault repair	fault cause	treatment measures	delivery date	maintenance signature
Note: Warranty service records can only record the performance of the fault range of maintenance. Testing, software updates, download services, product maintenance, decryption, etc. are not part of the maintenance model.					

FCC Warning Statement

Changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment. This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.