

美国区域
OR989牙刷说明书
材质工艺: 105G铜版纸
成品尺寸: 65*110mm

65mm

Electric toothbrush
OPERAN
User Manual



Model: OR989

110mm

First of all thank you choose our sonic electric toothbrush, gives us the opportunity to provide you with the best quality service, in order to make our service make you more satisfied, after the purchase please read this manual carefully.

Caution:

Check if there was any obvious damage each time before you use the toothbrush. If there is, please do not use it.

If you had a mouth or gum operation in the last two months, please consult your dentist before using this toothbrush.


If your teeth bleed after using this toothbrush please consult your dentist.

The toothbrush is used only for cleaning tooth and gum. Please don't use it for other purpose.

If you feel uncomfortable or even pain in using the toothbrush, stop using it immediately and ask for help from your doctor.

This product is not suitable for children, infants and those unable people, unless they are supervised by person in charge.

01 / Product Description

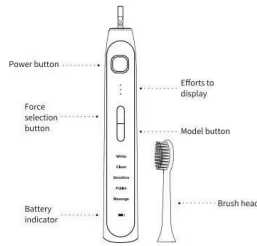

Your dentist recommends that you replace your electric toothbrush head every 2-3 months.



02 / Product specification

product name	Sonic Electric Toothbrush
Vibration mode	9 Modes
Rated voltage	3.7V ^{max}
Rated power	1.2W
Waterproof level	IPX7
Battery capacity	1300mAh

Plug one end of the USB charging cable into the charging adapter or COMPUTER USB port and the other end into the base, and it will take about 6 hours to fill the toothbrush.

Charging status: indicator blinks
Full: The indicator is steady on

Note: To protect the product, it is recommended to use the 5V1A adapter for charging.

03 / The warranty

OPERAN offers an extensive warranty for all products that extends through 1 YEAR of the delivery date. This warranty covers quality related issues and product defects. If you are past the 30-day return period on Amazon, contact us directly via the "Contact Seller" option on your Order Details. If the product is not working properly, we can send a replacement product after we receive the original defective item and our team has a chance to inspect and verify the issue(s). This warranty does not apply in cases of abuse, mishandling, or unauthorized repair of the product.

If you have any questions about this warranty or would like additional information, please contact us via the "Contact Seller" option. Please keep the original box, packing materials and all documentation in the event that service is required.

Wireless charging function introduction:
The principle of wireless charging is that the wireless charging device transmits the energy to the charging terminal device through near-field induction, and the terminal device converts the received energy into electric energy and stores it in the battery of the device.

The principle of energy transmission is inductive coupling, which can ensure that there is no exposed conductive interface. It can not only eliminate the messy entanglement lines between devices, but also be safer for electronic devices such as electric toothbrushes that often contact with conductive media such as liquids.

04 / Battery processing



Do not attempt to replace the battery by yourself, lest the waterproof performance or internal structure of the product be damaged. When discarding this product, please do not put it together with the general life garbage. The battery is not removable and should be given to an officially designated recycling center. The built-in rechargeable battery contains substances that may pollute the environment. When discarding products, please send the products to the officially designated recycling station to help the environment.

Even healthy teeth may start to tingle or bleed slightly when using a toothbrush. This is because the gums are not used to the irritation when first using the toothbrush. The symptoms usually stop within 1-2 weeks. More than two weeks may indicate that your gums are abnormal, in which case you should stop using this toothbrush and consult your dentist.

Email: support@operan.net

FCC ID: 2A152-OR989

FCC Radiation Exposure Statement:
This equipment complies with FCC radiation exposure limits set forth for an uncontrolled environment. This equipment should be installed and operated with minimum distance 20cm between the radiator & your body.

FCC Warning

This device complies with Part 18 of the FCC Rules.

Refer to FCC Part 18.213(c) and KDB 744748 A2 for A section A.6.

(a) Max. Electric Field Emissions 5.37V/m Max. Magnetic Field Emissions 0.43V/m

(b) The maintenance of the system is the responsibility of the manufacturer.

(c) To prevent interference, the product should be as far away from devices operating at the same frequency as possible.