

[Cover]

[Brand LOGO] **Emay Plus**

[产品名称] **Emay Plus ZERO-FAT 塑身仪**
Emay Plus ZERO-FAT 塑身儀
Emay Plus ZERO-FAT Body Slim

[产品型号] **EP-601**

[Icon or Photo]



User Manual

Please read this user manual carefully before use and keep it properly for future reference.

Contents

Safety Precautions.....	
Product Introduction.....	
Product Structure.....	
Product Specification.....	
How to Use.....	
Cleaning, Maintenance and Storage.....	
Causes and Troubleshooting.....	
Product Warranty.....	

Safety Precautions

⚠ Please read and follow the safety precautions carefully before use. Improper operation might pose a risk to personal safety and property damage. Do not use the product for other purposes other than as directed.

【DANGER】

1. Do not use the product on the head, eye area, laryngeal prominence or area of chest and back **area** near to the heart.
2. Stop using the product immediately if you feel any discomfort, and consult your physician for guidance.
3. The product is not intended for use by persons:
 - a) Under 14 years old;
 - b) Pregnant, lactating and menstrual women;
 - c) Have reduced physical, sensory or mental capabilities;
 - d) Feeling unwell or undergoing medical treatment;
 - e) Have or are undergoing medical electronic equipment (eg, cardiac pacemaker, artificial heart-lung or implantable medical devices);

- f) Have or are undergoing wearable medical electronic device (eg, ECG scanner);
 - g) Have chronic diseases, heart disease, cancer and tumor etc;
 - h) Have low immunity caused by immunosuppressive diseases, or taking drugs of immunodepressants;
 - i) Have hypersensitive skin, contact dermatitis or skin disease;
 - j) Light sensitive or taking prescription medicine of isoflavones;
 - k) Have undergone surgical operation, laser minimally invasive operation, cuticle removal on the treatment area for the last 3 months.
- * Please consult your health care physician directly before using the product if you are not sure of the conditions.
- 4. The electrode plate on the middle part of main unit will generate heating during working. Be careful to use the product if you are not sensitive to heat. Do not stay on the same treatment area for too long to avoid scald.
 - 5. Do not treat the same area for too long when using BURN mode. Make sure the treatment head is in close contact with your skin and moves slowly to avoid scald.
 - 6. Do not use the product while driving or while concentrating on work.
 - 7. Do not disassemble, repair or modify this product. Improper handling might lead to fire, electric shock, injury or product malfunction. Unauthorized disassembly or modification may result in exclusion from the warranty.

【WARNING】

- 1. The product is for personal care only, not medical device. Do not use the product for other medical purposes.
- 2. Do not use the product in conjunction with other medical devices.
- 3. Do not rinse the product in water while using.
- 4. Remove metal jewelry (eg. necklaces, rings, watches, earrings) before using the product to avoid possible strong electric shock from contact between metal part of jewelry and electrode of the main unit.
- 5. Do not look at LED light for a long time to avoid eye injury.
- 6. Do not use the product on below area:
 - a) Area with plastic surgery;
 - b) Injured and not fully healed area like abrasion, scratch and scald;
 - c) Area of skin inflammation caused by abscess, rash and sunburn;
 - d) Area of erythema caused by itching or physical stimulation;
 - e) Area with implanted metal, silicone, introducer in body, face and mouth;
 - f) Area with tattoo or artistic make-up.
- 7. The product can be used by person aged from 14 years and above, and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the product in a safe way and understand the hazards involved. Children shall not play with the product. Cleaning and user maintenance shall not be made by children without supervision.
- 8. Do not remove the built-in rechargeable battery. Do not charge or use the product with high temperature. Follow local government laws and regulations for proper disposal of battery. Please recycle battery if possible. Do not heat battery or device containing battery into fire.
- 9. The product should be charged with an approved safety isolation adapter with the output capacity of DC 5.0V/2.0A and complied with IEC/EN 60335-1.

10. Make sure your hands are dry before touching the plug of USB charging cable or power adapter to avoid electric shock.
11. Do not place any object on the surface of USB cable during charging to avoid malfunction or fire. Please switch off immediately and remove USB plug if any failure occurs to avoid electric shock or fire.
12. Make sure USB plug or power adapter is inserted properly to avoid possible electric shock and fire.

[CAUTION]


1. Keep the product out of reach of children and pets while using.
2. Start from low intensity level for initial use, and increase intensity level gradually if you are accustomed to the product.
3. Keep the product away from heat source, open flame nor place it in a room with high temperature. Do not expose the product to the sun or rain.
4. Do not use the product close to flammable material like alcohol, diluent, gas or spray.
5. Do not throw the product to avoid falling and collision.
6. Charge the product before initial use or when it is not used for a long time. Make sure button is working and indicator is normal.
7. Unplug and switch off the product when not in use.
8. Stop using the product immediately if it shows visible damage or product malfunction in case of mechanical or electronic failure. Contact our service center or local distributor in time if necessary.
9. Battery Safety Information: The product contains non-replaceable battery; the battery must be removed before product is scrapped; the product must be disconnected from the supply mains when removing the battery; the battery shall be disposed safely.

Notes: The latest user manual is available via www.emayplus.com

Product Introduction

The product contains BURN and SLIM two modes, through the perfect combination of ultrasound, EMS, RF and LED technology, to achieve the purpose of fat burning, slimming and skin elasticity restoration. It can also stimulate collagen regeneration, tighten skin and remove wrinkles for beauty care with three selectable light sources.

Product Structure

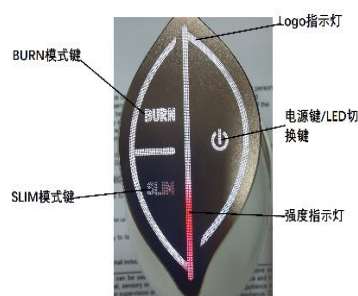

1. Main Unit
2. Base
3. Control Panel
4. Power / LED Switch Button
5. BURN Mode Button
6. SLIM Mode Button
7. Logo Indicator
8. Intensity Level Indicator
9. Type-C Charging Port
10. Charging Port Silicone Cover
11. Treatment Panel
12. Ultrasonic Electrode Plate
13. EMS/RF Electrode Plate
14. LED Light



Product Specification

Product Name	Emay Plus ZERO-FAT Body Slim
Product Model	EP-601
Input	DC 5V/2A
Rated Voltage	7.4V
Rated Power	MAX. 18W
Battery Capacity	7.4V/1000mAh
Functional Mode	2 Modes
Intensity Level	3 Levels
Charging Time	≤200 Minutes (The data is based on an input voltage of DC 5V/2A. Actual charging time will vary depending on usage and charging environment)
Working Time	≥50 Minutes (May vary from modes and levels)
Timer Function	10 Minutes
Waterproof Level	IPX6
Main Material	ABS, Stainless Steel
Product Size	115*80*97mm (Main Unit)
Product Net Weight	360g (Main Unit + Base)
Content List	Main Unit*1, Base*1, Storage Bag*1, USB Type-C Charging Cable*1, User Manual*1

How to Use

1. Functional Modes

Modes	Ultrasonic	EMS	RF	LED	Intensity Level	Timer
SLIM	•	•		•	3 levels	10 minutes
BURN	•	•	•	•		

2. Operation Instructions

- 1) Power on: Long press on "Power" button  around 1 second to start the device with two beeps. Power indicator, SLIM indicator and intensity

indicator are on respectively at this time, it enters into SLIM mode with low level, LED light is red on treatment panel.

- 2) Power off: Long press on "Power" button  around 1 second to turn off the device with two beeps. All indicators are off.
- 3) LED colour adjustment: Short press on "Power" button  under any operating state, LED colour will be adjusted in a cycle of <Red-Blue-Purple-Off>. Switch of modes or levels will not change the LED colour.
- 4) Mode switch: Short press on "BURN" button  when in SLIM mode, BURN indicator is on, and it enters into BURN mode with low level; Short press on "SLIM" button  when in BURN mode, SLIM indicator is on, and it enters into SLIM mode with low level.
- 5) Intensity adjustment: Short press on "SLIM" button  when in SLIM mode, or short press on "BURN" button  when in BURN mode, intensity level will be adjusted in a cycle from <Low-Mid-High> with corresponding red intensity indicator turned on. Each time the mode is switched, intensity will be reset to low level.
- 6) Timer function: For any type of mode, the device will shut down automatically with two beeps after working for 10 minutes. Countdown timing will restart if you switch mode, while it is not workable for intensity level adjustment.


3. Charging Instructions

- 1) Power adapter is not provided. Please charge the main unit with an adapter of DC 5V/2A.
- 2) Charge the main unit at room temperature.
- 3) When the battery is low, the indicator flashes with 5 beep tones before power off. Please charge in time for normal use.
- 4) Charge the main unit with supplied USB Type-C cable. Logo indicator lights on with a breathing light during charging and stays on after fully charged.
- 5) Unplug power adapter and remove USB charging cable after fully charged. Constantly doing a full recharge will shorten battery's lifespan.
- 6) The product cannot operate during charging.
- 7) The product does not support Type-C quick charge.

4. Treatment Method

- 1) Precautions for use
 - a) Cleanse the treatment area before treatment, make sure skin is clean;
 - b) Apply adequate amount of skin care products such as gels, toners, and slimming lotions to the treatment area; try to use less skin care product that contains much oil to avoid influence on ultrasonic performance;
 - c) Start from low intensity level for initial use;
 - d) Do not stay on the same treatment area for a long time to avoid scald;
 - e) Do not use the product in bathtub or hot spring with body wash, or in hot sauna with high temperature;
 - f) Electrode head of main unit shall be in contact with skin as vertically as possible and move slowly. Do not press or rub hard.

2) Body care method (recommend)

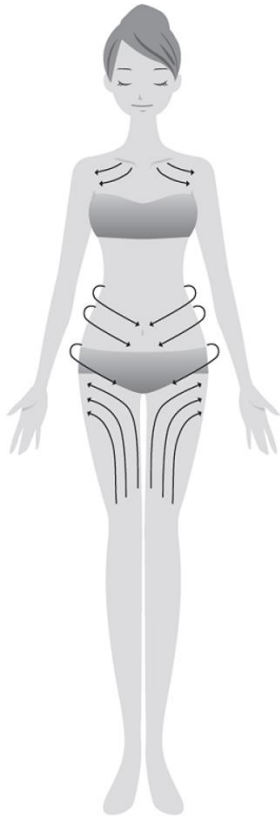

Move as indicated arrow (refer to picture) during treatment.

✧ **Waist and abdomen**

Move from back waist towards navel, or move slowly around navel in a circle way.

Recommend mode: SLIM mode + red LED

BURN mode + purple LED

Frequency of use: Recommend 10 minutes per time, 2-3 times per week.


✧ **Thigh**

Move from above the knee towards base of the thigh.

Recommend mode: SLIM mode + red LED

BURN mode + purple LED

Frequency of use: Recommend 10 minutes per time, 2-3 times per week.

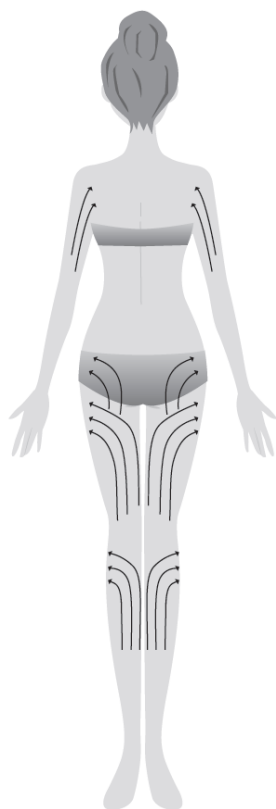

- ✧ **Upper arm**
Move from elbow towards underarm.
Recommend mode: SLIM mode + blue LED
Frequency of use: Recommend 5 minutes per time, 2-3 times per week.
- ✧ **Hip**
Move from lower part of hip towards outside.
Recommend mode: SLIM mode + red LED
BURN mode + purple LED
Frequency of use: Recommend 10 minutes per time, 2-3 times per week.
- ✧ **Calf**
Move from heel towards the back of knee.
Recommend mode: SLIM mode + red LED
Frequency of use: Recommend 5 minutes per time, 2-3 times per week.

Notes: Do not use the product on laryngeal prominence or area of chest and back right next to heart.

- 3) Please apply an appropriate amount of skin care product for moisture based on personal needs after treatment.

Cleaning, Maintenance and Storage

1. Switch off the product before cleaning.
2. Do not cleanse main unit with steel wire ball, alkaline detergent, abrasive chemicals, diluent, gas, nail polish remover or other corrosive liquids.
3. Clean the residue of cosmetics on treatment panel in time after treatment. Use a soft and wet cloth, wet tissue or cleaning cloth to wipe stains on treatment panel and surface of main unit. Wipe the treatment panel with dry cloth or tissue.
4. Do not immerse the product directly in water while cleaning.
5. Dry the product in a cool and ventilated place after cleaning, then store it in storage bag.
6. Do not keep the product in a place with high temperature or high humidity.
7. Charge the product every two months for not long-term use.
8. Do not put the product on uneven surface or high places to avoid injury from falling or damage to device itself.
9. Keep the product out of reach of children.

Causes and Troubleshooting

Problem	Possible Cause	Solution
Cannot turn on	Not properly press the power button	Long press on power button around 1 second
	Battery is low	Charge the main unit
Cannot charge	Use quick charge adapter	Charge with regular adapter
	Charging cable is not inserted properly	Make sure cable is fully connected to main unit
Shut down automatically during use	Wrong press on power button	Long press on power button around 1 second to start the device
	Battery is low	Charge the main unit
Cannot feel EMS	Skin is dry	Apply adequate amount of gel and toner
	Intensity level is low	Increase intensity level gradually
Cannot feel RF heating	Skin is dry	Apply an appropriate amount of gel and toner
	Intensity level is low	Decrease intensity level gradually
	Move too fast	Lower the speed of movement
Cannot feel ultrasound	Vibration of ultrasound goes too fast, almost cannot feel it on skin	Drop water in the middle part of metal electrode, if drip vibrates, then ultrasound is in normal use.

*This chapter summarizes the most common causes and solutions. If you have any inquiry about the product, please feel free to contact our customer service center or local distributor. It is strictly prohibited for unauthorized personnel to disassemble the product.

Product Warranty

We offer 12-month warranty for the product from the date of purchase (excluding any intentional damage). Please retain your receipt as proof of purchase.

FCC Warning

This device complies with part 15 of the FCC rules. Operation is subject to the following two conditions: (1) this device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

Changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

NOTE: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

The user manual or instruction manual for an intentional or unintentional radiator shall caution the user that changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment

Radiation Exposure Statement

The device has been evaluated to meet general RF exposure requirement. The device can be used in portable exposure condition without restriction.

Importer: XXXX

Address: XXXX

[Back Cover]



This marking indicates that this product should not be disposed with other household wastes throughout the EU. To prevent possible harm to the environment or human health from uncontrolled waste disposal, recycle it responsibly to promote the sustainable reuse of material resources. To return your used device, please use the return and collection systems or contact the retailer where the product was purchased. They can take this product for environmental safe recycling.

[Cert Logo] CE, FCC, TELEC, RoHS

[Emay Plus logo]

Distributor: Emay Global Limited

Address: 1/F., Centre 600, 82 King Lam Street, Cheung Sha Wan, Kowloon, Hong Kong

Tel: (852) 3595 3871

Email: cs@emayplus.com

Website: www.emayplus.com