Wireless Pressure Mat Doorbell And Alarm System

!!! Please visit our website: www.geemarc.com as there may be important updates and changes that you should be aware of !!!

Contents of the box



The box contains the following items:

- 1 Guide to use
- 1 Cousin drop detector with its transmitter
- 1 Receiver of alert
- batteries are not provided
- power adapter is provided as an option (6v 0.3a).


DESCRIPTION of receiver and transmitter

DESCRIPTION of receiver

1. Ringing level low (Lo)/strong (Hi)
- 2,3,4. Indicator lights of alarm
3. Low battery indicator light
5. Ringing of the bell
6. Alarm/Off/Chime switch

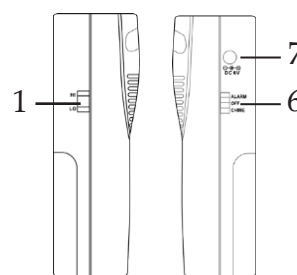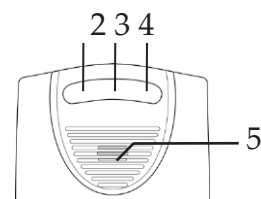
7. Power adapter connector 6V 0.3A
8. "Learn" pairing touch
9. Battery compartment cover
10. Wall mounting of the receiver

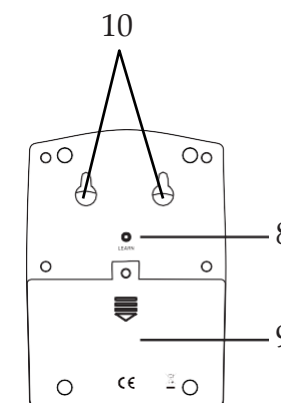Description of the issuer

1. Ringing of the bell (85dB)
2. Battery compartment cover (PP3 9V not included)
3. Interrupteur Alarme/Arrêt/bip


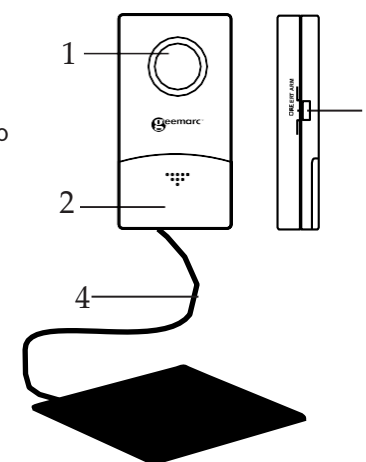
4. Presence detection pad


INSTALLATION

Operation of the transmitter

1. Open the lid of the battery compartment (Fig. 1) by sliding it down.
2. Connect the PP3 alkaline battery type (not included) to the connector located in the compartment respecting the polarity.
3. Close the lid of the battery compartment by sliding it To the top.
4. Connect the sensor pad cable to the transmitter (Fig. 2).
5. v You can fix your transmitter, make 2 vertical holes Separated by 44 mm.
6. Insert 2 wall fixing pins into the 2 holes.
7. Screw 2 screws into the pins and let them protrude 5mm, then slide down your transmitter to lock it.

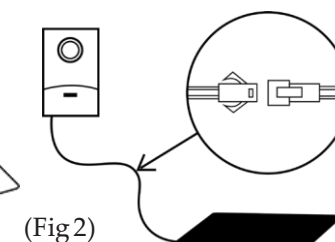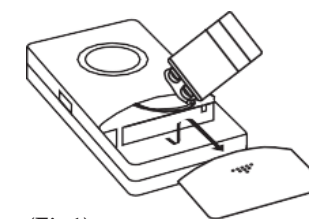
8. Remove the cushion from the chair and place the detection cushion

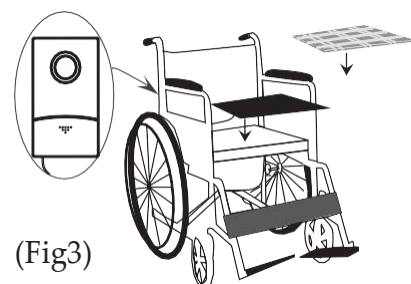9. Replace the chair cushion over it.

Note: a warning beep will alert you when the battery is to be replaced every 2 minutes.

Starting up your receiver

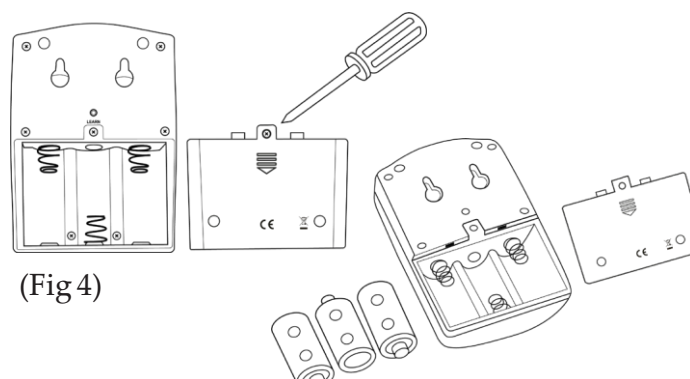(Fig 4)

(Fig 5)

1. Remove the marker screw (Fig. 3) with your cruciform and open the cover of the pile compartment (9) by sliding it down.
2. Insert the 3 1.5 V LR14 alkaline batteries (not included) respecting the polarity "+" and "-" (Fig. 4) Fig. 3 Fig. 4
3. Close the lid of the battery compartment by sliding it up.
4. Replace the reference screw (10) with your cruciform.
5. You can place your receiver in wall position, make 2 separate horizontal holes of 36 mm.
6. Insert 2 wall fixing pins into the 2 holes.
7. Screw 2 screws into the pins and let them protrude 5mm, then slide them into 11 marks by pressing down your receiver to lock it.


Note: you can connect a 6v 0.3a adapter which is sold as an option. You can leave the batteries in your device for operation in case of power failure.

The middle light (3) will start flashing to warn you that the battery will need to be replaced.

Using your device

The AmpliCall PSA 200 system you just purchased is very easy to use, requires no external wiring and can be used at home, in a hospital, retirement home or other. Your system is intended to warn a person or caregiver in the receiver room that the person no longer puts pressure on the sensing pad, such as in a wheelchair or chair.

Maximum Distance of 30 meters

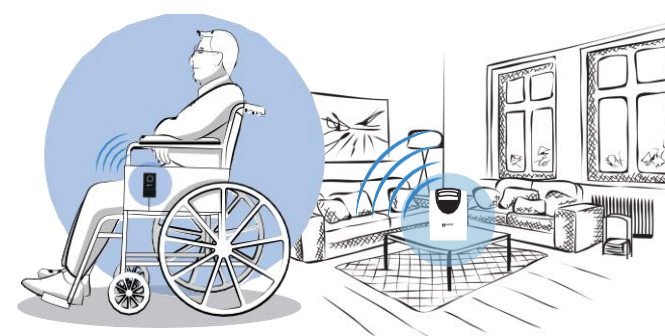
When the person is no longer sitting on the sensor cushion, the lack of pressure will trigger the sending of information to the Receiver which must be placed at a maximum distance of 30 meters from the transmitter. The receiver will start to ring or beep (in weak or strong position) depending on the chosen mode and he will have the 3 flashing lights.

1. The user is placed on the detector cousinr
2. Set the switch of your transmitter (4) to the "ALERT" position you will beep after 4 seconds to signal that your transmitter is ready. Your transmitter will just send the information to the receiver every 40 seconds as long as no pressure Is detected on the detection pad. "ARM" you will beep after 4 seconds to signal that your transmitter is ready. It will start to beep for 5 seconds if a drop is detected at the same time as your receiver. This procedure will be repeated every 40 seconds as long as no pressure is detected on the detection pad.

3. Placez Set the switch (6) on your receiver to the position "CHIME" "it will beep and blink for 5 seconds when it receives information from the transmitter. After the 5 seconds ring, you will have a 30seconds stop and it will resume for 5 seconds and so on as long as the sensor pad does not detect a pressure. "ALARM" it will ring and blink for 30 seconds when a signal is received from your transmitter. After the 30 seconds ring, you will have a 10seconds stop and it will resume for 30 seconds and so on as long as the sensor pad does not detect a pressure.

Note: place your switches on the OFF of your 2 systems to stop your PSA 200 amplifier ringing or to save your batteries.

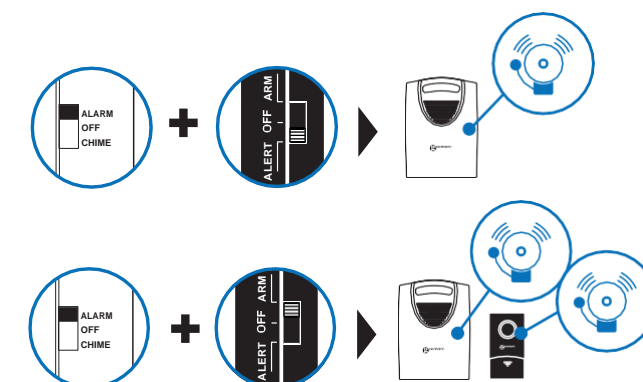Recording a new receiver on your sensor pad

Note: when you receive your Amplicall PSA200 package, the devices are already paired. So you don't need to follow the steps below

To pair a new receiver (up to 10), please follow the steps below :

Step 1: turn on your receiver in "Alarm" or "Chime" mode and enter recording mode on your receiver by remaining press the "LEARN" button with a pen for 3 seconds, release the button when you hear a beep.

Step 2: immediately afterwards, put your transmitter in position "Alert" or "Arm" without pressing the detection pad.

Step 3: if the recording is successful, your receiver will make 4 beeps and will start to ring according to the mode of your alarm.

Clear a receiver from the sensor pad

You can delete the recording of a receiver by pressing its "LEARN" button (8) for 12 seconds you will hear 2 long beeps sound To confirm the deletion.

Characteristic of your AMPLICALL PSA 200

Range: up to 30 meters
Frequency of use: 433.92Mhz

Level of ringtone:
Issuer of securities: 75dB
Receiver receiver : 85dB (fo) 70dB
Energy in Europe: Transmitter: 1 alkaline battery 9V type PP3
Receiver receiver: batteries Alkaline 1.5V type LR14
Number of receiver for a transmitter:10

Breakdown of work

- 1) **The receiver does not sound:**
please check that the receiver is either on function.
ALARM or either on CHIME.
check that the transmitter is operating on the.
ALERT or ARM position.
Check that the receiver is not more than 30 meters from The transmitter.
Check that there is no pressure on the sensor pad
Change the batteries.
- 2) **Transmitter does not ring:**
Check that the transmitter is operating on the
ALERT or ARM position
Your transmitter beep every 2 minutes, replace the battery.
Check that there is no pressure on the sensor pad
Change the batteries.

FCC Warning

This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions:
(1) This device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.


NOTE 1: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:
- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
-Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
• -Consult the dealer or an experienced radio/TV technician for help.

NOTE 2: Any changes or modifications to this unit not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment
The device has been evaluated to meet general RF exposure requirement. The device can be used in portable exposure condition without restriction

Directive IC

“This device complies with Industry Canada’s license-exempt RSSs. Operation is subject to the following two conditions:

(1) This device may not cause interference; and
(2) This device must accept any interference, including interference that may cause undesired operation of the device.”

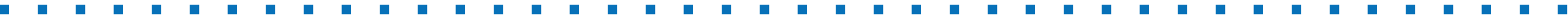Guarantee guarantee

Your Geemarc AMPLICALL PSA 200 product is guaranteed for 2 years, depending on the date of purchase. During this period, all repairs or replacements (at our sole discretion) will be made at no charge. If you have any problems, please contact our helpdesk or visit our website: www.geemarc.com. The warranty does not cover accidents, negligence or broken parts. The product must not be altered or disassembled by anyone who is not an approved representative Of Geemarc. Geemarc's warranty does not limit your legal rights.
.

IMPORTANT: YOYR INVOICE IS PART OF YOUR WARRANTY AND MUST BE KEPT AND PRODUCED DURING A WARRANTY CLAIM.

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eemarc Telecom SA hereby declares that its AMPLICALL PSA 200 Complies with the essential requirements as well as other relevant directives of standard 2014/53/UE.

*The declaration of conformity can be consulted on the following website **www.geemarc.com***

Guidelines for recycling

Weee (waste electrical and electronic equipment)
Was created for end-of-life products to facilitate their recycling. When this product is at the end of its life, please do not dispose of it in your household garbage bin.

Please use one of the following disposal options:
remove the batteries and put them in a suitable weee bin. Place the product in a suitable weee bin.

or, return the old product to the retailer. If you buy a new one, they should accept it. Thus, if you follow these instructions, you are ensuring the protection of human health and the environment

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