

*** USER'S MANUAL ***

FCC ID : 2AABU-WBT168

The Federal Communication Commission Statement

This equipment has been tested and found to comply with the limits for a Class B Digital Device, pursuant to Part 15 of the FCC rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instruction, may cause harmful interference to radio communication. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one of more of the following measures: -

- **Reorient or relocate the receiving antenna.**
- **Increase the separation between the equipment and receiver.**
- **Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.**
- **Consult the dealer or an experienced Radio/TV technician for help.**

Use only shielded cables to connect I/O devices to this equipment. You are cautioned that change or modifications not expressly approved by the party responsible for compliance could void your authority to operate the equipment.

THIS DEVICE COMPLIES WITH PART 15 OF FCC RULES. OPERATION IS SUBJECT TO THE FOLLOWING TWO CONDITIONS:

- 1 This device may not cause harmful interference and
- 2 This device must accept any interference received, including interference that may cause undesired operation.

The antenna used for this transmitter must not be collocated or operation in conjunction with any other antenna or transmitter.

Notice : The changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.


IMPORTANT NOTE: To comply with the FCC RF exposure compliance requirements, no change to the antenna or the device is permitted. Any change to the antenna or the device could result in the device exceeding the RF exposure requirements and void user's authority to operate the device.

WBT-168 Bluetooth Bathroom Scale User's Guide (Android system)

SPECIFICATIONS

- Capacity: 150 kg / 330 lb
- Graduation: 100 g / 0.2 lb
- Bluetooth connectivity – BT 2.1 + EDR
- Upload data each time automatically to smartphone (**Android system only**)
- Apps that calculate BMI, track trend, provide weight chart and allow multiple users
- Stainless steel platform (size: 256 x 267 mm)
- LCD display (size 68 x 33 mm) with backlight
- Tap on & auto off
- Low battery indication
- Over-load indication
- Error indication
- Power source: 4 x 1.5V AAA batteries

BEFORE USE

- 1) Turn over the scale, find the battery cover on the top, slide it to the right and take it out. Insert 4 pieces of AAA battery with correct polarity in battery compartment. Close the battery cover.
- 2) Put the scale on a flat table, ensure it is steady.
- 3) Download app "Satrue WBT-168" from Google Play to your smartphone. (<https://play.google.com/store/search?q=satrue&c=apps>) or scan the QR code on the right with your smartphone. Find the icon  on the screen of your smartphone.
- 4) Turn on Bluetooth of your smartphone.




OPERATIONS

- 1) Run the app you have previous downloaded from Google Play.
- 2) Establish your user data, please refer to "Advanced Operations – create new user".
- 3) Turn on scale - tap on the platform or press the white button on the bottom.
"8888" appears on the display, then 0.0 after 2~3 seconds. If no weight applied to the platform for more than 5 minutes, the scale will turn off automatically.




Important!

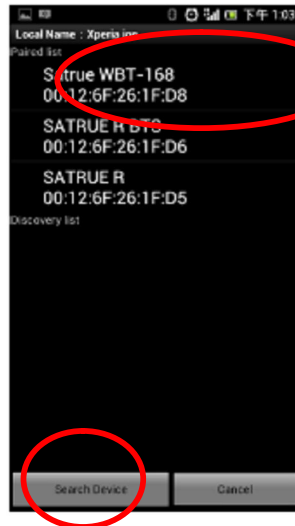
For first time use, press "BT connection" then "search device" to search Satrue WBT-168, click Satrue WBT-168 to pair your smartphone with the scale. If pair-wise key is required, enter 1234. After Bluetooth connection is enabled, the app will automatically memorize the scale, once the app is turned on, it will automatically search the scale it paired with and enable Bluetooth connection. If the automatic connection does not work, please go to setup page, disable Bluetooth and turn it on again.



- 4) Carefully step on to the platform and stand still. Your weight will display on scale LCD and transmits to your smartphone.
- 5) Press "Conform" beside weight column, the weight and weighing date will be recorded.



- 6) Press  to check your weight and BMI chart. Color of the curves is changeable, press   to select your prefer color.

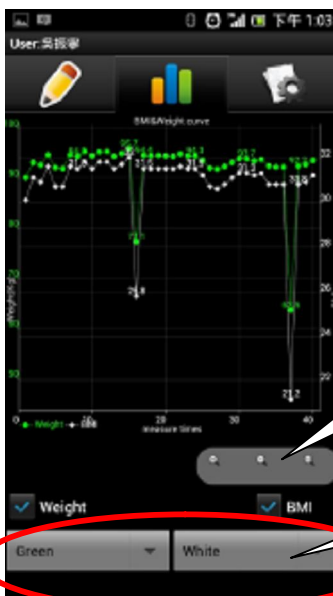


Zoom in and zoom out the screen of curves.


Select your prefer color for curve.

ADVANCED OPERATIONS

Create new user

Press , Press "Create New User", input your name, height, weight, birthday and gender, press "Confirm". If weight is left blank, the initial weight on the chart will be 0. Please input a correct date for birthday, a date in the future is not allowed.

Delete User

- 1) Delete a user unwanted. If all the users are deleted, the screen will go to Create New User page.
- 2) Press , Press "Delete User", press "Confirm" to delete current user.

Clear user REC.

Press to clear weight record of the current user.

Select Language

Select language in Traditional Chinese or English.

Send Record

Send weight record by e-mail. Input your e-mail address to the blank column above "Send Record", press "Send Record", select your e-mail software and send.



Press to turn on or turn off Bluetooth.



Press to escape APP.



ADVICE FOR USE AND CARE

- 1) The scale is made with sensitive electronic components; keep it away from any heat sources and environment with excessive moisture.
- 2) For accurate measurements, put the scale on a flat, solid and non-slip floor.
- 3) Only clean the scale with neutral detergent and a cloth.
- 4) Remove the batteries if you are not going to use the scale for a long period of time, and store it in a dry and clean place.
- 5) Do not attempt to disassemble or repair the scale.
- 6) Treat the scale with care. Never drop it or jump on it.

WARNING INDICATION

- 1) If weight exceeds maximum capacity of 150kg (or 330lb), "----" will be displayed. You should step off the scale immediately or else the scale may be damaged.
- 2) If "batt" is displayed, replace batteries.