You Brush.
We Care.
You Smile.

playbrush

Smart One
Smart OneX

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
INSTRUCTION FOR USE

Model number: A122

Manufacturer: Trisa AG

Distributor: Playbrush GmbH,
www.playbrush.com, support@playbrush.com

Electrical specifications:
see bottom of charging device

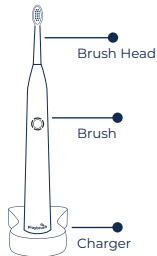IMPORTANT SAFEGUARDS

When using electrical products, especially when children are present, basic safety precautions should be followed at all times during

use, storage and disposal, including the following:

READ ALL THESE INSTRUCTIONS

DANGER – To reduce the risk of electrocution:

Children and people with reduced physical, sensory or mental capabilities, or lack of experience and knowledge

must be supervised during use of the device or are warned prior to use and understand potential dangers. Cleaning and user maintenance shall not be made by children without supervision. Children shall not play with the appliance. The product is not suitable for children less than 3 years old.

Check the power cord and charger regularly for cracks or other damage. Never touch the device with sharp or hot objects. Never unplug by pulling on the cable instead of the connector. The power cord cannot be replaced. If the power cable is damaged, dispose of the device. Do not continue to use a damaged charger or Playbrush.

Do not use charging base near water or outdoors. Never touch the plug with wet hands. Do not reach for a product that has fallen into water. Unplug immediately. Do not use while bathing. Do not place or store product where it can fall or be pulled into a tub or sink. Do not place in or drop into water or other liquid.

WARNING – To reduce the risk of burns, electrocution, fire, or injury to persons.

Use this product only for its intended use as described in this manual (tooth brushing). For the purposes of recharging the battery, only use the detachable supply

unit 4680 or 4656 provided with this appliance.

Symbol detachable supply unit

Only use with original brush heads. If the brush head breaks off, after fall, impact or improper use, there is danger of suffocation. Close supervision is necessary when this appliance is used by, on, or near children or

individuals with certain disabilities. Never use while sleeping or drowsy. Do not use outdoors or where oxygen is being administered. Do not repair or change the product. Do not use attachments not recommended by the manufacturer. Never operate this product if it has a damaged cord or plug, if it is not working properly, if it has been dropped or damaged,

or dropped into water. Return the product to a service center for examination and repair. Keep the cord away from heated surfaces. Batteries are non-replaceable.

SAVE THESE INSTRUCTIONS

WARNINGS

Use only at room temperature (0-40°C). This toothbrush is suitable for braces and dentures. Any person with a pacemaker or other electrical implant should consult a medical professional in advance.

HYGIENE

Replace brush heads every 3 months. Thoroughly rinse brush head under running water after each use.


BATTERY PACK

To extend the life of the battery, use the toothbrush at least every six months until fully discharged and fully charged again. The toothbrush should never be left uncharged and charged at least every 6 months. Place the toothbrush in the charger after each use to ensure optimal performance at all times. For the most energy-efficient use, it is recommended that the plug of the charger be removed from the socket after full charge.

REMOVING THE BATTERY CELL

Replacement of the battery is not possible. Only dispose of fully discharged batteries. Once opened the guarantee automatically expires.

1. At the end of the service life, open the cap with a coin.
2. Press the unit against a firm surface and push out the inner components.
3. Cut through connecting wires and remove battery. Do not short-circuit the positive and negative poles and mask off for protection. Take the battery to a designated collection point.

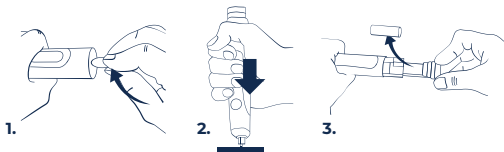RECYCLING

Electrical appliances must not be disposed of with household waste. Use the return options in the designated sales or collection points. This toothbrush contains a lithium ion rechargeable battery. For environmental reasons, it must be removed at the end of its life and disposed of in an environmental manner.

GUARANTEE INFORMATION

For this service, the manufacturer will provide you with a two-year guarantee from the date of purchase (purchase receipt). The guarantee covers the repair or replacement of a device with material or manufacturing defects. To maintain the warranty, use only original Playbrush brushes. Wear and tear in the context of normal use are not covered by the guarantee. The guarantee expires in the case of third-party entries or changes to the original condition, as well as in commercial use. The guarantee service requires that the defective device is sent with the purchase receipt.

The manufacturer disclaims any liability for damage caused by non-observance of these operating instructions, by improper use and repairs or incorrect operation. We accept no liability for damage to the smartphone or tablet that has arisen in connection with the use of Playbrush. Subject to change!

DECLARATION OF CONFORMITY

Hereby, Playbrush GmbH declares that Playbrush complies with the essential requirements and other relevant provisions of Directives 2014/30/EU, 2014/35/EU, 2014/53/EU, 2009/125/EC and 2011/65/EU. The compliance statement can be viewed at <https://www.playbrush.com/doc.pdf>.

FCC RULES

This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions:

- (1) this device may not cause harmful interference, and
- (2) this device must accept any interference received, including interference that may cause undesired operation. No changes shall be made to the equipment without the manufacturer's permission as this may void the user's authority to operate the equipment.

This device has been designed and complies with the safety requirements for portable RF exposure in accordance with FCC rule part At the end of the §2.1093 and KDB 447498 D01. FCC ID: 2AF47A122

Let's get you started...

1. Activate

The battery of your Smart One is pre-charged. To activate, place the brush on the supplied charging station, until a solid green light is shown.

Important: For maximum battery power, leave the brush on the charging station for the first 7 days. However, keep using the brush as usual.

2. Turn on Bluetooth

Turn on Bluetooth on your smartphone or tablet. Don't actively pair the Playbrush Smart One in the settings. This takes place in the app.

3. Download the Playbrush App

Download the Playbrush App from the Apple App Store or Google Play Store.

Android users: As Bluetooth falls under the location settings, please allow the app to access your location. Playbrush will not actually access your location.

A one-off registration is required when first launching the app. If you have purchased this product from the Playbrush website, please enter the email address you used for your order. Otherwise please use any email address.

4. Connect

Pick up your Playbrush and launch the app. If the Playbrush isn't showing a blinking blue light, press the start button. If your Playbrush shows a blinking blue light, it will automatically begin to connect within the app. Your Playbrush is successfully connected when it shows a solid blue light.

5. Explore & brush as you are

Cleaning performance & cleaning programs

Click here to change the cleaning program. Make sure that your brush is connected with the app.

Cleaning programs


You can easily adjust the intensity of the vibrations in the app.

Sensitive - for sensitive teeth & gums

Clean - for a particularly thorough deep cleaning

Massage - for gentle tooth polishing and massage of the gums

Ultra Sensitive - for Smart One X only

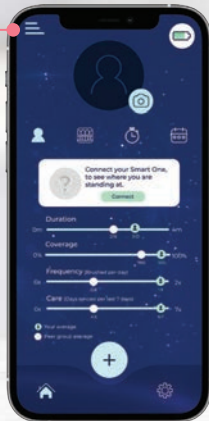Indicator bristles

If the gray bristles have faded by half, it is time to change the brush head, at the latest after 8-12 weeks.


The app also reminds you to swap brushes.


Hydro-clean vibrations

The smart micro-vibrations of your toothbrush can easily reach the spaces between the teeth without having to exert pressure. So slide slowly over your teeth for the best cleaning results.


Coverage Monitor

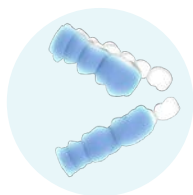
Left side
completed

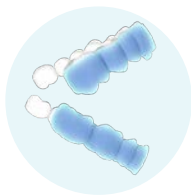Right side
completed

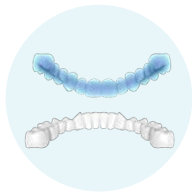Top teeth and
inner surfaces
completed

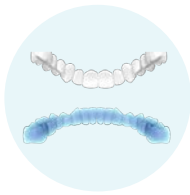Lower teeth and
inner surfaces
completed

Coverage Monitor

If all 4 lights are shining,
then all quadrants in the
mouth have been sufficiently
cleaned.

There is a small gap of 1mm
between the brush and brush
head, which is necessary to
support the Hydro-Cleaning
vibrations.

Cleaning Tipp: Clean the
metal rod underneath and
other open areas regularly.

Battery

-  Red flashing = Battery is low
-  Red continuous = Battery is empty
-  Green flashing = Battery is charging
-  Green continuous = Battery fully charged

Top Tip: To boost your battery, for the first 7 days of use return the brush to the charging station after each brushing session.

Light feedback during your cleaning session

Blinking blue = Ready to connect to the app

Continuous blue = Connected to the app

Flashing orange = You are brushing too fast or too hard.

Flashing green/red = Your Playbrush needs a performance upgrade. Connect to the app.


Short vibration (1 impulse):

A quarter of your cleaning time has elapsed.

Long vibration (3 impulses):

The cleaning time has completely elapsed.

The Playbrush App

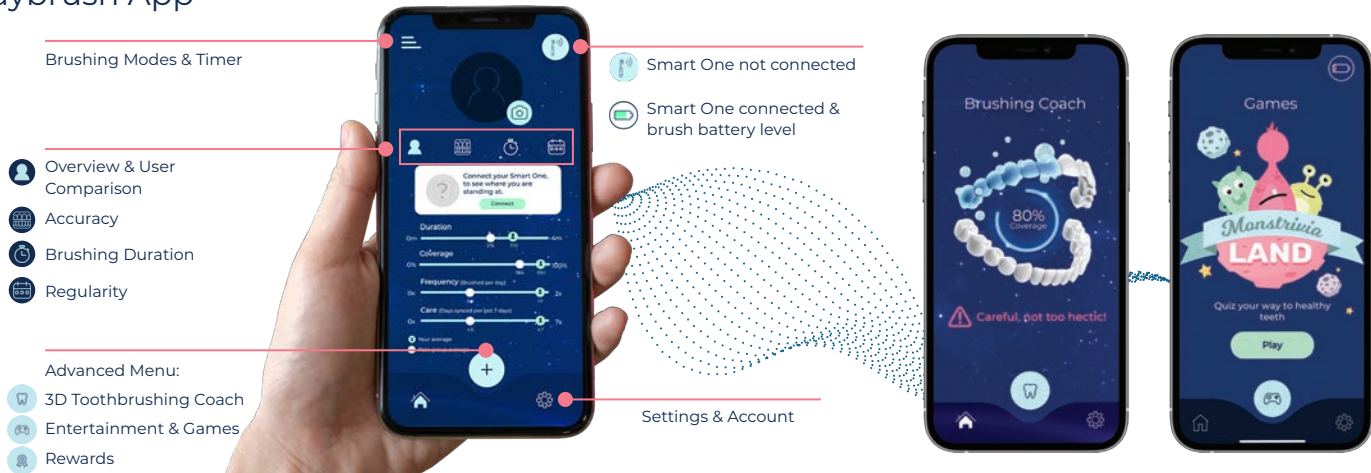All App Features in Detail

Personal Space

Monitor your status-quo as well as brushing progress over time and get highly personalised tips!

Personal Recommendation

This is your home base!

You control the Playbrush App from here. Also find a quick overview and key brushing tips.

Detailed Statistics

Want to monitor your progress?

Keep your brushing time, regularity, coverage and pressure in check and observe your scores improving over time.

Cash-Back

Pick up rewards

Collect cash with every cleaning session and have it paid out to your account for your next dentist bill. Premium feature for subscription customers*.

*in Germany and Austria

Training Area

Choose from a wide variety of oral hygiene tools to optimize your routine!

Brushing Coach

Perfect brushing your teeth

Get real-time feedback on your performance as you clean your teeth perfectly and 100% interactively with the help of a 3D simulation.

Brushing Entertainment

Bored while brushing your teeth?

Navigate through game adventures with your own brushing movements.

Settings

Customise your brushing preferences by adjusting the settings as you please!

Brushing Preferences

Want to customise your experience?

Feel free to set preferred brushing modes (clean / sensitive / massage), timer (2-4 minutes) or screen mood (dark/light).

App Connection & Use

To connect your Smart One sonic toothbrush to the app, follow the Quick Start Guide on page 2. For connection problems, consult the FAQs on page 27.

You don't want to constantly use the app?

No problem! To keep track of your progress, offline brushing sessions, i.e. those done without the app, can be stored on the brush for up to 6 months and transferred when synced.

Syncing your brushing progress

For necessary performance updates, connect the brush to the app regularly, even if you are not using the app.

Replacement Heads & Other Services

Discover our replacement brush heads deliveries by post to your home, and many other subscription services at **Playbrush.com**. **30 Days Free Trial!**

playbrush 

Attention parents: Playbrush Smart One also connects to the Playbrush Kids app. So you can easily clean together as a family.

Top Tip: Use the same login credentials for both apps so your progress can be synchronised!

Playbrush for Kids

We believe that good dental habits are a family affair. That's why we've also developed Playbrush for kids, smart toothbrushes for kids connected to gaming apps that make brushing not only more enjoyable but also help kids to interactively improve their brushing skills.

12 Game Adventures
(incl. multiplayer)

Statistics

Brushing Coach

Refill Services

Find out more on **Playbrush.com** and save 25% with the Code **FAMILYONE** when purchasing a Playbrush Kids toothbrush or plan.

Playbrush Smart
Attachment
for manual
toothbrushes for
kids 3+

Playbrush Smart Sonic
Electric toothbrush
for kids 6+

FAQ - Troubleshooting

Toothbrush does not work:

- Your Playbrush may need to be charged. Place the toothbrush on charge overnight.
- Try to check the socket with another electrical appliance.

Brushing performance is low or insufficient:

- The toothbrush head may need to be cleaned or replaced.
- Your Playbrush may need to be charged.
- Change cleaning mode.

Toothbrush is turning off during brushing:

- Playbrush will automatically turn off when the battery level is low. Place the Playbrush on charge for at least 1 hour or overnight if possible.
- During brushing Playbrush also vibrates for every quarter of the brushing duration (changing sides).

Toothbrush vibrates in unusual motions during brushing:

- A short interruption can indicate the termination of the brushing session.

- The vibrations will also quickly pause to instruct you to change sides.

App connection does not work:

- Make sure your Playbrush Smart One is charged and flashing blue.
- Activate Bluetooth on your mobile phone/tablet and connect your Playbrush directly via the app - do not actively couple it with the end device!
- Make sure the location settings have been activated:
 - Settings > Privacy > Location Services> Playbrush App> Location Shared
 - Settings > Playbrush app > Bluetooth active + location enabled

Note: We are not actually accessing the location.

Didn't find what you were looking for?

You can find more FAQs on **Playbrush.com** under the navigation point "Support".

Still have questions?

Ask us anytime **support@playbrush.com**

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