HOME

WAYSKIN

[care button] [care part] : manual selection : manual change (color change) 2 seconds : go back [screen] : manual selection 2 seconds : go back [measurement button] [skin measurement]



WAYSKIN HOME -Your Personal Skin Therapist

WAYSKIN HOME offers LED color therapy for your unique skin by using application interworking and the Wayskin app to analyse big data of your skin/environment, and it helps your skin improve effectively with micro-electricity to vitalize your skin's energy.

Effects of LED color therapy for each color







and Pigment

Rose Red Improvement of Elasticity

Increase of Skin Moisture Improvement of Skin Texture

Boost the Creation of Collagen and Ellastine Improvement of Blemishes Prevention of Skin Aging

Use it like this

Prevention of Trouble



1. Use LED Color Therapy after pasting the WAYSKIN SOOTHING CREAM or Moisturizing Product. (It's okay to use one product only if it is more convenient for you.)



Check the changes to your skin through the skin measurement before/after care.

1. Avoid looking directly at the LED light.



Direct visual exposure to the LED light may cause side effects such as loss of vision.


You could feel a slight stinging while undergoing LED care because of the micro-electricity.







The care part of the device could become hotter during the LED care.

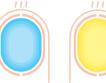



Don't worry. The heat creates a good effect in the vitalization of your skin cells.

4. The skin measurement is not working during charging.



Please use the skin measurement after releasing the charger.

How to use WAYSKIN HOME

1. Turn on the power

1) Press the [Care Button] or [Measurement Button] for two seconds.

2) If it does not turn on, please charge the device and try again.

Skin Measurement

1. Check the Skin

[Skin Measurement] Enter the skin measurement after pressing the button once.




[Skin Measurement] Start the skin measurement after pressing the button once.



Attach the part you want to measure to the [measurement part] of the device.



There will be three beeps after a while and the measurement ends with vibrations. (Don't remove the device from your face until you get the vibrations.)

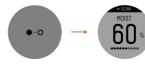

2. Finish the skin measurement

End the skin measurement after pressing the [Measurement Button] for two seconds.

LED Color Therapy

1. Start the care

Enter the LED Care after pressing the [Care Button] once.



Choose the LED care color after pressing the [Care Button] twice.



Start the care after pressing the [Care Button] once.



Attach your face to the [Care Part] of the device.



After a moment, the care will start. Rub the device on your face as if you were massaging your face with it.



The care will finish automatically after three minutes.



2. Finish the Care

Finish the care by pressing the [Care Button] for two seconds.

3. Pause and restart

- 1) While undergoing care, press the [Care Button] once to change the pause status.
- 2) If you press the [Care Button] once again, the care will resume.

Mobile Connection

Connect the device to your mobile device.

- 1) Search for Wayskin in the Apple app store or Google Playstore and download the application.
- 2) Press the home icon of Wayskin and enter the [WAYSKIN HOME] screen.
- 3) Press the [Connection] button and connect the device.
- 4) In some cases, the connection may fail
- When the distance between the device and the mobile device is too great.
- When the device is turned off or battery power is insufficient.
- When Bluetooth has not been enabled.

2. Using the device while the mobile device is connected

While your mobile device is connected to the device, you can use every function of WAYSKIN HOME through the inside buttons of the mobile device and device buttons.

ETC

When the screen or buttons of the device are not working, use a small, thin object like a paperclip to push the [RESET] button to reset the device.

2. Turning OFF (long-term hibernation)

- 1) Press the [Care Button] and [Measurement Button] of device together for at least five seconds.
- 2) Battery power is not consumed at this stage.
- So it is recommended that you turn off the device for long periods of time.

Homepage www.wayskin.com

E-mail help@wayskin.com

FCC Precautions

This device complies with part 15 of the FCC Rules for use with cable television service. This device complies with part 15 of the

FCC Rules. Operation is subject to the following two conditions: (1) This device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

The Grantee is not responsible for any cha nges or modifications not expressly approbed by the party responsible for compliance. such modifications could void the user's

no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be following measures: authority of operate the equipment.

This equipment has been tested and found to

comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules.

These limits are designed to provide

reasonable protection against harmful

interference in a residential installation.

—Connect the equipment into an outlet on a circuit different from that to which the recei

ver is connected.

/TV technician for help.

determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the —Reorient or relocate the receiving antenna. —Increase the separation between the equip ment and receiver. —Consult the dealer or an experienced radio

This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is

