NUFACE®

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NUFACE®

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DOWNLOAD THE APP

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Get more out of your device by downloading the NuFACE Smart App to unlock exclusive treatments, watch step-by-step tutorials, track your results over time and so much more.

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2. Activate + Smooth

2. Activer + Lisser

Apply a mask-like layer of FIX® Serum (Microcurrent Activator) in the section(s) you want to smooth. Push the ON/OFF button on your device. Perform quick, feathering and circular motions with your FIX+ device. Keep both spheres in contact with your skin during the treatment. Complete the same sequence in sections around the eyes, lips, and forehead.

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3. Finish

3. Terminer

Massage into skin any remaining FIX® Serum and finish by applying eye cream, moisturizer, and SPF.

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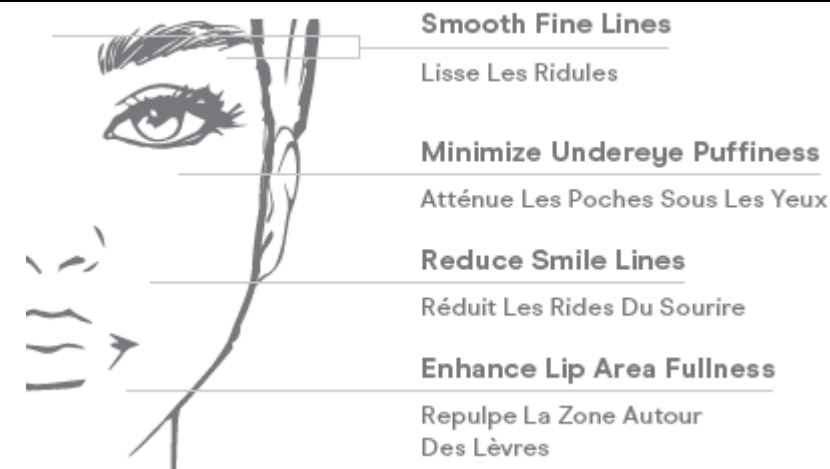
The FIX+ device produces low levels of microcurrent which is emitted through the two spherical probes as they both touch the skin. Prior to using the FIX+ device, the selected area of the skin to be treated is first prepped with a NuFACE® Microcurrent Activator to moisturize the skin and ensure optimal microcurrent conductivity. The two spherical probes are then gently glided over the skin to deliver the electrical impulses per the selected treatment setting. Microcurrent is an aesthetic modality providing electrical current in millionths of an ampere and has the ability to increase the firmness of the skin and muscles.

INTENDED USE

The NuFACE FIX+ Device is intended for facial stimulation and is indicated for over-the-counter cosmetic use on the face and neck as described. Any consequence resulting from misuse, application of the device to other parts of the body or face that are not approved for use, connection to improper voltage sources, dirty spheres, expired or non-NuFACE conductive alternative, or any other improper application is neither the responsibility of NuFACE nor its affiliates. **DO NOT** use on infants and children under the age of 18. Keep this device out of the reach of children and pets.

IMPORTANT INFORMATION

NUFACE®

- Use the FIX+ Device only as described within this User Manual.
- The risks and dangers of using the device in any way that does not follow this User Manual are unknown.
- Your NuFACE FIX+ Device has a unique serial number located on the back of the device. Record your serial number and register your device online at www.NuFACE.com or download the NuFACE Smart App.

The NuFACE device is designed for cosmetic use only.

The device should be operated, transported, and stored in the manner specified in this manual.

Failure to use and maintain the device and corresponding attachments in accordance with the instructions outlined in this manual will void the warranty and may result in product damage or physical injury.

The long-term effects of transcutaneous electrical stimulators are unknown.

Apply stimulation only to normal, intact, clean, healthy skin.

Avoid applying excessive pressure when using the device.

Skin irritation or hypersensitivity may occur due to the electrical stimulation and / or use of the NuFACE® Microcurrent Activator.

A microcurrent sensation may occur while using the device. Using a lower intensity and/or additional activator may be needed to reduce the sensation.

Skin may become dry during treatment in which Microcurrent Activator may need to be reapplied.

A perception of flashing lights may occur if device is used near or around the eye region due to stimulation of retina and/or optic nerve from the electrical current. The perception of flashing lights is not harmful.

DO NOT use the device if you see any signs of damage or leaks from the device.

DO NOT remove the battery of the device.

DO NOT store or leave the device under the sun or on a hot surface.

Follow all local regulations and waste management guidelines for safely disposing devices containing batteries.

The device has been designed to provide freedom from unacceptable risk when used under normal conditions and reasonably foreseeable hazards. In the event of complete loss of treatment function, the device does not pose a risk to the user.

CONTRAINDICATIONS FOR USE

Consult your healthcare provider if you have any pre-existing or potential health conditions that may be impacted by use of the device including those listed below.

- **DO NOT** use if you have an electronic implanted device such as a defibrillator, neurostimulator, pacemaker, or ECG Monitor. If you use a body-worn medical device, such as an insulin pump, consult your healthcare provider.
- **DO NOT** use this device if you have impaired cognition.
- **DO NOT** use this device if you have an infection, tuberculosis, or swelling of bone tissue.

N U F A C E®

• DO NOT use if you are pregnant.
• DO NOT use if you have any active cancer or over any suspicious or cancerous lesion.
• DO NOT use this device if you have a history of epilepsy or seizures.
• DO NOT use this device if you have a history of cardiac failure.
• DO NOT use in areas of impaired skin sensation
• DO NOT use if you have hemorrhagic conditions or bleeding disorders
• DO NOT use this device on any area of the body other than instructed in this manual including across the chest, the temples, upper back, and groin area.
• DO NOT use this device directly on the anterior neck, over the carotid artery/sinus.
• Consult your healthcare provider if any of the contraindications above apply to you.
PRECAUTIONS
• If you are under the care of a physician, consult your physician before using this device.
• Apply stimulation only to normal, intact, clean, healthy skin.
• Use caution when using this device near metal implants. locally over areas of swollen, infected, or inflamed areas or skin eruptions, e.g., phlebitis, thrombophlebitis, broken capillaries, varicose veins, etc.
• Use caution when using this device if you have deep-vein thrombosis, impaired circulation, chronic wounds, or skin disease (e.g. Eczema)
• Use caution to avoid applying electrical stimulation near the thyroid gland.
• DO NOT operate FIX+ Device in the presence of electronic monitoring equipment (e.g., heart monitors, heart tracing or in close proximity (e.g. 1 m) to shortwave, microwave or electromagnetic therapy equipment as it may cause the medical equipment to malfunction.
• Use of accessories and cables other than those sold by NuFACE® may result in increased electronic emissions, decreased electronic immunity of the system, damage to the device, or injury.
• Only use the external charging cable provided with the FIX+ Device. Optional 5VAC 1A USB to AC wall adaptors can purchased directly from NuFACE or authorized NuFACE sellers. The optional power adapter is marked NuFACE®. Always position the power adapter so it is accessible in case it is needed to be disconnected from the wall outlet.
WARNINGS
DO NOT apply stimulation or use the device:
• While bathing or showering.
• At extreme temperatures.

NUFACE®

- If you suspect the device is damaged or not functioning properly.
- While driving or operating machinery.
- **DO NOT** disassemble the device as this may cause damage, malfunction, electrical shock, fire, or injury. There are no user-serviceable parts inside. **Note: Disassembly voids the warranty**

Discontinue use of the device and contact your healthcare provider if you experience any of the following:

- Persistent pain or discomfort.
- Discomfort or any difficulty breathing or swallowing.
- Negative or harmful reaction(s) or injury.
- Symptoms of hypertension (dizziness or lightheadedness).
- Constant headaches and/or painful sensations.
- Skin irritation or burning.
- Unpleasant sensations on skin during use.
- Consistent flashing lights without the use of your device.
- If you experience any negative reactions while using the device and/or the microcurrent activator.

BEFORE USING YOUR NEW DEVICE

Follow these important steps before using your new device:


- 1) Remove from the box your FIX+ device, USB charging cable, and download the NuFACE Smart App.
- 2) Connect your FIX+ Device to a USB outlet to begin charging, safely away from any source of water, on a counter or tabletop. The lights will increasingly illuminate, indicating that your device is charging. All of the lights will stay illuminated when device is fully charged. **Note: The FIX+ Device will not turn on when it is charging.**
- 3) **Important: Charge your FIX+ Device for 4 hours before the first use.** The white LED light flashes to indicate that the beauty device is charging. When charging is complete, the white lights will light up continuously.
- 4) **Important: Test your skin to make sure you are not sensitive to the NuFACE Microcurrent Activator.** To do so, rub a dime-sized amount of Microcurrent Activator on your forearm and leave for 15-30 minutes before washing it off. If an allergic reaction (rash, redness, etc) develops, immediately contact NuFACE for recommendations on alternative NuFACE formulas. You should not initiate treatment with the FIX+ Device until your reaction to the Microcurrent Activator has been resolved.

N U F A C E®

- 5) When your FIX+ Device is ON it will emit an energizing hum and the white LEDs will illuminate.
- 6) **IMPORTANT: Always ensure that the two (2) spheres remain in direct contact with the skin.** Avoid pulling or stretching the skin during treatment.
- 7) The FIX+ Device energizing hum will emit 3 rapid pulses during treatment when the battery needs recharging. **Note:** The FIX+ Device will NOT contract the muscles. The current is soft, soothing, and mild when used correctly.

TREATMENT AREAS

The FIX+ Device is designed to smooth and tighten the appearance of fine lines & wrinkles around the lip and eye area.

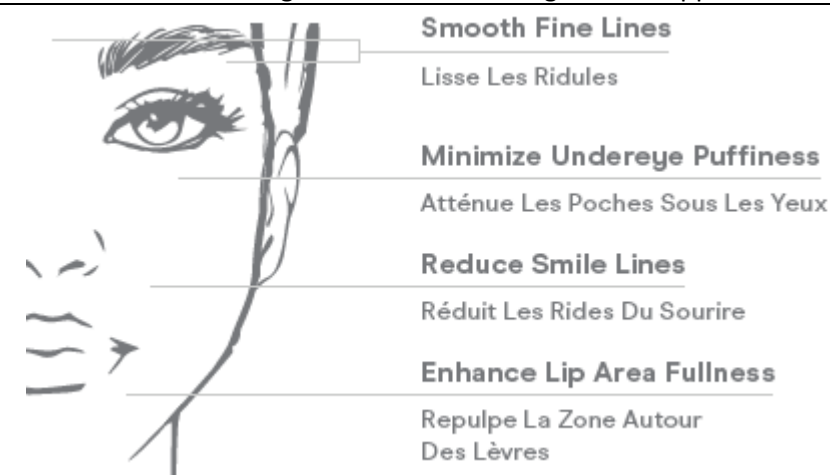
TREATMENT TECHNIQUES


Ensuring that both spheres are in contact with the skin, apply light to medium pressure while using the FIX+ Device feathering technique.

Circular Motion

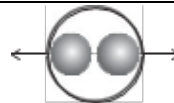Horizontal Gliding Motion

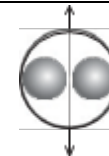Vertical Gliding Motion

NUFACE®

DEVICE INTENSITY INDICATOR		
LED LIGHT INDICATOR	INTENSITY LEVEL	
1 LED	 	Low Level
2 LEDs	 	Medium Level
3 LEDs		High Level
HOW TO USE YOUR FIX+ DEVICE		
Get ready to Firm, Smooth and Tighten.		
The instructions below are designed to provide the most effective results when followed correctly.		
1) CLEANSE Cleanser your face with an oil-free cleanser. Men should shave before treatment — facial hair can interfere with the conductivity of the spheres. Note: Any cream, lotion, self-tanner, and/or oil-based products on the skin can block the microcurrent.		
2) ACTIVATE IMPORTANT: Always apply a NuFACE Microcurrent Activator when using your device. • Apply a mask-like layer of the NuFACE Microcurrent Activator section-by-section as you treat. • Reapply activator if skin becomes dry during treatment in order to continue safely using device.		
3) PREPARE YOUR DEVICE Unplug the USB Cable from your FIX+ Device and turn it ON by pressing the Power Button; the white LEDs will illuminate, and the energizing hum will begin when the device is ON.		

NUFACE®

4) BEGIN TREATMENT

- Place the device in direct contact with the treatment area.
- Always keep both spheres in contact with the skin during treatment to keep microcurrent flowing.
- Using light to medium pressure, begin treatment.
- Begin treatment by using the recommended technique shown.
- You can treat one or more treatment areas for the full 3-minute cycle.
- The FIX+ plus device has three intensity levels. The intensity level can be changed with a quick press of the Power Button.

IMPORTANT: If you feel the electrical current and / or experience a tingling sensation on your skin, do the following:

- Add more NuFACE Microcurrent Activator to the area you are treating. The device will NOT contract the muscles when used as instructed.
- Lower the intensity on the device. Lowering the intensity may also reduce or eliminate this sensation.

5) FINISH TREATMENT

- Your FIX+ Device will automatically turn **OFF after a 3-minute** treatment.
- The device can also be turned OFF by holding the power button while the device is running until the indicator LEDs turn off, and the device stops humming.
- After treatment is complete and device is turned off, massage remaining NuFACE Microcurrent Activator into skin followed by eye cream, moisturizer, and SPF.

Go to the section '**CARING FOR YOUR DEVICE**' for instructions on cleaning and maintaining your device.

TREATMENT FREQUENCY

- The FIX+ Device may be used up to 2 times a day per treatment area.
- 5 days a week.
- Can be used as the finishing step to your NuFACE TRINITY+ or MINI+ treatment.

For more advanced treatment options and holds, download the **NuFACE Smart App**, or visit our website at **www.NuFACE.com**.

RESULTS

Results may vary based on skin type and skin condition. Optimal results are achieved with regular use.

NUFACE SMART APP

NUFACE®

Download the NuFACE SMART APP for additional how-to instructions, and customized treatment options.

- You can download the NuFACE Smart App from the App store or the Google Play Store.
- To pair your device, in the App, click on ‘Add Device’ and follow the on-screen instructions to connect.

FOR PRIVACY AND SAFETY, CONNECT YOUR DEVICE WITH THE NUFACE APP IN A SECURE LOCATION. ENSURE IT IS USED ONLY WITH YOUR PERSONAL SMART DEVICE FOR INDIVIDUAL APPLICATION.

CARING FOR YOUR FIX+ DEVICE

Storage

- The device is intended to be stored and used at room temperature. Use caution prior to using the device if it has been stored or not used for an extended period of time. Inspect the device and do not use it if there are any signs of internal battery leaks.

Cleaning

- Slightly dampen a soft, lint-free cloth with water or rubbing alcohol to clean your device after every treatment.
- **DO NOT submerge your device in water at any time.**

Charging

- Before first use, charge your FIX+ Device for 4 hours.
- The device will NOT turn on when it is charging.
- For continued optimal battery life, it is recommended to charge your device regularly.
- After being fully charged, the device shall operate for up to 40 full 3-minute treatment cycles.
- Keep the product out of the reach of children when charging.
- Keep the product away from high temperature, open flame and water sources when charging to prevent short circuiting or damaging the battery.
- Plug the charger into a wall outlet away from water, and within a location that can be accessed and disconnected easily.

Battery Care

For ideal battery life, we recommend charging your device when the low battery indicator **(1 LED light only)** is displayed on your device and the low battery indicator begins.

Check Battery

N U F A C E®

To check your device battery level, **hold the power button down on your device for one (1) second**. The battery level will be illuminated through the three (3) LEDs:

LED LIGHT INDICATOR		BATTERY POWER
NO LED LIGHTS		No Battery Power A full charge cycle of four (4) hours is recommended before use.
1 LED Flashing		Critical Battery. Not enough charge to start a treatment.
1 LED		Low Battery
2 LEDs		Medium Battery
3 LEDs		Full Battery

Note: The rechargeable battery of the FIX+ Device is not removable by the user. Use caution prior to using the device if it has been stored or not used for an extended period of time. Inspect the device and do not use it if there are any signs of internal battery leaks.

Disposal

IMPORTANT: If disposing the device, do not put into a fire. Dispose of the device and its batteries per local regulations.

TROUBLESHOOTING

Device Not Turning On:

Check to see if your device is fully charged.

- If the device does not turn ON, the battery may be drained.
- A full charge cycle of four (4) hours may be required for the device.
- For additional information go to the “Battery” section in this user guide.

Resetting Your Device:

- To reset your device, plug your device into the power supply.
- Hold down the power button for **10 seconds**.
- The 1st and 3rd LEDs will blink three **(3) times**.
- Once the LEDs stop blinking you may start using your device.

N U F A C E®

With a reset, you will need to reconnect your Bluetooth® pairing and settings located under the NuFACE App section.

DO NOT attempt to troubleshoot by disassembling the product as this will void the warranty and may cause injury.

CUSTOMER SUPPORT

If you are still experiencing issues with your device, contact NuFACE Client Services at info@nuface.com or call **888-360-9171** (US Only).

WARRANTY

NuFACE® warrants that all NuFACE® devices and attachments are free from manufacture defects in materials and workmanship for a period of two (2) years from the date of your original purchase. Our obligation under this Warranty is limited to replacement of the product which will not extend the period of your Warranty. This Warranty does not cover damage covered by abuse, improper care, application of harmful chemicals, product that may have been altered in any way, loss of the device and product purchased from an unauthorized seller.

NO OTHER WARRANTY EXPRESSED OR IMPLIED, WHETHER OF FITNESS FOR A PARTICULAR PURPOSE OR OF MERCHANTABILITY OR OF ANY OTHER KIND, SHALL EXIST IN RESPECT TO SUCH PRODUCT.

For Devices Purchased in the US: For more information about your warranty or to initiate a warranty claim and receive a return authorization number, contact our Client Services Team via our website at www.NuFACE.com or by calling **888-360-9171**. A receipt must be shown from an authorized seller as proof of purchase.

For Devices Purchased Outside of the US, International Returns, and International Warranties: For all product purchased outside of the United States, warranty and returns are per your place of purchase—spa, retailer, or distributor. Please consult your receipt or place of purchase for their return and warranty policy.

TECHNICAL INFORMATION

The NuFACE® FIX+ complies with the requirements of IEC 60601-1-2:2014 (EMC Collateral Standard) including the E-field susceptibility requirements at a 10 volts per meter, at frequencies from 80 MHz to 2.7 GHz. However, even at this level of device immunity, certain transmitting devices (cellular phones, two-way radios, cordless phones, paging transmitters, RFID devices, etc.) emit radio frequencies that could interrupt FIX+ operation if operated in a range too close to the device. Users should be aware of possible radio frequency interference if portable devices are operated in close proximity to the FIX+. It is recommended that transmitting devices such as cellular phones, two-way radios, cordless phones, etc. not be located within 10 feet of the FIX+ during treatment. The NuFACE® FIX+ complies with the requirements of IEC 60601-1-2:2014 for Radiated and Conducted Emissions in Group 1 Class B.

N U F A C E®

The NuFACE® FIX+ complies with the requirements of IEC 60601-1-2:2014 for Harmonic distortion in Harmonic classification A.
The NuFACE FIX+ incorporates Bluetooth® low energy communications at 2.4 GHz frequency single-ended RF with -20 to +4 dBm Effective Radiated power.
FCC Compliance FCC ID: 2AGNA-10800 This device complies with part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) This device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation. NOTE: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures: —Reorient or relocate the receiving antenna. —Increase the separation between the equipment and receiver. —Connect the equipment into an outlet on a circuit different from that to which the receiver is connected. —Consult the dealer or an experienced radio/TV technician for help.
Changes or modifications not expressly approved by NuFACE could void your authority to operate the equipment.
RF Exposure information (SAR) Cet équipement a été testé et répond aux limites applicables pour l'exposition aux radiofréquences (RF). Le taux d'absorption spécifique (SAR) fait référence au taux auquel le corps absorbe l'énergie RF. La limite SAR est de 1,6 watts par kilogramme dans les pays qui fixent la limite en moyenne sur 1 gramme de tissu. Lors des tests, les radios de l'appareil sont réglées sur leurs niveaux de transmission les plus élevés et placées dans des positions qui simulent une utilisation à proximité du corps, avec une séparation de 5 mm

Industry Canada Compliance

IC: 25861-10800

This device contains license-exempt transmitter(s)/receiver(s) that comply with Innovation, Science and Economic Development Canada's license-exempt RSS(s). Operation is subject to the following two conditions:

1. This device may not cause interference.
2. This device must accept any interference, including interference that may cause undesired operation of the device.

CAN ICES-3 (B) / NMB-3 (B)

This Class B digital apparatus complies with Canadian ICES-003.


Conforms to EN IEC/UL: Device - 60601-1, 60601-1-6, 62366-1, 60601-1-11, 60601-2-10, 60601-1-2, 60529, 62133-2, 60335-1, 60335-2-23
Power Supply - 60950-1, 60068-2-70

RF Exposure information (SAR)

This equipment has been tested and meets applicable limits for radio frequency (RF) exposure. Specific Absorption Rate (SAR) refers to the rate at which the body absorbs RF energy. The SAR limit is 1.6 watts per kilogram in countries that set the limit averaged over 1 gram of tissue. During testing, device radios are set to their highest transmission levels and placed in positions that simulate use near the body, with 5 mm separation.

SYMBOLS

One or more of the following markings may appear on the FIX+ Device, User Manual, or packaging.

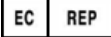	European Conformity
	UK Conformity Assessed
	Nemko Electrical Safety Certification Program
	European Representative
	Part Number
	Operator's Manual; Operating Instructions
	Refer to Instruction Manual/Booklet

N U F A C E®

	Manufacturer
	Complies with European Restriction of Hazardous Substances Directive 2002/95/EC
	Product must be disposed of in accordance with the WEEE directive
	Type BF applied part, which is equipment that provides a particular degree of protection against electric shock. Particularly regarding allowable leakage current and of the protective earth connection.
	FCC-approved equipment authorization mark
	Class II double insulated electrical device
IP54	Protection against dust: Protection against splashed water
	Detachable power supply part
FIX+ DEVICE GENERAL SPECIFICATIONS	
Device Measurements: 154mm/6.06" (L) x 27mm/1.06" (W) x 24mm/0.94" (D)	
Device Net Weight: 56 g / 2 oz	
System Kit Gross Measurements: 197mm/7.75" (H) x 89mm/3.5" (W) x 59mm/2.34" (D)	
System Kit Gross Weight: 108 g / 5.8 oz	
Antenna Peak gain: 2.21 dBi	
Antenna type: Chip Antenna	
Antenna impedance: 50Ω nominal	
Battery Rating: 3.7V, 350mAh	
Max Output Current: 250μA	
Max Output Voltage: 28VDC	
Pulse Duration: 10ms – 60ms	
Pulse Frequencies: 8.3 – 50 Hz	

NuFACE®

NuFACE AC Power Adapter DZ006THU050100H: 14030	
NuFACE USB Cable (5V) to Barrel KC-DZ081A: 10671	
ENVIRONMENTAL CONDITIONS	
Operating conditions: 10 to 40 ° C, 30 to 75%rH*, 700 to 1060hPa	
Storage conditions: -10 to 55 ° C, 10 to 90%rH*, 500 to 1060hPa	
Shipping condition: -10 to 55 ° C, 10 to 90%rH*, 500 to 1060hPa	
	
Manufactured by NuFACE® 1325 Sycamore Avenue, Suite A, Vista, CA 92081 1.888.360.9171 www.MyNuFACE.com	
©2024 NuFACE®	
All Rights Reserved.	
NuFACE® and FIX+® are registered trademarks of the Carol Cole Company DbA NuFACE.	
Device designed and Manufactured in California.	
Microcurrent Activator made in the USA from U.S. and Imported Materials.	
  	
Medical Device Safety Service GmbH	
Schiffgraben 41, 30175 Hannover, Germany	
Uk Responsible Affiliate: MDSS-UK RP	
6 Wilmslow Road, Rusholme, Manchester	
M14 5TP United Kingdom	
BIORIUS / EU: WAVRE, BE / UK: LONDON WC2H9JQ, GB	NuFACE®

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