PROJECT NAME GLOWCAP INSERT DATE 05/16/2016

GLOWCAP SETUP GUIDE

Material 40#, BRIGHT WHITE FINISH FRONT GLOSS/BACK GLOSS

PAGE SIZE 215.9 × 279.4 mm VIEW FRONT VIEW, BACK VIEW

FRONT VIEW

NantHealth®

GlowCap® Setup Guide

* RECEIVE EMAIL

Retrieve your welcome email from NantHealth. This email contains your username and password.

DOWNLOAD APP

Using your Bluetooth®-enabled mobile device, download the NantHealth DeviceConX Mobile™ app from either the Google Play™ store or Apple App Store®.



PAIR DEVICES

Once the application is installed, launch it and log in with your username and password.

- 1. Select Add Device
- 2. Select GlowCap
- 3. Open and close the GlowCap
- 4. Select your GlowCap once it appears in the available list
- 5. If desired, rename your GlowCap
 Note: Each GlowCap should have a unique name.






Access the URL provided in your welcome email by clicking it, or by copying and pasting it into your browser's navigation bar, to review your information, track your usage, and make changes to your schedule.

Technical Support Contact Information

Support Line Hours: Monday through Friday 8am to 7pm (Central) Phone: 1-844-626-8935 (NANTWELL)

*Before using your GlowCap, review the safety information and precautions on the back of this instruction sheet.

SUPPORT: 1-844-626-8935 WWW.NANTHEALTH.COM

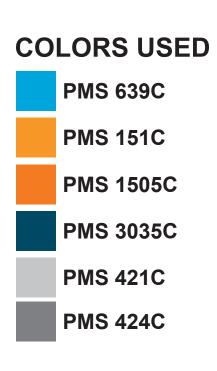

BACK VIEW

USER INFORMATION

Compliance Statement (Part 15.19)

The enclosed hardware device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) This device may not cause harmful interference, and (2) This device must accept any interference received including interference that may cause undesired operation.

Warning (Part 15.21

Changes or modifications not expressly approved by NantHealth® could void the user's authority to operate the equipment. Manufacturer is not responsible for any radio or TV interference caused by unauthorized modifications to this equipment.

Compliance Statement (Part 15.105(b))

This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following

- Reorient or relocate the receiving antenna
- Increase the separation between the equipment and receiver
 Connect the equipment into an outlet on a circuit different from that to which the
- Connect the equipment into an outlet on a circuit different from that to which the
- Consult the dealer or an experienced radio/TV technician for help

Warnings and Precautions

- It is important to take all medications (prescribed and over-the-counter) in accordance with your healthcare provider's exact instructions. Please consult your healthcare provider before taking any medications and consult the drug manufacturer's guidelines to ensure your safety and well being.
- Please consult your healthcare provider before using GlowCap.
- Please make sure that the information entered in your GlowCap account is accurate and up-to-date.
- Do not load the GlowCap with any medications not prescribed to the intended end user.
- Do not pour liquid/gel medications directly into the GlowCap.Please make sure the GlowCap is properly labeled.
- Please make sure that the GlowCap lid is properly secured.
- Please keep the GlowCap away from children and pets.
- Do not permit unauthorized or untrained individuals to program, load, or use the GlowCap.
- Do not misuse, abuse, alter, neglect, or attempt to repair the GlowCap. Any tampering or modification to the GlowCap will void the limited warranty.
- \bullet Do not reuse the GlowCap for other medications or prescriptions.
- $\bullet \ \, \text{Please report any malfunction of the GlowCap to NantHealth at 1-844-626-8935 or visit www.nanthealth.com}. \\$

Specifications

a. Device Power Requirements

- 3 VDC (Lithium Battery type CR2447)
- b. Environmental RequirementsOperating Temperature: 0° to 40°C
 - Humidity: 0% to 95% (non-condensing)
 - Indoor use only, dry locations only

Caution

Risk of explosion if battery is replaced by an incorrect type. For battery replacement, call NantHealth at 1-844-626-8935.

Additional Information:

For additional questions regarding your GlowCap, or for additional product information, visit www.nanthealth.com or call 1-844-626-8935.

6 NantHealth®

GlowCap® Setup Guide

*

RECEIVE EMAIL

Retrieve your welcome email from NantHealth. This email contains your username and password.

2

DOWNLOAD APP

Using your Bluetooth®-enabled mobile device, download the NantHealth DeviceConX Mobile™ app from either the Google Play™ store or Apple App Store®.



3

PAIR DEVICES

Once the application is installed, launch it and log in with your username and password.

- 1. Select Add Device
- 2. Select GlowCap
- 3. Open and close the GlowCap
- 4. Select your GlowCap once it appears in the available list
- 5. If desired, rename your GlowCap
 Note: Each GlowCap should have a unique name.



4

LOGIN

Access the URL provided in your welcome email by clicking it, or by copying and pasting it into your browser's navigation bar, to review your information, track your usage, and make changes to your schedule.

Technical Support Contact Information

Support Line Hours: Monday through Friday 8am to 7pm (Central)

Phone: 1-844-626-8935 (NANTWELL)

SUPPORT: 1-844-626-8935

^{*}Before using your GlowCap, review the safety information and precautions on the back of this instruction sheet.

USER INFORMATION

Compliance Statement (Part 15.19)

The enclosed hardware device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) This device may not cause harmful interference, and (2) This device must accept any interference received including interference that may cause undesired operation.

Warning (Part 15.21)

Changes or modifications not expressly approved by NantHealth® could void the user's authority to operate the equipment. Manufacturer is not responsible for any radio or TV interference caused by unauthorized modifications to this equipment.

Compliance Statement (Part 15.105(b))

This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- · Reorient or relocate the receiving antenna
- Increase the separation between the equipment and receiver
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected
- Consult the dealer or an experienced radio/TV technician for help

Warnings and Precautions

- It is important to take all medications (prescribed and over-the-counter) in accordance with your healthcare provider's exact instructions. Please consult your healthcare provider before taking any medications and consult the drug manufacturer's guidelines to ensure your safety and well being.
- Please consult your healthcare provider before using GlowCap.
- Please make sure that the information entered in your GlowCap account is accurate and up-to-date.
- Do not load the GlowCap with any medications not prescribed to the intended end user.
- · Do not pour liquid/gel medications directly into the GlowCap.
- Please make sure the GlowCap is properly labeled.
- · Please make sure that the GlowCap lid is properly secured.
- Please keep the GlowCap away from children and pets.
- Do not permit unauthorized or untrained individuals to program, load, or use the GlowCap.
- · Do not misuse, abuse, alter, neglect, or attempt to repair the GlowCap. Any tampering or modification to the GlowCap will void the limited warranty.
- Do not reuse the GlowCap for other medications or prescriptions.
- Please report any malfunction of the GlowCap to NantHealth at 1-844-626-8935 or visit www.nanthealth.com.

Specifications

a. Device Power Requirements

3 VDC (Lithium Battery type CR2447)

b. Environmental Requirements

- Operating Temperature: 0° to 40°C
- Humidity: 0% to 95% (non-condensing)
- · Indoor use only, dry locations only

Caution

Risk of explosion if battery is replaced by an incorrect type. For battery replacement, call NantHealth at 1-844-626-8935.

Additional Information:

For additional questions regarding your GlowCap, or for additional product information, visit www.nanthealth.com or call 1-844-626-8935.