

正



- Do not modify, dismantle, or repair its S1 toothbrush (hereinafter referred as "the Product") on your own.
- Do not soak the charging cradle in water or any other liquid.
- After cleaning the Product, make sure the plug of the charging cradle has been dried completely before connecting the power.
- The Product has a built-in rechargeable battery, so do not throw it into fire, heat it, or charge, use, or place it in high-temperature environment.
- Do not knock or drop the Product and do not use it any more if it has been damaged due to collision or fall.
- Keep the product far from electromagnetic field, such as induction cooker when it is being charged.

- Check for **visual damages** each time before using the product. Do not use it if it has been damaged because of drop or any other reason.
- Before connecting the product to power, check whether the local voltage is equal to the rated voltage of product.
- If the charging cradle has undergone any damage, do not use it.
- Do not use the charging cradle outdoors or close to high-temperature surfaces.
- If any breaking sign is found on the product including brush head and handle, do not use it any more.
- Do not use the product when bathing or showering.

If any operation was made on your mouth or gum in the past two months, contact your dentist for advice before using this product.

Your teeth are bound **to** keep using this Product or the **bleeding** lasts over one week, contact your dentist for help.

If a heart pacemaker or other device has been implanted into your body, your surgeon or the device manufacturer before using this Product. Do not use a brush head that is not recommended by the manufacturer.


This Product can be used only **to** clean the teeth, gum and tongue. Do not use it for other purpose. If you feel uncomfortable or painful, stop using it and go to see a doctor for help.

If you use the toothpaste containing hydrogen peroxide, sodium bicarbonate, or bicarbonates (usually "whitening toothpaste"), clean the brush head completely with soap and water after brushing tooth every time.

- Do not place any metal object, for example, coin or clip on the charging cradle.
- Do not twist the power lines around the adapter tightly for storage.
- Unplug to power the Product off if it will not be used for a long time, preventing fire or other hazard.
- Do not share one toothbrush with other family members or other people.

Note: Although your teeth are in good condition, tingling or bleeding may happen at the beginning when you use this Product, because your gum has not been accustomed to the stimulus for the first time. Generally, the tingling or bleeding will disappear in one or two weeks, otherwise your gum may have something wrong, and you should stop using this Product and contact a dentist for help.

The diagram shows three items: a protective case for travel, a toothbrush head, and a full toothbrush. The protective case is a simple rectangular box with a flap. The brush head is a rectangular head with bristles. The full toothbrush is a standard manual toothbrush with a handle and a head. The text 'General type (standard configuration)' is written below the full toothbrush.

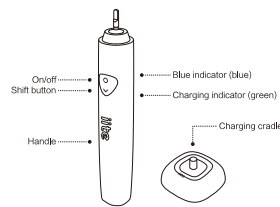

Perform the following steps before the first use:

1. Place the charging cord at a safe, dry, and shadow location to avoid **falling** into water or other liquid to prevent electric shock.
2. Connect the power cable.
3. Check the battery status: Refer to the charging cord with caution. If the charging indicator goes on, the Product is being charged.

Notes:


1. Change the Product for 24 hours before the first use.
2. It is recommended to use the Product in temperature ranging from 5°C to 35°C, and moisture ranging from 30% RH to 85% RH.
3. Charge the Product and use the following guidelines:
 - a. Before the first use.
 - a. To avoid not being used for at least three months.
 - a. The bottom can't be placed on a wet surface used.
 - a. To be used after affected by low battery during usage, please do the Product on the charging cord with caution.
4. With the increase of charging times, the Product can be used for longer times after each charging. If the Product is fully charged and out of power after being used only once, the battery is damaged.

	Over 20%	20%~1%	Less than 1%
State of Green LED	Off	Flash slowly (every 1s)	Flash quickly (every 0.2s)


5.5cm

反

For daily use, perform the following steps:

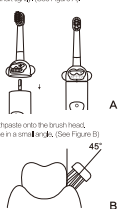

- A) Move the brush **along** the teeth **slowly** and **gently**, and do not use excessive force.
- B) **Starting** **below** the **hand**, **brush** **lightly** with **finger**, put the brush **into** mouth **now**, **switch** the brush **on**, keep the brush in a moderate **angle**, and **place** it **on** the **incisal** connecting teeth **first**. Rotate the brush **head** **slightly** on a fixed location for **several** seconds and then change to another location and rotate again. Perform the same operations on every part of teeth without any omission, including the internal surface and the chewing surface, as well as the new area. Make the brush **along** the gum line, no press or rub. This is the whole process of brushing teeth.
- C) Clean the brush head and hands. Wash the brush head and handle with clear water and knock the brush head **lightly** to remove food debris. (See Figure C)
- D) Put the brush head as shown in Figure D and wave the head to make it dry. (See Figure D)
- E) Put the hands back on the charging cord. If the charging indicator goes on, the Product is being charged.

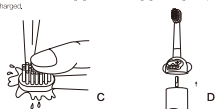

1. This Product is charged in inductive manner. Put the brush on the stereo-configuration charging cradle that has been connected to power and then the charging begins.
2. Charging indicating in the process of charging, the LED indicator flashes slowly, and the lighting when the battery has been fully charged. The charging takes 24 hours for the first time and 15 hours each time later.
3. When the charging begins, the brushing mode is shut off, the multi-function button is disabled (can not be pressed), and the motor stops rotating.
4. During charging, if the handle leaves the charging cradle, the Product will start standard operation. The LED indicator flashes rapidly (not rotate).
5. When the Product is switched on and the handle is placed on the charging cradle, the motor stops rotating and the charging mode begins.

	Charging	full charge
State of green LED	Flash slowly (every 1s)	Keep lighting

	Not bound APP	Bind APP	Bind APP success
State of the blue LED	Off	Flash quickly (every 1s)	Keep lighting

1. Charging cradle
 - 1) Wipe the charger with soft cloth.
 - 2) Wipe the metal part of the plug with dry cloth every 6 months to avoid poor contact caused by accumulated dust.
2. Brush head
 - 1) Take the brush head off the handle.
 - 2) Wash the brush head and brush with clear water.
 - 3) Take the brush head off the handle and wash the connection part at least once a week.
 - 4) Reassemble.
 - 1) Take the brush head off the handle.
 - 2) Wipe the outer surface of the handle with soft and wet cloth and put it back on the charging cradle.

Name: Ite S1 Kid's Super Toothbrush
 Model: Ite-BS1M
 Rated voltage: 200~240V/AC
 Rated frequency: 50 Hz/60 Hz
 Rated input power: 2W
 Product standards: GB4706.1-2005, GB4706.55-2008
 GB4706.56-2008

Around three years. If the number of times that the battery can be used after being fully charged has reduced significantly, the battery life has possibly come to an end. (The battery life varies according to usage and storage conditions.)

Faults	Possible causes	Corrective measures	
The indicator light is not lighting during charging	Product is made, bought or transported incorrectly by customer	At the beginning of charging, the indicator light is used to indicate whether the product is in a normal state. If the indicator light is not lit, it indicates that the product is not in a normal state. The indicator controller is faulty, contact the manufacturer or the supplier of the maintenance unit, and submit the problem to the manufacturer or the supplier of the maintenance unit, and submit the problem to the manufacturer or the supplier of the maintenance unit.	
The product battery can only stand for about 10 minutes after being fully charged	Insufficient charging time	Charge for over 10 hours	
The battery does not come to use	The battery is not connected	I need a repair, contact the store where you buy the product or the manufacturer, and submit the problem to the manufacturer or the supplier of the maintenance unit.	
The battery cannot support use	Accessories (tools) are used or not used	The battery is damaged or broken	Contact the store where you buy the product or the manufacturer, and submit the problem to the manufacturer or the supplier of the maintenance unit.
The battery does not work and the indicator does not light after being fully charged	The battery has been damaged	Contact the store where you buy the product or the manufacturer, and submit the problem to the manufacturer or the supplier of the maintenance unit.	



Warranty

Dear life toothbrush users:
Thanks for using our life moto™-driven toothbrush. It is our honor to provide the optimal service to you. To ensure satisfactory service, please read the following terms and keep the warranty (Valid only in PRC and Hong Kong, Macao and Taiwan are not applicable) carefully upon purchasing our product.

[illegible]


说明书2

正



iite S1 quick start guide
www.iite.cc

Components

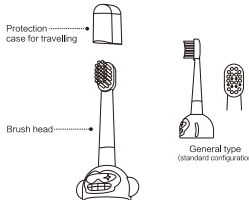


Protection case for travelling

Brush head

General type (standard configuration)

Body

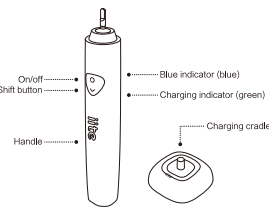

On/off Shift button

Blue indicator (blue)


Charging indicator (green)

Handle

Charging cradle

How to use


Step 1: Put the motor-driven brush on the charging cradle for 24 hours before the first use.

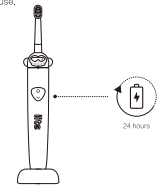


24 hours

Note: This product applies only to children above 4 years old.

Step 2: Cover the brush head on the shaft tightly.

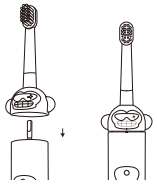


Step 3: Wet the brush head and squeeze some toothpaste onto the head.

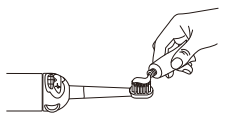


反

Step 4: Place the brush head on the line connecting teeth and the gum.

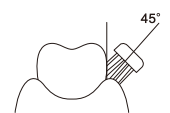


45°

Step 5: Switch on and begin to clean your teeth. Make sure every part is cleaned.

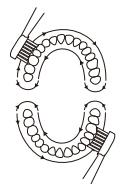


Step 6: To avoid sprinkling of the toothpaste, switch off and take the brush out of your mouth.

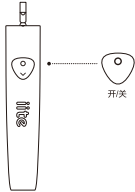


开关

Step 7: Wash the brush head with flowing clear water after brushing teeth, wave the brush head to make it dry, and wipe the handle with soft cloth.

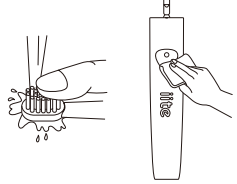


Step 8: Put the toothbrush back on the charging cradle.



Note: If the brush has not been used for a long time, unplug and clean the plug and then store it in a shadowy and dry location with no direct sunlight.

Download the iite smart toothbrush App as follows:

a. Swipe the following QR code.



b. Visit www.iite.cc and download the App.

c. Search for iite in Apple's App Store, SnapRes, or 91 Mobile Assistant and download the App.

FCC Statement

This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

This device complies with part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) This device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

Changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.