Instructions of intelligent falsies

Model: BP

Brand: Z private beauty



Brief Product Introduction

----Brief introduction of intelligent falsies


intelligent falsies realizes intelligent control with a mobile phone APP operation by connecting with mobile phone APP, utilizing phototherapy function of frequency conversion and microseism, and integrating the professional courses of "physiotherapy for a healthy breasts", "breasts enhancement care" and "breasts shaping" in APP, which has developed a comprehensive new area for intelligent breasts beauty.

----Function of intelligent falsies

By microwave vibration, it accelerates cell metabolism of breasts, activates the sleeping cells, and has the breasts cells recovered active status. Besides, it increases the elasticity of ligament cells, shapes the breasts and reveals a healthy and charming posture.

Physiotherapy to the breasts is conducted with red light, yellow light, blue light and other different colors, which is for the increase of metabolism and circulation of breasts cells, for valid activation of adipose cells, for stimulation of the generation of collagen, and for shaping of perfect breasts. Phototherapy has very strong penetrability, which can directly reach to corium layer, can rejuvenate the aging cells, and can have the breasts persistently maintained in a healthy states.

---Appearance of intelligent falsies (as picture below)---

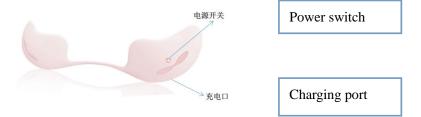

Note (picture above is only sectional drawing, not design drawing, which need a new drawing to show according to the picture)

Best ergonomics design, U-typed chest protector

Operation Description

-----Function introduction of intelligent falsies

Charging: Insert the contact pin of supporting battery charging line into the charging port of falsies, and connect the USD interface of the other side with a computer or a charger to fulfill charging. The indicator light will blink in the charging process, and the indicator light will keep light when fully charged.



Note (the above is only a picture of real products, not design drawing, which need a new drawing to show according to the picture)

Power on: Press the on-off button for three seconds. The indicator light shall quick flash three times as hints.

Power off: Press the on-off button for three seconds. The indicator light shall quick flash three times as hints.

Standby: Short press it for a one-time flash, which indicates the product is in a standby condition.

-----Downloading and Installation-----

Downloading and installation methods for Android mobile phones

Search "Simei(the meaning of "private beauty")" in Yingyongbao(application box), Baidu assistant, Peasecod, 91 assistant and other major application markets for downloading and installation.

Scan the two-dimension code to skip into the application market, and then download and install the APP according to the hints.

Downloading and installation methods for IOS iphones

Search "Simei(the meaning of "private beauty")" in the APP Store for downloading and installation.

Scan the two-dimension code to skip into the APP Store, and then download and install the APP according to the hints.

--- Registration---

Registration method: register---input your cell phone number---receive the verification code---input the verification code received by your cell phone---click the next step---set a password---complete the registration

Logging methods:

Logging in by account number: input your account number---input your password---click the logging

Logging in by a third party: Choose a direct logging in with the QQ, Sina Microblog and Wechat blow the interface.

----Bluetooth Connection---

Before the connection, please keep the intelligent falsies BU1.0 on a power-on condition, and the bluetooth of cell phone on a open condition. If the bluetooth of cell phone is off, the system will automatically inform to open the bluetooth of cell phones, which need to be opened according to the hints.

Operation methods of bluetooth connection: Click the "intelligence" in the navigation bar below the home page of dynamic condition to enter into the main

interface of products---click the bluetooth connection symbol on the top right corner of the main interface of products, and search the bluetooth name of intelligent falsies in APP---click the bluetooth name of intelligent falsies for connection. Successful connection: a $\sqrt{}$ hint before the bluetooth name of intelligent falsies. When bluetooth is successfully connected, the bluetooth symbol will show a picture as

. When bluetooth is disconnected, the bluetooth symbol will show a grey picture



----Control Mode----

Under the condition of successful connection betweenSimei(the meaning of "private

beauty") APP with intelligent falsies BU1.0, click the on-off button on the control interface of Simei(the meaning of "private beauty") APP to control the massage function and phototherapy function of intelligent falsies

Operational methods:

Enter into the course selection page, and click the "intelligence" in the navigation bar below the interface---enter into the course selection page and choose "physiotherapy for a healthy breasts", "breasts enhancement care" or "breasts shaping" to enter into relevant course.

Choose the massage or phototherapy mode on the top of interface after entering into course control interface. The on-off button remain grey before choosing this function, which shall turn purple after the selection.

Massage mode:

There are three-typed buttons for massage mode in the below of massage control interface. Click the button to start the massage, then the falsies will relatively vibrate-a second click will stop the massage, and the falsies will stop vibrating. Massage mode remains grey before the starting up, and turns into purple after the opening, besides, the button will become larger at the same time. The three-typed massage modes below the control interface can be switched at a free will by clicking.

Massage strength ajusting:

Drag the slide symbol of the strength set in the control interface to adjust the massage strength. The massage strength is the weakest when the slide symbol is in the far left "-" position, and the massage strength will become stronger when the slide symbol is dragged into the "+" position in the right direction.

Phototherapy mode:

Select phototherapy, and enter into the phototherapy control interface. There are three light buttons below the interface, which are "red light", "yellow light" and "blue light". Click the buttons to start phototherapy (falsies will deliver the relevant light), and a second click will turn off phototherapy(the light of falsies will disappear). Massage mode remains grey before the starting up, and turns into purple after the opening, besides, the button will become larger at the same time. The three-typed modes below the control interface can be switched at a free will by clicking.

Massage function and phototherapy can be used both in separate and at the same time. It can be switched to phototherapy mode while massaging, and to vibration mode in the process of phototherapy mode.

Click the return button on the left top of interface to get back to the main interface of products when stopping the exercise---click the bluetooth mark on the top right cornor---click the bluetooth name to break bluetooth connection.

--- Targeted Group of Intelligent Falsies

The intelligent falsies is for groups who have sub-healthy breasts, like hyperplasia of mammary glands, fibroid of mammary glands, galactoma, slack breasts, mastoptosis, external expansion of breasts and flat breasts.

-- Usage Methods of Intelligent Falsies

Usage steps

Clean the intelligent falsies with warm water or cleaning articles.

Keep the intelligent falsies on.

Open the installedSimei(the meaning of "private beauty") APP in the mobile phone, and start physiotherapy after selection of the necessary courses.

Turn off the control button on the control interface after completion of physiotherapy to stop the working of intelligent falsies Then break the bluetooth equipment, and quit the control interface.

Proposal use time

The massage mode and phototherapy mode can be used above thirty minutes respectively each day. It can be used in a long term without starting the function of intelligent falsies, which acts the role of shaping, promoting and locking the fat of breasts.

Maintenance notice of intelligent falsies

Please clean with cleaning articles or warm water before the initial usage. The effect will become better with a synchronized combination usage of breast care products. When using the intelligent falsies with breast care products, you need first take appropriate breast care products, smear it evenly on the breasts, and then massage the breasts in circle until total absorption, with the combination usage of intelligent falsies, you will have a comprehensive experience and care.

Put it in a dark and dry location after cleaning up and charging.

Please do not use damaged charging lines, USB and socket. Please do not place it into fire or heating equipment(like microwave oven or oven etc.)

In order to guarantee the service life of product, please do not soak the products into water while washing, and try to avoid the charging port during the process of cleaning.

In order to better ensure the service life of intelligent falsies, please charge the intelligent falsies to the full regularly(at least one time each week) if the intelligent falsies has not been use for a long time).


In order to give you a better experience, please charge the intelligent falsies to the full before use. The intelligent falsies will not work during the charging process.

In order to provide you a better effect, a timely switch of intelligent falsies after a 6-month usage is advisable.

Wearing methods:

Tuck the concave surface of intelligent falsies facing the direction of skin into underwear. The intelligent falsies will be next to the lower limit outside breasts. Adjust the intelligent falsies into a comfortable position.

Pictures of wearing methods for intelligent falsies

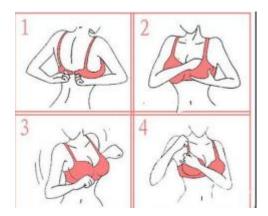

(note: all the pictures like these shall be drawn with lines, which need a communication in person with Sun Xianling)

Precautions for use

Forbidden for seriously injured breasts.

Forbidden for breast cancer.

Forbidden for cardiac patients.

Adaptive mobile phones:

Apple: iPhone 4S, iPhone 5, iPhone 5C, iPhone 5S, iPhone 6, iPhone 6 Plus

GALAXY: GALAXY S3, GALAXY S4, GALAXY S5, GALAXY Note 2, GALAXY

Note 3

Google: Nexus 4, Nexus 5, Nexus 6

Meizu: MX3, MX4, MX4 pro Mi: Mi 3, Mi 4, Mi note

And most mobile phones with Android system which are issued in 2014 and

afterwards.

the company: Shenzhen Beautifying Technology Co., Ltd.

address: 5A-066, NO. 2 Trading Plaza, South China City, Pinghu Street, Longgang District, Shenzhen, China

CE Maintenance

- 1. Risk of explosion if battery is replaced by an incorrect Type . Dispose of used batteries according to the instructions.
- 2. The product shall only be connected to a USB interface of version USB2.0.
- 3. EUT working temperature: -20 °C to 40 °C

4. Adapter

Input:100V-240V,50/60Hz 0.2A (EU Plug) Output: DC5V,500mA (USB Port)

Declaration of Conformity

Shenzhen Beautifying Technology Co., Ltd. hereby declares that this intelligent falsies is in compliance with the essential requirements and other relevant provisions of Directive 1999/5/EC. a copy of the original declaration of conformity may be found or obtained at http://www.beautifying.com/

FCC Caution.

This device complies with part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) This device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

Any Changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

Note: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- -Reorient or relocate the receiving antenna.
- -Increase the separation between the equipment and receiver.
- -Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- -Consult the dealer or an experienced radio/TV technician for help.