



OcClean Smart Sonic Electric Toothbrush

User Manual

Thank you for using OcClean X Smart Sonic Electric Toothbrush.

Important Safety Tips


- It is recommended to use this product by users who are at least 12-year-old or adults.
- Due to the high-speed vibrations, it is normal to experience a tingling sensation while using the toothbrush for the first time. This feeling will gradually go away by regularly using the toothbrush.
- It is a normal situation for the toothbrush to heat up while charging. If the toothbrush becomes extremely hot, stop charging immediately and power off, and then contact OcClean customer service center for checking.
- The toothbrush can only be used with OcClean special head.
- The toothbrush is for oral use only, like brushing the teeth, cleaning the tongue or massaging gums. Do not use it for other purposes.
- The toothbrush is a personal item. It is not intended for use in the dental clinic and should not be used for multiple patients.
- The toothbrush can be used indoors only. It should not be placed under direct sunlight for a long period of time. Also, it should be stored in a cool and dry environment.
- Please cut off the power supply when the charging base is not in use or when it is cleaned. Clean the power plugs regularly to prevent dirt as well as in order to avoid reducing the efficiency of charging.

People with Sensitive Oral Cavity

The following people should use this product with great caution.

- People who have undergone oral periodontal surgery or are being treated.
- People who have mouth diseases, loose false teeth or gums with hemorrhage.
- People with sensitive and susceptible oral cavity, they should consult a professional dentist before using this product.

The toothbrush is not suitable for people with reduced mental and physical capabilities. Not suitable for infant and children under the age of 12. Children who use this product should be strictly supervised by an adult.

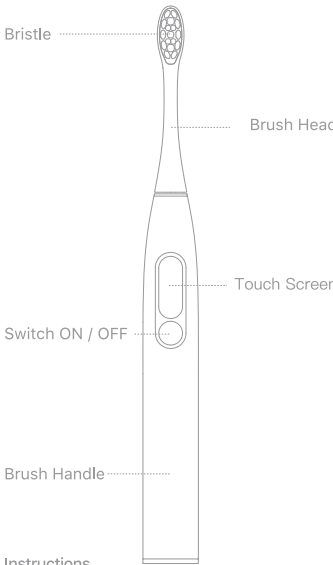

Instructions

- First Time Use**
The bluetooth is on for a long time after the first time use with pressing the Switch ON/OFF button for the brand new toothbrush. After downloading the OcClean App, please enter into App to complete the connection and settings of toothbrush.
- Install App & Toothbrush Settings**
Scan QR code to install App and customize your own special brushing plans through App. After that it will complete the whole process of tooth cleaning, tooth whitening and gum massaging automatically in one brushing. No need to open OcClean App in daily brushing process but it is available to open App in every few days to sync the brushing data to get the bonus points for accessories exchange.
- Direct Use without App Settings**
It is available to brush directly without connection with toothbrush. When the toothbrush pops up the tip of "Connect Toothbrush", please choose "Skip". It is supportive to touch control by touch screen and Switch button of X1 smart toothbrush.
- Automatic Awake Mode & Automatic Sleep Mode**
When the toothbrush is removed from the holder, it will feel the action and wake up automatically. When there is no any operation for 8 seconds after brushing, it will automatically enter into sleep mode.
- Default Brushing Plan & Direct Brushing**
Press Switch button directly when the toothbrush is woken up, it starts working with default mode" Standard Cleaning Plan "with soft intensity and lasts for 2 mins. You can also change the intensity and mode in the settings.


- Record Time & Change Brushing Area**
Toothbrush will start default mode with Standard Cleaning Plan with 2 mins if no new customized brushing plan setting in App.

There will be a reminder with vibration to remind user to change the brushing area in every 30s. Toothbrush will stop operating automatically after 2 mins.

- Over-voltage and Frequency Reduction**
In order to prevent excessive pressure in the process of brushing teeth, this product will automatically reduce the frequency of operation to protect enamel and gums and prevent excessive brushing. This frequency change can be sensed by the sound of motor operation, and the display screen will display too large icons. This is a protection mechanism instead of motor fault, but also can turn on or off this function in the settings.

- Brushing Report**
After brushing, the toothbrush's display will show the effective brushing time, the brushing points, and the tips of the brushing area. You can intuitively see the insufficient brushing time and intensity through the schematic diagram.

- Brushing Mode Setting**
Enter advanced setting by sliding up and down the display. In Brushing Mode Setting, there are 4 modes with 32-level intensity can be chosen, such as Cleaning Mode, Whitening Mode, Massaging Mode. After the new setting is completed, it will be automatically remembered and the toothbrush will run it for the next cycle.

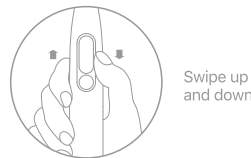

Swipe up and down



Swipe left and right



Click



- Weather Icon**
There will no weather information on the display if the brushing data of toothbrush doesn't sync to App over 7 days.



- Bluetooth Icon**
The Bluetooth icon will show up on the display when the toothbrush connects to App



- New Message Icon**
The new message icon will show up on the display when there is new activity information, software update notification, or other prompt information.



- Cleaning Intensity Reference Value**
The cleaning effect is in the range of high, medium and low. This range is an evaluation of the selected or customized brushing plan. It may be given a low value if the brushing time is too short. The cleaning intensity value will change accordingly by adjusting the brushing time or brushing strength.



- Brushing Points**
It is available to convert the scores into points after each time brushing. It also supports to transfer the points to App by synchronizing the data. There is no points limit in App but it only can be stored 999 points in the toothbrush. The new points will become invalid if the previous points up to 999 and no sync to App.



- Low Battery Warning**
Please charge the toothbrush if a low battery pop-up warning message appears on the display. The display will show the percentage of power when the toothbrush is charging. It will show up the 100% battery power on the display once toothbrush fully charged.

Brushing Plan & App Settings


You can either choose or customize different brushing plans through the OcClean App. You can also check the exact battery power level and you can upgrade the firmware.

Charging

Please take the bottom of toothbrush on the top of magnet charging base. The display will show the percentage of power when the toothbrush is charging. Fully charged usually around 2.5 hours. It will take over 3 hours to get fully charged if you use a low electric current or connect to computer for charging. Recommended adapter: 5V/1A – 5V/2A

Important Charging Tips

The toothbrush runs with a fast charge mechanism. It is recommended to put the toothbrush on the holder when is not charging. Please don't put the toothbrush on the charging base for a long time. It may cause the built-in battery to decay prematurely and affect the user experience.

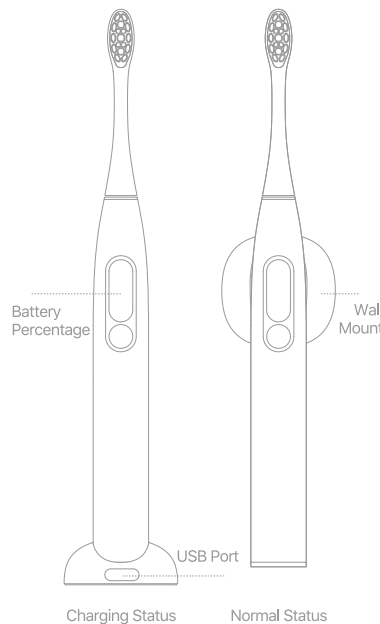



Charging Status Normal Status




- No disassembly, impact, extrusion or throwing into fire
- Please stop using the toothbrush if the battery has expanded
- Please don't put the battery in high temperature environment
- Do not use the battery if water entered the device.

Recommended Tooth Brushing Method


- Tooth Brushing Time**
Oral experts (Dentists) recommended that each brushing should be 2 to 3 minutes long.

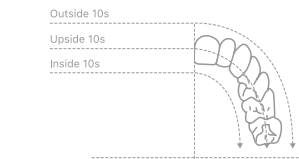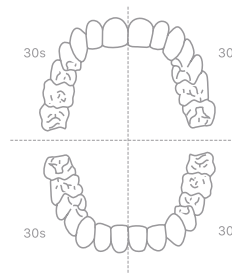

Recommended Tooth Brushing Method

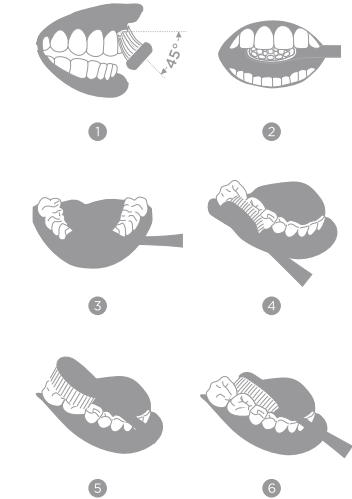

Direction

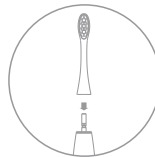


Insert the brush head carefully into the body of the toothbrush. A small gap should be left between the head and the main body for work space so the brush head can work normally.



To prevent toothpaste from splashing, smear a small amount of toothpaste on the head before using it. (A sonic electric toothbrush saves more toothpaste than a manual toothbrush.)



Soak the toothbrush head and put it into your mouth, then turn on the toothbrush to start brushing.

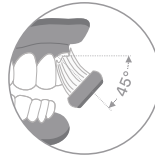

Slowly move the bristles at 45° to the gum line to brush teeth properly.




Hold the toothbrush handle and slowly move it in circular motion. Do not press too hard.



Please rinse with clean water after brushing.


Maintenance

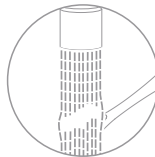


Cleaning Brush Head
Clean the brush head after use and air dry it. Do not pull the brush head for this act as it may result in deformation of the toothbrush.



Cleaning Body
Use clean water to wash the surface of handle and scrub it with soft cloth and air dry it. Do not press the rubber sealing ring with sharp item or metal object, otherwise it will be damaged.

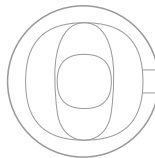

Cleaning Charger
Scrub the dirty parts of the charger with a soft cloth. Pull out the battery plug from the socket prior to cleaning the charger. Scrub the surface of charger with wet cloth.

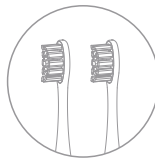

Change the Head
It is suggested by dentists that the head should be changed every 2-3 months to maintain optimal brushing.

Specifications	
Product Name	OcClean Smart Sonic Electric Toothbrush
Type	OcClean X
Magnetic Power	220gf.cm
Power of Subject	5W
Charger	OcClean USB Cradle
Input	5V / 1A
Service Time	20 Days (*)
Battery Capacity	800 mAh
Charging Time	2 Hours
Bluetooth	Bluetooth 4.0 above

* Validated by 2 minutes brushing period, twice a day.

FCC Statement

This device complies with Part 15 of the FCC rules. Operation is subject to the following two conditions: 1) this device may not cause harmful interference, and 2) this device must accept any interference received, including interference that may cause undesired operation.

Note: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
 - Increase the separation between the equipment and receiver.
 - Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
 - Consult the dealer or an experienced radio/TV technician for help.
- Changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

FCC Radiation Exposure Statement

This device complies with FCC RF radiation exposure limits set forth for an uncontrolled environment. This transmitter must not be co-located or operating in conjunction with any other antenna or transmitter.

Remark

- The content of toxic and harmful substance in all homogeneous materials of this component is less than the limited quantity as specified in GB/T 26572-2011;
- The content of toxic and harmful substance in all homogeneous materials of this component is less than the limited quantity as specified in GB/T 26572-2011;

Environment Protection		10 Years		5 Years	
Toxic and harmful substance or element	PBDE	○	○	○	○
	PBB	○	○	×	○
	Cr6+	○	○	○	○
	Cd	○	○	○	×
	Hg	○	○	○	○
Name of parts	Pb	×	×	×	×
	Vibrating motor Electroacoustic device Accessories Battery	○	○	○	○

Environment Protection		10 Years		5 Years	
Toxic and harmful substance or element	PBDE	○	○	○	○
	PBB	○	○	○	○
	Cr6+	○	○	○	○
	Cd	○	○	○	○
	Hg	○	○	○	○
Name of parts	Pb	×	○	×	×
	PCBA Plastic parts Metal Hardware FPC	○	○	○	○

Please ask for a maintenance report from the maintenance center and properly keep it with the warranty certificate.
Shenzhen Yunding Information Technology Co., Ltd., in strict accordance with the Responsibility Regulations on Maintenance, Replacement and Sales Return of Household Appliances, offers guarantee for customers who purchase OcClean Smart Sonic Electric Toothbrush. In case of inconsistency between the service articles and the provisions of guarantee, the provisions of guarantee and the contents in the guarantee certificate shall govern.

Shenzhen Yunding Information Technology Co., Ltd.
#32G, Tower E, CR Land Building, Tonggu RD 5#, Nanshan District, Shenzhen, GD, China
Zip Code: 518057
If there is fault occurring in the product you purchase, please contact email: service@occlean.com

OcClean Facebook Page



www.facebook.com/occleanofficial

Follow the OcClean Official Facebook Page

Download App



www.occlean.com/download



www.occlean.com