Model: CM10B

MASSAGER GUN

User Manual

Arboleaf Corporation
5465 Legacy Drive, Suite 650 Plano, TX 75024 USA
www.arboleaf.com

EN	5
DE	15
IT	24
ES	33
FR	42

IMPORTANT SAFEGUARDS

1. Read all instructions.
2. This device is not intended to diagnose, treat, cure, or prevent any disease.
3. Consult your physician before use if you experience any pain or have any following health conditions:
 - (1) You are pregnant or have diabetes, retinal problems, or neuropathy.
 - (2) You have internal medical devices in your body.
 - (3) You have other medical conditions or health concerns.
 - (4) You are injured, after surgery, or fracture. Or you are infants, or elderly people with poor physical fitness.
4. This device is charged by a 5V / 2A charger or some fast charging power adapter.
5. Fully charge the device before its first use. Do not turn on this device during charging.
6. Do not replace the attachment head while it is working.
7. Do not immerse this device in the water or other liquid.
8. Do not use this device to massage your head, chest, genital, or any area with wounds. Avoid massaging the surgical area during the healing process.
9. Do not attempt to repair the device yourself.

CAUTION

This device has a built-in rechargeable Li-ion battery.

1. Do not dispose of the device with other ordinary trash.
2. Properly handle it by following your local recycling guidelines.

SPECIFICATIONS

Model: CM10B

Dimensions: 3.74 x 1.89 x 5.63 inches (95 x 48 x 143 mm)

Net weight: 0.82 lb / 370 g / 13 oz

Speed settings: 5

Speed: 1500 - 3000 rpm

Auto-off: Every 15 mins

Power supply: Rechargeable Li-ion battery

Battery: 2000mAh*2

Charging time: about 120 min.

Working time: about 500 min.


Max. voltage: 8.4 V

Max. current: 7.2 A

Max. power: 60 W

Rated input: 5 V \Rightarrow 2 A


Working noise: < 45 dB

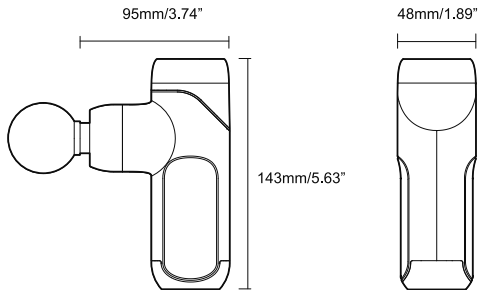PRODUCT OVERVIEW

Device (Massager gun)

- ① Replaceable massage head
- ② Handle
- ③ Speed indicator
- ④ Battery indicator
- ⑤ Power on-off / Speed switch button
- ⑥ Charging port

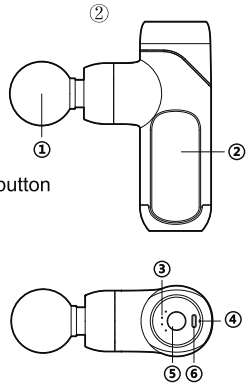Accessories

1. Attachment heads

- (1) Ball head
- (2) Fork head
- (3) Bullet head
- (4) Flat head
- (5) Air cushion head

2. Type-C charging cable

QUICK START GUIDE

1. Massage head installation

(1) Before installing the massage head, please turn off the machine first.

(2) Select an attachment head according to your needs:

Ball head: for large muscle groups, like back, leg

Fork head: for both sides of the spine and Achilles tendon

Bullet head: for point massages such as palms, feet, and soles

Flat head: for relaxing all body muscle parts

Air cushion head: for first users or sensitive muscles

(3) Insert the selected attachment head into the device.

Note: Never change an attachment head when the device is on.

2. Turn on / off the device

(1) Press and hold the Power on-off / Speed switch button

⑤ for 1.5 seconds, the speed indicator light will flash on and off in turn, and then one blue light is left to enter normal working mode, at this time you can open the App for Bluetooth settings and control. If you do not want to connect the App, you can directly adjust the speeds to use.

(2) Press and hold the on-off / Speed switch button ⑤ for 1.5 seconds to turn off the power and the device stops working.

3. Set your desired speed

(After the device is on, tap the on-off / Speed switch button ⑤, the speed will be shifted up and the speed indicator light will also change. Keep tapping the button until your

desired speed indicator light is on. The speed setting is done.

Note: The higher the speed, the higher the intensity of the massage. Below is the general guideline for your reference.

(1) Speed 1: One blue light on, suitable for muscle arousal.

(2) Speed 2: Two blue lights on, suitable for fascia relaxation.

(3) Speed 3: Three blue lights on, suitable for decomposing lactic acid.

(4) Speed 4: Four blue lights on, suitable for deep massage.

(5) Pressure-Adaptive Speed and Percussion: This innovative system utilizes pressure-responsive technology to automatically adjust speed and percussion based on sensed pressure. By incorporating pressure-sensing mechanisms, the system adapts to the user's input, providing an optimal experience with precise and tailored speed and percussion levels.

4. Charging instructions

(1) When the power is low, the power indicator flashes red.

(2) When charging, the power indicator is red and always on.

(3) When charging is completed, the power indicator will be green and always on.

5. Auto power off protection

(1) System protection (low voltage prompt): When the voltage is lower than 7.3V during normal operation, the red light flashes; when the voltage is lower than 7.0V, the

system automatically shut down.

(2) Blocking protection: when the blocking current reaches the protection value, the motor stops rotating; after entering the blocking state, short-press the function button to enter the next speed to work again.

(3) Maximum run-time protection: This device will auto-off every 15 mins. To re-start, press the Power on-off / Speed switch button for 1.5 seconds to turn on the massager.

6. "arboleaf" App Download and Use Instructions

a. To download the "arboleaf" App, Scan the QR code below or search "arboleaf" App in the App Store or Google Play

b. Create an arboleaf account. Sign up and complete your profile.

c. Turn on the Bluetooth of your phone

(a)

(b)


(c)

- d. Press and hold the Power on-off / Speed switch button ⑤ for 1,5 seconds to turn on the massager, and the green light comes on to enter the power-on state.
- e. Tap "My Device" and select "Add Device" in the App to connect the device.
- f. Follow the App instructions to complete adding the device.
- g. You can use the App to operate and control the device.


(d)

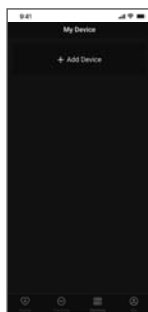(e)

(f)

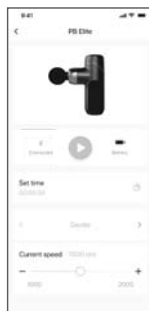(g)

WARRANTY

Customer Service:

Arboleaf Corporation warrants this product purchased from an authorized seller to be free from defects in workmanship and materials, under normal use and conditions, for a period of one (1) year following the date of purchase. Proof of purchase is required.

This warranty does not cover failures due to misuse, abuse, accidental damage, or unauthorized repair. Charging cable is excluded from the warranty.

1-800-658-1148 (US)

support@arboleaf.com (US)

support.eu@arboleaf.com (EU)

www.arboleaf.com

Arboleaf Corporation

5465 Legacy Drive, Suite 650 Plano, TX 75024 USA

FCC Statement

This device complies with part 15 of the FCC Rules.

Operation is subject to the condition that this device does not cause harmful interference (1) this device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

Changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

NOTE: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation.

If this equipment does cause harmful interference to radio or television reception,

which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.

-- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.

-- Consult the dealer or an experienced radio/TV technician for help.

The device has been evaluated to meet general RF exposure requirement. The device can be used in portable exposure condition without restriction.

FCC ID: 2ANDX-CM10B

IC Caution:

- English:

This device contains licence-exempt transmitter(s)/receiver(s) that comply with Innovation, Science and Economic Development Canada's licence-exempt RSS(s). Operation is subject to the following two conditions:

- 1.This device may not cause interference.
- 2.This device must accept any interference, including interference that may cause undesired operation of the device.

The device has been evaluated to meet general RF exposure requirement. The device can be used in portable exposure condition without restriction

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