

Anything else we can help with?

Call us!
1-844-633-8729



Hello!
Let's get started.





Any changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate this equipment."

This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) This device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.


From our lives to yours

Back in 2012, a group of friends were given the news that their friend had been diagnosed with a rare form of cancer. Their friend, Nick, was put on a new regimen of medications to manage on a daily basis. Seeing Nick struggle, the group of friends set out to create what would eventually become the TowerView Health pillbox.

Since then, Nick and those friends have worked hard to be able to provide that same pillbox to others. We're so happy to be able to provide this pillbox to you now, and we hope that having this pillbox will help you stay in good health.

We are always here to help. If you ever need to talk to us, please call **1-844-633-8729**, even if it's just to say hi!

What you will receive

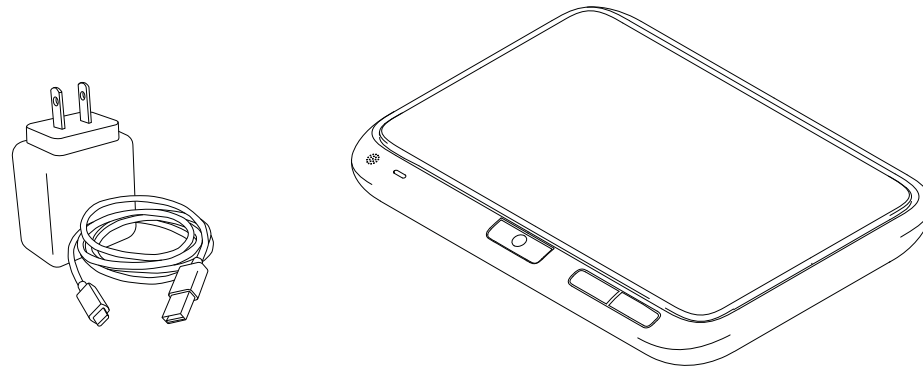

USB Charging Cable

Pillbox

Box 1

Find an uncluttered surface where you can place your pillbox. It should be at table-height. Then plug in your pillbox and it will set itself up!

Do's and Don'ts

- Do keep your pillbox in a dry environment and close to a window.
- Don't keep your pillbox in a closed cabinet or drawer.
- Don't keep your pillbox in extremely hot or cold environments.



Safety

- Interactions between the pillbox and other medical devices are not likely, but can be possible.
- The best prevention for adverse interactions is separation of devices: Three feet or about one arm's length for an adult is sufficient.
- When first using the pillbox, watch any other medical devices for strange operations. If you notice anything unusual, unplug and move your pillbox and contact us at **1-844-633-8729** for instructions.

Customize your reminders

Beyond just the default reminders, we love to set up tailored reminders that work best for you. Some options that our members find useful are:

- Sending text message, call, and email reminders to your phone.
- These can also be sent to family, friends and loved ones.
- Turning off reminders completely.

You can set these up and change them at any time and as often as you like by calling **1-844-633-8729**.



Power Light


The power light has two modes based on whether the pillbox is plugged in

Pilbox plugged:

-  Charging
-  Fully Charged

Pillbox unplugged:

-  Good Battery
-  Low Battery

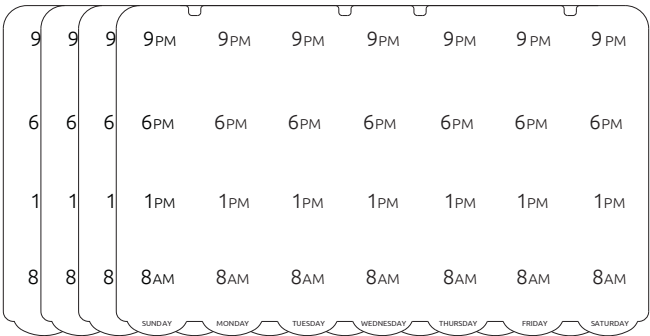


Medication Trays

Box 2

Your trays will come in a separate shipment and the scheduled times will be customized for you!

Get to know your pillbox

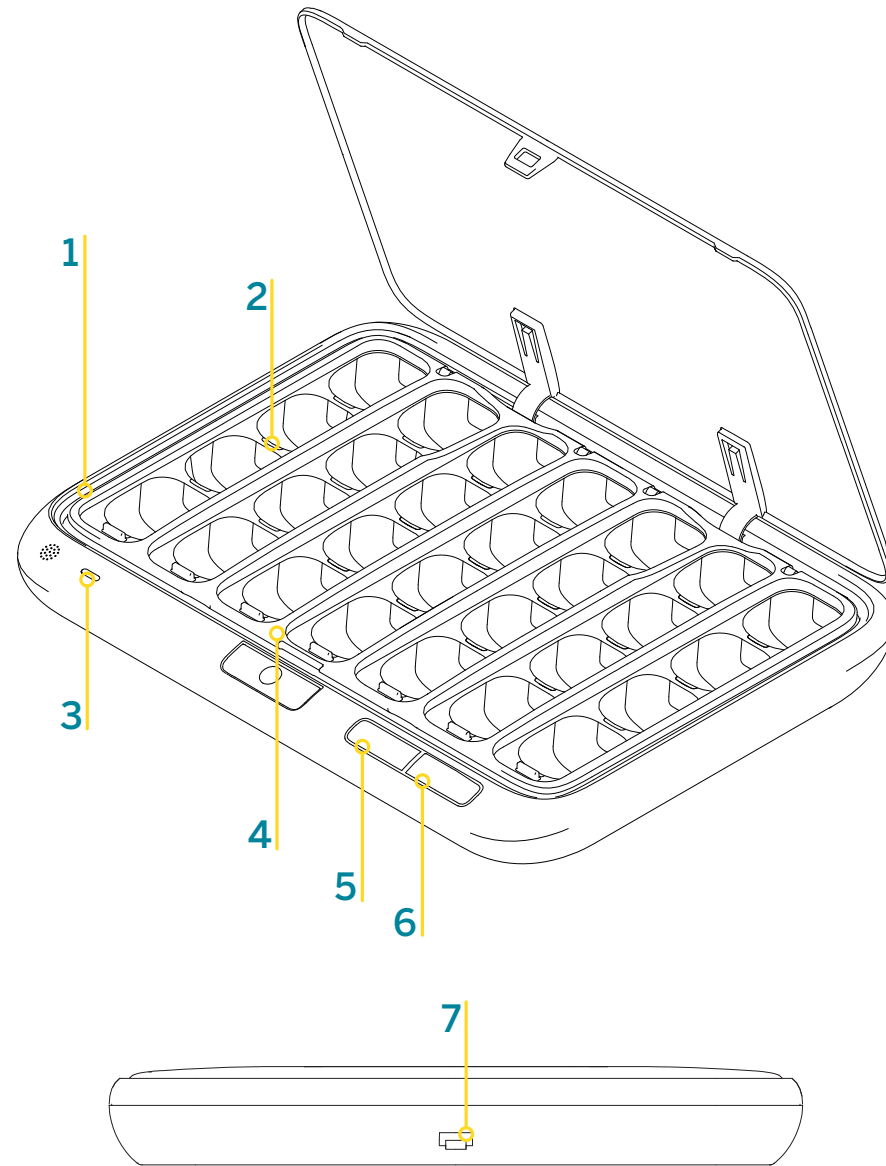

How and when you will receive reminders

- Once your pillbox is activated, by default, the light ring will light up and play a reminder chime at:
 - Your scheduled dose time
 - 30 minutes after your dose time
 - 90 minutes after your dose time
- If you take your medication at any point during this period, the subsequent 30 minute and 90 minute reminders will not turn on.

For example:

- If your morning dose is at 8:00 AM:
 - Your pillbox will light up and chime first at 8:00 AM.
 - Only if you do not take your medication before 8:30 AM, your pillbox will light up and chime again.
 - Last, only if you do not take your medication before 9:30 AM, your pillbox will light up and chime a last time.

How your pillbox works

A few notes about your pillbox:

- Your pillbox can sense when medication is removed.
- Your pillbox will send you reminders to take your dose if you've forgotten. Call **1-844-633-8729** to customize your notifications. You can choose from:



Sound



Light Ring



Call



Text Message



E-Mail

- No need to connect your pillbox to the internet! The pillbox automatically uses a cellular connection to always download your schedule and stay up-to-date.
- The pillbox has no on/off switch. Keep your pillbox plugged in at all times so you never have to worry about the battery dying.

1 Light Ring

Lights up at scheduled dosage time. Turn to page **9** for more information.

2 Well Light

Lights up to indicate the scheduled dose.

3 Power Light

Pillbox battery power indicator. Turn to page **10** for more information.

4 Tray Cage

Locks medication tray into pillbox. Lock and unlock the cage by pushing until you hear a click. Cage should sit on top of your medication tray.

5 Multi Button

Multipurpose button. Can silence reminder chime while chime is playing.


6 Volume Toggle

Turns sound volume down, \ominus , or up, \oplus . To mute pillbox, press \ominus repeatedly until pillbox stops beeping.

7 Power Port

Connects your pillbox to the USB charging cable for charging.

Get to know your medication tray

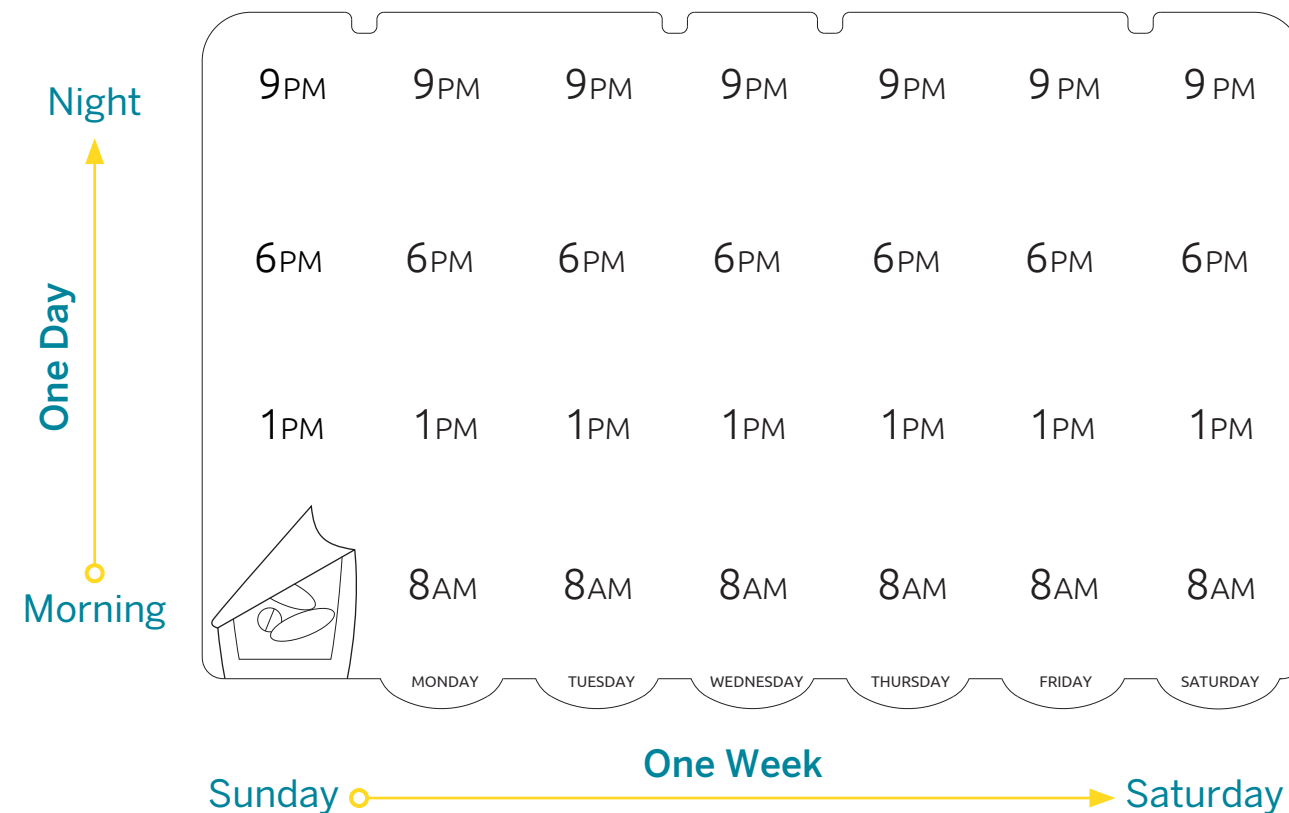

Your tray will be customized to your own medication schedule.

How to set up your medication tray with the pillbox:

1. Open the pillbox lid by pressing the central lid button.
2. Release the cage by pressing the central “push” button located on the cage.
3. Remove the medication card on the left of the tray by tearing along perforated edge.
4. Insert the tray so that the Sunday column is to your left and the morning row is nearest to you.
5. Make sure the tray sits level in the pillbox — you should be able to close the cage easily!
6. Lock the cage by clicking it closed.
7. Close the pillbox lid.

How to remove your medication:

1. Find the tab for your dose. It will be labeled with the current date and time. It will also be illuminated with a blue well light.
2. Pull up and peel back the tab for the appropriate well.
3. Remove your medications.
4. Once you remove the tab sticker, throw it away!*

* Leaving the label on or in the tray may affect the pillbox's ability to sense your medications.