

Table of Contents

1	Functions and effects of the SOQI Bed	34
2	Diagram of the SOQI Bed	36
3	Parts List	37
4	Installation Steps	38
5	Connection/Diagram of the Electrical Devices	52
6	Electric Installation	53
7	Operation Instructions	56
8	Technical Specifications	62
9	Warranty Card	63
10	Service Record	64
11	Warranty	65
12	Troubleshooting	66

Functions and effects of the SOQI Bed

Functions and Effects of the SOQI Bed HotHouse™

- This product provides FAR Infrared Heat, it has been reported that it has had a positive effect over arthritis of the shoulder and pain in the lower back area.
- The Infrared Rays are a type of electromagnetic wave, which is good for the health and the growth of the organism. The person's organisms absorb the Infrared Rays, regulating the physiologic functions.
- The Infrared rays can promote the blood circulation, relieving fatigue and high blood pressure; also, regulating our physiological functions, which will give us a healthy life.
- Through the warm and efficient FIR Rays, we can get thermotherapy benefits in a natural and secure way.

Arthritis of the Shoulder

This condition usually presents the following symptoms: Pains of long duration in the shoulders, spasms and atrophy in the muscular area of the shoulder, tendonitis and restriction of the movement of the joints in every direction. The pain can increase if movements of turning the back are made. Sometimes the pain can extend to the neck and back.

Pain in the Lower Back

It usually involves long lasting pain in the waist, which sometimes can be unnoticed, sometimes others can even be intolerable. The pain can extend from the buttocks to the calf of the leg and it can intensify if you are sitting or standing for long periods of time. The pain can be relieved by standing up with a straight back or by putting your hands on your waist. (Usually a cold sensation is felt in the lower back and the movement of your waist is limited. The intensity of the pain can increase when your waist does not have a warm temperature.



Chi Machine®


The Chi Machine aids in the therapy to assist with the pain at the lower back area caused by a lesion. The design of the Chi Machine can be described in this way: the feet of the patient rise and the lower part of the extremity is turned, the force can spread from the joint of the damage area to the muscle. The inflammation and pain can be relieved and restored to the normal function of the joint.

When you use The Chi Machine, the person must be in a hard and flat surface and the area shouldn't be pressured. When you put your feet in the resting area of The Chi Machine, the feet will elevate in a natural way, which will promote the blood circulation. Also, in addition to the frequency of the movement is a determinant factor for the benefit of the blood circulation and the metabolism. It also has some effects in the health maintenance and prevention of diseases.


- The oscillatory movement is designed accordingly with the needs of the physiological system.
- The oscillatory movement can promote the blood circulation, relieving fatigue and pressure.
- Relieves pain and the sensation of tiredness, improves the functions of the body and release stress from the body and mind.



Diagram of the SOQI Bed

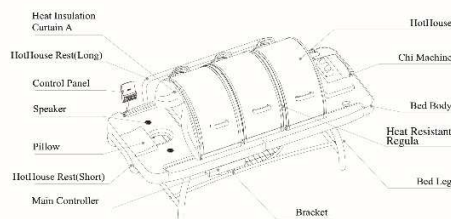

Parts List

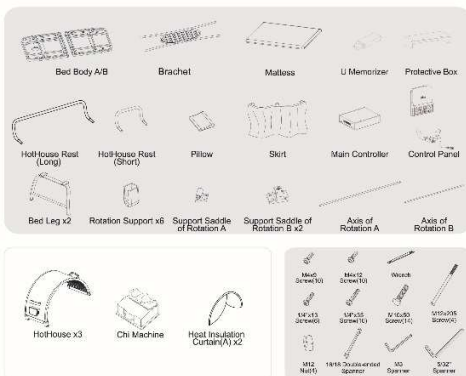

Installation Manual


Step 1

Open the package, take out the bed structure and put it on the floor with the interior part facing up. Insert the bed supports to the bed frame and secure them with the M12x20S screw (medium) like is shown in the illustration.



Important

The bed structure must be placed in a clean and flat surface in order to avoid damage and dirt.

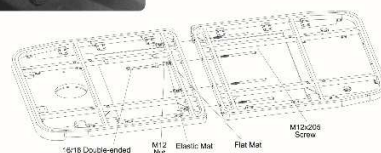


Step 2

Install the second support of the bed the same way you did in step one.



Important

When you install the supports of the bed, make sure that you align them in a way that the holes of the bed and the supports are in the same place. Make sure the screws are firmly tightened.

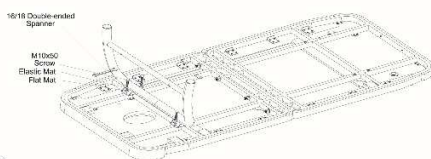

39


Step 3

Take out the bracket and insert it into the bed frame with the M10 x 50 Screw as shown in the illustration.



Important

When you install the bracket, make sure the holes of the bracket are perfectly aligned with the holes of the bed support. When the parts are aligned, firmly tighten the screws.

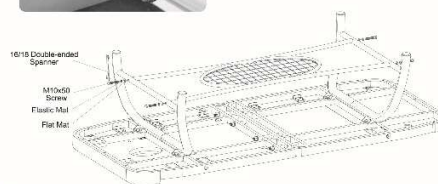

40


Step 4

Take the MP3 player control box and put it under the surface formed like an "I." of the bed frame connect it from the side with the M4 x 9 screws.



Important

Take special care in handling with the MP3 player, the side with the switch must face the front.
After you finish this step and make sure that everything is in place, you can now turn the bed.

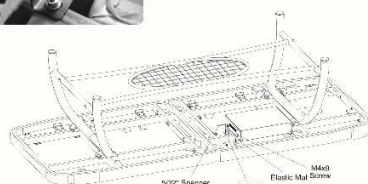

41


Step 5

Once the bed is facing up, take the support saddle of rotation A/B, fix them with the M4x9 screw and M4x12 screw at the right side corresponds to the location of bed board like shown in the picture.



Important

After the bed body is assembled, if you find the bed is not stable, please use the wrench from the tool bag to adjust the height of the legs. There should be no gap between the legs and the ground. Adjusting the legs will ensure the bed body is even.

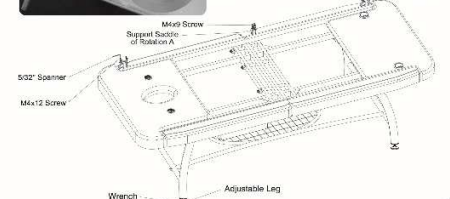

42


Step 6

Remove the rotation axis A/B, insert them from the both sides of the support saddle of rotation B to the support saddle of rotation A, and then firmly tighten the screws like shown in the picture.



Important


1. Before the rotation axis A/B are installed to the support saddle of rotation, the rotation support must be installed at first.
2. Pay attention to the relative position of the rotation support while installation.

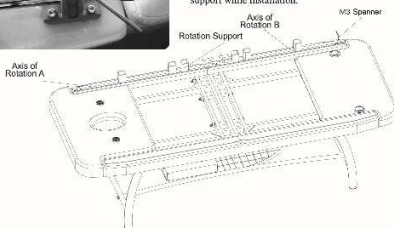

43

Step 7


Take the Hot House Rest (short) and insert it in to the groove located at the head of the bed. Fix it with the 1/4"x35 screw (short) like shown in the picture.

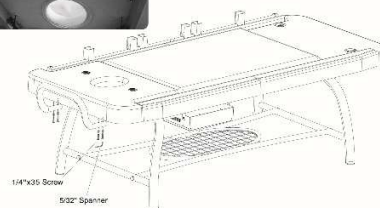

44

Step 8


Take the HotHouse Rest(long), fix it with the 1/4" x 35 screws from underneath the bed like it is shown in the picture.

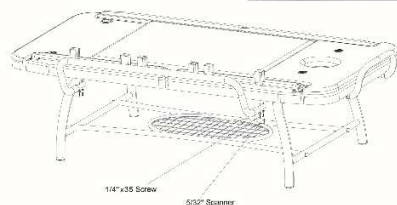

45

Step 9

Take the Control Panel, and slide the cables through the hole design for the cables and fasten it with the 1/4" x 13 screw (short).

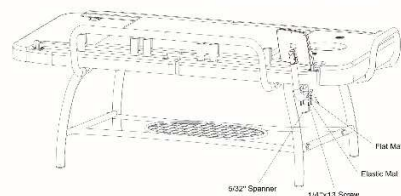


46

Step 12




Place the Chi Machine at the end of the bed. On one side is a gray cap, press the edge and pull up the cap in order to pass the cables to the electrical socket. Afterwards, replace the gray cap.

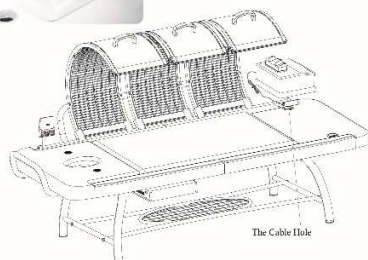

49

Step 13

After all the devices are installed and aligned properly, plug the cables of the machines to the back of the MP3 player/Control Box.

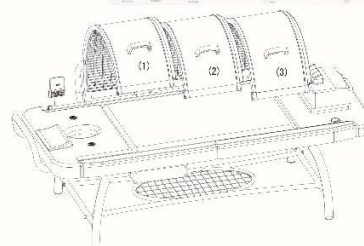


50

Step 10



Take the HotHouse (1) and slide the edges in that go inside the rail and have them coincide with the interior edges of the HotHouse as shown on the illustration.

Important
Please install the HotHouse in order and sequence.

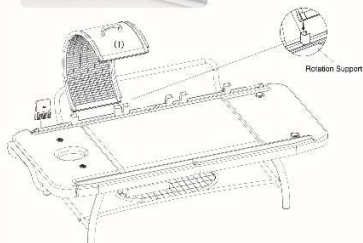


47

Step 11



Install the HotHouse (2-3) using the same procedure mentioned before.

Important
After all the devices are installed and aligned properly, plug the cables to the electrical socket which is located at the back of the MP3 player.

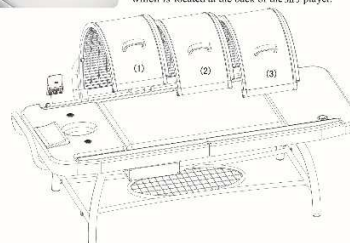

48

Step 14

1. Take out the bed skirts from the packaging and spread them. The two long pieces are for the left and right sides while the two short ones are for the front and the end.
2. First, install the longer piece with the cover on the side of the main controller by hanging the 7 rings of the bed skirt onto the 7 hooks of the bed (as shown in picture no. 1). Follow the same steps for the other 2 pieces.
3. Tie the ribbons on the corresponding tubes (as shown in picture no.2)
4. Fasten the sticking nubs.
5. Overall view (as shown in picture no.3)



Step 15


1. Take the Heat Insulation curtains A from the package and spread them out, as shown in Picture 1.
2. Install curtain A on the HotHouse from both the front and backside, as shown in Picture 2.
3. Organize Curtain A.



51

Installation Manual

Connection /Diagram of the Electrical Devices

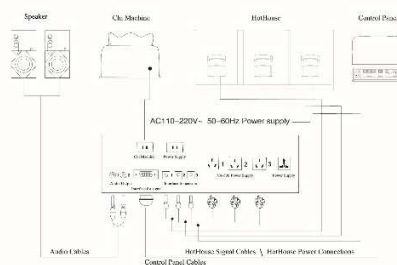

52

Electrical Installation

1. Assembling of Cables.

Connect the speaker cables into the audio socket located at back of the MP3 player/Control Box.

2. Assembling of the Control Panel

According to the diagram, plug the cable to back of the MP3 player/Control Box and tighten the screws of the plug.

3. Assembling of the HotHouse

Insert the cables of each HotHouse into each correspondent socket located on the back of the MP3 player/Control Box while following the diagram, "Connection of the Electrical Devices". For example: if the HotHouse is labeled with a no.1, the electrical cord should be plugged in socket no. 1.

4. Assembling of the Chi Machine


Insert the plug into the socket "SDM" according to the connection diagram.

Important :

Make sure that all the cables are plugged properly. Secure the cables with ties.

53

5. Take out the protective box and insert it into the bed frame with the screws M4x9 screw as shown in the illustration.

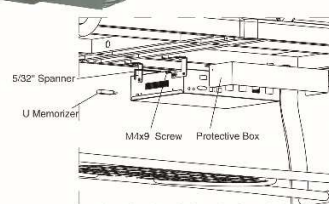


54

- (6. When all of the connections are made and have no problems, release the main electrical cord, which is located at the back of the MP3 player/Control Box, and plug it to the room outlet.

Important :


The sockets of the holes 2-3 are designed to plug in other devices.

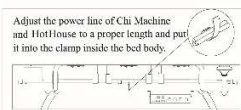

55

Operation Instructions

1. Operation may start after the multifunctional physiotherapy couch has been installed and modulated properly.
2. Centralized operation mode is applied to control the couch. It is consisted of three parts: Far infrared apparatus, physiotherapy apparatus and music player. They can work either simultaneously or separately.
3. Bluetooth-controlled couch is an upgraded version of the touch-controlled couch, mainly for the user's convenience. Table/cellphone can connect the couch via Bluetooth and hence realize the control. Operation procedures are listed below:

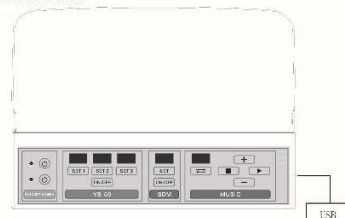

I. Controller:

1. The controller is redesigned and added with the Bluetooth system. The Bluetooth beacon is set inside the master control box.
2. The controller is also added with USB ports that are available for charging the table/cellphone.

56

3. User can operate the couch directly via the controller apart from tablet/cellphone Bluetooth. The functions of buttons are introduced below:

- ◆ “ ” U disk music play/next. Touch once to play, and touch again for next song;
- ◆ “ ” U disk music stop. Touch once to stop music;
- ◆ “+” “Touch once to increase the volume;
- ◆ “-” “Touch once to reduce the volume;
- ◆ “ ” Single cycle playback. Touch this button to play the song shown on the digital display during the music play of U disk. Touch it again to cancel the cycle and go to next song after current play;
- ◆ “ ” The nixie tube shows the track number of music that is playing in the U disk;
- ◆ “SET 1” “SET 2” “SET 3” “YS-50 work time setting. Time can be set for 15/30/45/60min;
- ◆ “ON/OFF” “Control YS-50 to start/stop work. After the work time is set, press the button to start work and press again to stop. Time returns to zero;
- ◆ “SET” “Set the SDM work time. Time can be set for 1/2/3/4/5/6/7/8/9/10/15/20/25/30min;
- ◆ “ON/OFF” “Control SDM to start/stop work. After the work time is set, press the button to start work and press again to stop. Time returns to zero;

Button settings are detailed below:

- ◆ “ ” The upper power button is to control the backup power on/off from the 2-core socket behind the master control box. The left red indicator turns on after touch and the power is on. It turns off after touch again and the power is off;

- ◆ “ ” The lower power button is to control the backup power on/off from the 3-core socket behind the master control box. The left red indicator turns on after touch and the power is on. It turns off after touch again and the power is off;

II. Tablet/cellphone operation:

1. Install the “iSOQI” App:

- a. Open the tablet/cellphone;
- b. Download the “iSOQI” install package;
- c. Click the package to install the “iSOQI” App onto the tablet/cellphone;

2. Connect the Bluetooth of the master control box:

- a. Turn on power of the master control box;
- b. Open the tablet/cellphone and turn to setup“ ”;
- c. Select the Bluetooth and enable it“ ”;
- d. Click the search;
- e. Among the search results, find out the Bluetooth device number soqibed***** (e.g.: soqibed00016), and then connect it (if connect successfully, it will

show: soqibed00016);

- Note: 1. After Bluetooth is matched, click soqibed***** directly;
2. Installation is required for the first time. User can begin with the third step later. (If a text of “please enable your Bluetooth first” appears, user needs to enable the Bluetooth manually).

3. Start the “iSOQI” App to enter the control interface;

- a. Click the “iSOQI” App to enter into the master interface“ ”;

b1. Click the Visitor Login:

- Select the matched Bluetooth device (e.g.)

b2. Click the Store Owner Login (tablet/cellphone must be connected to the network before login):

- Enter the owner's account and password, and then click

confirm (e.g.)

- Select the matched Bluetooth device (e.g.)

4. Functions of the controller interface:

- ◆ “ ” “User can select the VIP mode for individual customization (not available to visitor);

Notes:

Individual customization may not be provided only after the SOQI apparatus and its dedicated Bluetooth box are equipped and tested by professionals.



- ◆ “ ” 7 standard modes, 1 customized mode and monthly care mode are available;
- ◆ “ ” After selecting the mode, press the start button and then the system tips pop up. It will start work after waiting for 60s or by pressing the immediate start;
- ◆ “ ” By pressing the stop button, the couch will stop work and all settings return to zero;
- ◆ “ ” By pressing the VIP fresh button, the software will update the VIP data automatically;
- ◆ “ ” Music control buttons from left to right are: last; play; stop; next; single cycle. Notes: these control buttons are to control the music from SOQIBed App only. They cannot control the music from the U disk;
- ◆ “ ” Volume control button (it can control the volume of music from “iSOQI” App only. The track number is the music number in the “iSOQI” App, totaling 01-06. The music varies with different work modes);
- ◆ “ ” After the work is over, press the exit button to quit the App.



Additional attentions


1. The positions of far infrared apparatus and physiotherapy apparatus may be adjusted according to actual conditions. Please pay attention to the position of conductor wire and extend it if necessary.
2. The far infrared apparatus cannot be independent of the couch controller for separate use.
3. Please refer to "Step 14" for use and installation of thermal insulation curtain.
4. Please refer to "Step 15" for use and installation of enclosing curtain.
5. Please clean the couch with leather detergent and soft cloth to maintain its daily sanitation.
6. Changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.
7. This equipment has been tested and found to comply with the limits for a Class B digital device pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:
 - Reorient or relocate the receiving antenna.
 - Increase the separation between the equipment and receiver.
 - Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
 - Consult the dealer or an experienced radio/TV technician for help.

61



Technical Specifications

Voltage of Power Supply : <input type="checkbox"/> AC100-120V <input type="checkbox"/> AC220-240V	Rated Power : approx. 1000W
Frequency of Power Supply: 50-60Hz	Applicable Environment : 0-40°C
Weight : 63 kg (Bed Body) 92 kg (Bed Body and Devices)	Exterior Dimensions(mm) : Shown as Fig

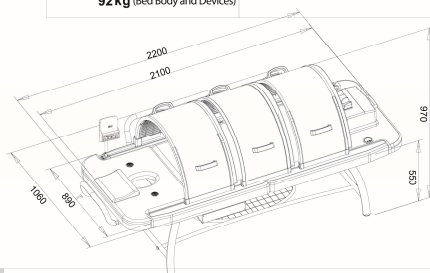

62



Warranty Card

Warranty Card

NO. _____

Name	Address
Phone Number	City : State/Province : Zip : Country :
Purchase Date	1. The product can be exchanged in 14 days after purchase, get free warrant for one year since the date of purchase (warranty card should be shown).
Product Name	2. If warrant card can not be shown, material cost will be charged as per the condition.
Model #	3. Service telephone number: _____
Reference #	Free warranty policy
Stamp of the Dealer	1. If the problem occurs within warranty period and it is correctly used according to the notes in the Operation manual, the maintenance will be free.
	2. Problem caused by user dismantling or contrived problem will not get free maintenance or exchange.
	3. Free warranty of exchange will be refused under following conditions even within warranty period
	* The warranty card can not be shown;
	* Contents in the warranty card are altered;
	* Missing, incorrect repaling of relic;
	* Problems by falling or mis-operating after purchase;
	* Problems or damages by fire, earthquake, or other disaster.
	4. For other questions, contact dealer.
	5. *Please keep this card carefully. It will not be reissue again if it is lost.


63



Service Record

Date	Damages	Examine and Repair Records	Signature of Repairer	User's Signature


64



Warranty

- ◆ This product has one year warranty since the date of purchase.
- ◆ The customer shall not dismantle the devices to avoid injury.
- ◆ If problem occurs during the usage, contact the dealer.
- ◆ If the following condition happen, certain maintenance cost will be charged as per the condition even within the warranty period,
 - ◆ Problem caused by incorrect use;
 - ◆ Problem cause by modification by costumer or incorrect repair;
 - ◆ Problem or damaged by fire, earthquake, flood or abnormal voltage;
 - ◆ Problem by abnormal condition of the power supply system;
 - ◆ Parts consumption;
 - ◆ Without warranty card;
 - ◆ The warranty card is not complete if there is no stamp of the dealer or its modified.
- ◆ The warranty card will not be issue again if it is lost.

65




Troubleshooting

Before returning the unit to the company for service, please read the following problem check list carefully.

Problem	Possible Reason	Item to check
<ul style="list-style-type: none"> • The unit does not work normally after the power is turned on. • Indicator light doesn't show functionally. • The HotHouse and the Chi Machine won't work at the same time. 	Power outlet or plug is not connected correctly.	<ul style="list-style-type: none"> • Is the plug broken? • Is the fuse broken? • Is it completely plugged into the wall? • Is the power well connected? • Is the timer set to "0"? • Is the power switch turned "On"?
<ul style="list-style-type: none"> • The Unit overheats. 	<ul style="list-style-type: none"> • Power outlet or plug is not connected correctly. • Cables are plugged in the wrong order or plugged in improperly. 	<ul style="list-style-type: none"> • Is the plug broken? • Is the fuse broken? • Is it plugged in completely and correctly? • Is the power well connected?
<ul style="list-style-type: none"> • Speakers won't work. • The volume adjust function won't work. 	Power outlet or plug is not connected correctly.	<ul style="list-style-type: none"> • Is the plug broken? • Is the fuse broken? • Is it plugged in completely and correctly? • Is the power well connected?

66

<div>  <div> <div>SAFETY CHECK LIST</div> </div> </div> <div> <p>For your safety and to avoid fire, if any of the following conditions occur, stop using the unit, turn the Unit off, and pull out the plug. Return it to the closest HTE Branch for inspection and repair.</p> <ul style="list-style-type: none"> • The radiant plate is sometimes hot and sometimes cold while the unit is "On." • The power cord is broken or frayed. • The cord or plug becomes excessively hot. • The Unit emits smoke or an unpleasant odor. • The unit is not working properly. </div> <div> <p>To avoid damaging the unit or causing injury DO NOT take the Unit apart.</p> </div> <div> <ul style="list-style-type: none"> • This unit is installed with a fuse. If the fuse breaks, please purchase one from the company and replace it as soon as possible. • When first using the unit, you may see a little liquid. This is normal and you may safely keep using the unit. </div> <div>67</div>	
--	--



3. Start the “iSOQI” App to enter the control interface;

a. Click the “iSOQI” App to enter into the master interface“



b1. Click the Visitor Login:

-- Select the matched Bluetooth device (e.g.



scqibet00176
00:11:11:11:11:11

)

b2. Click the Store Owner Login (tablet/cellphone must be connected to the network before login):

-- Enter the owner's account and password, and then click

confirm (e.g.



)

-- Select the matched Bluetooth device (e.g.



scqibet00176
00:11:11:11:11:11

)

4. Functions of the controller interface:



“



”User can select the VIP mode for individual customization

(not available to visitor);