

## Replacing Brush Head

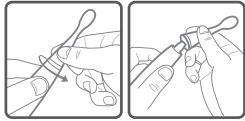Turn the used head to the right and pull out it.  
Put the new head to the top of the head and turn left to lock in.

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## Replacing Battery

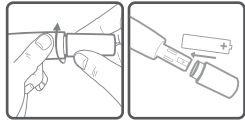When the vibration power decreases, separate the bottom battery cap by twisting the cap to the right. Insert the new battery on the right pole as shown below. Then, turn the battery cap to the right to lock in.

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## How to solve problems

• **Before inquiry** Before asking for help or AS, follow the instruction for solving the below lists of problem. If the problem is not fixed with them, please call to the service center(+82-70-7620-0405) • **Product is not working** It may discharge completely. Replace to the new battery. Checking the battery is inserted in right way. (check the positive(+) and negative (-) terminals.) • **Vibration works less than usual** It may discharge almost. Replace to the new battery. Vibration power is controlled in 4 levels on the App. • **Fail to find the Brush Monster on Bluetooth lists** There is no process for connecting to the App. Just turn on the Bluetooth and run the App. After following the instruction, push the button on the body. It will be connected automatically. • **Bluetooth connection is not available or disconnected** Turn on the Bluetooth on the smartphone. Check the version of application. If it is not updated version, update it. Using both of the smartphone and product in 10 meters. Distance is affected by the usage environment. Reboot the smartphone and open the application again. • **Product is over automatically** It is intended to powered off automatically for guide of correct brushing.(160 Secs with App. 120 Secs without App.)If the product is off earlier, ask to the service center.

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• **Inquiry about product** Before inquiry, let us know about Model name, Symptom, Contact Kitten Planet service center : +82 - 70 - 7620 - 0405

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## WARRANTY

### One Year Limited Warranty

Name of product	BRUSH MONSTER Smart Toothbrush
Model	BMT100
Date of Purchasing	
Place of Purchasing	
Your Name	

1-year warranty applies to this product from the date of purchase. According to the condition of the product, defects caused by the device itself or during production will be replaced free of charge for 1-year. For further inquiries, please contact Kitten Planet E-mail: contact@kittenplanet.net, Tel: +82-70-7620-0405

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Country of origin Made in Korea

Copyright

Copyright©2018 Kitten Planet Co., Ltd. This manual is protected under copyright laws. No part of this manual may be reproduced, distributed, translated, or transmitted in any information storage and retrieval system, without the prior written permission of Kitten Planet Co., Ltd.

## Caution

• This product is for children ages 3+. Avoiding use for children ages 3-. • Do not bite or chew on the product. It may damage to brush head • Do not wash the toothbrush with hot water. It may cause the brush head to deform. • Excessive use of UV tooth-brush sanitizers may deform the brush head. • It is recommended for replacing the brush head monthly for the dental hygiene and effectiveness of brushing. But if it is transformed or worn out, replace it right away. • Do not use unauthorized batteries. Their usage may reduce the lifespan or cause damage to the product. Only use authentic batteries that have been authorized. • If you treat for gums or teeth, consult a dentist before using the product. • If you experience excessive bleeding that persists for more than a week, consult a dentist. • With domestic safety standards, the products certified with the electromagnetic interference test. If your body is implanted with pacemakers and other electronic medical devices, consult a doctor or the device manufacturer before usage. • Children or adults with the physical or mental illness have to use under the guardian or a person who is responsible for their safety.

**FCC Warning Statement** FCC Part 15.19 This device complies with part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) this device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation. FCC Part 15.21 Any changes or modifications (including the antennas) to this device that are not expressly approved by the manufacturer may void the user's authority to operate the equipment. **FCC Radiation Exposure Statement** This equipment complies with FCC Radiation exposure limits set forth for an uncontrolled environment. This device and its antenna must not be co-located or operating in conjunction with any other antenna or transmitter. This equipment should be installed and operated with a minimum distance of 10cm between the radiator and your body.

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## User Guide

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## BMT100

Thank you for purchasing BRUSH MONSTER

We hope your child will enjoy using BRUSH MONSTER toothbrush and application.  
Before guiding your child, please read the instructions thoroughly.  
Until your child becomes familiar with the product, parental supervision is recommended.

www.brushmon.com