


BRUSH MONSTER SMART BRUSH(BMT200_Y)

Thank you for purchasing the BRUSH MONSTER product.

Please read the user manual before using the product.

- Components

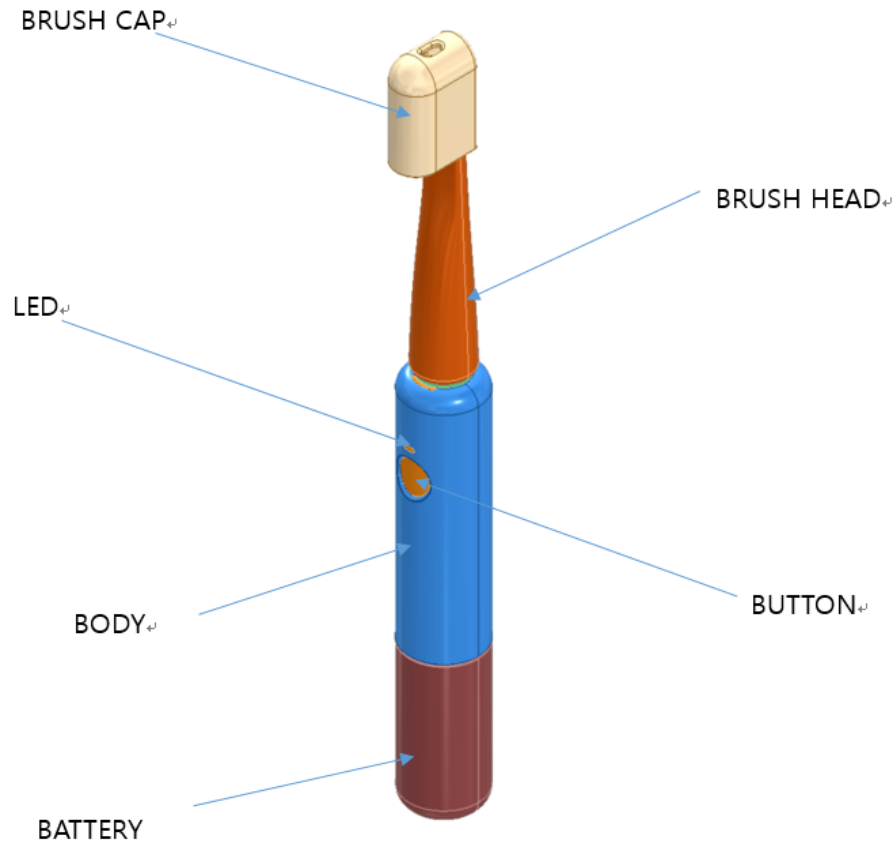

- Body: The mainbody of the product.
- Battery Cap: A cap the can be separated when the battery is replaced.
- Brush Head: Replaceable toothbrush head.
- LED: A light that shows the current status of a smart toothbrush.
- Button: Button for product operation.
- Head Cap: Cap to protect the brush head.
- Battery: Battery for smart toothbrush operation. (AAA)

- Specification
 - Product name: BRUSH MONSTER Smart toothbrush
 - Model: BMT200_Y
 - Size: 177 x 24 x 21 (mm)
 - Weight: 32 g(excluding battery)
 - Sonic Vibration: 20,000 RPM
 - Power: AAA Dry battery 2 ea
 - Bluetooth: 5.0, 2.4GHz
 - Material
 - Body: ABS, Silicone
 - Button: TPU
 - Brush Head: PP
 - Brush: PBT
 - Water Proof: IPx7
 - Age: Over 36 Months old

- Brief Instructions on how to use it.
 - With App: Push the button once before starting brushing, then it is connected to the app and starts to vibrate.
 - Without App: Push the button once, then it starts to vibrate (Duration is 2 minutes, and it alerts every 30 seconds.)
- How to use
 - Install the application and create account
 - Search For 'Brush Monster' on the Google Play or App Store.



- Run the app after installation.
- Register your profile after create account.




- Starting with the product


- Remove the battery cap of the product by turning it counterclockwise.



- Insert the enclosed battery (AAAx2) to match the polarity (+/-) inside the product.

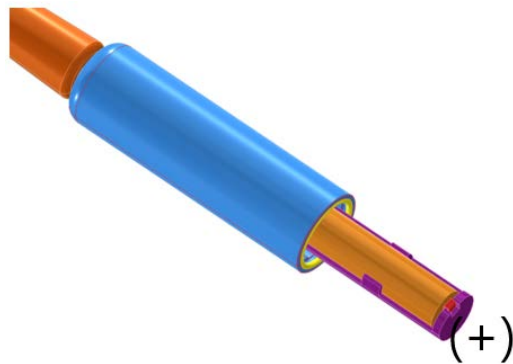


- Turn the battery cap of the product clockwise and insert it.

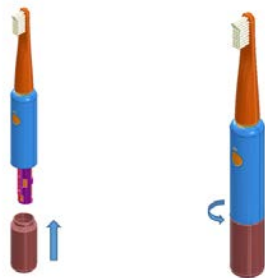


- Connect to device.

- Run the Brush Monster app.
- Activate Bluetooth on mobile devices such as smartphones or tablets.
- Press the button on the toothbrush once to connect to the device.

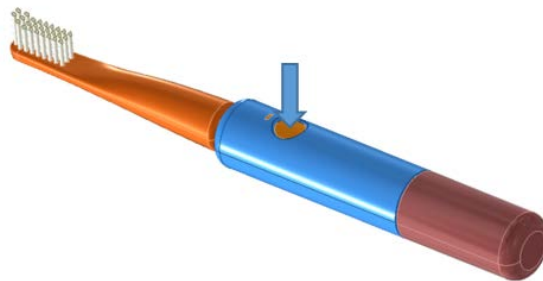

- Brushing teeth

- The product and application are connected, and brush your teeth according to the tooth brushing guide on the screen. (The sensor of the smart toothbrush directly guides the direction and location of the tooth brushing.)

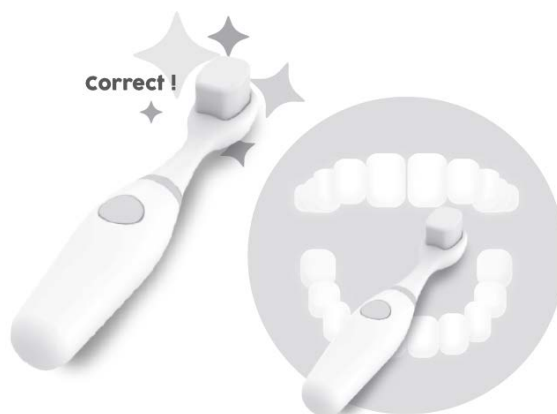


- Check the dental report

- Check the brushing report when the entire process of brushing is completed. (You can check where to brush well and where to brush more carefully in the brushing report.)



- Without app

- Press the button of the product once.

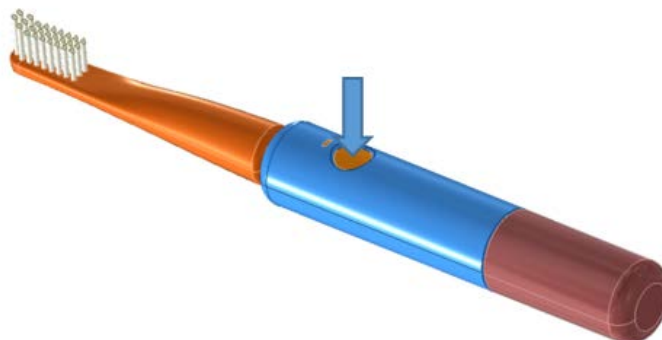

- The vibration starts from the toothbrush. (Vibration is maintained for 2 minutes, and haptic alarm goes off every 30 seconds.)

- Replacing brush head

- Turn the toothbrush head counterclockwise and pull it out upward.



- Turn the new toothbrush cap clockwise and insert it.




- Replacing battery


- Remove the battery cap of the product by turning it counterclockwise.



- Insert the enclosed battery (AAA) to match the polarity (+/-) inside the product.

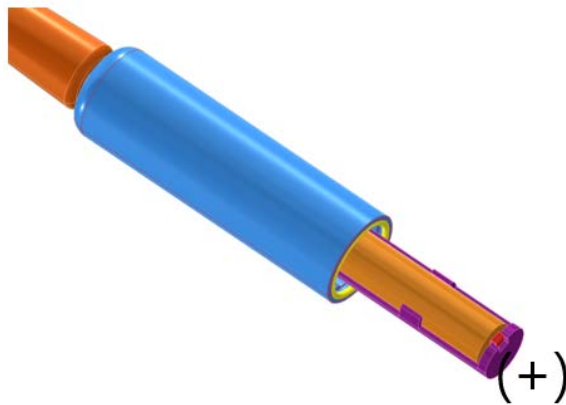

- Turn the battery cap clockwise and insert it.

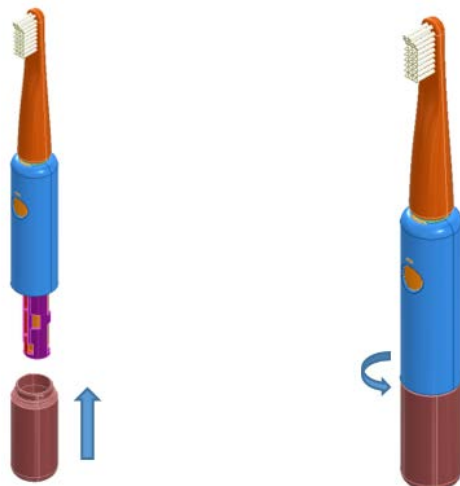

- Failure self-test before service request
 - The smart toothbrush is not working.
 - Check if the battery's polarity is correct, and if it doesn't work even though the polarity is correct, try putting in a new battery.
 - The vibration of the smart toothbrush is weak.
 - Check the vibration intensity in the app. The vibration intensity may be set low.
 - Replace the battery with a new battery (AAA) and try operating it.
 - Smart toothbrushes are not searched when searching for Bluetooth devices.
 - Without a separate interlocking process, once you press the button on the smart toothbrush, the app recognizes and connects the smart toothbrush.
 - The product is not connected or disconnected from the mobile device.
 - Change the Bluetooth function of your mobile device to active.
 - Make sure you have the latest version of the Brush Monster app installed on your mobile device. (If you're not up to date, try updating it to the latest version.)
 - Use it within 10 meters without a wall between the product and the mobile device.
 - Try turning off and on your mobile device and running the Brush Monster app again.
 - The product turns off automatically.
 - For the correct brushing time guide, the product will be automatically turned off after about 160 seconds when interlocked with the mobile device and about 120 seconds when using the smart toothbrush alone.

- Inquiry about product
 - Kitten planet Service center: +82-70-7620-0405
 - Kitten planet Service e-mail: contact@kittenpla.net

- Warranty

- Product name: BRUSH MONSTER Smart toothbrush
- Model: BMT200

This product provides a one-year warranty service from the date of initial purchase. During the warranty period, defective products resulting from the product itself or production process will be repaired or replaced depending on the condition.

Damage due to normal wear and tear, accident or abuse cannot be guaranteed.

If you want to apply for an exchange or remuneration service, please contact the service center of Kitten Planet Co., Ltd. or e-mail.

- Kitten planet Service center: +82-70-7620-0405
- Kitten planet Service e-mail: contact@kittenpla.net

- Caution

- This product is made for children over the age of 3. Infants under the age of 3 need to be careful.
- Biting or chewing the toothbrush cap can cause damage to the toothbrush cap.
- Don't wash or rinse your toothbrush in hot water. It causes toothbrush hair deformation.
- Excessive sterilizer use can deform toothbrush hair.
- Please replace toothbrush hats for one month or for deformation/damage/wearing considering the effect and hygiene of brushing teeth.
- Only use certified AAA batteries.
- If you are being treated for an oral disease, consult a dentist before using the product.
- If excessive bleeding occurs after using this product or if bleeding continues a week later, consult a dentist.
- This product has been tested for electromagnetic disturbance according to domestic safety standards. If you have received a transplant of aids such as heartbeat aids, contact your doctor or manufacturer of aids to fully understand and use the product.
- Use it under the guidance of a guardian when using children or adults with physical or mental disabilities.

- FCC :

This device complies with part 15 of the FCC Rules, Operation is subject to the following two conditions; (1)This device may not cause harmful interference, and (2)This device must accept any interference received, including interference that may cause undesired operation

Any changes or modifications (including the antennas) to the device that are not expressly approved by the manufacturer may void the user's authority to operate the equipment

This equipment complies with FCC RF Radiation exposure limits set forth for an uncontrolled environment. This device and its antenna must not be co-located or operating in conjunction with any other antenna or transmitter.