



METAWELL
变啦

METAWELL Intelligent Body Fat Scale

Instruction Manual

Product Introduction

In addition to measuring body weight, body fat scale can also measure body fat, water, muscle, bone and other components after Bluetooth connection with Metawell app and send the measured data to Metawell.





Test Tips

1 Body fat scales are only suitable for pregnant women in mother and child mode. Other measurement modes are forbidden for pregnant women.

2 In order to ensure the accuracy of the measurement, the intelligent calibration method is needed after the first use

of the body fat scale or the moving body fat scale: light step on the screen, and the calibration can be completed after 5 seconds of rest.

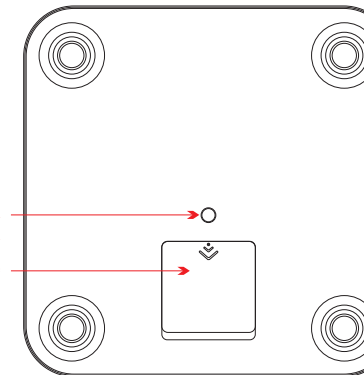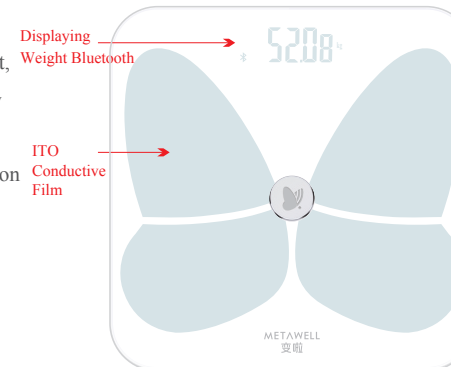

Weight Measurement

Body fat scale can be directly used as weight scale without the Metawell app open, and the weight can be displayed when standing up.

- 1 Stand on the scales as shown in the picture and the LED window displays 0.0 kg. After about 1 second, there are data showing, and constantly changing. After the weight is stable, the data flickers two times and stays fixed which is the final measured weight.

- 2 You can get down the scale when completing the measurement.



Weight Display

*Open the battery lid on the back of the scales before use, and install 4 AA batteries according to the right direction. If "Lo" is prompted when power on, replace the batteries following the above operation.

Body Component Measurement

Body component measurement is required to use Metawell app. Please download and install the app according to the following methods and read the help in the app.



Cautions:

- 1 Stand on the scale after the Bluetooth is turned on.
- 2 Make sure your feet stand on the scales barefoot and have good contact with the conductive film (butterfly part)
- 3 If your feet are too dry to measure, use a wet wipe to wipe the scales surface before measurement.
- 4 Please place the scale on a flat and hard ground before measurement.
- 5 If you need to move the body fat scales for measurement, please carry out intelligent calibration before measuring.



Forbidden for those who has a pacemaker or other metal devices inside the body. Pregnant women do not use.
Age of application: 7-99 years old

APP Testing

After downloading Metawell app, open Bluetooth, ensure network connection, and enter "register" interface.

- 1 New User Registration---Register through mobile phone, after obtaining authentication, setting password and confirming, enter "Set Information" interface.
- 2 Set Information---Set gender, birthday, height to complete your information.
- 3 Measure body component.

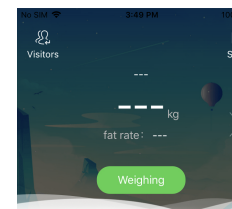


Figure 1

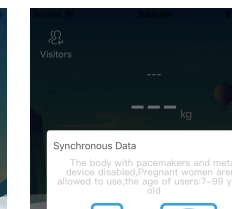


Figure 2

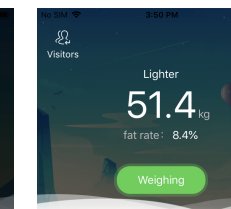


Figure 3

When you enter the main interface (as shown in Figure 1), click on the "Stand on" button in the middle. After displaying the "synchronize data" window (as shown in Figure 2), stand on the scale barefoot and then the Bluetooth sign on the scale LED screen lights up and the weight numbers appear, which is automatically synchronized to the app. Body fat scale LED screen displays from 0.0kg to weight test figures. After twice of flickering, figures fixed and the test is completed. All data can be viewed on the app. The body fat data shows BF when measuring and shows oL when overweight.



Warranty Description

Dear users,
Thank you for using Metawell intelligent body fat scale. In accordance with the Consumer Rights and Interests Protection Law of the People's Republic of China and relevant regulations of the State Administration of Technical Supervision, we will provide you with the following services with this specification and valid machine purchase invoices:

This product is a quality-tested product.
Keep the instructions properly as your certificate of warranty service.
If you need more information or have any questions, please visit our website or call our customer service center at 400-685-8885.

If the machine can not be troubleshooting in the process of use, please follow the following methods to contact for maintenance: Within the warranty period: please contact our after-sales service center for maintenance with the instruction manual and valid machine purchase invoices for repairing products.

Beyond warranty: please contact our after-sales service center for maintenance and pay according to regulations. Do not disassemble the machine yourself. It should be maintained by a professional or a maintenance center. Otherwise, danger may happen.

Hangzhou Metawell Network Technology Co., Ltd.

Warranty Year Limit

The whole machine is free of charge for one year warranty (from the date of purchase, subject to the date of invoice). The following circumstances are not within the scope of free service:

- ❶ Failure to provide initial warranty or invoice
- ❷ Invoice altered.
- ❸ Damage caused by accident or improper use.
- ❹ Damage caused by self repair without the permission of our company.
- ❺ Exceeds valid period of three guarantees

Product Parameters

Name: Metawell intelligent body fat scales

Model No.: m1801

Size: 300mm*300mm*30.5mm

Weight: 2.5kg

Power supplies: 4 AA batteries

Boot mode: weight boot (>5kg)

Maximum weight: 180kg

Display Instruction

Lo : Low power. Please replace the batteries.

oL : Over weight

BF : Body fat measuring

⚙ : Mother & child mode

KG: Kilogram

Download Metawell APP

Scan the QR Code to download Metawell APP



For IOS system users: Please go to the App Store to search for "Metawell" and download.

For Android system users: Please search for "Metawell" in the Android market and download. Login our official website <https://www.bianla.cn> to download the appropriate model.

iOS version: suitable for iOS 8.0 or above.

Android version: suitable for Android4.3 or above.

Bluetooth 4.0 model.

Packing List

Packing List		
No.	Name	Quantity
1	Body fat scale	1
2	AA size battery	4
3	Instruction manual	1

Hangzhou Metawell Network Technology Co., Ltd. Level 5, Building 1, Meidaxiaocheng Tiandi, Jianggan District, Hangzhou City, Zhejiang Province.
 变啦 METAWELL
<https://www.bianla.cn>
After-sale Hotline: 400-685-8885

This device complies with part 15 of the FCC Rules. Operation is subject to the following two conditions:

(1) this device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation. Any changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

NOTE: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation.

If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

The device has been evaluated to meet general RF exposure requirement. The device can be used in portable exposure condition without restriction