



DoseSmart

DS UNIT USER GUIDE

Please read the following before proceeding

Note: This device complies with part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) This device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation. Changes or modifications not expressly approved by DoseSmart could void the user's authority to operate the equipment.


Note: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:


- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

1

Download DoseSmart

Download the “DoseSmart” app from the App Store on your iOS device.

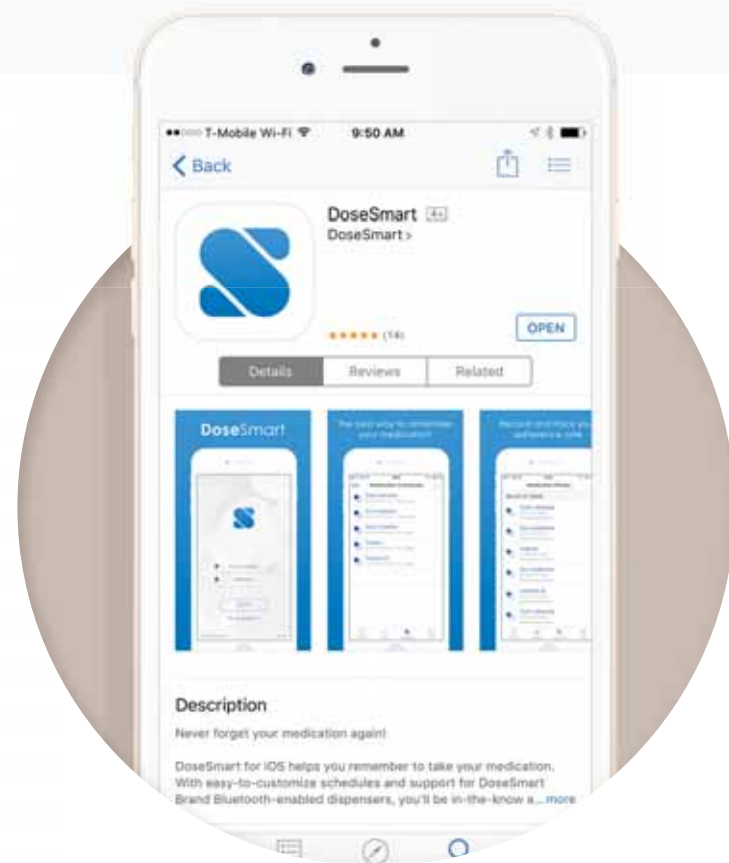

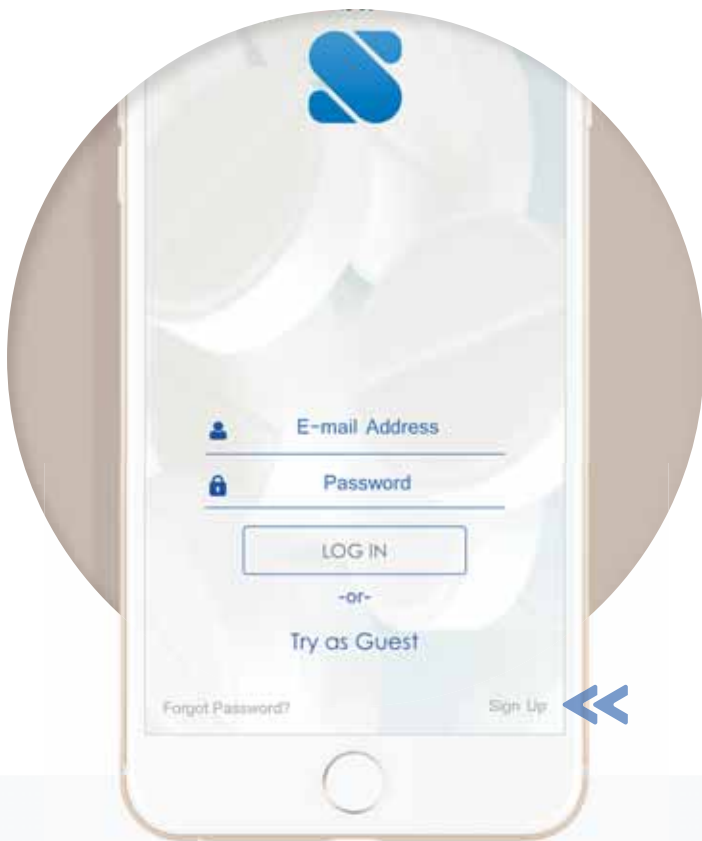

Sign up to DoseSmart


To register to use the DoseSmart app, tap the '**SIGN UP**' button in the lower right hand corner of the screen.


3

Create an account

Tap '**New Account Email Address**' and enter your email address, then tap '**New Account Password**' and create a password of your choosing.

Press the **SIGN UP** button to complete your account creation.

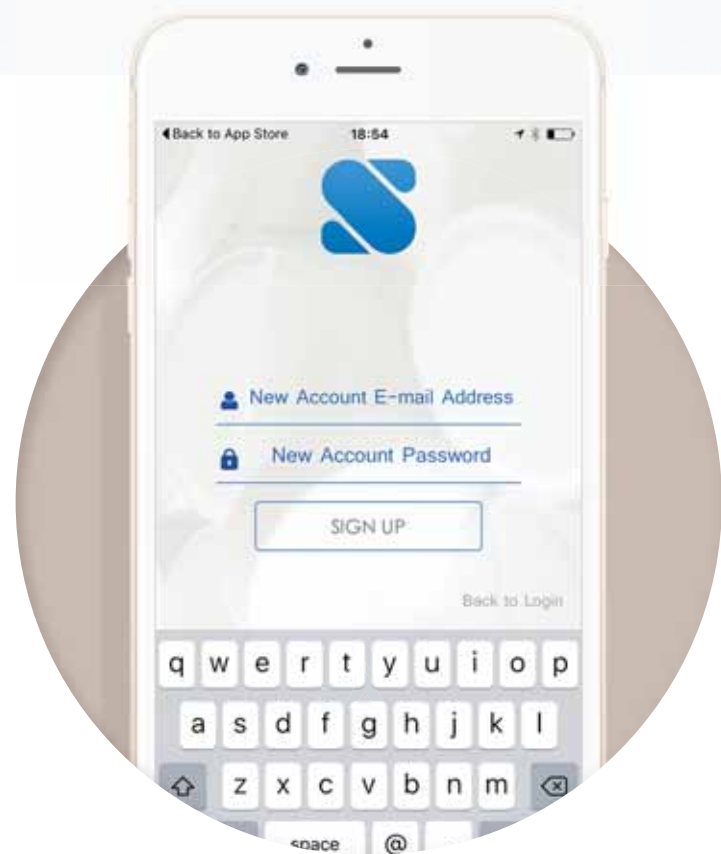

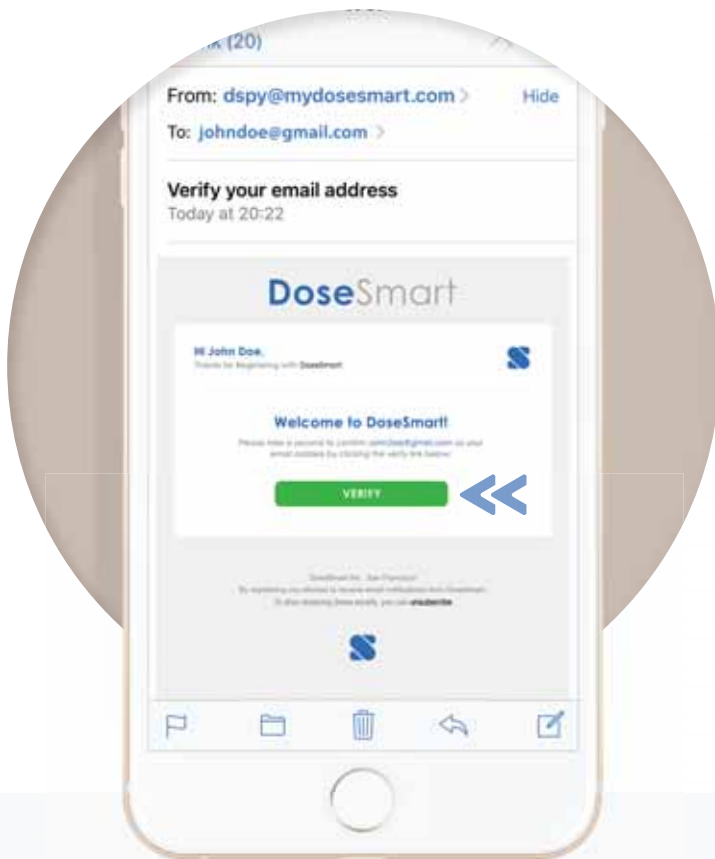


Verify Your Email

You will receive an email from DoseSmart. Please click on the link and verify your email address.

5

Login to your Account

Login with the user name and password
you just created

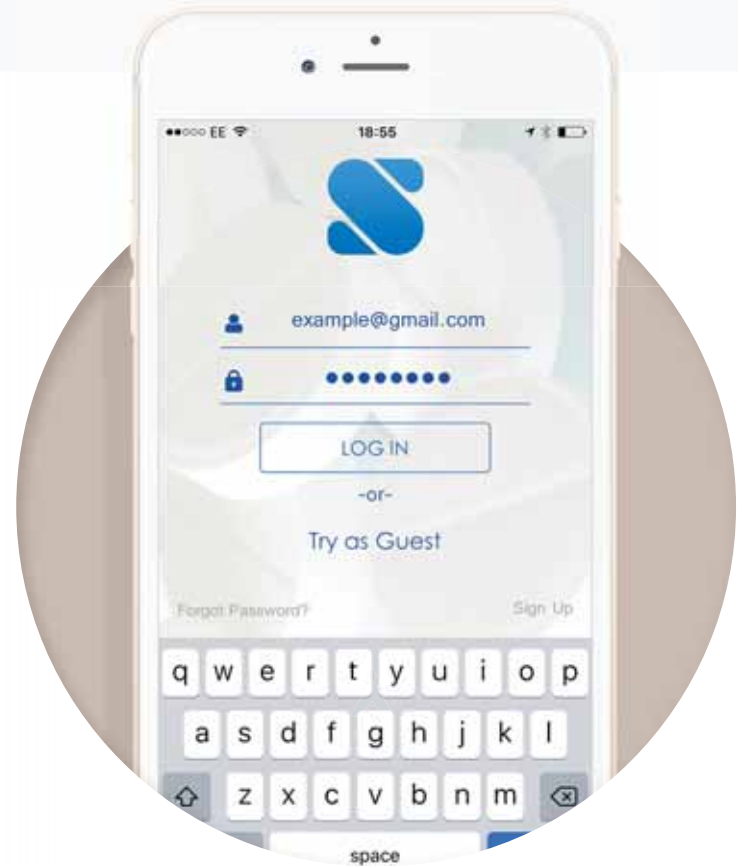



First Launch

When you first launch the DoseSmart App you will begin on the Activity screen.


Initially this page will be empty as you will not yet have any scheduled medications.


Tap on the 'medications' button at the bottom to set up your medications.

7

Add a Medication

When on the medication tab, press the '+' icon in the top right hand corner of the screen to add a new medication.

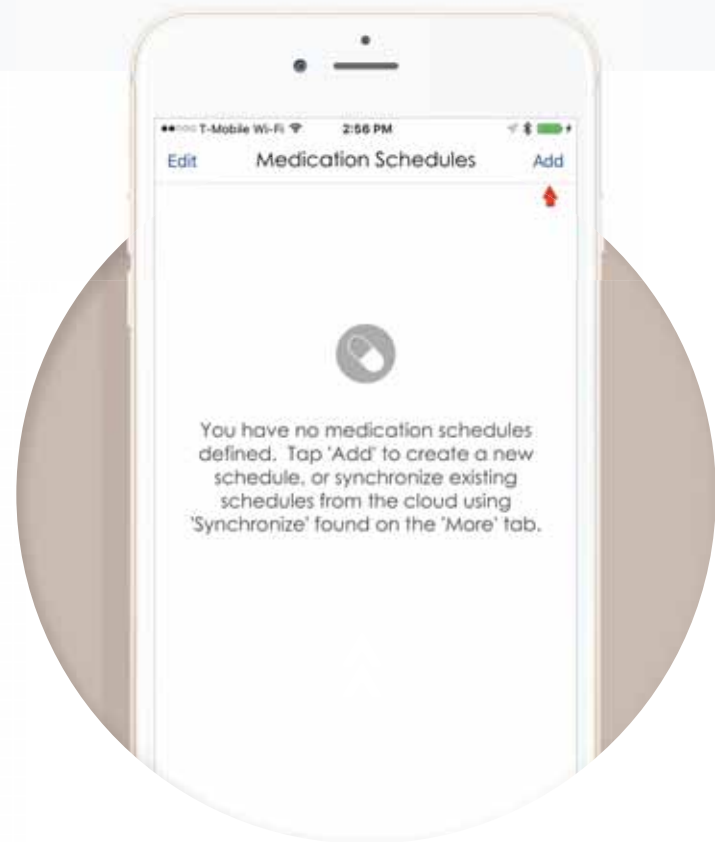

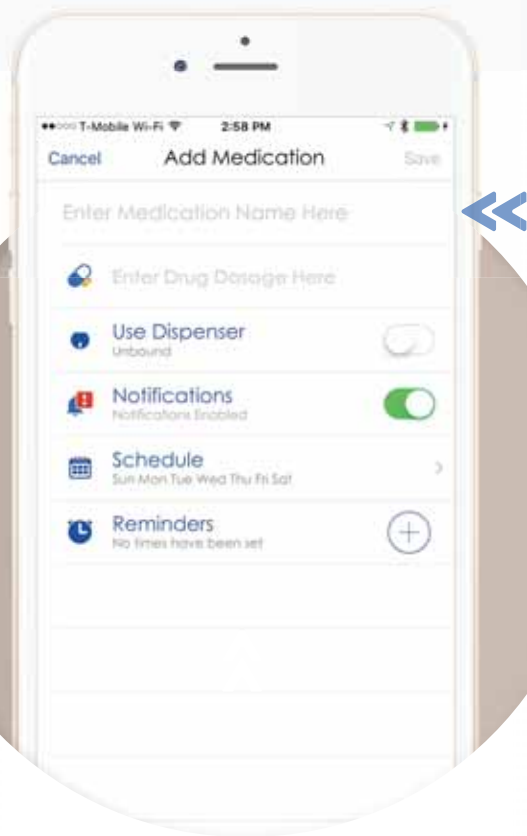


Name your Medication


Tap '**Enter drug name**' and type the name of the medication you would like to add.

9

Enter Drug Dosage

Tap '**Enter Drug Dosage**' and type the quantity of the dosage.

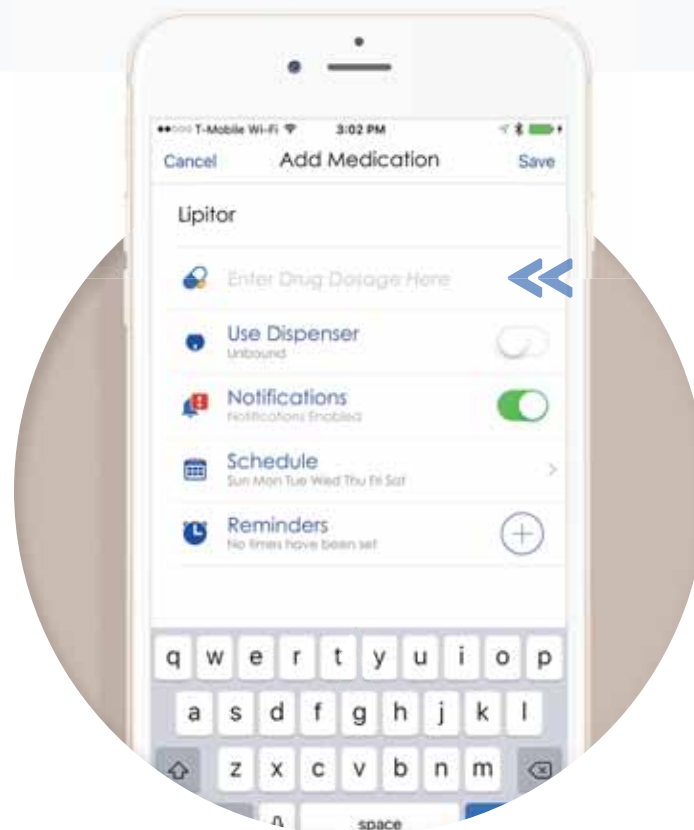

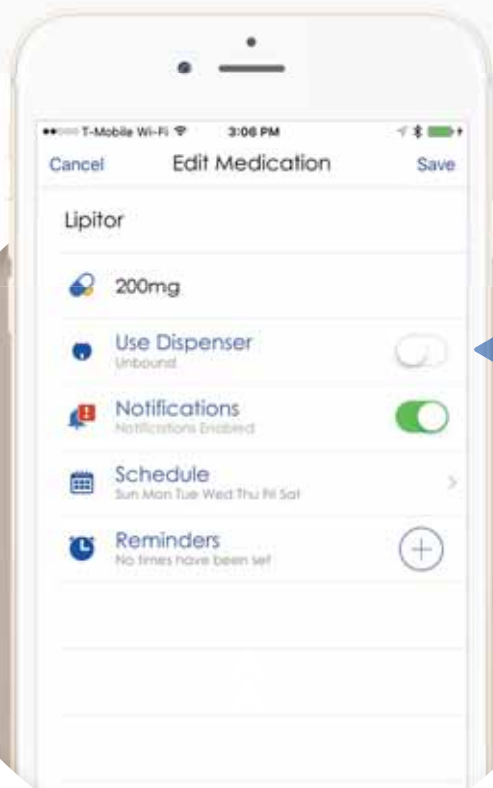


Use Dispenser

The DoseSmart App can be used with or without a DoseSmart dispenser.

To use with a dispenser slide the “Use Dispenser” toggle button to the right.

Start Pairing

When prompted tap the **'Start Pairing'** button to begin syncing your device with the DoseSmart Unit.

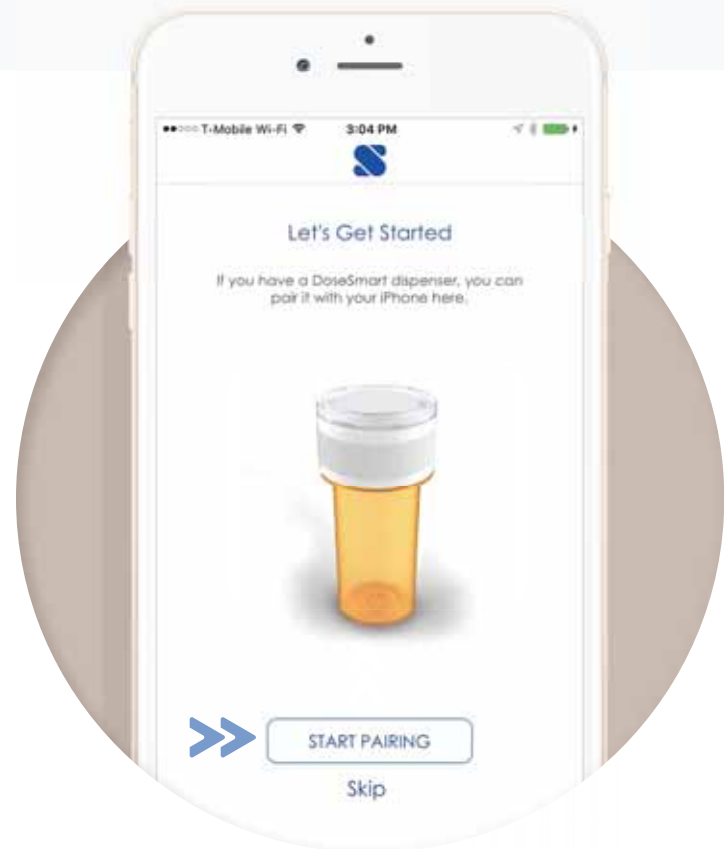



Wake Up Dispenser

Remove the Smart Cap from the base of the bottle

Press the Scan button in the app when the device has been woken up.


Note: Please make sure the DoseSmart unit is charged first.

Device Detected

Once the device is found, press the **'Link Device'** button at the bottom of the screen.

If the device is not found for some reason, repeat the steps and double check to make sure your bluetooth is activated.



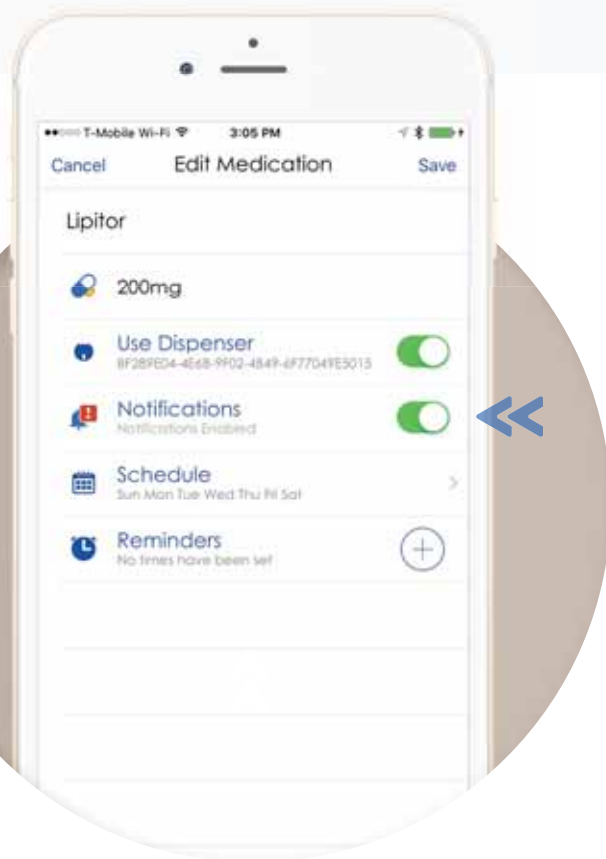

Notifications

Notifications are automatically turned on.

- Sms
- Email
- and phone call reminders


To adjust notification settings go to the **+More** tab and select **Notification Settings**


15

Set Schedule

Your Schedule is set for everyday of the week.

If you only take your medication select days, than you can press Schedule and adjust as needed.

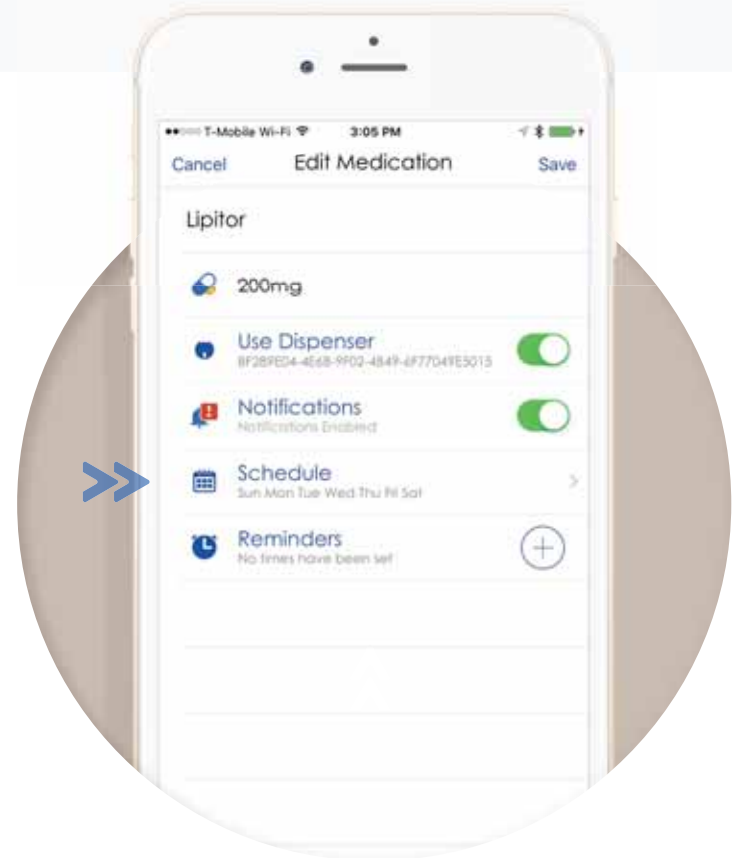

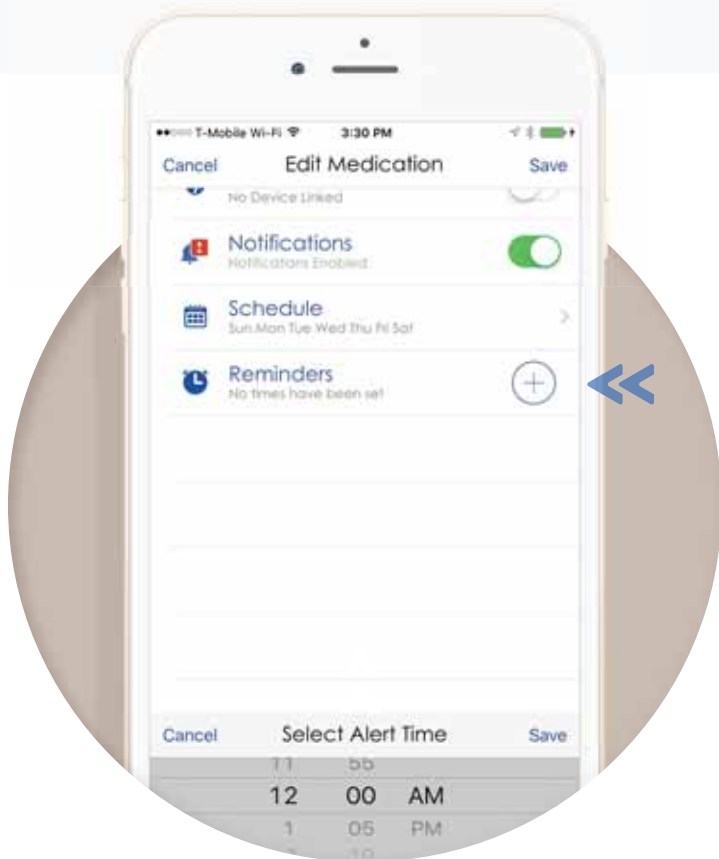



Set Reminders

Simply press the + then select a time and press 'Save' to add a reminder time. You can add as many reminders as needed.

Save Medication

Once you have completed setting up your medication details and schedule, tap **'Save'** in the top right of the screen to complete the process.

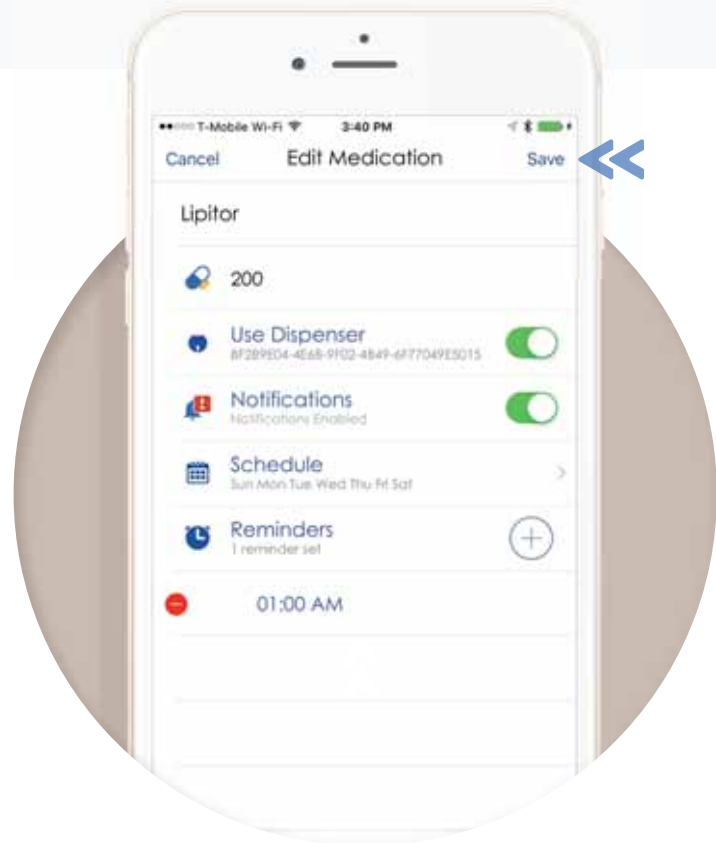

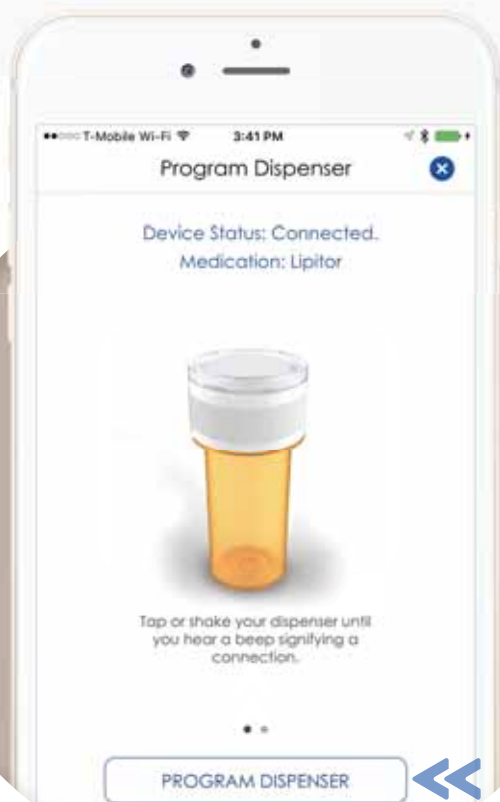

Program Dispenser


You may hear a beep automatically on your device. Otherwise, Tap or shake your dispenser until you hear a beep.

Once you hear a beep then press 'Program Dispenser' button.

Setup Complete

Congratulations. Now you're all set up!

When you take your medication from the device the pill icon will turn green registering that you have taken your medication

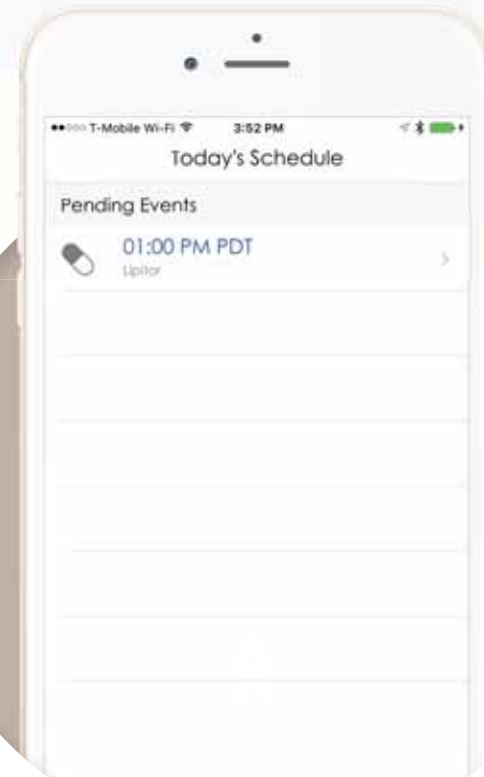

Contact information:

David MacVittie

8186186375

info@mydosesmart.com



DoseSmart

DS UNIT USER GUIDE