


# INSTRUCTION MANUAL

## MASSAGE CHAIR

**YH-7260L-HY**

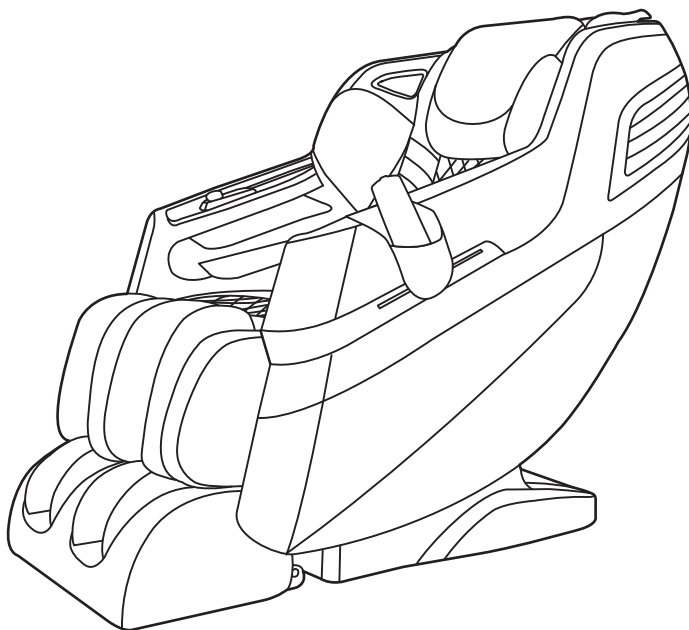

PLEASE BE SURE TO READ THIS MANUAL CAREFULLY TO LEARN THE OPERATIONS OF THE THIS PRODUCT. DO NOT DISPOSE OF THE MANUAL AND KEEP FOR FUTURE REFERENCE.

# IMPORTANT SAFETY INSTRUCTIONS

## OPERATING SURROUNDINGS

- OPERATING ENVIRONMENT
- Do not use or store the massage chair under damp conditions:
  - i.e. bathroom, sauna, spa, to avoid any mold development, electric shock, fire or mechanical failure.
- Don't use the massage chair outdoors
- Don't use the massage chair under dusty and temperature varying conditions.
- Don't place the massage chair near heater(s), stove or direct sunlight to avoid inner components overheating which may lead to electrical fires and shock.
- Make sure to place the massage chair on flat level ground surface
- The chair should be placed on carpet, wood or hard surface flooring
- Make sure the power plug is properly grounded to avoid shock and injury.
- Don't use the machine while smoking



## FOLLOWING CONDITIONS REQUIRING CONSULT FROM A MEDICAL DOCTORS BEFORE USE

- Suffered from malignant tumor
- Serious heart disease
- Serious skin disease
- Bone disease
- While pregnant or menstruating
- On rest by the doctor
- Spinal condition including injuries associated to the back and spine
- Hypertension and cramping
- Arthritic bone disease
- Any implanted medical device such as pacemaker
- Suffering from internal conditions (i.e. gastritis, enteritis, diarrhea, Zepatitis)
- Serious blood circulation trouble
- Consciousness issues such as vertigo and other equilibrium matters

## WARNING

Children under the age of 18 years old or with physical or mental condition are not permitted to use this product unattended by an able adult and with the permission of a medical doctor. These precautions are taken to prevent injury and damage to the product

# IMPORTANT SAFETY INSTRUCTIONS

For safety reasons and in order to clearly indicate potential dangers should the product be used improperly, instructions have been marked in the following manners and should be strictly noted

## ⊘ CAUTIONS

**Non-authorized technicians are not permitted to dismantle or repair this product to avoid fire, electrical shock or serious injury**



**Do not clean or wipe down this product with anything wet or moist when this symbol is shown.**



**When storing or during maintenance, be sure to disconnect the power plug from socket. Do not detach the plug with wet hands**



**Make sure to turn the power switch to the "off" position and detach the power plug when cleaning or storing the product**



**During cleaning or when moving the chair, make sure to detach power cord from the wall and chair to avoid fire or electrical shock  
Do not place the machine on a unstable, uneven and unlevelled surface to avoid damage, insulation failure or electrical shock**



**When storing the machine, make sure to disconnect power plug from the outlet and chair to avoid fire caused by electrical surges.**



**Clean power plug periodically to avoid buildup leading possible fire**



**Do not use it at any unsafe area or for any other purposes than personal use. Not for commercial use.**



**Do not store the chair in areas that are not climate controlled.  
Must be in room temperature with very low humidity.**



**Keep the power cord unobstructed at all time, i.e. not twisted, knotted, over bent or pulled and do not place anything heavy on it**



# IMPORTANT SAFETY INSTRUCTIONS

## CAUTIONS

**Do not make any modifications or change components by an unauthorized person**

---



**Do not dismantle the machine by unauthorized person to avoid electrical shock, short-circuit or fire**

---



**Do not store it near a fire, heater or flammable substances It is one of the major reasons of causing product damage.**

---



**Should any abnormality occur(i.e. burning smell), stop using immediately and detach power plug from outlet.**

---



**please contact your local seller or authorized service agent for maintenance to avoid electrical fire and injury.**

---



**Only proper rated voltage must be used . Stop using immediately if power cord or power plug is damaged or not tightly fitted into the outlet or the base of the chair**

---



**Be sure to keep children and pets away while operating to avoid any injury  
Do not place the chair in direct sunlight at all times**

---



**If the supply cord is damaged,it must be replaced by the manufacturer or an authorized technician**

---



# IMPORTANT SAFETY INSTRUCTIONS

## **RF exposure statement**

This equipment complies with the FCC radiation exposure limits set forth for an uncontrolled environment. This transmitter must not be co-located or operating in conjunction with any other antenna or transmitter.

## **FCC Warning**

This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions:

(1) This device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

Note: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

?Reorient or relocate the receiving antenna.


?Increase the separation between the equipment and receiver.


?Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.

?Consult the dealer or an experienced radio/TV technician for help.


NOTE 2: Any changes or modifications to this unit not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

# NAME OF EACH PART

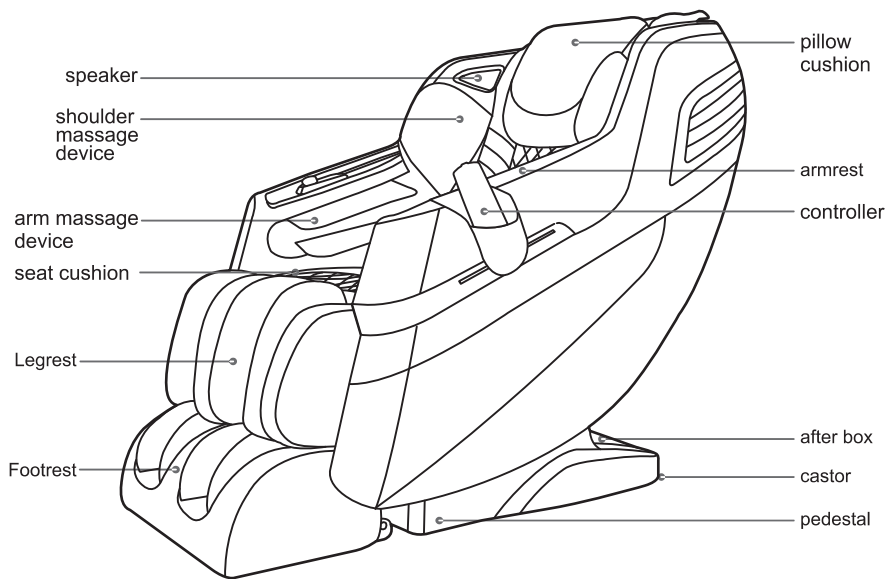

Voice Control Enter			
Serial No	Automatic Program	Voice Control Entry	Answer Entry
1		Hello Message	I Am Here
2		Massage On	Powered On
3		Switch Off	Shut Down
4		Pause Massage	Paused
5		Resume Massage	Massage Resumed
6	Z:G	Zero Gravity On	Zero Gravity On
7		Stretch Mode	Stretch Mode Is On
8	DEMO	Vitality Mode	Vitality Awakening Has Been Turned On
9	GENTLE	Soothing Mode	Soothing Mode Is On
10	NECK&SHOULDER	Shoulder Neck Mode	Shoulder Neck Mode Is On
11	VIGOROUS	Motion Mode	Sport Mode Is On
12	RELAX	Relaxation Mode	Relaxation Is On
13	LOWER BACK	Hip Mode	Beauty Mode Is On
14		Turn On Heating	Heating On
15		Stop Heating	Heating Off
16		Sole Massage	The Sole Roller Is Turned On
17		Stop Foot Massage	Sole Roller Closed
18		Back Up	Backrest Raised
19		Lie Down A Little	Backrest Lowered
20		Backrest Stop	Backrest Stopped
21		Increase Efforts	Efforts Have Been Intensified
22		Reduce Force	The Force Has Been Reduced

# INSTRUCTIONS ON CONTROL PANEL

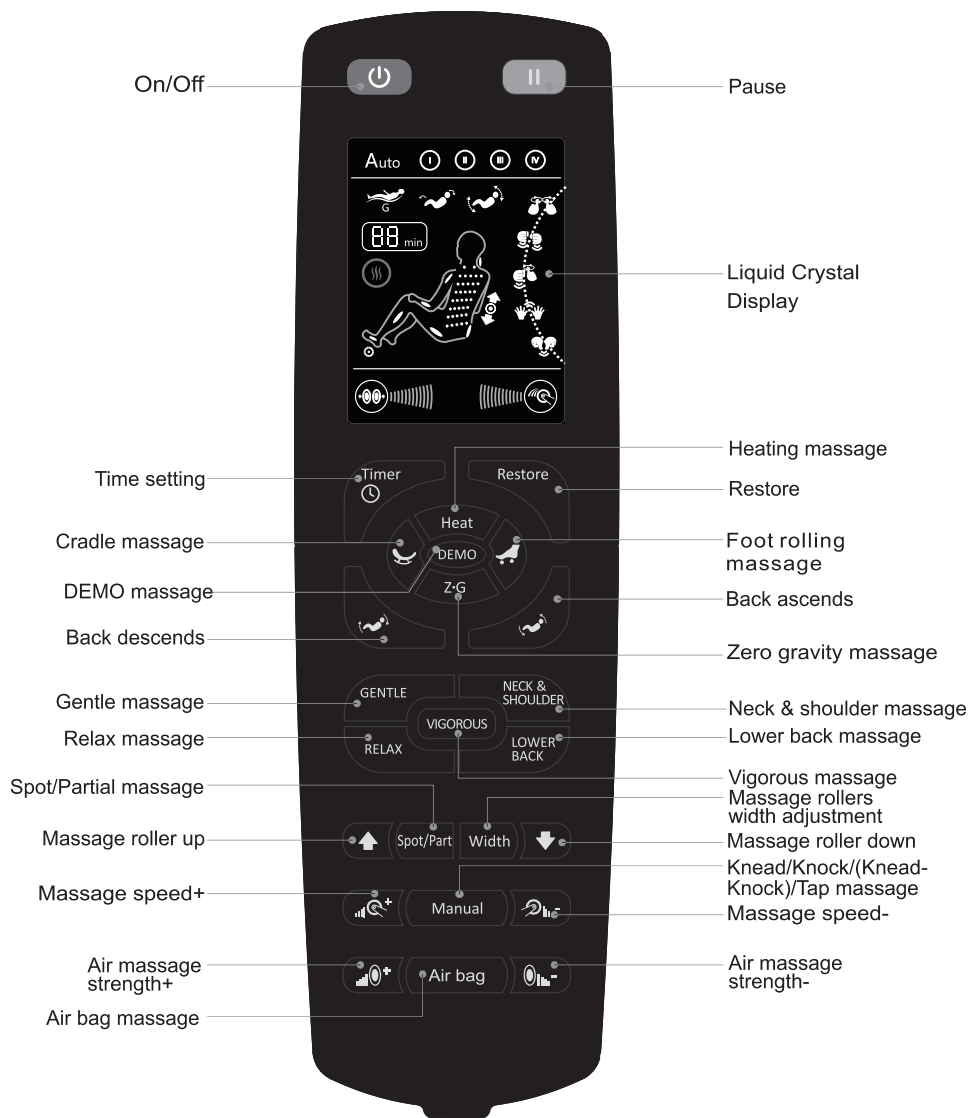


**Warning:** Only rated supply 110V~ must be used to avoid damage to electronic components.

# OPERATION GUIDE

## 1. Powering up the massage chair

first be sure the power switch at the back of the massage chair is in the off “O” position. then insert the power cable into the port to the left of the switch. Once firmly plugged in, proceed to plug the power cord into the outlet. Then put the switch to the on “I” and the remote should light up.

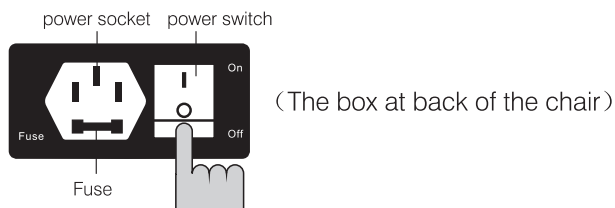

## 2. Starting your massage

Sit in the chair with your buttocks firm pressed against the backrest to ensure that the chair will provide a consistent massage pressure from top to bottom. The press “on” button and select a massage.



On/Off



Time setting

Press the "time setting" button to adjust massage time. You can adjust time to maximum 30 minutes. The unit of increment is 5 minutes.

## 3. Adjusting the angel of the recline

Please adjust the backrest and legrest angle prior to selecting a massage function. Getting into position before selecting the massage function allows for the chair to deliver a more accurate massage. The chair and legrest can be adjusted simultaneously or the independently, depending on the button selected on the remote as illustrated below.



## 4. Stop or finish massage

Normal stop : the timing function automatically fix time is 15 minutes,time out,all the program will stop and keep the setting posture.

Force stop: when press the restore during massage, it will stop massage immediately, press again to continue.



Pause










# OPERATION GUIDE

## 5. electing a preset programmed massage

**DEMO, Gentle, Neck & shoulder, Relax, Lower back, Vigorous, Cradle massage.**

Automatic programs combine shiatsu, kneading, tapping and air massage.

Button	Description	LOWER BACK Main function
	DEMO massage program	Press for a 7 minute demonstration of all available massage techniques.
	Gentle massage program	Kneading massage, air massage and tapping massage. The recline of the backrest can be adjusted during the massage.
	Neck & shoulder massage program	Kneading and beat massage program. The angle of backrest cushion and leg rest cushion could be adjusted during the process.
	Relax massage program	This course combines mainly with quick kneading, tapping and air massage. The angle of backrest cushion and leg rest cushion could be adjusted during the process.
	Lower back massage program	Kneading and scraping massage program. The angle of backrest cushion and leg rest cushion could be adjusted during the process.
	Vigorous massage program	A slower pace massage , delivering a soothing intensity with full airbag massage.
	Cradle massage program	This course combines mainly with cradle massage. The angle of backrest cushion and leg rest cushion could be adjusted during the process.

※ For above programs, the fix time is 15 minutes and the default of airbag intensity is low.







※ The above automatic full-body massage courses can not be combined with other programs.

# OPERATION GUIDE

## Manual setting course

### Introduction of manual setting course






Partial air massage is the course which targets different needs and areas of body.

Button	Description	Main function and adjustment guide
	 Kneading	Press button one time to start kneading massage. The speed is adjustable.
	 Knocking	Press two times to start the knead-knocking massage. The speed is adjustable.
	 Knead-knocking	Press three time to start the knocking massage. The speed and width is adjustable.
	 Tapping	Press four times to start the tapping massage. The speed and width is adjustable.
	 Shiatsu	Shiatsu massage.

※ The above manual setting course could not be combined with other courses.

※ For the above manual setting course, the default setting for the massage is spot-specific with medium width. Speed and width is adjustable by pressing relative button.





### Function introduction for manual setting course

Button	Description	Main function
	Speed adjustment	There are five levels to choose.
	Width adjustment	There are three levels to choose from narrow, medium and wide.
	Spot /Partial massage	Press the button one time the rollers will stay in one location. Press two times to start the partial massage. The third time to start full course massage .
 	Message roller move up Message roller move down	Press down the button to move message roller, release the button to stop.

# OPERATION GUIDE

## Functions on buttons

### Introductions on buttons used for manual setting course

Button	Description	Main function
	Air pressure intensity adjustment	Air pressure intensity is adjustable It could be adjusted only in air pressure massage. There are 5 levels of strength can be adjusted.
	Air massage	Press one time to start air massage for upper body. Press two times to start air massage for lower body. Press three times to start air massage for full body. Close by four times.
	Rolling massage	Foot rolling control key. Press the button to start foot rolling massage, press again to stop.
	Heating function	Waist heating function. Press the button to start waist heating massage, press again to stop

## Music Function

**Connecting your mobile device:** Turn on bluetooth of mobile phone or tablet, then select bluetooth device labeled (code 7260-HY, password “9999” ), after it syncs with massage chair speakers, you can now play your favorite tunes or audio books from your massage chair. Enhance your massage experience!

**⚠ Warning:** Supervision and guidance are necessary when people with disabilities and mental illness. Please advise with your doctor before using the chair if you have ANY medical conditions. Children under 18 should not use the chair and consult with a doctor for approval if attempted. Do not operate the chair if there are any children or pets near by. Serious injuries may occur if these instructions are not followed.

# SAFETY INSTRUCTIONS

1. Should any abnormality occur, shut off and detach power plug from socket.
2. When power cord damaged, it must be changed from special power cord or available from other manufacturer or service department.



**Warning:** Do not dismantle, modify or change components by unauthorized person to avoid electric shock or fire.

## DAILY CARE AND MAINTENANCE

### Storage guide

- ⦿ Do not dismantle any part the chair or make modifications. Serious injury may occur.
- ⦿ Do not attempt to repair the machine by unauthorized person to avoid electric shock or fire.
- ⦿ Make sure to shut off and detach the power plug when cleaning to avoid electric shock or hazard.
- ⦿ Keep the massage chair clean and free dust. Do not place the massage chair in a humid area.
- ⦿ If not in use for a long period, please cover it with a cloth to protect from dust.
- ⦿ Do not place the massage chair under sunlight or heating source like a heater or directly by a window that has direct sunlight.

### Maintenance

- ⦿ Use a neutral detergent only that is free of acid and alcohol.
- ⦿ Use a neutral detergent to clean the back cushion, head cushion and leg cushion. Let it dry.
- ⦿ With a damp cloth wipe down the plastic parts, pipe and Synthetic leather material of the chair.
- ⦿ Make sure not to use any gasoline, amyl acetate, benzene based cleaners to avoid leather color fading or discoloring.
- ⦿ Do not spray pesticide or aerosols near or on the chair to avoid electric shock or fire.
- ⦿ Please clean with dry cloth to avoid mechanical failure.

### Controller maintenance

- ⦿ Do not place the massage chair near heater or fire.
- ⦿ Keep away from humidity to avoid electric shock and damage to the chair.
- ⦿ Please clean dust with soft cloth or small banister brush.
- ⦿ Make sure not to use detergents containing benzene, gasoline, amyl acetate.

# TROUBLE SHOOTING

Checking is necessary in frequent usage.

Make sure the machine work normally before usage, after not in used for long time.

Abnormality situation	Solution guide
Failure to start and indicator light is off	<ul style="list-style-type: none"><li>⊙ Check the plug is well connected</li><li>⊙ Check the switch on the chair is turned on in the “I” position</li></ul>
Sound in operation	<p><b>The following sounds are normal</b></p> <ul style="list-style-type: none"><li>⊙ The mechanism will produce a sound for kneading or knocking</li><li>⊙ Sound changes after speed adjustment</li><li>⊙ Sound from when the rollers go up or down</li><li>⊙ When the chair does a tapping massage it will make the chair vibrate making a sound</li><li>⊙ When the air inflates and deflates it will make a swooshing sound</li><li>⊙ Sound for airbag inflating and venting. Airbags may crackle</li><li>⊙ Sound for leg adjustment</li><li>⊙ Sound for motor when the massage chair is loaded</li></ul>

# SPECIFICATION

Product name	Massage Chair
Model	YH-7260L-HY
Power supply	110V~ 60Hz
Rated power	150W
Rated working time	15min