



Meet delights for your home here



Home of Delights

Massage Chair  
Fauteuil de Massage

**JL10053**

USA office: Fontana

AUS office: Truganina

FRA office: Saint Vigor d'Ymonville

GBR office: FDS Corporation Limited, Unit 4, Blackacre Road, Great Blakenham, Ipswich, Ip6 0FL, United Kingdom

THIS INSTRUCTION BOOKLET CONTAINS **IMPORTANT** SAFETY INFORMATION.  
PLEASE READ AND KEEP FOR FUTURE REFERENCE.



### Before You Start

- ⚠ Please read all instructions carefully.
- ⚠ Retain instructions for future reference.
- ⚠ Separate and count all parts and hardware.
- ⚠ Read through each step carefully and follow the proper order.
- ⚠ We recommend that, where possible, all items are assembled near to the area in which they will be placed in use, to avoid moving the product unnecessarily once assembled.
- ⚠ Always place the product on a flat, steady and stable surface.
- ⚠ Keep all small parts and packaging materials for this product away from babies and children as they potentially pose a serious choking hazard.

## IMPORTANT SAFETY INSTRUCTIONS

### OPERATING SURROUNDINGS

- OPERATING ENVIRONMENT
- Do not use or store the massage chair under damp conditions; i.e. bathroom, sauna, spa, to avoid any mold development, electric shock, fire or mechanical failure.
- Don't use the massage chair outdoors.
- Don't use the massage chair under dusty and temperature-varying conditions.
- Don't place the massage chair near heater(s) stove or direct sunlight to avoid inner components overheating which may lead electrical fires and shock.
- Make sure to place the massage chair on flat level ground surface.
- The chair should be placed on carpet, wood or hard surface flooring.
- Make sure the power plug is properly grounded to avoid shock and injury.
- Don't use the machine while smoking



### FOLLOWING CONDITIONS REQUIRING CONSULT FROM A MEDICAL DOCTOR BEFORE USE

- Suffered from malignant tumor
- Serious heart disease
- Serious skin disease
- Bone disease
- While pregnant or menstruating
- On rest by the doctor
- Spinal conditions including injuries associated to the back and spine
- Hypertension and cramping
- Arthritic bone disease
- Any implanted medical device such as pacemaker
- Suffering from internal conditions(i.e.gastritis, enteritis, diarrhea, hepatitis)
- Serious blood circulation trouble
- Consciousness issues such as vertigo and other equilibrium matters

### WARNING

Children under the age of 18 years old or with physical or mental conditions are not permitted to use this product unattended by an able adult and with the permission of a medical doctor. These precautions are taken to prevent injury and damage.

## IMPORTANT SAFETY INSTRUCTIONS

For safety reasons and in order to clearly indicate potential dangers should the product be used improperly, instructions have been marked in the following should be strictly noted.

### CAUTIONS

Non-qualified technicians are not permitted to dismantle or repair this product to avoid fire, electrical shock or serious injury.



Do not clean or wipe down this product with anything wet or moist when this symbol is shown.



When storing or during maintenance, be sure to disconnect the power plug from socket. Do not detach the plug with wet hands.



Make sure to turn the power switch to the "off" position and detach the power plug when cleaning or storing the product.



During cleaning or when moving the chair, make sure to detach power cord from the wall and chair to avoid fire or electrical shock Do not place the machine on a unstable, uneven surface to avoid damage, insulation, failure or electrical shock.



When storing the machine, make sure to disconnect power plug from the outlet and chair to avoid fire caused by electrical surges.



Clean power plug periodically to avoid buildup leading to possible fire.



Do not use it at any unsafe area or for any other purposes than personal use. Not for commercial use.



Do not store the chair in areas that are not climate controlled. Must be in room temperature with very low humidity.



Keep the power cord unobstructed at all time, i.e. not twisted, knotted, over bent or pulled and do not place anything heavy on it.



## IMPORTANT SAFETY INSTRUCTIONS

### CAUTIONS

Do not make any modifications or change components by an unqualified person.



Do not dismantle the machine by unauthorized person to avoid electrical shock, short-circuit or fire.



Do not store it near a fire, heater or flammable substances it is one of the major reasons of causing product damage.



Should any abnormality occur(i.e. burning smell), stop using immediately and detach power plug from outlet.



Please contact your local seller or authorized service agent for maintenance to avoid electrical fire and injury.



Only proper rated voltage must be used. Stop using immediately if power cord or power plug is damaged or not tightly fitted into the outlet or the base of the chair.



Be sure to keep children and pets away while operating to avoid any injury. Do not place the chair in direct sunlight at all times.



If the supply cord is damaged, it must be replaced by a unqualified technician.



## IMPORTANT SAFETY INSTRUCTIONS

### RF exposure statement

This equipment complies with the FCC radiation exposure limits set forth for an uncontrolled environment. This transmitter must not be co-located or operating in conjunction with any other antenna or transmitter.

### FCC Warning

This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions:

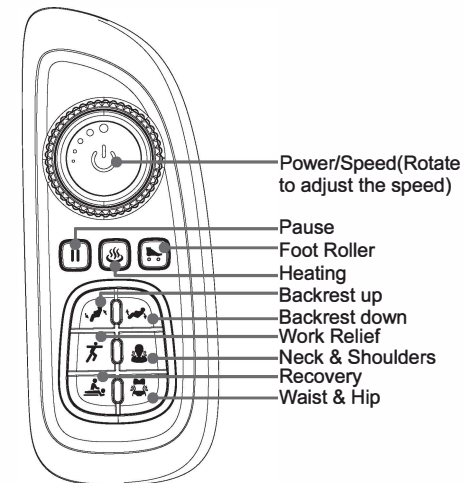
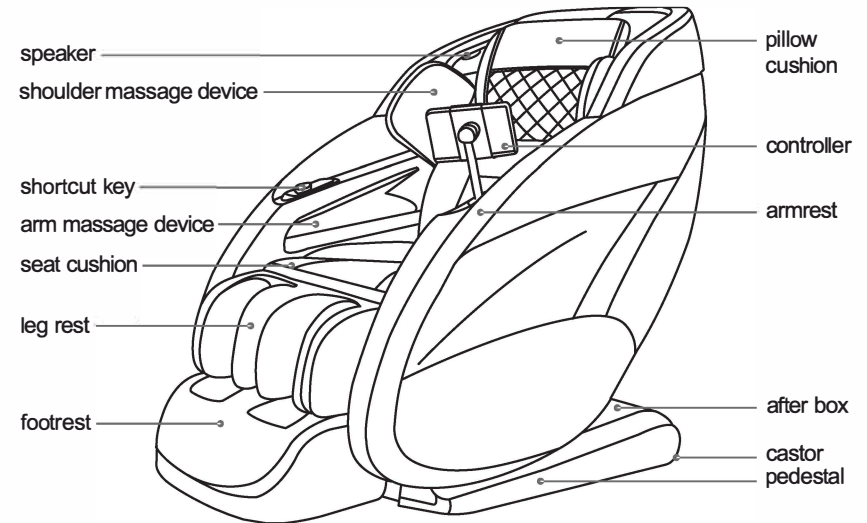
(1) This device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

**NOTE 1:** This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

**NOTE 2:** Any changes or modifications to this unit not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

## NAME OF EACH PART



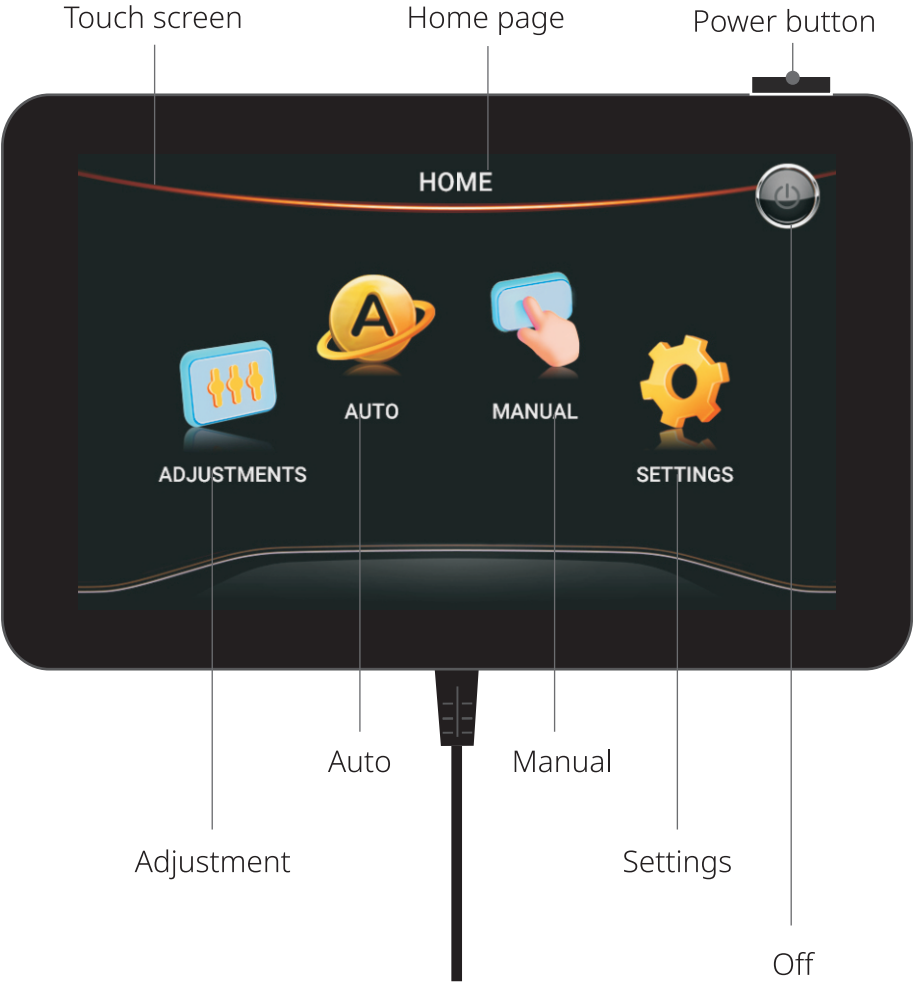
**SHORTCUT KEYS**  
Inside the right handrail

**Warning:** Only rated supply 110V must be used to avoid damage to electronic components.

OPERATION GUIDE





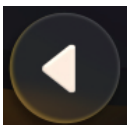


CONTROLLER MENU

Always press the Power Button to start massage.




TABLET CONTROLLER INSTRUCTIONS

1. BASIC KEYS



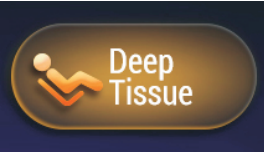
KEYS	DESCRIPTION	GUIDELINES
	Off	Touch to stop massage and shut off the screen.
	Home	Touch to return to the home page.
	PAUSE	Touch once to halt current massage and countdown timer.
	RESTART	Touch to resume previous massage functions, as well as countdown timer.
	BACK	Touch to return to the previous page.
	ADJUSTMENTS	Touch to <b>ADJUSTMENTS</b> page.
	Air Massager	Touch to <b>AIR MASSAGER</b> page.

## 2. HOME PAGE MENU

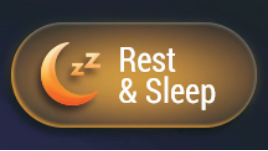
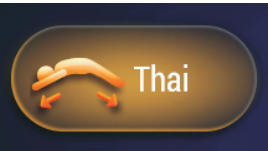
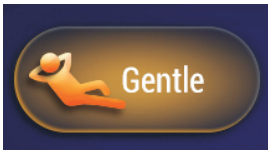
KEYS	DESCRIPTION	GUIDELINES
 <b>AUTO</b>	Auto	Touch to choose from 12 <b>AUTO</b> massage programs. <b>Work Relief, Recovery, Deep Tissue, Rest&amp;Sleep, Thai, Gentle, Athletic, Vertebral, For Her, For Him, Neck &amp; Shoulders, Waist&amp;Hip.</b>
 <b>MANUAL</b>	Manual	Touch to go to the <b>MANUAL</b> massage programs page, choose your preferred massage functions at desired intensities, fully customize your massage.
 <b>ADJUSTMENTS</b>	Adjustments	Touch to go to the <b>ADJUSTMENTS</b> page for zero gravity and massage intensities.
 <b>SETTINGS</b>	Settings	Touch to go to the system <b>SETTINGS</b> page for timer and other preferences.


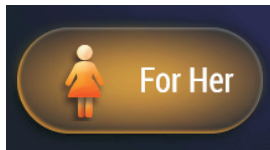
## 3. AUTO MASSAGE PROGRAMS




KEYS		GUIDELINES
	Work Relief	This is an exclusive massage program for the workplace, according to the features of long time desk work and business trip, it can ease muscles and recuperate body.
	Recovery	A full body deep massage includes Shiatsu, Knocking and Tapping to ease and relax muscles after surgeries.
	Deep Massage	A full body deep massage includes Shiatsu, Knocking and Tapping to ease and relax muscles after intense exercise, especially promote the body circulation system to improve the muscle recovery after exercise.






	Rest & Sleep	<p>A program designed to provide a mild massage that is great for deep relaxation. This is the program to help improve the quality of lunch break and sleep. The mild massage kneads body parts from heavy strength to light one, from light to no strength to help soothe tired body gradually.</p> <p>Product features front and rear swing function that takes you into the fantasy space swinging state and a thorough relaxation, best used before going to sleep.</p>
	Thai Stretch	<p>This massage will stretch your calf muscles and hamstrings with a combination of gentle tapping, kneading, rolling, and air compression.</p> <p>This massage is inspired from Thai Massage, the strong kneading and stretch of the body can effectively stretch legs, it's the best choice to relax tired muscles and increase physical vigor.</p>
	Gentle	A full body gentle massage includes Kneading and Heating, designed for daily use.

	For Athletics	An auto programme generated by the system, focus on the tension points collected from athletics, especially promote the body circulation system to improve the muscle recovery after exercise.
	Full Back Massage	A full body massage, mainly focus on your full back with a deep massage that includes air compression to relieve the symptoms of lumbar spine pain.
	For Her	A full body gentle massage includes Kneading and Heating, specially designed for female.
	For Him	A full body strong massage includes Kneading, specially designed for the male.
	Neck & Shoulders	A full body massage, mainly focus on your neck and shoulders with a slow massage that incorporates soothing heat to relieve neck and back pain symptom. The curve design of the chair backrest makes the massage rollers completely fit the shoulders and the neck, which can effectively relieve fatigue and stiffness.

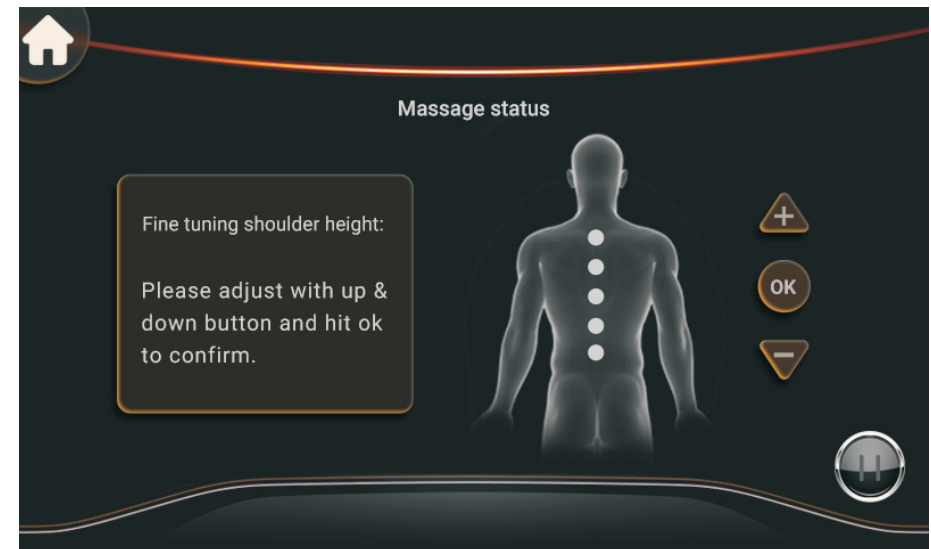
	<p>Lower Back, Waist and Hips</p>	<p>A full body massage, mainly focus on your lower back, waist and hips with a deep massage that includes air compression to relieve full body pain.</p>
---	-----------------------------------	--

## Backrest Scanning Instructions

A body scan is performed every time the power is on and the user selects an Auto Mode, it is designed to ensure an accurate and personalized massage for every individual. After the power is on, sit on the chair slowly with your back against the chair back, your head against the pillow so the massage sensors can detect your Shoulder Height and the whole back position.

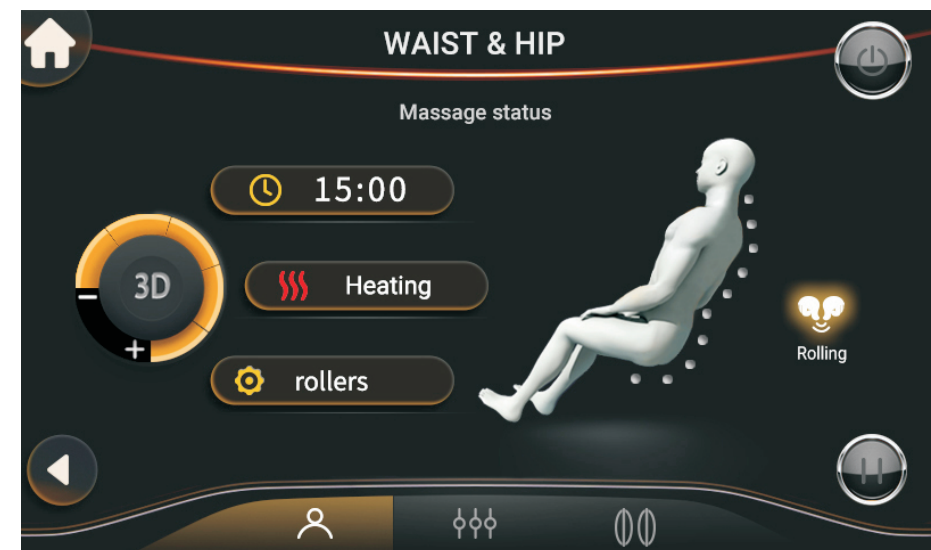
The Tablet Control gives beep sounds when completing a body scan. When micro-adjustment is needed, just touch Roller up/down   , then  to start massage. Make sure the body scan is finished correctly for the best massage.

Once body scan is completed, it will skip to the massage status page without further actions.



### Note:

Heating, foot rollers can be turned on/off manually during any auto mode, see below.



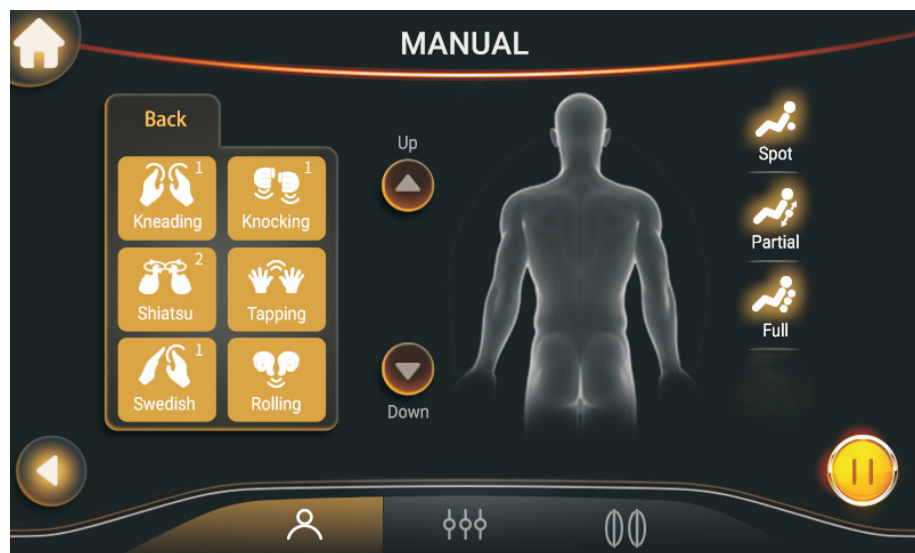


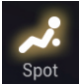
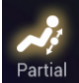
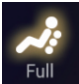

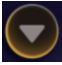
## 4. MANUAL MASSAGE PROGRAMS

6\*2 massage functions on demand.

### Kneading, Shiatsu, Swedish, Knocking, Tapping, Rolling.

- Touch to choose a manual massage function for a full back massage at medium strength.
- Manual setting modes can not be combined with other modes.
- The default setting is spot-specific massage with medium width.



KEYS	DESCRIPTION	GUIDELINES
 	Roller Spot/ Partial Massage	<p>Touch Spot to activate localized Spot Massage while the original massage mode continues to run at the latest position.</p> <p>Touch Partial to activate Partial Massage, rollers take the latest massage spot position as the center, and will travel approximately 3 inches up or down. Only effective under manual massage mode.</p>
	Full Back Roller	Full Back Massage.
 	Roller Up or Down	Up or down adjustment.