


THERAGUN

Multi-Device Wireless Charger

Theragun Multi-Device Wireless Charger


The Theragun Multi-Device Wireless Charger is compatible with Theragun PRO™, Theragun Elite™, and Wave Roller™. Only ONE device can be charged at a time.

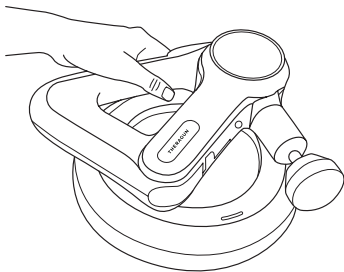

Charging your Theragun PRO or Theragun Elite

1. Plug the power adapter into the Theragun Multi-Device Wireless Charger's charging port.
2. Turn off your Theragun PRO or Theragun Elite.
3. When holding the Theragun PRO or Theragun Elite by the top handle, attachment away from you, place the left side of the device on the plate of the Theragun Multi-Device Wireless Charger, following the shape on the top surface.
4. The LED light of the Theragun Multi-Device Wireless Charger will light up when properly connected. The Theragun PRO and Theragun Elite are fully charged when the battery icon on the device's OLED screen is full. The Theragun PRO and Theragun Elite can remain on the Theragun Multi-Device Wireless Charger even when fully charged.

Charging your Wave Roller

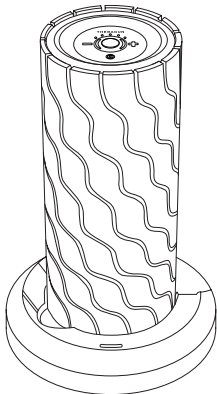

1. Plug the power adapter into the Theragun Multi-Device Wireless Charger's charging port.
2. Turn off your Wave Roller.
3. Place the base of your Wave Roller on the plate of the Theragun Multi-Device Wireless Charger, with the traditional charging port and functional buttons facing upwards.
4. The LED light of the Theragun Multi-Device Wireless Charger will light up when properly connected. The Wave Roller is fully charged when the battery indicator light is solid green. The Wave Roller can remain on the Theragun Multi-Device Wireless Charger even when fully charged.

THERAGUN

Born in Los Angeles, CA.
Engineered for every**body**.