



























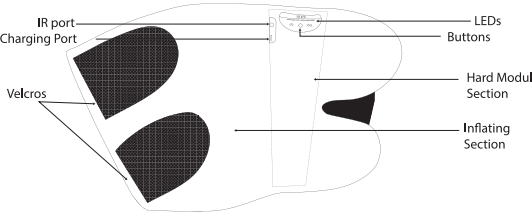















<div></div> <div>Spryng™ User Guide</div> <div></div> <div><p>This device complies with Part 15 of the Federal Communications Commission Rules (47 C.F.R. Part 15). Operation is subject to the following two conditions:</p><p>(1) This device may not cause harmful interference, and</p><p>(2) this device must accept any interference received, including interference that may cause undesired operation.</p><p>That the product is manufactured at a facility registered with the Food and Drug Administration.</p></div>	<div>TABLE OF CONTENTS</div> <div><ul style="list-style-type: none"><li>• Warnings</li><li>• Labels</li><li>• Know your SPRYNG™ device</li><li>• Operating Instructions</li><li>• Know The Led Indications</li><li>• Caring for Your SPRYNG™</li><li>• Technical Information</li><li>• Troubleshooting</li></ul></div>	<div>WARNINGS</div> <div><p>IMPORTANT SAFETY INSTRUCTIONS: Read the entire instruction manual before using the SPRYNG™ devices. These warnings apply to the appliance and, where applicable, to all accessories.</p><p>WARNING: Keep away from water. As with most electrical appliances, electrical parts in this product are electrically live even when the switch is off.</p><p>1.To reduce the risk of electric shock, do not use the product near or place in water.</p><p>2.Do not place or store where it can fall into a bathtub, sink, or washing machine.</p><p>WARNING! If you experience severe pain, any unusual symptoms or want to remove excessive pressure during use:</p><p>1.Stop the device by pressing the Power Button.</p><p>2.Remove the SPRYNG™ device(s) from your limb(s).</p><p>3.Promptly consult your licensed healthcare practitioner.</p><p>Do not attempt to take apart the system. If repairs are required, it must be done by the manufacturer. Please contact customer service at support@spryngme.com for more information. Do not use the product if it is damaged or requires repairs.</p><p>Use the product only as recommended and with the charger and power cord provided. The use of other chargers and power chords may harm the operation of the system.</p><p>Do not place the system, charger, charging cable, or any accessories where they could present a fall hazard. ☐</p><p>Keep the USB charging ports of the SPRYNG™ system free of debris.</p><p>Do not puncture or otherwise damage the inflatable parts of the SPRYNG™ device. The system should not be used if damaged in any way.</p><p>DO NOT USE this product if you are experiencing inflammation, an infection, pain of unknown origin, bleeding (internal or external) at or near the site of application, or if you have a wound at or near the site of application.</p></div> <div>1</div>	<div>DO NOT USE this product if your doctor advises against the use of Spryng.</div> <div><p>If you have any existing medical condition consult your doctor before using Spryng.</p><p>DO NOT USE this product if you are allergic to any of the materials, and it should be used only on intact/undamaged skin.</p><p>DO NOT USE this product if you have any of the following conditions:</p><ul style="list-style-type: none"><li>• Acute pulmonary edema</li><li>• Acute thrombophlebitis</li><li>• Acute congestive cardiac failure</li><li>• Acute infections</li><li>• Deep vein thrombosis (DVT)</li><li>• Episodes of pulmonary embolism</li><li>• Wounds, lesions, or tumors at or near the site of application</li><li>• Where increased venous and lymphatic return is considered unsafe</li><li>• Bone fractures or dislocations at or near the site of application</li><li>• Skeleton injuries and diseases</li><li>• Bone fractures in area of usage</li><li>• Currently being treated for Cancer</li><li>• High/abnormal blood pressure</li><li>• Severe atherosclerosis or other cardiac diseases</li><li>• Congestive heart failure</li><li>• Gangrene, untreated or infected wounds, recent skin graft, and dermatitis in the legs</li><li>• Known presence of malignancy in the legs</li><li>• Active limb infections, including cellulitis</li><li>• Presence of Lymphangiosarcoma</li></ul><p>Do not use this product:</p><ul style="list-style-type: none"><li>• During pregnancy or for at least six months after giving birth</li><li>• If you have an implantable device such as a cardiac pacemaker</li></ul><p>Children should not use this product. Do not use on anyone who is unconscious or incapacitated as this may present an increased risk of harm.</p><p>Each SPRYNG™ unit contains a rechargeable Lithium-ion Battery. It should be disposed in accordance with any applicable local regulations. It should be disposed with the unit once the battery fails to hold a charge.</p></div> <div>2</div>	<div>LABELS</div> <div><p>The following labels and symbols appear on the SPRYNG™ recovery wrap and removable innerLayer.</p><table><tr><th>Symbol</th><th>Description</th><th>Location</th></tr><tr><td></td><td>Read instructions before use</td><td>On SPRYNG™ calf wrap</td></tr><tr><td></td><td>Fragile, handle with care</td><td>On SPRYNG™ calf wrap</td></tr><tr><td></td><td>Keep dry</td><td>On SPRYNG™ calf wrap</td></tr><tr><td></td><td>Do not wash</td><td>On SPRYNG™ recovery wrap</td></tr><tr><td></td><td>Separate collection for waste electrical and electronic equipment</td><td>On SPRYNG™ calf wrap</td></tr><tr><td></td><td>Do not dry-clean</td><td>On SPRYNG™ calf wrap and removable inner layer</td></tr><tr><td></td><td>Serial number of the console</td><td>On SPRYNG™ calf wrap</td></tr><tr><td></td><td>Power On/Off Pause</td><td>On SPRYNG™ calf wrap – user interface</td></tr><tr><td></td><td>Change the intensity</td><td>On SPRYNG™ calf wrap – user interface</td></tr><tr><td></td><td>Change the mode</td><td>On SPRYNG™ calf wrap – user interface</td></tr><tr><td></td><td>Wash at 30° C</td><td>Washable inner layer</td></tr><tr><td></td><td>Do not tumble dry</td><td>Washable inner layer</td></tr><tr><td></td><td>Do not iron</td><td>Washable inner layer</td></tr></table></div> <div>3</div>	Symbol	Description	Location		Read instructions before use	On SPRYNG™ calf wrap		Fragile, handle with care	On SPRYNG™ calf wrap		Keep dry	On SPRYNG™ calf wrap		Do not wash	On SPRYNG™ recovery wrap		Separate collection for waste electrical and electronic equipment	On SPRYNG™ calf wrap		Do not dry-clean	On SPRYNG™ calf wrap and removable inner layer		Serial number of the console	On SPRYNG™ calf wrap		Power On/Off Pause	On SPRYNG™ calf wrap – user interface		Change the intensity	On SPRYNG™ calf wrap – user interface		Change the mode	On SPRYNG™ calf wrap – user interface		Wash at 30° C	Washable inner layer		Do not tumble dry	Washable inner layer		Do not iron	Washable inner layer	<div></div> <div><p>Charging Port – This is used to recharge your Spryng device using a USB charging cable</p><p>Velcro – These are used to fasten your Spryng devices once wrapped around your legs</p><p>IR Port – These are for the Infra-Red emitter and receiver</p><p>Button – These are used to interact with the device, power on/off, to change modes and increase/decrease level of intensity</p><p>Hard Module Section – This contains the electronics of the device. It is meant to sit on your shin when using Spryng</p><p>Inflating Section – This contains an air bag which provides compression. This is meant to wrap around your leg when using</p><p>SpryngLED – Will provide light indication when using Spryng</p></div> <div>4</div>
Symbol	Description	Location																																													
	Read instructions before use	On SPRYNG™ calf wrap																																													
	Fragile, handle with care	On SPRYNG™ calf wrap																																													
	Keep dry	On SPRYNG™ calf wrap																																													
	Do not wash	On SPRYNG™ recovery wrap																																													
	Separate collection for waste electrical and electronic equipment	On SPRYNG™ calf wrap																																													
	Do not dry-clean	On SPRYNG™ calf wrap and removable inner layer																																													
	Serial number of the console	On SPRYNG™ calf wrap																																													
	Power On/Off Pause	On SPRYNG™ calf wrap – user interface																																													
	Change the intensity	On SPRYNG™ calf wrap – user interface																																													
	Change the mode	On SPRYNG™ calf wrap – user interface																																													
	Wash at 30° C	Washable inner layer																																													
	Do not tumble dry	Washable inner layer																																													
	Do not iron	Washable inner layer																																													