


planck[®] O1

Quick Start Guide

v1.0 2020.7

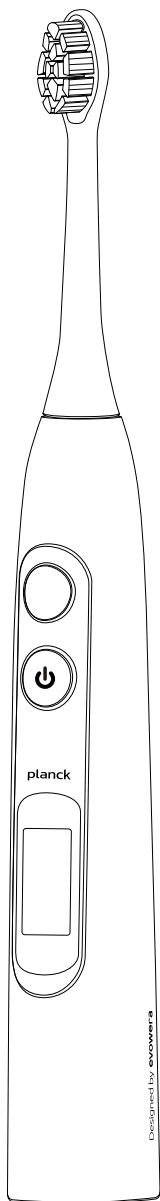

Disclaimer and Warning

Congratulations on purchasing your new **evowera**[®]**planck**[®]O1 adaptive control sonic electric toothbrush (hereinafter referred to as 'product'). Read this entire document carefully to ensure proper configuration before use and it is recommended to keep this document for future reference. This product is a personal care device. By using this product, you hereby agree to this disclaimer and signify that you have read it fully. Please use this product in strict accordance with the document, user manual, related regulations, terms of use and be sure to pay attention to the Important.

Shenzhen EVOWERA Technology Co., LTD assume no liability for damage(s) or injuries incurred directly or indirectly from improper use of this product. Misuse includes, but is not limited to, using unapproved chargers, or modifying the product.

evowera and **planck** are trademarks of Shenzhen EVOWERA Technology Co., LTD (abbreviated as "evowera"). Names of products, brands, etc., appearing in this document are trademarks or registered trademarks of their respective owner companies. This product and document are copyrighted by evowera with all rights reserved. No part of this product or document shall be reproduced in any form without the prior written consent or authorization of evowera.

This disclaimer is produced in various languages. In the event of divergence among different versions, the Chinese version shall prevail when the product in question is purchased in China, and English version shall prevail when the product in question is purchased in any other countries or regions.

Important

Danger


- Keep the power adapter away from water. DO NOT place the power adapter over or near any container with water. DO NOT drop the power adapter into water or any kind of liquid. After cleaning, make sure the power adapter is completely dry before connecting it to a power outlet.

Warning

- Always use officially approved power adapter to charge to avoid electric hazards.

- DO NOT use the power adapter outdoors or near heated surfaces.
- If any part of this product gets damage in any form (e.g:brush head, handle, charging base or power adapter), stop using it immediately.
- DO NOT let children use, clean or maintain this product without supervision.
- DO NOT let children play with this product.
- DO NOT clean the brush head, the handle, the charging base or the power adapter in the dishwasher.
- Consult your dentist before using this product if you have had oral or gum surgery in the previous two months.
- Consult your dentist if excessive bleeding occurs after using this product or if bleeding continues to occur after one week of use or if you feel discomfort or pain when using planck O1.
- planck O1 complies with the safety standards for electromagnetic devices. If you have a pacemaker or other implanted device, contact your physician or the device manufacturer of the implanted device before use.
- If you have medical concerns, consult your doctor prior to use.
- This product is only designed for cleaning teeth, gums, and tongue. DO NOT use it for any other purpose. If you experience any discomfort or pain, stop using this product and contact your doctor.
- planck O1 is a personal care device, and is not intended to used by dental clinics or institutions for multiple patients.
- Stop using a brush head with crushed or bent bristles or when planck O1 reminds you the brush head has exceeded its lifespan. DO NOT use other brush heads than the ones recommended by the manufacturer.

Overview

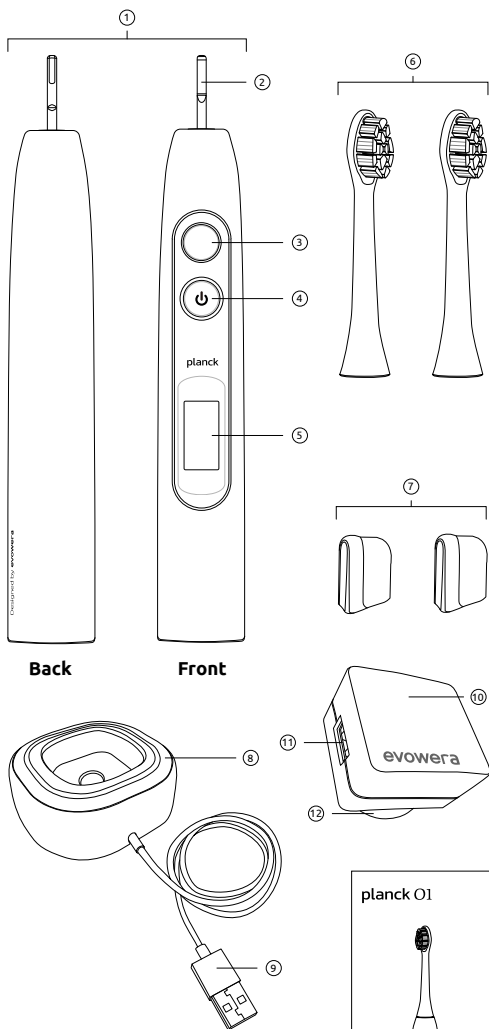


- | | |
|-------------------------------|---------------------|
| ① Handle | |
| ② Metal Shaft | |
| ③ Function Button(Mode/Force) | |
| ④ On/Off Button | |
| ⑤ LCD Screen | |
| ⑥ Brush Head | |
| ⑦ Hygienic Brush Head Cap | |
| ⑧ Charging Dock | ⑪ USB Port |
| ⑨ USB Connector | ⑫ Power Connector |
| ⑩ Power Adapter | ⑬ Quick Start Guide |

Charging


Charge the toothbrush before first use. While charging, please note:

①

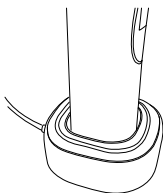

Gently place the toothbrush on the charging dock to avoid the toothbrush from being damaged.

②

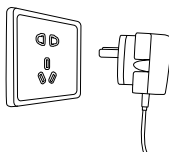

Make sure to connect the charging base to a AC power outlet of a voltage between 100 to 240 V, frequency between 50-60Hz.

Downloading and Installing evowera App

Scan the QR code below or go to app store to download evowera app. Please follow the instruction, register, and log in evowera account.



QR Code

Activating

To ensure optimal usage experience and to enjoy better after-sales service, please follow the steps below to activate the toothbrush using evowera app before first use.


Wake up the toothbrush

- Press **function button once or press on/off button once or slightly shake** the toothbrush to wake it up.


Connection and Activation

- Launch evowera app and log in. Enter device connection page, and then place your toothbrush near your phone. Follow the prompts to complete activation.

Note: Make sure your mobile device is connecting to the internet during activation. Retry if activation fails. Contact evowera if activation fails for several times.

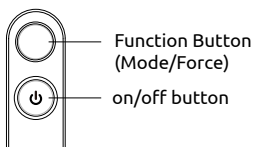

Button configuration

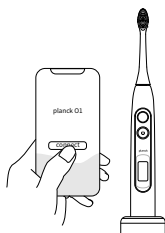

Device connection configuration

Basic Operations

Brushing Teeth

Press on/off button to switch on the toothbrush after it has been woken up.

Selecting Brushing Mode

After waking up your toothbrush, when the LCD screen is displaying home page and the toothbrush has not started working:

- ① Press **function button once** to enter mode setting page. Follow the prompts and **short press function button** to select the mode.
- ② **Press and hold function button** to confirm mode and exit to home page. **Press on/off button once** to start teeth brushing.


Adjusting Brushing Force

Note: You can only adjust the brushing force while using the toothbrush.

- Level by level: the brushing force will increase to a higher level every time you **press function button**.

Note that when the brushing force reaches the highest safe level, it will go down to the lowest if you press function button again.

- Skip through levels: **Press and hold function button**, and the brushing force will increase continuously **until you release function button**. Note that when the brushing force reaches the highest safe level, it will go down to the lowest and then increase again.

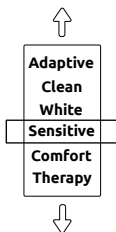

Checking Tooth-brushing Data

When the brush head is vibrating:

- ① Press on/off button once to stop the brush head vibration. The LCD screen will enter "Coverage" page.
- ② Press function button once to enter "Rating" page.
- ③ Press function button again to return to home page.

When the toothbrush is waking up, but not vibrating:


- ① **Press function button twice** to check record about your previous tooth brushing.
- ② **Press function button again** to return to home page.

Recording Your Periods

When the LCD screen is displaying home page and the brush head has not started vibrating yet:

- ① **Press and hold function button for 3 seconds** to enter "Female Period Tracker" page.
- ② Following the prompts, **press and hold function button** to record your period. The notice of "record with success" will appear on the LCD screen.
- ③ Check the LCD screen to see if you successfully record your period. If not, follow the instructions to retry.



*The LCD screen will return to home page if the toothbrush is idle for a while.

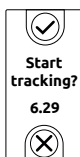



Using the Toothbrush

Follow the steps below to use the toothbrush:



- ① Attach the brush head
- ② Wet the bristle

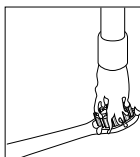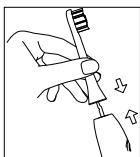

- ③ Apply a small amount of toothpaste
- ④ Place the toothbrush bristles against the teeth at a slight angle (around 45 degrees)

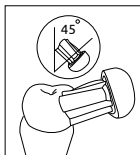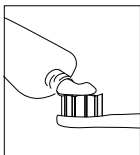

- ⑤ Switch on the vibration
- ⑥ Gently move the brush head slowly across the teeth.


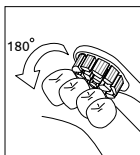
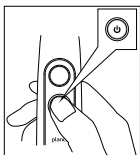


planck O1 can be safely used with:

- orthodontic appliances (it will accelerate bristle wear)
- restored tooth (dental filling, teeth crown, porcelain teeth)

Cleaning

- ① Remove the brush head , rinse the mental shaft area. Be sure to remove any remaining toothpaste.
- ② Use wet wipe to clean the handle.
- ③ After each use, rinse the brush head and bristle.
- ④ Remove brush head from handle, rinse the joint area at least once a week.
- ⑤ Use wet wipe to clean wireless charging dock and adapter. Be sure to remove any remaining toothpaste or other residues.

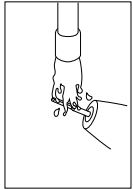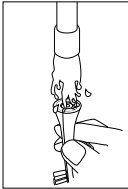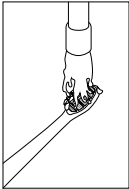

Storage

If you are not going to use the toothbrush for an extended period, unplug the charging base from the wall socket. Clean the toothbrush thoroughly and store it in a cool and dry place away from direct sunlight.

Specifications

Product Name: planck O1 Adaptive Sonic Electric Toothbrush

Product Model: planck O1

IP Rating: IPX7

Rated Voltage: 7.4V

Rated Power: 4.0 W

Operating Frequency: 2.402– 2.48GHz

EIRP: ≤2dBm

SAR: 000

Battery Capacity: more than 1040 mAh

Charger Input: 100 ~ 240V , 50 ~ 60 Hz , ≤ 350mA

Manufacture: SHEN ZHEN EVOWERA TECHNOLOGY CO.,LTD

Warranty and Support

If you encounter any problem in using the toothbrush while having no time to go to local evowera authorized dealer for support, feel free to request mail-in repair service provided by evowera.

How to apply for mail-in repair service

Via evowera app:

- ① Launch evowera app and tap to select "Mail-in Repair Service". Fill in and submit the form.
- ② If you have already submitted your request for mail-in repair service, tap "Status" to check the repair progress.

Warranty and Support

The terms of the international warranty do not include the following:

- ① Brush Heads
- ② Damage caused by use of unauthorized replacement parts.
- ③ Damage caused by misuse, abuse, neglect, alterations or unauthorized repair.
- ④ Normal wear and tear, including but not limited to chips, scratches, abrasions, discoloration or fading. Visit evowera official website to see detailed information on warranty policies.

This content is subject to change.

Download the latest version from
<http://www.evowera.com/plancko1>

Compliance Information

FCC Compliance Notice

This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

This device complies with part 15 of the FCC Rules.

Operation is subject to the following two conditions:

- (1) This device may not cause harmful interference, and
- (2) this device must accept any interference received, including interference that may cause undesired operation.

Radiation Exposure Statement

This equipment complies with FCC radiation exposure limits set forth for an uncontrolled environment.

This transmitter must not be co - located or operating in conjunction with any other antenna or transmitter.

Any changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

EU Compliance Statement

Shenzhen EVOWERA Technology Co.,Ltd. hereby declares that this device is in compliance with the essential requirements and other relevant provisions of the RED Directive.

A copy of the EU Declaration of Conformity is available online at www.evowera.com.

RoHS

