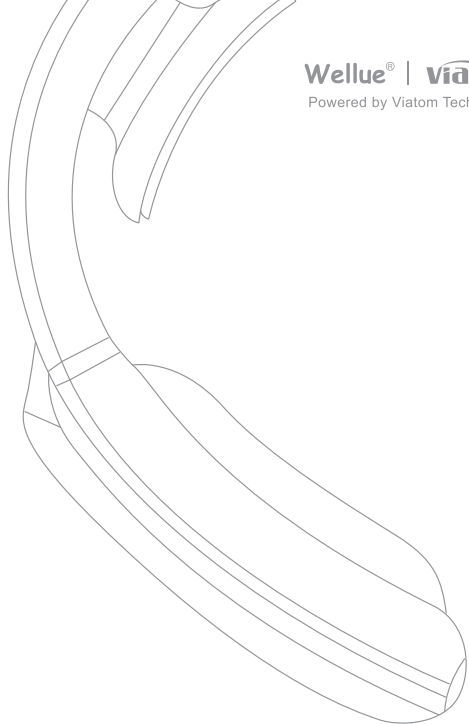
Intelligent Neck Massager

CONTENTS

English 01-14

Deutsch..... 15-29

Italiano 30-44


Français..... 45-59

Español 60-74


How to connect the Vihealth app

- 1 Scan the QR code or Search "ViHealth" in App Store and Google Play.

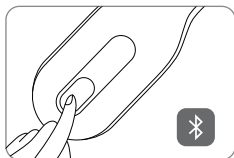- 2 Long press the power on button of the massager for 2 seconds and release it after hearing the "beep" and Indicator light is blue and always on.

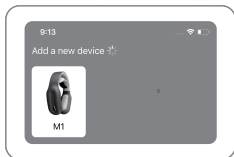
- 3 Open ViHealth App, it will search your device automatically. Choose and add your device to ViHealth.

- 4 Select the mode and function you need, adjust the intensity to start the massage.

Product diagram

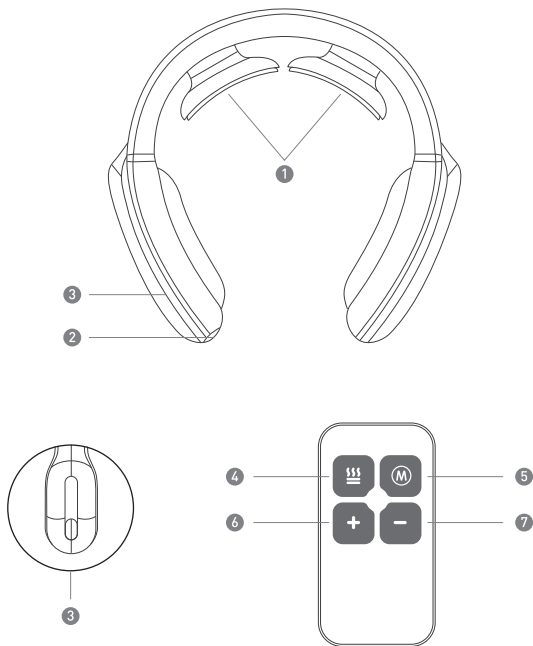❶ Electrode

❷ Charging interface

❸ Switch

❹ Heating button

❺ Mode button

❻ Intensity
increase button

❼ Intensity
decrease button

Operation method

1. Before use, first wet the skin of the cervical spine with a wet towel, and then put the product around the neck, make the electrode close to the skin of the cervical spine, adjust the position of the product, host and adjust the appropriate massage intensity and mode.then put the product around the neck, make the electrode close to the skin of the cervical spine, adjust the position of the product, host and adjust the appropriate massage intensity and mode.
2. Long press the host “” button for two seconds, turn on the host (with the buzzer “ beep ” sound and prompt light), enter the default mode with the default level (0). Long press host “”for two seconds when the machine is on, hear two"beep" sounds and turn off the products.
3. Mode selection (vitality mode、 beat mode、 dynamic mode、 relieve mode and auto mode.)“”on the remote control is the mode button. press the mode button to switch the mode.
There are five modes. After each press, it will send out “ beep ” sound and switch to higher mode. The five modes can be switched in cycle(Note: M1 mode, M2 mode, M3 mode M4 mode M5 mode). Each time a moda is switched, the intensity returns to the default o level.
4. Gear selection
 - 4.1 “” on the remote control is the intensity increase button.

After each press, it will send out “beep” sound and switch to higher mode. The intensity is divided into 15 levels (0, 1, 2...15). When the intensity is at level 15, When pressing “+”, it will send out two “beep” sound with same intensity at level 15.

- 4.2 “-” on the remote control is the intensity decrease button. After each press, it will send out “beep” sound and switch to lower mode. The intensity is divided into 15 levels (0, 1, 2 ... 15). When the intensity is at level 0, it will send out two “beep” sound with same intensity at level 15 after pressing “-”.

5. Heater selection

- 5.1 “” on the remote control is the heating button. Press the heating button and start host heating function with a “beep” sound. The two electrodes will be not at the temperature of 42°C.

- 5.2 “” on the remote control is to heating off button. Press the heating off button and stop the heating function with two “beep” sounds.

6. Timing function

- 6.1. The host is scheduled for 15 minutes. When the host is turned on, it will automatically shut down after 15 minutes.
- 6.2. When the host is turned on, the product will automatically shut down in 3 minutes when the products is not put on normally.

7. Low voltage alarm & charging function

- 7.1 Low voltage alarm: when the host is short of power (battery voltage is 3.5v), the blue light of the indicator flashes quickly(frequency: 3 times per second), you can insert one end of micro-usb cable into the host and the other end into the mobile phone Charger (output is 5V 1A/2A), green light breathing flash when charging, green light always on when full.
- 7.2 Battery model of remote control: CR2032. When replacing, please note that the battery positive pole is facing up.

Tips

- The product will be automatically turned off after 15 minutes of operation. It is recommended to use it for no more than 30 minutes per day.
- The product should be stored in a dry and ventilated plane after use.
- Do not place other weights on the product.
- Don't pull on the rubber plug at the USB charging port, as it may cause rupture and falling off.

Precautions

❗ Danger

- Please never use the product with the following medical electronic devices:
 1. Heart pacemakers and other implantable medical electronic devices;
 2. Artificial heart and lungs and other life-supporting medical electronic devices;
 3. Electrocardiographs and other medical electronic devices;
- * These medical electronic devices may cause malfunctions and can be life-threatening.

❗ Warning

- If the patient applies high-frequency surgical device and the product on the body at the same time, the electrode of the product may cause burns and damage to the product.
- If the product is used near (1M)a short-wave or microwave therapy device, the output of the product may be unstable.
- The electrode should prevent contact with wounds and scars during treatment.
- When you use the product, don't place the electrode near the chest heart, otherwise the risk of heart fibrillation will increase.
- The product is recommended for use by one user. When it's used by different users, clean the electrode on the product before using it.

- When you use the product, stick the electrode tightly and evenly to the skin.
- The electrode of this product has a heat-generating function. When the hot compress function is applied on the same part for a long time (more than 15 minutes), it may cause low temperature burns.
- If you feel unwell during or after use, stop using it and seek medical attention.
- In the following cases, consult your doctor before use, otherwise it may cause accidents or cause physical discomfort.
 1. Patients under medical treatment and feel discomfort.
 2. Patients with abnormal heart and brain nerves
 3. Patients whose body temperature exceeds 38 degrees (fever period).
 4. Patients with infectious diseases.
 5. Patients with skin perceptual disorders or skin abnormalities.
 6. Patients under medication.
 7. Patients with abnormal blood pressure.
 8. Patients who are not free to express their wishes.
 9. Patients with sensory impairment due to high peripheral circulatory disorders such as diabetes.
 10. Patients with skin allergies to stainless steel electrodes.
- Usage under the following conditions is prohibited:
 1. Malignant tumor patients
 2. Patients with bleeding tendency
 3. Pregnant women, lying-in women

4. Acute suppurative inflammation patients
 5. Heart pacemaker users
 6. Patients with metal implant
 7. Near the heart
- Don't apply it near the heart, above the neck, head, mouth, genitals, or skin disorders and don't treat the soles of the feet simultaneously. In addition, don't compress the internal organs when using. Otherwise, an accident or physical discomfort may occur.
 - Don't use with other medical devices or more than 2 products. Otherwise, nausea or physical discomfort may occur.
 - Don't use for purposes other than the ones described in the instruction manual. Otherwise, accident or malfunction may occur.
 - Don't use the electrode with a chipped or damaged surface. Otherwise, accidents or malfunctions may occur.
 - Don't disassemble, repair or modify this product without authorization. Otherwise, malfunction or electric shock may occur.

i Note

- If you feel unwell due to an abnormality in this product, stop using it immediately.
- If your body or skin feel strange, consult your doctor and follow the doctor's advice.
- If you want to move the pad or place it on another person

durine massage, be sure to turn off the power before moving it. Otherwise you may receive a strong electric shock.

- Don't treat while having electronic devices. Otherwise, the electronic device may malfunction(for example, the electronic watch disorder).
- If abnormal operation occurs, immediately switch off the power supply and remove it. Otherwise, high temperatures or short circuits may occur and lead to an accident. At the same time contact customer . service for after-sales advice.
- Don't plug or unplug the USB charging cable with wet hands. Doing so may result in electric shock or accident.
- Don't use it for children or people who have no expressiveness, otherwise it may cause an accident or cause discomfort.
- Don't use it in a ;oace with high humidity, such as a bathroom, or use it while bathing.Otherwise,you may receive a strong electric shock and damaged the product.
- Don't use it when you are sleeping. Otherwise, it may cause the main body to malfunction. In addition, the pad may move to other parts and cause discomfort.
- Don't use while driving a car, otherwise it may cause physical discomfort due to excessive stimulation, resulting in an accident.
- Don't stick the electrode to the affected area for a long time otherwise it may cause skin inflammation.
- When massaging, don't use the metal parts of belts, necklaces, etc to touch the electrode, otherwise you may receive a strong

electric shock.

- Don't apply the product on the same place for long time(more than 30 minutes), otherwise it may cause muscle fatigue at the massage site to cause physical discomfort.
- Please wear this product correctly and use it on the cervical spine. Wear it symmetrically. Otherwise, accidents or physical discomfort may occur.
- When unplugging the USB charging cable from the body do not pull on the cable, but rather, pull the cable plug and pull it out. otherwise, it may result in electric shock or damage to the product.
- The product electrode has a surface heating function for users who are not sensitive to heat or subject to heat allergies.
- Please consult a doctor before use.
- When the product is working for a long time, the surface will be hot, Which is normal.
- When the remote control battery is exhausted, replace it with a new battery of the same model as soon as possible. Battery model: CR2032
- When replacing the remote control battery, do not install the pole reversibly.
- After the remote control battery is exhausted dispose of it according to the relevant regulations in your place of residence

① Suggestion

- Please make sure that the host button is working properly.
- When you use it for the first time, or if the battery loses power after the product has been left idle for a long time, please perform the charging operation according to the instruction manual to ensure that it can work normally.
- When disposing of this product, please dispose of it in accordance with the relevant environmental protection regulations of the city, otherwise it may cause environmental pollution.
- It is recommended for people aged 16 and above.
- Don't use or store this product in the following locations:
 - Places exposed to direct sunlight
 - Places subject to high temperature and humidity
 - Places with water
 - Dusty places
 - Places near the fireworks
 - Places susceptible to vibration shocks
 - Strong electromagnetic places

(i) Routine maintenance

- When cleaning and maintaining this product, please turn off the host, use a wet paper towel to clean the electrode. When it is particularly dirty, use a little medical alcohol (alcohol concentration 75%) to wipe and disinfect.
- Please wipe with neutral detergent and soft cloth during daily maintenance, do not use corrosive detergent and don't place it into water for soaking or washing.
- During the cleaning process, water or other liquids should be prevented from entering the body to avoid short circuit and electric shock during use.

(i) Troubleshooting

- If you have problems using the product, please refer to the following tips and solutions.
- If the product is still in a fault condition or cannot be used, turn off the power.
- If the product is still in a fault condition or cannot be used, turn off the power.
- Users should not disassemble or repair the product by themselves. Dismantling or repairing by themselves will cause damage to the product and cause the invalidity of warranty. The company will not be responsible for any consequences caused thereby.

Problems	Possible cause/treatment
Failure to start	• The internal rechargeable battery is running low • The remote control battery is dead. • Please confirm whether the host product switch has been pressed.
Failure to charge or no charge light	Please check if the USB cable and host and phone power adapter are connected.
Normal display but no stimulation	• Does the pad contact the skin directly? Adjust the position so that it is close to the skin. • Is there hair or clothing, jewelry, etc. between the pads and the skin? Move away or remove items. • Is the intensity level at 0? Increase the massage intensity.
Tingling sensation	• Is the massager's electrode and neck skin closely contacted? adjust the position so that it's stick to the skin. • Is the massage time too long? It is recommended to massage for 15 minutes at a time, twice a day. • Is the massage intensity too high? Adjust to low level . • Is the skin allergic to the electrode? Those with severe allergies need to be desensitized first. • Is the climate cold and the skin dry? Use a wet towel to moisten the skin and use it.
Weak simulation	• Is the massager's electrode and neck skin closely contacted? adjust the position so that it's stick to the skin and adjust intensity. • Is the electrode dirty or oily? Clean the electrode surface. • Is the battery running low? Please charge if the battery is low. • Is massage intensity level too low? Increase massage intensity.

Product parameters

Weight	160g
Battery cell	Lithium cell
Battery capacity	900mAh
Charging voltage	DC 5V
Rated power	5W
Executive standard	IEC/EN60335-2-32
Manufacturer	Lepu Zhixin Medical Technology (TianJin) co., Ltd.

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FCC Warning

This device complies with part 15 of the FCC rules. Operation is subject to the following two

conditions: (1) this device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

Changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

NOTE: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

Radiation Exposure Statement

The device has been evaluated to meet general RF exposure requirement. The device can be used in portable exposure condition without restriction.