

DR.BEI 聲波電動牙刷 S7
使用說明書

DR.BEI Sonic Electric Toothbrush S7
User Manual

DR.BEI 音波式電動齒ブラシ S7
取扱説明書

DR.BEI Звуковая электрическая зубная
щетка S7
Руководство пользователя



使用產品前請仔細閱讀本說明書，並妥善保存

Read the user manual carefully before using and keep it with care for future reference.

ご使用の前に、この取扱説明書をよくお読みの上、正しくご使用ください。後も大切に保管してください。

Внимательно ознакомьтесь с руководством пользователя, прежде чем использовать устройство, и сохраните его для справки.

感謝使用DR.BEI口腔護理產品

請在開始使用前仔細閱讀本手冊，掌握正確使用本產品方法，以充分發揮其出色性能，獲得最佳口腔護理效果。

Thank you for using DR.BEI oral care product

Please read this instruction carefully before using it to master the correct method of using the product, so as to maximize its performance and help you enjoy the best oral care.

この度、DR.BEIの口内ケア用品をご利用いただき、誠にありがとうございます。

製品の素晴らしい性能を発揮させて、最良の口内ケア効果を得るために、ご使用前に、本取扱説明書をよくご覧になって本製品の正しい使用方法を了解してください。

Благодарим за использование продукта по уходу за полостью рта DR.BEI

Пожалуйста, внимательно прочитайте данное руководство перед использованием, чтобы управлять правильным методом использования данного продукта, полностью выявить выдающиеся свойства, получить наилучший эффект ухода за полостью рта.

Contents

I. IMPORTANT SAFEGUARDS	1
II. Instructions and recommendations for the handling of a used battery ..	2
III. Introduction	3
IV. How to Use	5
V. Overpressure Induction and Intelligent Frequency Conversion	8
VI. Stepless Power Adjustment	9
VII. Charging	11
VIII. Specifications	12
IX. After-sales Service	13
X. Troubleshooting	14

1. IMPORTANT SAFEGUARDS

Statement

First time use may accompany itch or slight bleeding due to the fact that your gum just starts getting used to the high frequency cleaning of electric toothbrush. Kindly suggest begin with Soft Mode or Sensitive Brush Head. Generally speaking, this phenomenon will disappear after 1–2 weeks. If you do not feel any significant improvement after more than 2 weeks, please stop using this product and promptly consult your dentist. This product is not suitable for users that under the age of 6.

* Tips

1. Do not share brush head with others.
2. Do not wrap the charging cable around the charging base while charging.
3. Do not put the travel case/bag in water or any liquid.

* Notes

1. This product is designed for oral cleansing. Do not use it for any other purpose. If you experience any discomfort or pain, please stop using the product and consult your dentist.
2. Make sure to check this product for any apparent damage before using it. If damaged or dropped, the product should not be used.
3. Before connecting to the power supply, please check whether the local power supply voltage meets the specification of the product.
4. Do not use the product if any damage is found on the charging base.
5. Replace brush head every 3 months.


* Warnings

1. Consult your dental professional before use if you once had orthodontic treatment, dental prosthesis or implantation.
2. This product complies with the safety standards for electromagnetic devices. If you have a pacemaker or other implanted device, contact your physician or the device manufacturer before use.
3. This product contains a built-in rechargeable battery. Please keep the product away from flames, and avoid charging, using, or placing it in high temperature.
4. Don't put the charging base in water or other liquid.
5. Make sure the power plug and interface are completely dry before connecting to the power supply.
6. Do not bump or drop the product. If it is damaged due a collision or being dropped, please stop using the product.
7. Do not modify, disassemble or repair the product by yourself.
8. Do not use the product in shower or bath.
9. Do not place metal objects on the charging base.

II. Instructions and Suggestions Regarding Battery Waste Disposal

1. Please take out the battery before disuse and send the battery to the battery collecting center designated by the government.
2. When taking out the battery, the product must be powered off.

Method of removing the battery from a disused product:

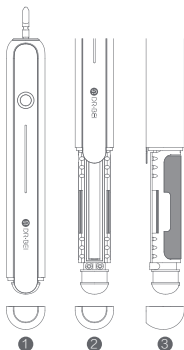


- ① Remove the bottom cover;
- ② Remove the internal bracket;
- ③ Cut off the battery nickel bar behind the support and remove the battery;

It is not recommended that the battery is replaced by any user. Forcing the battery out will damage the product's functionality. Please exercise with caution.

*The illustrations of the product, accessories, user interface, etc., in the manual are schematic diagrams for reference only. Due to the updating and upgrading of products, there may be slight differences between the actual product and the schematic diagram. The actual product may be a more up-to-date version.


III. Introduction

Thank you for choosing DR.BEI to assist on your oral health.

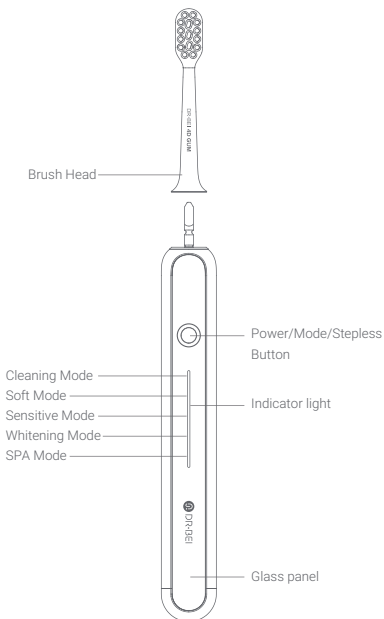

Product List

Thank you for choosing DR.BEI to assist on your oral health. Please see the list below to check if you have all the accessories.

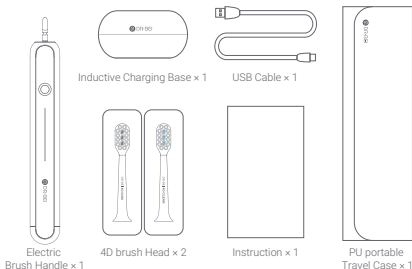

Bush Heads



DE.BEI provides you two brush heads for different oral conditions.

Soft 4D Sensitive Brush Head for people with oral sensitivity .

Standard 4D Cleaning Brush Head is suitable for people with more teeth loss, such as cigarette stains and tea stains.

Soft 4D Sensitive Brush Head
DR.BEI **4D GUM**

Standard 4D Cleaning Brush Head
DR.BEI **4D CLEAN**

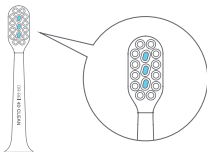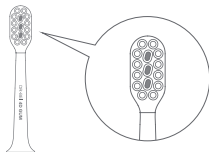

IV. How to Use

Power On

Press down the power button to switch on the product. For first use, the product will work in the cleaning mode by default.

Power Off

After the product is switched on, run a complete one-way cycle of the five modes from cleaning mode, soft mode, sensitive mode, whitening mode to SPA mode, and then turn off the toothbrush.

After the product is switched on, use any mode for more than 5 seconds, the toothbrush can be turned off by pressing the power button.


After booting, use any mode for more than 2 minutes, the toothbrush will be turned off automatically.

Switch Between Modes

After the product is switched on, it will work in cleaning mode by default.

Within 5 seconds, press the power button again to switch to the next mode.

The sequence shall follow the one-way cycle among 6 modes, namely "soft mode-sensitive mode-whitening mode-SPA mode-power off".

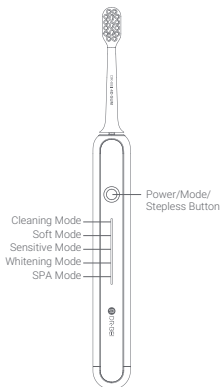

Power/Mode/Stepless Button


When certain mode is used for over 5 seconds, the product will store the mode as the memory mode. It will start in the same mode at the next boot.

Brushing time limit

Either mode has reminders for four parts of the mouth, each part has 30 seconds. There is an interval of 0.1 second after each reminder. Each time of brushing lasts for 2 minutes, with 3 reminders. After 2 minutes, the toothbrush will turn it off automatically.

Install brush head

Please hold the toothbrush bristles to face the same direction of the glass panel as shown in the picture below, and gently push the brush head down along the motor shaft until it is fixed.

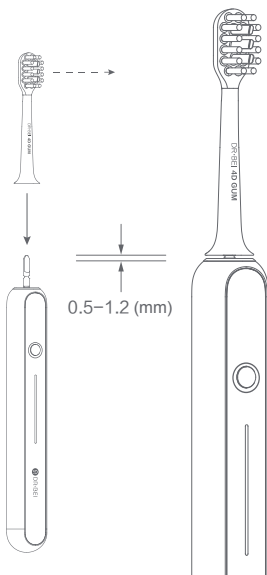

How to Brush Teeth

DR.BEI's Sonic Electric Toothbrush enables a better support for Bass method recommended by American Dental Association (ADA) and effective cleaning of plaque in your mouth.

The following suggestions can help to deliver superior cleaning results :



Bristles and teeth forms a 45°angle can better clean the plaque in gingival sulcus.



Dual frequency and low noise Brushless Magnetic Levitation Sonic Motor.



Smart reminder triggered every 30 seconds, please change to another area after a pause.



Gradual increase of vibration, no splashes exerted Improved comfort and cleanness.

Kindly advise that you press the brush against your teeth with a force less than 2N and move at every 2 teeth in your first use. With the 30s smart reminder, you will get a thorough cleaning of both inside and outside of your teeth.


Cleaning



Clean the brush head and keep it ventilated.



Whole body washable
Able to wash under running water.

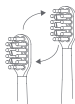


Change brush head every 2-3 months.

V. Overpressure Induction and Intelligent Frequency Conversion

The product can monitor the load force of the brush head when brushing teeth at all times, so as to sound a warning when the brushing force exceeds a certain range, and intelligently convert frequency to avoid damaging the oral cavity.

1. When the load force of the brush head is too high, the warning light at the end of the brush handle will emit red light, indicating that your brushing force may be too high at that time.
2. When warning light is on, the brushing power will be automatically reduced. As long as the force exerted on teeth by the brush head is appropriately reduced, the toothbrush can automatically resume its original state.
3. Frequency conversion here refers to changing the output power of the brush head to alleviate the possible damage to the oral cavity caused by excessive brushing force.

Illustration of the end of electric brush handle:

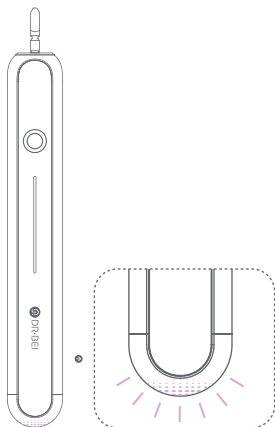

Warning Light On

VI. Stepless power adjustment

What is "stepless power adjustment"?


According to different oral care needs, the product can smartly adjust the swing scope of the brush head and provide different brushing force. How to use stepless power adjustment:

1. At power-on state, press the power button and hold it for a longer time (more than 1 second). The toothbrush enters the stepless power adjustment mode. After selecting the power that you feel appropriate, release the power button and the toothbrush will automatically exit the stepless power adjustment mode.

2. Each time the power is adjusted, the toothbrush will memorize the selected power as the default power. The next time you are adjusting power with this mode, the toothbrush will adopt this default power.


3. If you want to adjust the power again, press and hold the power button for a long time according as Step 1, and the power can be adjusted again based on the previous power stage.

Different power stage is indicated by different lengths and colors of light bar:

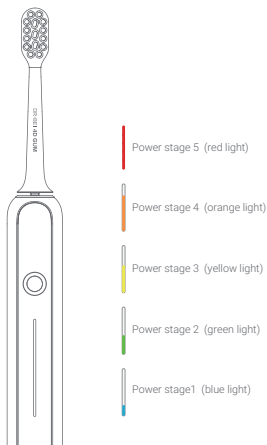

Factory Reset

If you want to factory reset power of the toothbrush, please put the electric brush handle onto the inductive charging base to begin charging after the inductive charging base is powered on. Press the power button for more than 2 seconds until all light turns white, which means that the toothbrush has already been factory reset.

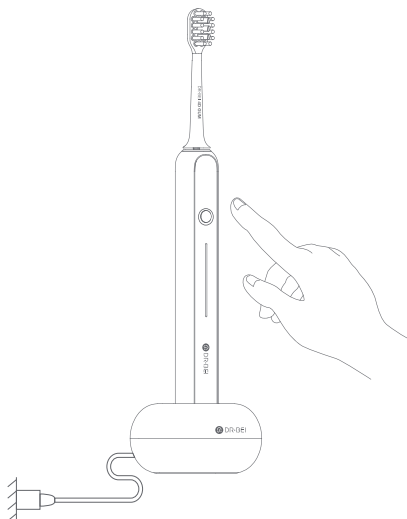

*Precautions:

1. The fifth mode, SPA mode, has no stepless power adjustment function for the time being in order to ensure the stable effect of gum massage.
2. Stepless power adjustment mode can give you more choices of brushing force. It is recommended that you adjust the power for the best experience and protect gums while cleaning teeth.

VII. Charging

Place the electric brush handle onto the inductive charging base, connect the base and the power supply through the power cord. Mode and light display of electric quantity at operation state

Blue light bar: 10%–100% electric quantity; the toothbrush can sustain 3 months' daily use.

Orange light bar: 5%–10% electric quantity; the toothbrush can sustain 7.5 days' daily use.

Red light bar: 0%–5% electric quantity; the toothbrush can sustain 3 days' daily use.

During operation, the light standing for mode is displayed as pink and white, while the colors of the other four lights represent the electricity quantity.

Light display at charging state

When charging in not completed: the light gradually changes and flows;

when charging in completed: the light is stable.


Stable green light: 100% charged

Green flowing light: 20%–100% charged

Orange flowing light: 10%–20% charged

Red flowing light: 0%–10% charged

When beginning charging, it displays flowing light; then entering energy-saving charging mode, it displays breathing light.

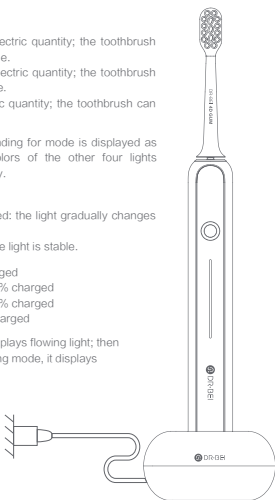

*Precautions for charging:

1. Do not immerse the base in water or any liquid;
2. Do not place the base at high temperature and humidity when charging and using.
3. Do not use the electric toothbrush when charging.

*Wake-up the charging display light: When the electric brush handle is placed on the charging base and begins charging, the light indicates current electric quantity, and flowing light will last for 30 seconds. Then, it enters the current electric quantity breathing light mode. During the charging, you can press the power button to wake up the breathing light, it will display the flowing light for 10 seconds, and then resume the breathing light.

VIII. Specifications

DR.BEI Sonic Electric Toothbrush S7		
Electric Brush Handle	Product Model	DR.BEI S7
	Rated Voltage	3.7V $\overline{\text{---}}$
	Rated Power	4.2W
	Waterproof Level	IPX7
	Charging Time	At least 12 hours from empty to fully charged.
Charging Base	Rated Input	5V $\overline{\text{---}}$ 1A
Date of Production		Printed on the package.
Place of Production		China

VI. Regulatory Compliance Information

WEEE Information
WEEE Disposal and Recycling
WEEE Disposal and Recycling Information



All products bearing this symbol are waste electrical and electronic equipment (WEEE as in directive 2012/19/EU) which should not be mixed with unsorted household waste. Instead, you should protect human health and the environment by handing over your waste equipment to a designated collection point for the recycling of waste electrical and electronic equipment, appointed by the government or local authorities. Correct disposal and recycling will help prevent potential negative consequences to the environment and human health. Please contact the installer or local authorities for more information about the location as well as terms and conditions of such collection points.



EU Declaration of Conformity

EU Declaration of Conformity

We, Wuxi QingHeXiaoBei Technology Co.,Ltd hereby, declare that this equipment is in compliance with the applicable Directives and European Norms, and amendments.

FCC WARNING

This device complies with part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) this device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

Any changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

NOTE: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

To maintain compliance with FCC's RF Exposure guidelines, This equipment should be installed and operated with minimum distance between 20cm the radiator your body: Use only the supplied antenna.

FCC ID: 2AWCX-DRBEIS7

IX. After-sales Service

If malfunction in DR.BEI Sonic Electronic Toothbrush Performance Failure List occurs, contact local authorized distributor or retailers.

DR.BEI Sonic Electronic Toothbrush S7 Performance Failure List

Name	Performance Failure
DR.BEI Sonic Electronic Toothbrush	Cannot turn on or off.
	Brush head cannot connect to power shaft when there is no foreign matter.
	Abnormal vibration of brush head.
	Power indicator does not work.
	Cannot charge the appliance.

WARRANTY EXCLUSIONS

1. Damage caused by unauthorized repair, misuse, collision, neglect, abuse, liquid intake, accidents, alterations or incorrect use of product accessories that is not recommended by the manufacturer, or tear and alter labels or anti-counterfeiting marks.
2. No longer in the validity period of warranty service (repair, refund, replacement).

3. Damage caused by force majeure.
4. Failure to comply with the malfunction listed in DR.BEI Sonic Electronic Toothbrush S7 Performance Failure List.
5. Failure of the product and its accessories in the DR.BEI Sonic Electronic Toothbrush S7 Performance Failure List but caused by users.

X. Troubleshooting

If malfunction in DR.BEI Sonic Electronic Toothbrush Performance Failure List occurs, contact local authorized distributor or retailers.

Problems	Reasons	Solutions
Brush head cannot attach to power shaft correctly.	Stuff clog between the brush head and power shaft; Brush head is not genuine or not matched; Incorrect orienting.	Remove the blocking stuff; Purchase genuine brush head recommended by manufacturer; Align the parts correctly according to the instructions.
Cannot power on.	Low power; Product in charging.	Charge the product; Remove brush from the charging base.
Weak brush head vibration.	Incorrect or loose attachment of brush head; Stuff clog between the brush head and power shaft.	Attach correctly; Remove the blocking stuff.
Handle cannot attach to charging base.	Stuff stuck at the bottom of handle;	Remove the blocking stuff.
The brush handle not showing charging status while attaching to the charging base.	The charging base failed to connect to the power supply.	Check the power supply.
Power indicator does not work.	Low power or used up.	Charging correctly will soon restore it correctly.
Low vibration.	Inappropriate stepless power adjustment.	Factory reset.

製造商：無錫青禾小貝科技有限公司
地址：無錫市濱湖區五三零大廈2號樓十一層1106室
官方網站：www.doctorbei.com
台灣代理商：瑞旭科技股份有限公司
台灣代理商地址：新北市中和區連城路268號17樓
台灣官方網站：www.doctorbei.com.tw
台灣客服專線：40555088（行動加02）

Manufacturer: Wuxi Qinghe Xiaobei Technology Co., Ltd.
Address: Room 1106, 11th floor, NO.2 of Building 530,
Binhu District, Wuxi City, Jiangsu, P.R. China
After-sales Service: service@doctorbei.com

メーカー：Wuxi Qinghe Xiaobei Technology Co., Ltd.
アドレス：Room 1106, 11th Floor, NO.2 of Building 530,
Binhu District, Wuxi City, Jiangsu, P.R. China
カスタマーサービス：service@doctorbei.com

Производители: Wuxi Qinghe Xiaobei Technology Co., Ltd.
Адрес: Room 1106, 11th Floor, NO.2 of Building 530,
Binhu District, Wuxi City, Jiangsu, Китай
почтовый ящик: service@doctorbei.com



Scan to connect and
discover more on Facebook.

 DR.BEI

了解更多DR.BEI信息請查詢官方網站

For more information about DR.BEI products,
please visit the official website

詳細については、をご覧ください

Для получения дополнительной
информации перейдите на
официальный сайт DR.BEI

www.doctorbei.com