Hyperice X 2

EN Operating Instructions
ES Instrucciones de uso
FR Instructions d'utilisation
PT-BR Instruções de uso
TH คู่มือการใช้งาน

IMPORTANT SAFETY INSTRUCTIONS

Read the entire instruction manual before using the Hyperice X 2.

⚠ WARNING! No modification of this equipment is allowed.

⚠ WARNING! If you experience severe pain, any unusual symptoms, or want to remove the device in an emergency during use:

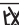
- Stop the temperature control by pressing the power button.
- Remove the device from your body.
- Promptly consult your licensed healthcare practitioner, as required.
- Use heated surfaces carefully. Always start at lowest heat setting to determine how your skin reacts to elevated temperatures.
- Do not use on an infant or invalid or on a sleeping or unconscious person. Do not use on insensitive skin or on a person with poor blood circulation. Check the skin in contact with the heated area of the device frequently during first or on new temperature level. Immediately remove from body and turn off at first sign of discomfort.
- This device must not be used by persons insensitive to heat and other very vulnerable persons who are unable to react to overheating.
- Consult your physician if you have a pre-existing medical condition that might prevent you from using this device. Please do not use the Hyperice X 2, or any heat or massage device without first obtaining approval from your doctor if any of the following apply: Pregnancy, diabetes with complications such as neuropathy or retinal damage, wear of pacemakers, recent surgery, epilepsy or migraines, herniated disks, spondylolisthesis, spondylolysis, or spondylosis, recent joint replacements or IUDs, metal pins or plates or any concerns about your physical health. Frail individuals and children should be accompanied by an adult when using any heat or massage device. These contraindications do not mean that you are not able to use this device, but we do advise you to consult a doctor first.

⚠ CAUTION!

- Do not attempt to take apart the device. The device has no user-serviceable parts. When service or repairs are required, please contact customer service at +1.949.565.4994.
- Only use the power supply provided with the device. Using a different power supply may cause the device to operate incorrectly.
- To avoid risk of electric shock, do not use the device near water, such as near a bathtub, kitchen sink, laundry tub, or swimming pool.
- To avoid damage and risk of electric shock, never spill liquid of any kind on the device.
- Do not place the device, power supply, or any accessories where they could be damaged, present a fall hazard, or become an obstruction to others.
- Keep the module vents and power inlet free of debris.
- To avoid damage and decreased performance, do not block or cover module vents.
- If the power supply is damaged, the device is dropped or damaged, liquid is spilled on the device, or the device does not operate normally when the operating instructions are followed, turn the device off by pushing the device's power button and then unplugging the device from the wall outlet. Contact customer service at +1.949.565.4994 for assistance.
- Do not puncture or otherwise damage the device. This may cause the device to operate incorrectly.
- To avoid risk of strangulation, do not leave a baby or child unattended with the power supply.
- Choking hazard, small parts. Keep away from small children.
- Do not leave the device, power supply, or any accessories where they could be damaged by children, pets, pests, or liquids. If you suspect your electronics box is damaged, contact customer service at +1.949.565.4994 for assistance.
- Do not allow lint or dust to accumulate on the modules. If lint or dust accumulates, wipe down the device with a dry cloth and clear hard-to-reach areas with compressed air before use.
- The IPX2 classification means the electronics box is protected against dripping water.
- The expected service life of the device and the integrated battery is two years.
- Do not hold the unit by the electronics box or modules.
- Please do not use the Hyperice X 2, or any contrast therapy device, without first obtaining approval from your physician if any of the following apply: are under the care of a physician or have a contraindication requiring the use of any medical device, if you are experiencing inflammation, an infection, pain of unknown origin, bleeding (internal or external) at or near the site of application, if you have a wound at or near the site of application, have sensitive skin, have bone fractures or dislocations at or near the site of application.
- Do not use the Hyperice X 2 device near eyes, nose, mouth, or ears. Keep away from hair.
- Use by unconscious or incapacitated persons may be dangerous without supervision.
- Make sure the charging port on the electronics box is always easily accessible in order to disconnect power if required.
- The device has a heated surface. Persons sensitive to heat must be careful when using the device.
- Use only on a dry clean surface of the body as shown on instructional material.
- Do not use Hyperice X 2 on wet skin or over topical ointment, especially ones that contain heat producing ingredients as skin burns could result.
- Use a thin fabric, clothing or thin towel between product and skin to reduce chance of burns.
- Do not place any metallic objects between Hyperice X2 and body, including pins or other metallic fasteners to fasten the Hyperice X 2.

The Hyperice X2 contains a lithium-ion battery which must be complied with safely at an appropriate e-waste disposal or recycling facility.

Save this manual for future reference.

Symbol	Description	Location
IP22	Degree of protection against ingress of water	On product
	Power	On product
	Cold	On product
	Heat	On product
	Pressure	On product
	Time	On product
	Start/stop	On product
	Level of protection type BF equipment	On product
	Double insulation	On charger
	Direct current	On product
	Alternating current	In manual
	Manufacturer's name and address	On product
	Date of manufacture	On product
	Separate collection for waste electrical and electronic equipment	On product
 XXXXX	Serial number of the console	On product
	Fragile, handle with care	On package
	Keep dry	On package
	This side up	On package
	Keep away from sunlight	On package
	Transportation & storage humidity limitation	On package
	Transportation & storage atmospheric pressure limitation	On package
	Transportation & storage temperature limitation	On package
	Do not wash	On product tag
	Do not dry clean	On product tag
	Do not tumble dry	On product tag
	Do not bleach	On product tag
	Do not iron	On product tag
	EU RF transmitter symbol	In manual
	Warning symbol to identify a hazard that may lead to death or serious injury	In manual and on product
	Caution symbol to indicate the need for the user to consult the instructions for use for important cautionary information such as warnings and precautions that cannot, for a variety of reasons, be presented on the medical electronics box itself	In manual
	Tip to provide guidance to make use easier. Risk to user is considered negligible.	In manual

	Refer to instruction manual	On product
	Manufactured In China	On product tag
	Indicates the manufacturer's model number so the device can be identified.	On control unit

Indications for use

Hyperice X 2 device is a portable electric hot, cold, and contrast device intended to reduce pain and swelling and improve range of motion in the treated areas. The device will monitor the cold/heat temperature and pressure to the selected level.


Risks and benefits

The risks and benefits of using the Hyperice X 2 device are the same as using traditional hot/cold therapy mediums such as an ice or heat compress. If the temperature of the treatment area feels uncomfortable, you can reduce the intensity or stop the session. Like using an ice or heat compress, the benefits include temporary pain and inflammation relief, as well as soothing joints and relaxing muscles in the treatment area.

Illustrations

Hyperice X 2 - Knee device

- 1 Electronics box
- 2 Modules
- 3 Kneecap ring
- 4 Straps
- 5 Strap extenders

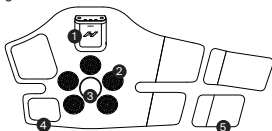Hyperice X 2 - Knee & Shoulder electronics box

- 6 Start/stop
- 7 Settings
- 8 On/off
- 9 Battery level indicator
- 10 Status light

Hyperice X 2 - Shoulder device

- 11 Loops
- 12 Modules
- 13 Electronics box

Operating Instructions

⚠ WARNING!

Before operating this device: Read all warnings at the beginning of this manual. If you do not understand the operating instructions, contact Hyperice at +1.949.565.4994

Set up the device

- To charge the device before use: Plug the power supply into an electrical outlet and then into the Hyperice X 2 device.

This device is equipped with a lithium ion battery. The battery automatically charges when the power supply is connected to the Hyperice X 2 device and an electrical outlet.

- Step 1: Place the device on body part. For the knee, align the colored ring to the kneecap.
- Step 2: Fasten straps so the device is snug against the body. For the knee, it is easier to put on the device if your leg is straight.
- Step 3: Press and hold the power button for 0.5 seconds to turn on the device. The screen and logo will illuminate.
- Step 4: Select your desired cold, hot, and pressure level.
- Step 5: Press the play/pause button to start the desired treatment.
- Step 6: Once the session begins, you can adjust the desired session time in 5 minute increments.
- Step 7: Press and hold the power button for 0.5 seconds to turn off the device.


Connecting to the Hyperice App

To connect your device to the Hyperice App via Bluetooth® make sure that your device is turned on, Bluetooth® is turned on in your phone, and your device is within close proximity. Select a routine within the Hyperice App and if prompted tap “Scan for Devices.” Select your device when it pops up on the screen. HyperSmart™ will automatically start your device. If you need to adjust the temperature or duration of your routine, simply open the Hyperice App and change the settings to your preference.


Changing Hyperice X 2 Shoulder strap from one shoulder to the other

Remove the large strap from its current loop and attach it to the other loop, leaving the long end of the strap loose.

Left shoulder strap placement

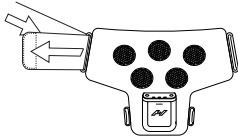Right shoulder strap placement

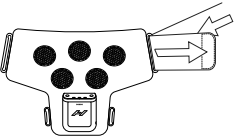Caring for the device

Cleaning the device

- Wipe down the pad with a damp, clean cloth.
- Dry thoroughly with a clean cloth.
- Do not machine wash or dry.
- Do not dry clean.

Maintaining the device

The device and power supply require no routine maintenance or service except for the care in this section.

Storing the device

Make sure the device is powered off and battery charger is not attached. Store in a clean, cool, dry place, not in direct sunlight, when not in use.

Replacement parts

Please call customer service at +1.949.565.4994 or visit our website at hyperice.com for information regarding available replacement parts and accessories.

Technical information

Do not attempt to take apart the device. The device has no user-serviceable parts. There are no user-replaceable fuses.

Bluetooth® Wireless Technology

The Bluetooth word mark and logos are owned by Bluetooth SIG, Inc., and any use of such marks by Hyperice is under license. In the unlikely event of loss of a stable Bluetooth connection, the device will attempt to re-establish its connection automatically. The Hyperice X 2 device is completely autonomous, and will continue operating normally, even during a loss of connectivity. If this device does cause interference, which can be determined by turning the device off and on, the user is encouraged to try to correct the interference by reorienting or relocating the device, increasing the separation between equipment and the device, or connecting the device to a different outlet on a circuit if it is plugged in.

The Hyperice X 2 device uses Bluetooth 5.0 wireless technology with the following radio specifications:

Frequency	2.4GHz ISM (2.402 – 2.480 GHz)
Modulations	GFSK at 1 Mbps, 2 Mbps data rates
Transmit Power	0 dBm
Security	128-bit AES ECB/CCM/AAR co-processor

FCC ID: 2AY3Y-HX20001

IC: 23655-HX20001

This device complies with part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) This device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

This device complies with Industry Canada license-exempt RSS standard(s). Operation is subject to the following two conditions: (1) this device may not cause interference, and (2) this device must accept any interference, including interference that may cause undesired operation of the device.

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Internal battery information

This device is equipped with a rechargeable lithium ion battery. The internal battery is designed to allow use of Hyperice X 2 anywhere—even when power outlets are not available. The Hyperice X 2 device may need to be plugged in before first use. The battery will provide power for approximately 90 minutes of continuous use (of cold therapy). It takes approximately 4 hours to fully charge the battery when the device is plugged in and not in use.

The rechargeable lithium ion battery is intended to be changed only by authorized service personnel with the use of a special service tool.

Product specifications

- Model: Hyperice X 2 Knee
- Model: Hyperice X 2 Shoulder
- Hyperice X 2 Knee dimensions: 25.75" x 16.5" x 1.5" [65.4 x 41.9 x 3.81 cm]
- Hyperice X 2 Knee weight: 2.5 lbs [1.13 kg]
- Hyperice X 2 Shoulder dimensions: 16" x 13" x 1.5" [40.64 x 33.02 x 3.81 cm]
- Hyperice X 2 Shoulder weight: 2.5 lbs [1.13 kg]
- Hyperice X 2 electrical requirement: 12V--8.5A
- Temperature (operating): +32° F to +113° F [+0° C to +45° C]
- Temperature (storage): -13° F to +95° F [-25° C to +35° C]
- Relative Humidity (operating): 15% to 90%, non-condensing
- Relative Humidity (storage): 15% to 90%, non-condensing
- Atmospheric pressure (storage and transportation): 190hPa to 1060hPa
- Atmospheric pressure (operating): 700hPa to 1060hPa
- Temperature of applied part: cooling +40° to +60° F [+4° to +16° C]; heating +105° to +121° F [+41° to +49° C]

AC/DC Adapter

⚠ WARNING! Only use the AC-DC adapter model number RSS1018-1200850-T2 provided with the device.

RSS1018-1200850-T2

- Input: 100-240 VAC, 50~60Hz, 2.0 A
- Output: 12V DC 8.5A, 102W

Using a different adapter may cause the device to not operate correctly.

Electromagnetic compatibility

This device is designed to operate reliably in intended for use in the electromagnetic environments of clinics, hospitals, athlete training, or home environments and is tested to meet the required immunity levels according to IEC 60601-1-2 standard for electromagnetic immunity, ensuring reliable operation even when exposed to potential sources of electromagnetic interference such as other medical devices, nearby electrical equipment or wireless communication signals and power lines. During the immunity testing the equipment continued to provide therapy normally.

General notes

Medical electrical equipment requires special precautions regarding electromagnetic compatibility (EMC). Portable and mobile RF communications equipment can affect medical electrical equipment. Cables and accessories not specified within the instructions for use are not authorized. Using other cables and/or accessories may adversely impact safety, performance, and electromagnetic compatibility (increased emission and decreased immunity).

Equipment classification

- Protection against electric shock: Class II/internally powered equipment
- Degree of protection against electric shock: Type BF applied part (Hyperice X 2 device)
- Ingress protection: IP22
- Equipment not suitable for use in the presence of a flammable anesthetic mixture with air or with oxygen or nitrous oxide
- Continuous operation

Call Hyperice customer service at +1.949.565.4994 if further assistance is needed.

Warranty

This product is covered by a limited warranty from Hyperice. Please visit hyperice.com/warranty to review the warranty in your country.

FOR MEXICO ONLY

- This device is not intended for use by persons (including children) whose physical, sensory or mental capabilities are different or reduced, or who lack experience or knowledge, unless such persons have been given supervision or training in the operation of the device by a person responsible for their safety.
- Children should be supervised to ensure that they do not use the devices as toys.
- The device must only be powered at the safety extra-low voltage indicated on the device.
- The device must only be used with the charger unit provided with the device.
- The battery must be removed from the device before disposal. The device must be disconnected from the charger when the battery is removed. The battery is removed for safe disposal.

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