


80 mm

120 mm

SeeBrush-SBH01

Video Electric Toothbrush

User Manual



ALWAYS HERE TO HELP YOU

For any questions or help, please send an email to support@opolar.com

HOW TO USE

Brushing technique

Wet brush head and apply foamless toothpaste on the both ends of the brush head. Wet your teeth and do not put too much water in your mouth, guide the brush head to your teeth before switching on the appliance. Long press to turn the video function on, and short press to turn the brushing mode on. Then move the brush head slowly and gently across the teeth in a small back and forth motion so the longer bristles reach between your teeth. Do not scrub like a manual toothbrush. Apply light pressure to maximize the toothbrush's effectiveness and let it do the brushing for you. The Quadpacer is an interval timer that beeps and briefly pauses to remind you to brush the 4 sections of your mouth. You may also consult your dentist or dental hygienist about the right technique for you.

In the first days of using any electric toothbrush, your gums may bleed slightly in general, bleeding should stop after a few days. Should it persist after 2 weeks, please consult your dentist or dental hygienist. If you have sensitive teeth and gums, it is recommended to use the 'Sensitive' mode.

Attentions

1. The foam toothpaste doesn't work with our toothbrush, which will stop the vision of the micro camera.
2. Put the foamless toothpaste on the ends of the brush head in case it falls into the hole to affect the vision of the camera.

Brushing Mode

1.Daily Clean Mode- Standard mode for daily mouth cleaning (the mode indicator turns blue)

2.Sensitive Mode- Gentle mode for sensitive teeth and gums (the mode indicator turns pink)

3.Gum Care Mode- Gently stimulate and massage your gums to improve gum health (the mode indicator turns green)

Note: Your toothbrush automatically starts in the "Daily Clean" mode. To switch to the other modes, successively press the button, it will cycle among "Daily Clean", "Sensitive", "Gum Care" and "Off" .

APP Connection

1.Please search "SeeBrush" in your APP store.



2.Long press the on/off button until the light of the mini camera turn on, then it will automatically release the Hotspot.

3.Active your mobile phone's WiFi feature and find the Hotspot of "SEEBRUSH" to connect it.

4.Click the APP "SeeBrush" and see your tooth condition.

5.Menu on the APP:

BACK
ROTATE
FILES
VIDEO
PHOTO




Note: if the APP is not working, please reset the WIFI APP and the micro camera (turn off the Hotspot and the App and re-connect them again)



Long press the button to turn the light of the camera on

Connect the hotspot "SEEBRUSH" (See "APP connection")

Click the APP "SEEBRUSH"



Apply foamless toothpaste on both ends of the brush head

Press the button shortly to choose a brushing mode

Start tooth brushing and view your oral condition

IMPORTANT SAFEGUARDS

SeeBrush has been carefully designed to offer you and your family a unique brushing experience that is both safe and effective. Read all instructions before using.

DANGERS

To reduce the risk of electrocution:

- Always keep the product away from water when charging.
- Unplug the charging cable immediately if it falls into water.
- Never charge the product with a damaged cord or plug.

WARNINGS

To reduce the risk of burns, electrocution, fire or injury:

1. Close supervision is necessary when this product is used by, on, or near children or infants.
2. Use this product only for its intended use as described in this manual. Do not use attachments not recommended by the manufacturer.
3. Keep cord away from heated surfaces.
4. Never drop the appliance or insert any object into any opening of the appliance.
5. Never operate this product if it has a damaged cord or plug, if it is not working properly, if it has been dropped or damaged or if it has been dropped into water. If any of these occur, please contact a service center and return the product for examination and repair.
6. Never force the plug into an outlet, if the plug does not easily fit into the outlet, discontinue use.

7. This product contains no user-serviceable parts. Do not modify or repair the product. This may cause fire, electric shock or injury. Do not disassemble the product except when disposing of the battery.

8. When taking out the battery for disposal, use caution not to short the positive (+) and negative (-) terminals.

9. Use this appliance only for its intended use as described in this document or as recommended by your dental professional.

10. The toothbrush is a personal care device and is not intended for use on multiple patients in a dental office or institution.

11. Do not touch the power plug with wet hands. This can cause an electric shock.

12. This product is not intended for use by children under age three.

13. Replace brush (head) every 3 months or sooner if brush head becomes worn.

MEDICAL WARNINGS

1. Consult your dental professional before you use this product if you have had oral or gum surgery in the previous 2 months.

2. Contact your dental professional if excessive bleeding occurs after using this product or bleeding continues to occur after 1 week of use.

3. Discontinue use of this product and contact a physician/dentist if discomfort or pain is experienced.

4. This product complies with the safety standards for electromagnetic devices. If you have a pacemaker or other implanted device, contact your physician prior to use.

5. Consult your physician if you have other medical concerns.

MEDICAL WARNINGS

IMPORTANT INFORMATION

• Your video toothbrush can be used with your Smart Phone (Find details under < "APP Connection" >).


• To avoid electromagnetic interference and/or compatibility conflicts, deactivate the radio transmission before using it in restricted surroundings, such as aircrafts or specially marked areas in hospitals. Deactivate radio transmission when your toothbrush is turned off.

• Long press the on/off button to reactivate transmission again.


• Persons with pacemakers should always keep the toothbrush more than six inches (15 centimeters) from the pacemaker while turned on.

• Anyone who suspect interference is taking place, long press to turn off the radio transmission of your toothbrush.

PRODUCT DETAILS

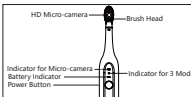

HD Micro-camera
Brush Head
Indicator for Micro-camera
Battery Indicator
Power Button
Indicator for 3 Modes



Charging Port

PACKAGE LIST

If there is anything missing, please contact support@opolar.com for help.

- 1* Video electric Toothbrush
- 3* Brush Heads
- 1* Foamless Toothpaste
- 1* Phone Holder
- 1* Toothbrush Holder
- 1* USB-C Charging Cable
- 1* User Manual

HOW TO CHARGE

Your toothbrush has a waterproof handle, which is electrically safe and designed for use in the bathroom.

1. Remove the bottom cover to plug the usb-c cable in, connect it to a power source.
2. The gradient red light indicate that the toothbrush is charging, once it is fully charged, the light stays in red.

3. A full charge can take up to 3-4 hours and enables up to 45 days of regular brushing (twice a day 2 minutes).

Note: If the rechargeable battery is low, the battery indicator light is flashing in red.

INSTALL INSTRUCTIONS

If you prefer to mount the video toothbrush on a wall, mirror or other surfaces, use the provided toothbrush holder and adhere it to a clean and dry surface. It will not work on dirt-repellent surfaces.

Please stick the phone holder on a clean and dry surface like a wall or somewhere you could easily see your tooth condition while brushing.

CLEANING AND STORE

After brushing, rinse your brush head under running water with the handle switched on. Switch off and remove the brush head. Clean handle and brush head separately, wipe them dry before assembling. Do not clean the brush head, handle, charging base or charging travel case in the dishwasher. The bottom part of toothbrush should only be wiped clean with a damp cloth.


If you are not going to use the toothbrush for an extended period of time, unplug it, clean all parts and store them in a cool, dry place away from direct sunlight.

REMOVING BATTERY

The rechargeable battery inside cannot be replaced, but it could be removed for recycling. At the end of product life and prior to disposal, please remove the rechargeable battery from the handle by following the instructions below.

Do not throw away the appliance with the normal household waste at the end of its life, but hand it in at an official collection point for recycling. By doing this, you help to preserve the environment.

This process is not reversible! To remove the rechargeable battery you need a flat head (standard) screwdriver. Observe basic safety precautions when you follow the procedure below. Be sure to protect your eyes, hands, fingers and the surface on which you work.



WARRANTY

To ensure that we meet your needs, we offer one year warranty from the purchase date for defective products. The incorrect human operation, not limited to disassembling, burning, soaking, and dropping. For any customer product question or honor your warranty, please write to support@opolar.com. Our usual response time is within 24 hours.

To become OPOLAR Member:

1. Visit the below link:
2. Click "Register" on the top right corner.
3. Type your email address and password.
4. Click Register.

The benefits you can get:

1. Extended 6 months warranty on Amazon order, which means 18 months warranty instead of 12 months.
2. Get our latest deal notifications and more via email.

WE NEED YOUR SUPPORT

We welcome you to POST AN UNBIASED REVIEW on our Amazon product page so that others can be benefited from your experience. As a growing business, we know how precious your time is, so we especially appreciate your time and input. Your support is very important to us! Thanks again!

For patents details, please visit the below website address:

www.opolar.com/page/patents

Warning:

This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) this device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

NOTE: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

NOTE: This device and its antenna(s) must not be co-located or operation in conjunction with any other antenna or transmitter

RF Exposure Statement

This device and its antenna(s) must not be co-located or operation in conjunction with any other antenna or transmitter.