

正面



EYE MASSAGER
CH-101

Instruction manual

- Thank you for purchasing our products
- Please read the instruction manual carefully before use
- This manual contains service guide, please keep it properly

Catalog

Product introduction	01
Scope of application	01
Product features	02
specification	02
Accessories list	03
Product structure	03
Preparation before use and charging	04
Method of use	05
Transportation and storage	06
Cleaning and storage	07
Matters needing attention	07
Common problems and solutions	08
After-sales service	09

Product introduction

The eye massager adopts low-frequency vibration motor, environmentally friendly PU airbag and PI heating film, graphene carbon fiber far-infrared heating wire and microcomputer chip control technology. It helps the eye muscles relax and keep away from eye discomfort by applying air pressure/vibration/hot compress to the eye. Constant temperature hot compress can effectively improve the ambient temperature around the eyes, and can drive the eye muscles to move together with the micro-motor vibration massage and the air bag simulated pressing point, bringing a full range of comfortable experience.

Product features

- Quick charging, automatically stop after working for 15 minutes
- Constant temperature hot compress, have no harm to the skin
- Folding and portable, super long endurance
- 5 modes switching, giving users a different experience
- Intelligent voice, meditation music

Scope of application


This product is suitable for people who use computers, mobile phones and other electronic devices for a long time. It can effectively avoid eye fatigue, blurred vision and other symptoms caused by eye problems such as using eyes for a long time, wearing glasses, insomnia, myopia, farsightedness, dark circles and eye bags.

Warm reminder: This product is personal health care massage product, cannot replace medicine or medical apparatus and instruments.

Specifications

Product name: Eye Massager
Model: CH-101
Rated voltage: 5V
Rated power: 3W
Product material: ABS shell + Lined with nano leather
Working time: 15 minutes by default
Battery capacity: 1200mAh
Size: 80 x 110x 220mm
Normal working environment/Ambient temperature: -20°C~40°C, Humidity: RH≤85%, Atmospheric Pressure: 70kPa~106kPa
Transportation and storage/Ambient temperature: -20°C~40°C, Humidity: RH≤85%,

Product structure

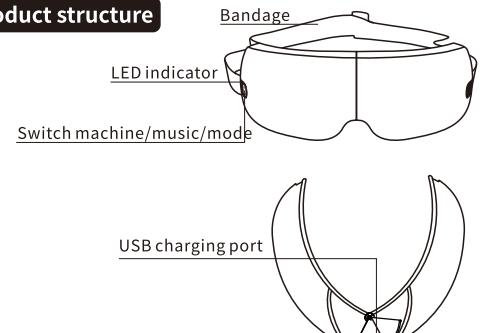

Accessories list

If you find the following items are inconsistent with this information, please contact the manufacturer.

Number	Accessory	Quantity	Remark
1	Eye massager host	1	
2	USB data cable	1	
3	Instruction	1	Includes certificate and maintenance card
4	Velvet bag	1	

Preparation before use and charging

Preparation before use

Align the eye massager with your eyes, wear the host, adjust the elastic band to the most comfortable size, and adjust the host to the most comfortable position.

Note:

1. Take off your glasses, contact lenses, hair clips and hair accessories before wearing them to avoid discomfort.
2. Do not use this product when your hands are wet.
3. Do not drop or step on the product.
4. After wearing and using the eye massager, please close your eyes for 1-3 minutes and then slowly open your eyes to avoid glare from the outside and discomfort.

Charging

1. Plug the USB data cable into the USB socket of the power adapter.
2. Plug the power adapter into the power socket.
3. The Type-C USB plug of the USB data cable is inserted into the type-C USB socket of the product.

Warm tips

1. During the charging process, the indicator light flashes: After the charging is completed, the indicator light is always on.
2. During the charging process, the product cannot be turned on.
3. If the product is charged while it is working, the product will automatically stop working.
4. The charging time is 3-5 hours. After charging is completed, please disconnect the power adapter.

背面

Method of use

Turn on the product:
Finger keep pressing the "O" button for two seconds to turn on the product. After power on, the indicator light is always on and the host starts to work. The automatic mode is defaulted, and the technique changes every 30 seconds. The working time is 15 minutes and music is played.

Warm tips:
The product has timing function, The default time is 15 minutes. After reaching the timing time, the product will stop working automatically. In addition, the timing time is the time period of each mode.

Switch mode:
Press the "O" button for less than one second to switch mode. Operate once to switch to clear mode; operate twice to switch to vitality mode; operate three times to switch to beauty mode; operate four times to switch to sleep mode; operate five times to switch back to automatic mode and cycle/loop in this way.

Mode introduction:
Automatic mode = change technique every 30 seconds + music
Clear mode = air pressure + hot compress massage + music
Vitality mode = air pressure + vibration massage + music
Beauty mode = hot compress massage + music
Sleep mode = air pressure massage + music

Music control:
Press the "O" button two consecutive times for less than one second to turn off the music.
Press the "O" button two consecutive times for less than one second again to turn on the music.

Shutdown product:
When the product is working, press the "O" button continuously for 2 seconds and the host is off.

Transportation and storage

Storage environment:
Room temperature, non-corrosive gas, well ventilated, dry indoor.
Transportation:
1. The product is equipped with simple earthquake-proof facilities, suitable for aviation, railway, highway and ship transportation. It should avoid rain, exposure, inversion and collision and other harsh conditions.
2. When the storage time exceeds three months, you should open the box for inspection and power on for half an hour.

Cleaning and storage

Cleaning:
Before cleaning, turn off the power and unplug the power adapter from the socket. Use a dry soft cloth to clean the surface of the product.
Warning:
Be sure to cut off the power when cleaning. Do not use corrosive cleaners such as benzene or diluents to clean the product, otherwise it will cause discoloration, surface aging, structural damage and other conditions.
Storage method:
If you do not use it for a long time, please cover it with a dust cover (provided by yourself), and place the product in a dry and ventilated room to avoid sun exposure.
Warning:
This product contains batteries, do not approach high temperature or put in the fire.

Matters needing attention

Please read carefully and understand all the contents of the manual, and operate according to the instructions.

1. Do not use when you have eye disease symptoms such as glaucoma, cataracts, and retinal detachment.
2. Do not use it if you have undergone eye surgery.
3. Do not use this product for stroke, hypertension and hypotension patients.
4. Do not use this product if you are slow of movement, have physical or mental dysfunction or intellectual disability.

Common problems and solutions

Problems encountered	Possible Causes	Processing method
Automatically shutdown after boot	Low battery	Please charge
Unable to boot	Low battery	Please charge
Hot compress temperature is not enough	Low battery	Please charge
Poor indicator light	Poor connection of USB socket/Poor connection of power supply	Unplug USB charging cable/power cord

If you have carried out the above treatment and there are abnormal phenomena, please turn off the power and contact the dealer you purchased.

Please do not disassemble and repair the product without authorization. We are not responsible for any product damage, personnel damage or other losses caused thereby.

After-sales service

Thank you for purchasing our electronic health products! While ensuring the high quality of this product, if you find that the product has a performance failure caused by the manufacturing process or raw materials, we will provide you with free maintenance services. Please fill out the warranty card when purchasing, and ask the purchaser for the proof of purchase as a proof of warranty.

Warranty description

1. The starting date of the warranty period shall be subject to the sales invoice of the product. Please show the invoice when repairing.
2. During the warranty period, free repairs will not be provided under the following circumstances:
 - 1) The machine fails to operate according to the requirements of the Operation Manual
 - 2) Defective machines caused by self-repair, modification and disassembly by users
 - 3) Due to natural disasters, falls, collisions, improper voltage caused by the failure of the machine
 - 4) Machines whose appearance is naturally dirty due to use

Certificate of conformity

Product name: Eye Massager
Product model: CH-101
Production date: see the body barcode
Inspectors:
★ This product is approved to leave the factory after inspection

Warranty Card

product name: _____
Product number: _____
Purchase date: _____
Customer name: _____
Customer phone: _____
Customer address: _____
Fault description: _____
Repair result: _____

尺寸:75X105mm

材质:120g双胶纸

工艺:单色双面印刷+风琴折

This device complies with part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) This device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

Any Changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

Note: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.