

# Instruction



## Reminder

The product is designed for household use only, excessive use in public is not covered by warranty

Make sure you read through instruction before use to avoid any unnecessary damage to the product or user

## **USB Wireless Back & Waist Massager X1**

1. The product is not designed for people (including children) lack of experience or knowledge, use this product under supervision
2. This product has heating function, use with care for people who are insensitive of heat
3. Sitting or stepping is not allowed, max safe load is 20kg

### **Directory**

1. Safety Instruction . . . . . 1
2. Packing List . . . . . 5
3. Operation Instruction . . . . . 6
4. Daily care / Troubleshooting / Emergency / Specificatio . 7

### **Warning**

1. Company is responsible for repairing products placed under special condition of corrosion such as outside islands, sulphur zone and seashore
2. Repair service uses OEM parts, company is not responsible for any damage if changed without using OEM parts

### **Emergency Handling**

When anything strange happens, please unplug and shut down the power

## **Safety Instruction**

### **Warning**

1. This product is not recommended for people under 16
2. The product is not designed for people (including children) lack of experience or knowledge, use this product under supervision
3. This product is not suitable for diagnosis, medical treatment or any business practices
4. Please do not use for goal other than massaging
5. Do not use with other massaging devices
6. Do not use by two or more people at the same time
7. Do not use when you are drunk or feel uncomfortable
8. Do not use for bones with protrusion, fasciitis or injured part of body
9. Check if power and function buttons work
10. Stop using immediately if you feel any discomfort when using
11. Use the product at least 1 hour after meal, exercising or showering
12. This product has heating function, please do not use for more than 10 minutes in case of inflammation or burn
13. This product has heating surface, be careful for people who are insensitive about heat
14. Take 15-20 minutes break before next program, we suggest 20-30 minutes every day
15. Please take off your jewelries before using to avoid
16. This product contains precision massaging parts, please do not step on or use upside down, company is not responsible for any damage caused above actions
17. Do not rinse with water, this will cause electronic components shortage or components failure
18. Do not spray anything like oil, gas, coating, pesticide or

any other corrosive liquid, it will cause electronic components shortage or damage of product shell

19. Please leave the product in safe cool ventilated place when not using, do not place it outdoors or under sunshine with high temperature
20. Please leave the product for 1 hour before using when moving from low temperature area to relatively high temperature. This can avoid circuit shortage caused by water due to sudden temperature change
21. Please power off the product before moving it
22. Please do not drop, collide or place heavy things on top of the product
23. Do not let kids play with the plastic packing bag to avoid suffocation
24. The product has been double checked before leaving the facility, please do not disassemble or repair the product on your own. Use OEM parts when repairing, company is not responsible for any damage caused using other unauthorized parts (Please contact the company for repair service if needed)
25. Please stop using the product when there is damage to cloth cover, and contact customer service, please do not take off or sew up the cloth cover and report to customer service immediately for exchange. (Remark: Cloth cover is consumable item, company will charge based on real time quote)
26. Please charge for 4 hours before first use. Please leave the product in ventilated area when not using, charge every two months
27. This product is personal hygiene, no exchange is allowed once the package is unsealed, except for defects

## **Safety Instruction Note**

**Please avoid using or consult a doctor when you have following**

1. Just recover from sickness, pregnancy, after delivery or surgery
2. Patients with acute disease
3. Patients with malignancy
4. Patients with malignancy
5. Patients with paralysis or disorder of perception
6. Patients with severe blood circulation problem
7. Wound at waist or back
8. People with dermatitis or skin issue
9. Long term treatment, severe wound near back or waist
10. People with osteoporosis, cataclasis, acute pain or big surgery
11. People with thrombus, aneurysm, varicosity, heart or blood disease, skin disease, wound or cellulitis
12. People with fever
13. People use pacemaker or any other implantable medical device
14. People with autonomic disorders
15. People who are above 70 years old
16. People under 16
17. People who cannot operate alone
18. For other situations and questions, please consult a doctor


## **Safety Instruction**

### **Note**

### **Electrical Safety**

1. Do not damage, twist, alter, pull or any other inappropriate action on power line. Do not put heavy things on top of power line, especially when moving around to avoid fire hazard. If there is any damage to lines, do not change on your own, please call customer service for help.
2. Stop using immediately when power line is damaged to avoid power shortage, sparks or electric shock, do not change on your own and contact customer service for help.
3. Do not unplug the power when massaging to avoid any injuries.
4. Please do not hold plug itself or pull power line when unplugging.
5. Please insert the plug in completely, take off pin, dust or any other things to avoid power shortage, electric shock or fire hazard.
6. Please do not use this product when there is thunder outside
7. Please do not touch with wet hands to avoid electric shock or getting water into the product
8. Please use proper power according to instructions, otherwise might cause damage to product, electric shock or fire hazard
9. Please unplug and power off after using or sudden power outage
10. This is rechargeable product, please do not charge when no one is around, and watch closely for charging status. Do not power on when charging, and please use proper power line and adapter for charging. Do not charge for more than 4 hours, and unplug after charging.

## Parts Name

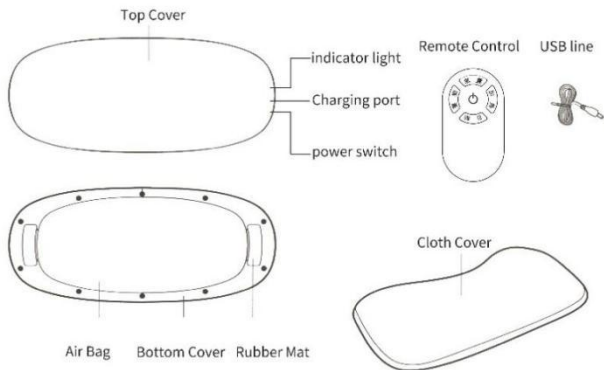

## Suitable for



**Whole Back**



**Shank**

Air Bag/Vibration/Heat/Auto

※The product is only suitable for above parts of body

# Instructions

## Power On/Off

Press this button to power on or off

## Air Bag

Press this button to have air bag filled, default is weak mode. Press again for strong mode. Press one more time to turn off. (Two Modes)

## Vibration

Press this button to turn on vibration mode, default is weak mode. Press again for strong mode. Press one more time to turn off. (Two modes)


## Heat

Under air bag mode or vibration mode, press this button to turn on heating function. Heating cannot be turned on alone for safety concerns, press again to turn off.

Heating Temperature is set to  $45^{\circ}\text{C} \pm 5^{\circ}\text{C}$

## Auto:

Press this button to turn on auto massaging, default is weak mode. Press again for strong mode. Press one more time to turn off. (Two Modes)

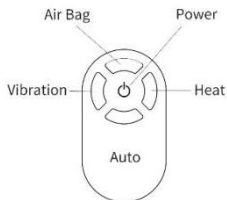

## Daily Care

- ※Unplug when not or finish using
- ※Do not pull or drag the product with power line
- ※Use soft cloth to wipe exterior
- ※Do not let water or liquid into product

## Emergency Handling

Symptom	Cause	Solution
Unable to power on	Insufficient Battery	Charge
Sudden stop	Use for too long, cause high temperature inside	Restart after on for 30 minutes
	Use for more than 10 minutes, the default setting is 10 minutes	Restart

In case anything happens, and you still cannot find solution:

1. Power off and stop using immediately
2. Contact the dealer

## Specification

Model	X1	Input	DC 5V-2A
Output	24W	Net Weight	1kg
Dimension	401mm x 173mm x 58mm	Material	Rubber, Metal, Cloth

This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) this device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation

NOTE: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

Warning: changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment