


成品尺寸：92*92mm，展开尺寸：366*92mm，材质：100g铜版纸，4折页，风琴折，双面单黑印刷_20200423




English

SMART BODY SCALE

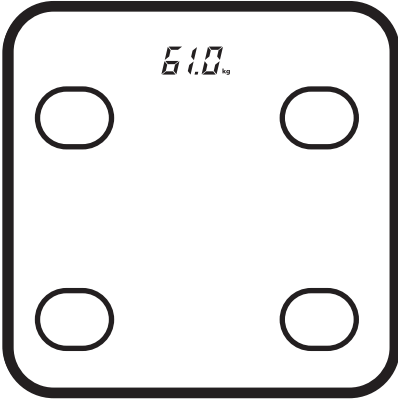


HOW TO INSERT BATTERIES

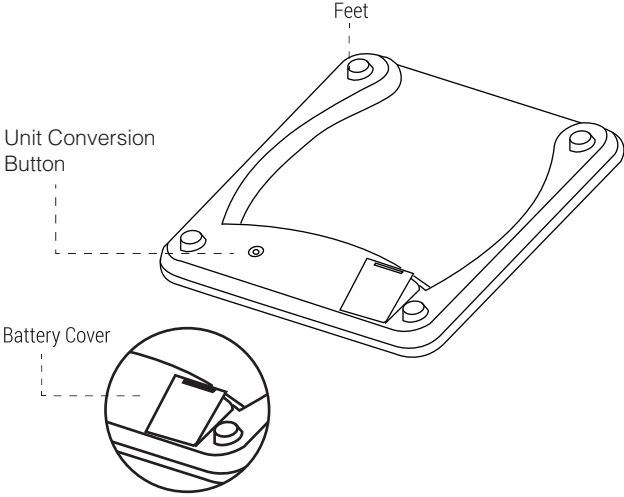

INSERTING THE BATTERIES:


- Turn the scale over and place it on a flat surface.
- Open the battery cover and insert four AAA batteries.
- Make sure the positive and negative ends are facing the correct polarity directions marked in the battery compartment.
- Close the battery cover and wait until the digits "0.0" are shown on the LED display.

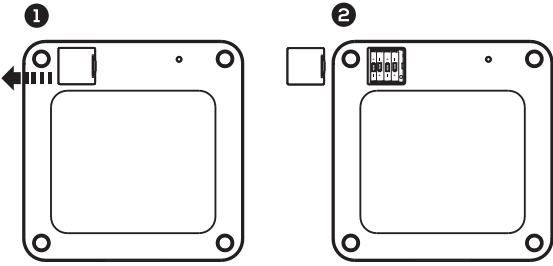

Questions?

Contact Health
855-633-5600

HOW TO WEIGH ON SCALE

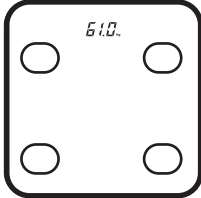

HOW TO WEIGH ON SCALE

- Make sure you stand bare ,clean and dry feet. Place one foot on the device and remove it when the display turns on.
- The device cannot measure body fat percentage if you wearing shoes and socks
- Stand up straight on the device to begin measurement .When the weight reading stops changing , it will blink three times and then stay on the display.

Note -- Please keep the scale on hard surface for accurate readings.

NOTIFICATION ON SCALE ALONG WITH THE SYMBOLS

 Low power (Need to insert new battery)	 Measuring body fat
 Measure Abnormally (Step down and try to weigh again)	 Overweight
 No user information is issued, and no user is identified (Turn on the Bluetooth of the mobile phone, use the mobile APP to establish a connection with the scale, and then measure again.)	 Initialize
 Multiple users identified in offline testing (During the offline test, the scale recognizes multiple users by weight, and cannot recognize a specific user parameter. Please connect to APP directly for online test)	

FCC Statement:

This device complies with part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) this device may not cause harmful interference , and (2) this device must accept any interference received,including interference that may cause undesired operation.

Note: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception,which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

The users manual or instruction manual for an intentional or unintentional radiator shall caution the user that changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment. In cases where the manual is provided only in a form other than paper, such as on a computer disk or over the Internet, the information required by this section may be included in the manual in that alternative form, provided the user can reasonably be expected to have the capability to access information in that form.

CARE AND MAINTENANCE

Care and Maintenance :
Please clean the surface of the scale with damp cloth.
Don't use any detergent and cleansers to clean the device.
Scale is not recommended for the following groups:

- 1) Underage minors
- 2) Pregnant women
- 3) Patient wearing pacemaker or any other internal medical devices