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
感谢您购买本产品,使用前请仔细阅读本说明书
Thank you for purchasing our product,please read
this manual carefully before use it.

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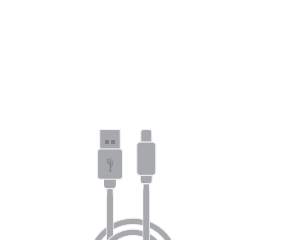
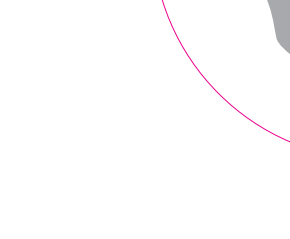
The manual operation method is as follows:
Tap the switch button to turn on the device,
and continue to press the button to adjust
different modes.

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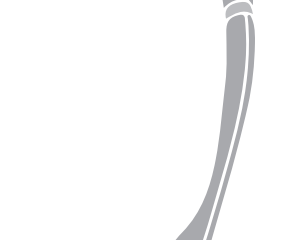
Indicator light description:
The flashing indicator light indicates that the
battery is low.
The indicator light is always on during charging,
The indicator light will go out when fully
charged.

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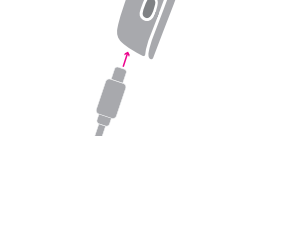The APP control method is as follows:
1. Scan the QR code of the APP on the right
2. Register and log in: Enter your phone number
and get the SMS verification code to complete
the registration
3. Add device: Make sure that the Bluetooth of
the mobile phone is turned on, turn on the
device, and click Add Device on the APP page

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Fading improvement:
Crow's feet (corner lines), nasolabial folds, dark
circles.

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technical parameter

Model: CY-1003A
Rated input power: 5w
Battery capacity: 700mAh
Charging input: 5V 1A
Charging interface: TYPE-C

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It is recommended to use it for 20-25 minutes a day, the next day, and temporarily stop using it during menstruation.
*This product should be stored in a ventilated and dry place after use
*Do not place other heavy objects on it

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Care and maintenance

*When cleaning this product, please turn off the power first, use wet paper to clean the surface of the product, do not use alcohol to clean, do not wash in water.
*Please use neutral detergent and soft cloth to wipe during maintenance.
*Prevent liquid from entering the product during cleaning

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Precautions

*If you feel unwell during or after use, it is recommended to stop using it and consult a doctor.
*This product contains lithium batteries, please do not disassemble, modify, immerse in water

or throw into fire by yourself, otherwise the product may be damaged and exploded.
*Please close your eyes; do not open your eyes during use.
*Do not use for purposes other than the functions described in the manual, otherwise it may cause discomfort or damage the product.
*Do not disassemble, repair or modify this product without authorization, otherwise it may cause malfunction.
*Do not plug or unplug the USB charging cable with wet hands, otherwise it may cause electric shock or accidents.
*Do not use for children.
*Do not use the product in a place with high humidity such as a bathroom, or use the product while taking a bath, otherwise it may cause discomfort or damage the product.
*Do not use it while sleeping, otherwise it may cause the product to malfunction.
*Please wear this product correctly and use it on your face, symmetrically.
*When the product is working for a long time, the surface will heat up, which is normal.
*Please confirm that the main unit buttons are working properly before use.
*When disposing of this product, please follow the city protection regulations for disposal, otherwise it may cause environmental pollution.
*Do not use or store this product in the following places: Places exposed to direct sunlight, places with high temperature and humidity, places with water, places with a lot of dust, places susceptible to vibration and shock, near fireworks, places with strong electromagnetic.
*Do not use this product during exercise.

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Troubleshooting

1. If the product fails during use, please shut it down immediately and contact customer service.
2. Users should not disassemble and repair the product by themselves, otherwise the company will not bear any consequences caused by this situation.

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Product and parts recycling instructions

In order to protect the environment, after this product or its components are scrapped, you, as a consumer, are responsible for separating it from the biochemical waste and sending it to a specially set recycling site, which will be classified, disassembled and recycled in accordance with relevant national regulations. Use etc.

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Cautions when disposing of the product

*The battery used in this product is a special rechargeable battery, which cannot be used in other products. Please do not charge the rechargeable battery that has been taken out.
*Do not throw in the fire

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Manufacturer: Beijing Truwin Optoelectronic Medical Co., Ltd