

PO* · 颈部美容仪

Thank you for purchasing our product. Please read this instructions carefully before using it.

①



①


① PO Light Neck Beauty Instrument



②

② TYPE-C charging cable

②

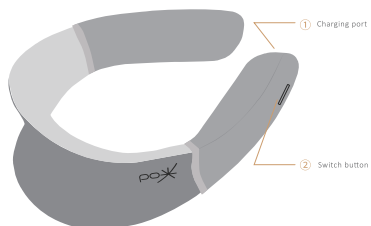


① Charging port

② Switch button

The manual operation method is as follows:

For manual operation, please tap the switch button to turn on the device, and switch to other modes by continuing to press the switch button.

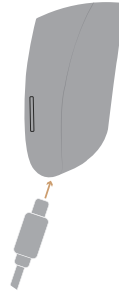

Instructions for mode selection:

1. Press the switch button to power the device.
2. Press the switch button to shift modes:
 - ① Professional mode: red light, the unit run time of the device is 20 mins;
 - ② Daily mode: red light and left-right alternating vibration, the unit run time of the device is 20 mins;
 - ③ Stress-free mode: red light and left-right simultaneous vibration, the unit run time of the device is 20 mins;(Different types of devices have slightly different modes.)
3. Hold on the switch button to turn off the device.

Indicator light description:


Indicator light - Red light flashes means low battery.

During charging - white light always on, The indicator goes off when product is fully charged.

3

Product Features:

Reduce and improve neck wrinkles + neck care with vibrative massage

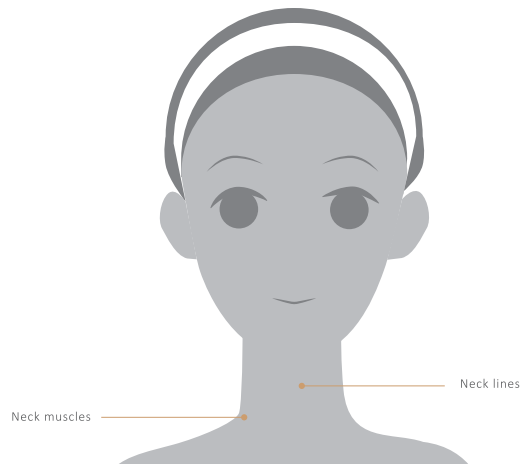

4

Wearing method

1. Suitable for neck (front or reverse wearing)
2. Turn on the device and choose a suitable mode.





**IMPORTANT SAFETY INSTRUCTIONS
BEFORE USING THIS APPLIANCE READ
ALL INSTRUCTIONS AND CAUTIONARY
MARKINGS IN THIS MANUAL AND**

ON THE APPLIANCE

When Using an electrical appliance, basic precautions should always be followed, including the following:

WARNING

**THESE WARNINGS APPLY TO THE APPLIANCE,
AND ALSO WHERE APPLICABLE, TO ALL
TOOLS, ACCESSORIES, CHARGERS OR MAINS
ADAPTORS.**

**TO REDUCE THE RISK OF FIRE, ELECTRIC
SHOCK, OR INJURY:**

1. This appliance is not intended to be used by persons (including children) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction by a responsible person concerning use of the appliance to ensure that they can use it safely.
2. Do not allow to be used as a toy. Close attention is necessary when used by or near children. Children should be supervised to ensure that they do not play with the appliance.
3. **WARNING:** Do not use this appliance near bathtubs, showers, basins or other vessels containing water.
4. If the supply cable is damaged, it must be replaced or repaired by the manufacturer, its service agent or similarly qualified persons in order to avoid a hazard.
5. Do not use the appliance for any other purpose.
6. Do not use any lubricants, cleaning agents, polishes or air fresheners on any part of the appliance.
7. Contact the us when a service or repair is required. Do not disassemble the appliance.
8. If the appliance is not working as it should, has received a sharp blow, has been dropped, damaged, left outdoors, or dropped into water, do not use and contact with us.

Technical parameters:

Model: CY-1004B

Battery capacity: 1500mAh
Charging Input: 5V 1A
Charging Interface: TYPE-C
Weight: 109g
Dimensions: 130mm *150mm *55 mm

Tips

- * It is recommended to use it every other day, and do not use it during menstruation.
- * This device should be placed in a ventilated and dry place after use.
- * Do not put other heavy objects on the device.

Care and maintenance

- * When cleaning this device, please turn off the power first, then use wet tissue paper to clean its surface. Do not clean it with alcohol or put it in water.
- * Please use neutral detergent and soft cloth for cleaning.
- * Please do not let liquid flow into the device during cleaning.

Notes:

- * If you feel uncomfortable after using the device, stop using it and consult a doctor.
- * This device contains a lithium battery, so please do not disassemble or modify the device, or immerse the device in water or place the device in fire; otherwise it may be damaged or

Contraindication:

- * Users who are allergic to light or suffering from porphyria, autoimmune diseases (such as systemic lupus erythematosus) should not be exposed to light therapy
- * Those who have taken photosensitizers within 2 weeks should not be exposed to light therapy
- * Ladies during pregnancy, lactation and menstruation;
- * Those with hypertension or uncontrolled blood pressure. People with severe cardiovascular and cerebrovascular diseases (such as multiple sclerosis, atherosclerosis) should use the product with caution;
- * Patients with high fever;
- * Those using pacemakers;
- * Those with malignant tumors in skin or other parts of body, or those with precancerous lesions
- * Patients with severe chronic diseases, endocrine diseases, major infectious diseases and acute inflammation (such as fever and skin suppurative inflammation);
- * Patients within 2 weeks after invasive cosmetic surgery and facial filling surgery.

Troubleshooting

- 1.If the device fails to function during use, please turn it off immediately and contact our customerservice.
2. Users shall not disassemble and repair the device by yourselves; other wise our company will not be responsible for any consequences by thesituation.

Instructions for device and components recycling:

For the sake of environment protection, after the device or its components are scrapped, you are responsible for separating them from

even explode

- * Do not use the device for purposes other than those described in the instructions; otherwise it may cause discomfort or the device may be damaged.
- * Do not disassemble, repair or modify this device without authorization; otherwise it may fail to function normally.
- * Do not plug and unplug the USB charging cable with wet hands; otherwise it may cause electric shock or accident.
- * Do not let children use the device
- * Do not use the device in the bathroom and other places with high humidity, or use the device while taking a bath; otherwise it may cause discomfort or it may be damaged.
- * Do not use the device while sleeping; otherwise the device may fail to function normally.
- * Please use this device in a correct way.
- * When the device works for a long time, the surface will be heated, which is a normal phenomenon.
- * Please confirm that the main switch button-works properly before use.
- * When discharging the device, please dispose of it according to the urban regulations on environmental protection; otherwise it may cause environmental pollution.
- * Please do not use or store this device in the following places: Places with direct sunlight, high temperature and humidity; places with water and a lot of dust; places vulnerable to vibration or strong magnetic field; near fireworks; places in strong electromagnetic fields
- * Do not use this device while exercising

biochemical wastes and send them to designated recycling stations for classified disassembly and recycling according to relevant national regulations.

Matters needing attention when discarding the device:

- * The battery used in this product is a special rechargeable battery, which shall not be used in other devices. Please do not charge the rechargeable batteries that have been taken out from the device.
- * Do not place the device in fire.
- * Do not pierce the device with nails, beat it with force, or dismantle or modify it.
- * Do not connect the positive and negative electrodes of the battery directly with metal.
- * Do not charge, use or place the device by fire or in the sun.
- * Please keep the removed parts, battery and other parts out of the reach of children to prevent them from swallowing by mistake.
- * Do not touch the battery liquid

Name of parts	hazardous substance					
	(Pb)	(Hg)	(Cd)	(Cr6+)	(PBBs)	(PBDEs)
Circuit board module	○	○	○	○	○	○
Plastics	○	○	○	○	○	○
Metal	○	○	○	○	○	○
Battery	○	○	○	○	○	○
Accessory parts	○	○	○	○	○	○

This form is prepared according to the standard of SJ/T 11364-2014

○: indicating that the content of the harmful substance in all homogeneous materials of the part is below the limit specified in GB/T 26572-2011.

X: indicating that the content of the harmful substance in at least one homogeneous material of the component has exceeded the limit specified in GB/T 26572-2011.

Note: The reason why some substance is marked with an "X" is that there are no alternative parts or technologies available at this stage.

This device complies with part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) this device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

Any changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

NOTE: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation.

If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

The device has been evaluated to meet general RF exposure requirement. The device can be used in portable exposure condition without restriction