


Parts

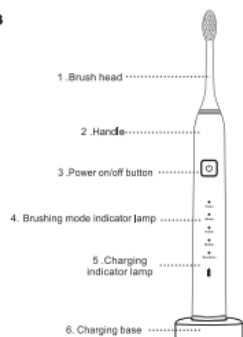

(Figure 1)

Introduction

Intelligent sonic-wave electric toothbrush (Figure 1)

1. Brush head
2. Handle
3. Power on/off button
4. Brushing mode indicator lamp
5. Charging indicator lamp
6. USB charger base

Note: items in the packaging box may vary depending on the model purchased.

Product Instructions


Input Voltage:DC 5V	Rated power: 2W
Product material: PC / ABS	Brush head: Nylon

Preparation for use

Brush head installation steps:

- 1.Align the brush heads so that the direction of the bristles is the same as the front of the toothbrush handle (Figure 2).
- 2.Press the brush head down hard to the metal axle until it will not go down anymore.

Instructions

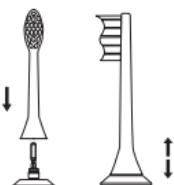

(Figure 2)

Charging the product

- Plug the charger into the power socket(Figure 3).
- Please make sure that the machine body and the charging base are perpendicular to the horizontal plane during charging, otherwise the machine body will fall down and cause safety hazards.
- Place the body on the charger base, and the charging indicator lamp will flicker until it is fully charged when the indicator is on.
- The body is pre-charged, so you don't need to fully charge it before you use it for the first time. This charger can only charge one toothbrush handle at a time normally.

Charging base



(Figure 3)

Use this product

Use steps:

- 1.Insert the brush head carefully into the oscillating axle, and leave a certain gap between the brush head and the main body to reserve a working space for the normal work of the brush head.
- 2.To prevent toothpaste splashing, take a small amount of toothpaste on the brush head before use (generally speaking, sonic-wave electric toothbrush saves more toothpaste than manual toothbrush).
- 3.First soak the brush head and put it into the mouth, then start the toothbrush.
- 4.Move the toothbrush slowly with the bristles at 45 ° to the gum line.
- 5.Hold the toothbrush handle gently and move slowly in a circular motion. Do not use press it hard.
- 6.Please rinse your mouth with water after brushing.

This product can be used safely with/for:

Dental braces (when used with the dental braces, the brush head will wear faster).
Restoration of teeth (filling, corona, porcelain teeth).

Brushing mode

Cleaning mode

Effectively remove periodontal plaque (default mode). The duration of cleaning mode is 2 minutes, including four 30s intervals.

Whitening mode

2 minutes to deep whitening teeth, effectively improving tooth surface color.

Polishing mode

It can brighten and polish the front teeth, and users can press and hold the power button to stop polishing in advance according to their own needs.

Milder mode (beginner mode)

It is suggested that this mode can be used for the beginner which can effectively clean teeth and gently massage gum to improve gum health.

Sensitive mode

It is a super gentle mode for sensitive teeth and gums which lasts 2 minutes.

Toothbrush handle

1.Remove the brush head and wash the metal axle with warm water. Make sure to remove all remaining toothpaste.

Caution: do not press sharp objects on the rubber sealing ring on the metal axle, otherwise it may cause damage

2. Wipe the entire surface of the toothbrush handle with a wet cloth.

Caution: do not use rubbing alcohol of isopropyl, vinegar or bleach to clean the handle, otherwise it may cause discoloration of the handle.

Cleaning, service and storage

Clean the brush head: wash the brush head and bristles after each use, and dry them with natural air. Do not pull the bristles, as this will deform them and cause them to fall off.

Clean the body: clean the handle surface with clean water, and scour it with soft cloth, and dry it with natural air. Do not press sharp objects on the rubber sealing ring on the metal axle, otherwise it may cause damage.

Clean the charger: clean the stains on the charger with a soft cloth. Be sure to unplug the power plug from the socket before cleaning the charger. Wipe the surface of the charger with a wet cloth and use it after it is completely dry.

Replace brush head: professional dentists recommend replacing the brush head every 2-3 months for the best results.

Clean the bottom of the brush handle: before charging, be sure to clean the bottom of the brush handle, wipe the bottom surface with a wet cloth, and then charge after the cleaning.

Storage

If you do not use this product for a long time, please fully charge the toothbrush handle, unplug the charger from the power socket, clean the two parts and store them in a cool and dry place away from direct sunlight.

Recovery

- When discarding this product, please do not put it together with the general domestic waste. Instead, hand them to the recycling center designated by the government. This will be beneficial to the environment.
- This product contains a built-in rechargeable battery; do not discard the rechargeable battery together with general domestic waste
- Follow your country's regulations on classifying and recycling electronic and electrical products and rechargeable batteries. Proper discarding of products will help to avoid adverse consequences on the environment and human health. Warranty period of intelligent sonic-wave electric toothbrush.

FCC WARNING

This device complies with part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) This device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

Note: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

Note: The Grantee is not responsible for any changes or modifications not expressly approved by the party responsible for compliance, such modifications could void the user's authority to operate the equipment.

The device has been evaluated to meet general RF exposure requirement. This equipment complies with FCC's RF radiation exposure limits set forth for an uncontrolled environment. This device and its antenna(s) must not be co-located or conjunction with any other antenna or transmitter.

To maintain compliance with FCC's RF exposure guidelines, the distance must be at least 20 cm between the radiator and your body, and fully supported by the operating and installation configurations of the transmitter and its antenna(s).



Description of charging time: it can last for 30-40 day (2 times/day, 2 minutes/time) after continuous charging for 8-10 hours. It takes about 20-22 hours to fully charge, which can be used for 50-60 days (2 times/day, 2 minutes/time).

