
---

SC-500

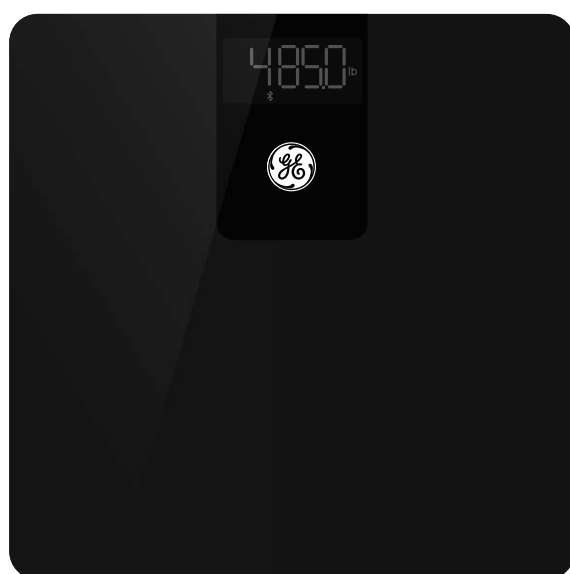**User Manual**

**EN**

**Digital Body Weight Scale**

**1- 4**

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**FR**

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Precautions


Read carefully before use:

- ⚠ Do not use the scale on a wet surface. This could lead to slipping.
- ⚠ Dispose of the scale properly following your local regulations.
- ⚠ Do not jump on the scale. This scale has a weight limit of 500lbs.

Scale Overview


GE SC-500 Digital Body Weight Scale measures weight and BMI. This digital scale is intended for home use only.

BMI (Body mass index) is a person's weight in kilograms divided by the square of height in meters. This indicator will be calculated and shown in Fit Profile App.

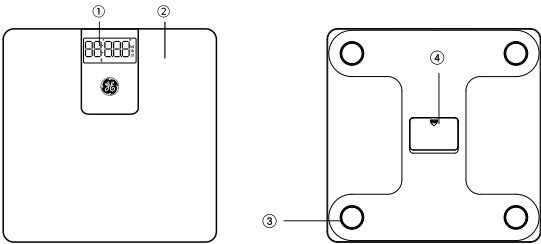Scale Layout

- ① LED Display
- ② Tempered Glass Surface
- ③ Scale Feet
- ④ Battery Box/Compartment

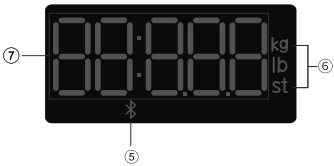LED Display

- ⑤ Bluetooth Connection
- ⑥ Unit
- ⑦ Weight Measurement

Scale Overview

|                 |                                             |
|-----------------|---------------------------------------------|
| Model:          | SC-500                                      |
| Dimensions:     | 11.81" × 11.81" × 1.06" / 300 × 300 × 27 mm |
| Weight:         | 3.7lb/1.68kg                                |
| Weighing Range: | 0.4 - 500 lb/0.2 - 226.8kg (±0.2kg)         |
| Power:          | 3*AAA battery                               |

## Get Started

### 1. Download the app and create an account

- (1) Download the app. Go to the App Store for iOS or Google Play Store for Android, and search for "Fit Profile", or scan the QR code below.
- (2) Sign up and complete your profile.

Scan QR code to download the app

Sign up

### 2. Setup your scale

- (1) Before putting your digital scale into operation, please open the battery compartment on the underside of the device and insert the three 1.5 V batteries (type AAA) included.
- (2) Place the scale on a hard flat floor surface.
- (3) Open the "Fit Profile" app and turn on the Bluetooth.
- (4) Turn on the scale by briefly stepping on and off it.



Your scale will be paired with your smartphone automatically within seconds. If not,

- (1) Tap "+ Pair the Device".
- (2) Turn on the scale by briefly stepping on and off it. On the app, tap the digital scale icon, then select "SC-500" and click "Connect."

*Your scale has now been successfully paired with your smartphone and is ready to start taking your measurements.*

### 3. Take measurements

- (1) Open the "Fit Profile" app, and then keep the app open on the "Home" screen.
- (2) Place the scale on a hard flat floor surface, do not place it on carpet or any other soft surfaces.

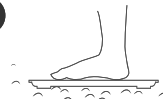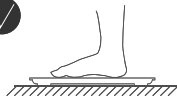- (3) Turn on the scale by briefly stepping on and off it. Then, wait for the scale to set to 0.
- (4) Step on the scale while barefoot till the numbers display on the scale. Your measurement results will also display on the app.

## Recommendations

**To have accurate, consistent, and comparable measurement results, please follow the recommendations below.**


#### 1. DO's

- Ensure the personal information in your profile is correct.
- Place the scale on a hard flat floor surface when you weigh.
- Make sure the scale displays a 0 before stepping on to measure. While barefoot, stand still on the scale until the measurement is completed.
- Weigh yourself under similar conditions and around the same time every day.
- Empty your bladder before weighing.

## 2. DON' Ts

- Do not place the scale on a carpet or a soft surface when you weigh.
- Do not measure with socks or shoes on.
- Do not measure when your body is wet. For example, immediately after a shower, bath, sauna, or strenuous exercise.
- Do not measure after drinking a lot of water or alcohol.
- Do not measure if you have a cold or fever.

## Scale Display
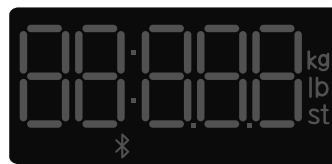
### Measurements

|                |                       |                        |                                  |
|----------------|-----------------------|------------------------|----------------------------------|
|                | No display            |                        |                                  |
| Normal working | Auto-off, or no power | No wireless connection | Low power                        |
|                |                       |                        |                                  |
| Overload       | Weight is not stable. | Unable to calibrate    | The scale needs to be calibrated |

### Messages Displayed

## Troubleshooting

| Symptom    | Cause                 | Solution                                                                                                               |
|------------|-----------------------|------------------------------------------------------------------------------------------------------------------------|
|            | Low power             | Replace new batteries.                                                                                                 |
|            | Overload              | Use the scale within its limits                                                                                        |
|            | Weight is not stable. | Please step off and measure again.                                                                                     |
|            | Unable to calibrate   | Place the scale on a flat hard floor surface gently and then turn it on. Let the scale calibrate.                      |
| No display | Auto-off, or no power | Auto-off: Normal condition. Turn on the scale (briefly step on then off)<br>No power: change and use the new batteries |

| Symptom | Cause                            | Solution                                                                                                                                                                                                                                |
|---------|----------------------------------|-----------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|
| 1360 lb | No wireless connection           | Reconnect the scale. On the app tap the device icon located in the top-right-hand corner, and then tap "+Add device". Next, turn on the scale then tap the smart scale icon on the app screen. Then, select "SC-500" and tap "Connect." |
| CAL     | The scale needs to be calibrated | Step on the scale with one foot to turn it on, and step on it again after it shows "0.00"                                                                                                                                               |

## Warranty

This product is warranted by the manufacturer against defects in materials and workmanship for one (1) year from the original date of purchase. Proof of purchase is required.

This warranty does not cover accidental damage, mistreatment, or unauthorized repair. Detailed warranty policy: [www.fit-profile.com/warranty](http://www.fit-profile.com/warranty)

If you have any issues or questions, please contact Global eFit Inc.

**Contact:** Global eFit Inc.

[cs@fit-profile.com](mailto:cs@fit-profile.com)

(866) 786-5988 (Mon – Fri. 9:00 am – 5:00 pm, US Central Time)

## FCC Compliance Statements

This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) this device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

Note: This equipment has been tested and found to comply with the limits for a Class B digital device, under Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses, and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation.

If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

Caution: Changes or modifications not expressly approved by the party responsible for compliance could void the user authority to operate the equipment.

The device has been evaluated to meet general RF exposure requirement. The device can be used in portable exposure condition without restriction

## ISED Compliance Statements

This device contains license-exempt transmitter(s)/receiver(s) that comply with Innovation, Science, and Economic Development Canada's license-exempt RSS(s). Operation is subject to the following two conditions:

- (1) This device may not cause interference.
- (2) This device must accept any interference, including interference that may cause undesired operation of the device.

This equipment complies with IC RSS-102 radiation exposure limits set forth for an uncontrolled environment.

Global eFit Inc.  
17304 Preston Road, Suite 800, Dallas, TX 75252, USA  
[www.fit-profile.com](http://www.fit-profile.com)  
[cs@fit-profile.com](mailto:cs@fit-profile.com)  
1-866-786-5988  
(toll-free US & Canada; M-F 9:00am - 5:00pm Central Time)

FCC ID: 2BAI6-SC500  
IC: 30163-SC500

GE is a registered trademark of General Electric Company.  
Manufactured under trademark license.

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Si tiene algún problema o pregunta, comuníquese con Global eFit Inc.