

Product and parts recovery instructions

Name and content list of harmful substances in products

| Component name | Hazardous substance | | | | | | |
|-------------------|---------------------|--------------------------------------|---------|----------------|--------------------------------|----------------------------------------|--|
| | Lead Pb | Marcury Hig | Cadmium | Cr (VI) | Patrominano lephende PEB | Folderstrand Springs offers PEDE | |
| hardware | ٥ | 0 | ٥ | 0 | 0 | 0 | |
| PCB board | × | 0 | ۰ | 0 | 0 | 0 | |
| Plastic parts | 0 | ۰ | 0 | ٥ | 0 | 0 | |
| | | content of the h ents appoiled in | | nce is alterno | erosus ristorial s | The post exceeds | |
| ® | | | | | enterprise scoot | ing to the product | |



to case of accidental inguistics, such medical intention immediately. Therefore y Each Sales, she not taken the Equid with years have branch and has it in the Salesing ways. 8 In the National value, there is a risk of counting of intrinse, places abort so, pure wheating it to eye, there is a risk of counting of intrinse, places abort so your eyes of other in value of the transit of the counting of other counting of other counting of other counting, when the interest pould no the other or address, may lost if you interest one of your places from with a south of counting of the co Marring The product was only be decomposed when the character, Pakin to comply you may receip in the extent once, years and cover accesses, if you need to specify receip to tany replacement, please seem. I the steep when you begin the product.

- The battery must be removed from the appliance before discarding it.
 The appliance must be powered off when the bartery is removed;
 Batteries should be disposed of safely.

- Daily maintenance

- When manifoling the place of the manifoline to make an additional dates the discrete with a well assembled and dates the discrete with a well assembled as the discrete with a well assembled as the discrete dates and discrete dates and a powerful dates and and a powerful dates and a discrete dates and a powerful dates and a discrete dates and a powerful dates and a discrete dates and a



Product warranty card






Bluetooth headsetcervicalmassager

Instruction manual

Stereo care instant comfortable experience

Please read this manual carefully before using the product and keep it properly for later use

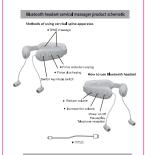

| Product name | Resolt teabelonical teasopr | Model number | KH 900(Rechargeable | | |
|----------------------------|-------------------------------------------------------|--------------|-------------------------------|--|--|
| Charging power supply | DOSV == | Rated power | 2.6W | | |
| Implementation standard | GB4706,1-2005 GB4706,10-2008 GB4943,1-2022 | Pover scotly | 3.7V tolk in little or sarroy | | |
| Sindian | Tigral by Inquency Pulse Massage - (Australia Headorn | | | | |
| Product size | 140X35X130mm | Colorbax si | as 152004200164mm | | |

Operation method

Operation method:

Methods of the first order of policy appearance.

In this was, for seal or self-order order ord

Correct values department in the counts connected continuously present every to excoust Product changing or before them between performance of the country of the country

Troubleshooting

If you encounter pull-eas when wing the polace, place even to the following sign and studions.
If they product is self-leadly or convention counts, have of the power.
Users are not abbout to dissease this or much the power.
Users are not abbout to dissease this or much the product by the medium. If the product is disseased for membranes, the preduct between your bank and the Corporary will not be taken up consequents and legislation as the product in dissease produces and legislation as the product in disease.

Condition Possible cause Handling method

| Unable to start | The battery is low. Check whether the switch button of the host product is pressed. | |
|---------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|---------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|--|
| No charging or no charging light display | Check whether the USB cable is properly connected the host and mobile phone power adapters. | |
| Display normally But without feeling | Are the electrodes in direct contact with the skin? Adjust the position for it against the skin. Is from any hin; dolling, accessedies between the electrode and the skin? To remove or remove objects. Is the power goor 1? Adjust and increase message strangth. | |
| There is a trigling sensition | Bit the electron is content and the most what it is easier to be adjusted to select the select | |
| The debather interest in susceptibilities and the second of the sec | | |

Matters needing attention

A Danger

OPLease never use with the following medical electronic devices:

- 1. Implanted medical destroits devices such as cording processation;
 2. Medical delicterion instruments to the apport such as articular hard lang;
 3. Bedrockendingstein and other medical declarate instruments.

 1. In all course mograture of the executive delicterion instruments.

 1. In all course mograture of these resoluted delicterion devices, which may endurgue fact, and processing supplication workings, give and other inflammable and enduring delicterion problems.

- - ◆ The following circumstances are prohibited.
 1. Basels with neligonst tamor 2. Potents with Identity teachers.
 2. Pregnant occurs and partners in 4. Potents with scale supposative inflammation 5. Progressive course. Or Protests with mobil englands.
 7. Protest the Journal Intel[®]
 7. Protest the J

- ne probled.

 **Timing

 **Diming

 **D

- S. Patients with side processors of the control of

- Do not use it near the heart, above the neck, on the head, in the mouth, in the public area,
 or on skin discribers. Do not use it on the soles of both feet. Also, do not compress the
 internal groups when using. Otherwise it may occur accident or love I othersial discorder.

- Do not with other medical devices or more than two (2) the products are used together, otherwise it may bow | nausea or physical discomfort.
- Die des Er purposes der train massage as describés.
 Die des Er purposes der train massage as describés in the liebución manual, or traing case document and fallam.
 Die des demanges de massage describés colores, citaveries, accidente and fallame may coox.
 Die des describés page or massage describés describés, describés and fallame may coox.
 Die de déscribés page or massage describés de product without authoritation, otherwise it may cause let laive or élabris shock.
 Pay attention

- Pay statistics
 If you feel (in do to be abromally of this product, please stop using II immediately.
 If you he follow feel of firerent, coroust your cooper and fallow your doubtr's darking.
 During the monages process, you want to more the decreade or over people, level as the cut of the poems before before increase; otherwise your my produced a strong electric strong.
 Do not before glacetions during so the instruction. Celebration says make electronic excipement images and the substitute during the artists.
- misspendin (such as deductive subt them dement discrete).

 In cose of aboutment subtem that content in places and the previous supply missistant that and the places and the size of the provide supply and table if it. Otherwise, high temperature or and discrete may occur and accordent may core. All the some time consiste contractions.

 Do not senter or remore the full fill impregnents with what have, Otherwise, electric shock, or accordant may come and contraction or accordant may come and contraction.

 Did not give the time them of proposite without the subtly to express consociatorsess, otherwise furning cause accordant or clause playsial discontrols.
- I may give an adolestic or claim is highly active plant of the price of the product of the price of the price
- exposition, oversimen you may receive a solid granulus of local.

 Pleases do not find a seare place for a long time (30 minutes) to use, otherwise may make back muscle folgate [bowlo massage the body unwest.

 Pleases wear this product convectly for cardial vertakon left and right side symmetry, otherwise address or playead descendining occur.
- scaderis or physiol-disconfoling, cour.

 When nerrowing the USB charging coller from the looks do not pall the coller instead into the play of the coller hand the coller instead in the coller instea
- If the product is not used for a long time (more than 3 months), please charge the product once every 3 months.

FCC WarningThis device complies with part 15 of the FCC rules. Operation is subject to the following two conditions: (1) this device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

Changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

NOTE: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

Reorient or relocate the receiving antenna.

Increase the separation between the equipment and receiver.

Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.

Consult the dealer or an experienced radio/TV technician for help.

Radiation Exposure Statement
The device has been evaluated to meet general RF exposure requirement. The device can be used in portable exposure condition without restriction.